# Midmark® Steam Sterilizers

For Models:

M9® (-050/-053)

M11® (-050/-051/-053)

[English](#)

[Español](#)

[Français](#)

## Setup and User Guide

Style U

003-10534-99 Rev. AA9 (3/25/24)  
TP202 20-42-FO-00014 Rev A3 ECO088522

# Setup

## **Attention Canadian Users - Action Required!**


Thank you for purchasing a M9® or M11® Sterilizer.

To comply with changes in Canadian regulatory requirements you must now provide documentation, a Manufacturer's Data Report (MDR), to your appropriate provincial regulatory agency when purchasing and using a steam sterilizer.

To obtain the MDR from Midmark®, record the sterilizer Serial Number and the American Society of Mechanical Engineers (ASME) National Board Number of your sterilizer. See the illustration below for the location of these numbers. Then, contact Midmark customer service at 1.800.MIDMARK and provide the recorded numbers.

After completing the above steps, Midmark will complete the MDR and send it to you. Once you have the MDR, you must submit it to the appropriate provincial regulatory agency.

We apologize for any inconvenience and appreciate your assistance in helping maintain compliance to Canadian regulations for sterilizers. If you have any questions or concerns, please contact 1.800.MIDMARK.

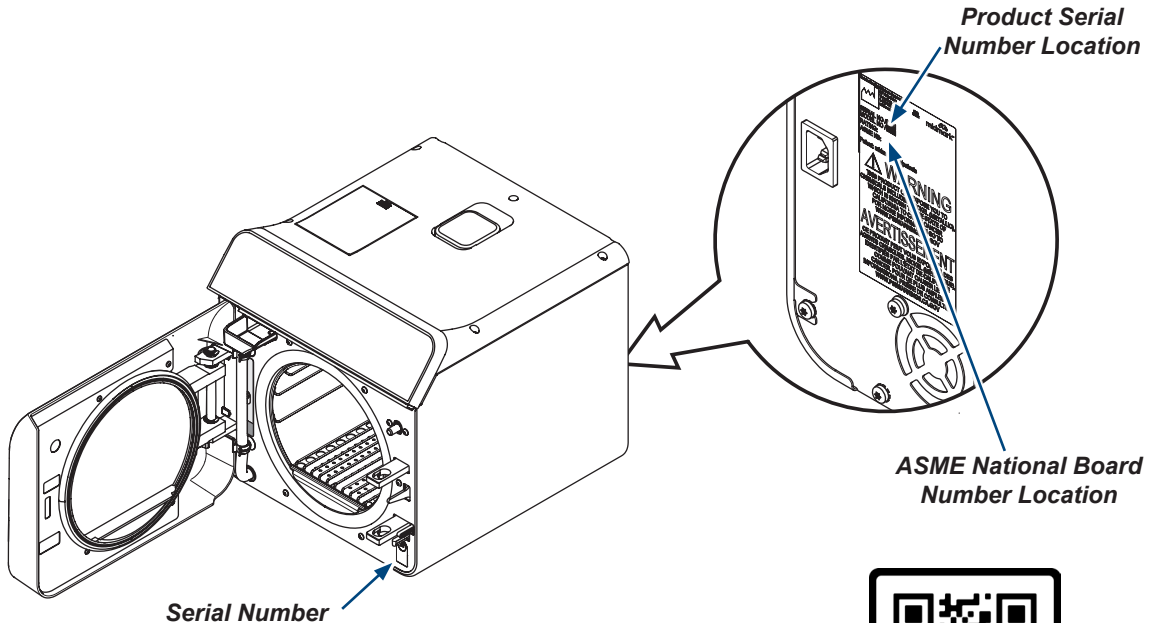Scan the QR code for access to the Canadian Install Checklist per standard: (CSA Z314:23 section 16.4.2.)

# Table of Contents

## Setup

Canadian Regulatory Requirements.....	2
Safety Notification .....	4
Shipping Symbols .....	4
Transportation / Storage Conditions .....	5
Operating Environment.....	5
Electrical Requirements .....	6
Location Requirements.....	7
Unpacking M11® Steam Sterilizer.....	8
Unpacking Steam M9® Sterilizer .....	9
Leveling Procedure .....	10
Connecting the Power Cord.....	10
User Interface Setup.....	11

## Important Information

Intended Use .....	12
Sterilization Technology .....	12
FCC Compliance Statement.....	12
Electromagnetic Interference .....	13
Wi-Fi and Ethernet Connections .....	13
myMidmark.com Cloud Service Information ..	14
Unregistering a Sterilizer.....	14
Cybersecurity .....	14
Software Update .....	14
Authorized Representatives.....	15
Accessories, Tools and Common Parts List..	16
Component Location .....	17
Sterilization Monitoring Guidelines .....	18

## Operation

Before Operating the Sterilizer.....	19
Controls and Interface .....	20
User Interface Icon Key .....	20
Additional Menus .....	21
Biological Indicators (BI) .....	21
Cycle History .....	22
Settings .....	23
Filling the Reservoir.....	26
Qualification Testing .....	26
Guidelines for Loading.....	27
Standard Cycle Parameters.....	35

## Operation Continued...

Custom Cycle Parameters .....	36
Standard Cycle Operation .....	37
Post-Sterilization Handling .....	40
Custom Cycle Operation .....	40
Bowie-Dick Cycle .....	42

## Routine Care

Daily Care .....	44
Weekly Care.....	45
Monthly Care .....	47
Extended Use Care .....	50

## Accessory and Tools Details

Using the Optional Tray and Cassette Tool...	51
VistaCool™ Direct-to-Drain System .....	52
USB Kit.....	52
Auto-Fill .....	52
Vertical Cassette Rack (M11® Only).....	52

## Troubleshooting

Troubleshooting Charts.....	53
Troubleshooting Menus.....	54
Error Codes.....	55
Calling for Service .....	57

## Product Compliance and Specifications

Specifications Chart: M9® Steam Sterilizer ..	58
Specifications Chart: M11® Steam Sterilizer	59
Certification + Standards Compliance.....	60
Electromagnetic Disturbances Guidance .....	60
Connectivity Options.....	61
Water Purity Specifications.....	61
Disposal of Equipment.....	61

## Warranty Information

Scope of Warranty.....	62
------------------------	----

# Safety Notification

The primary concern of Midmark® is that this equipment is operated and maintained with the safety of the patient and staff in mind. To ensure safe and reliable operation:

- Before attempting to setup or operate the sterilizer, read and understand this user guide and document this activity following your standard operating procedures.
- Ensure that the appropriate personnel are informed of the contents in this guide (this is the responsibility of the purchaser).
- Ensure that this guide is located near the sterilizer and that appropriate personnel know where it is located.

The M9® -053 and M11® -053 models are for veterinary use only and not FDA cleared for human use.

# Safety Symbols

## WARNING

*Indicates a potentially hazardous situation which could result in serious injury.*

## Caution

*Indicates a potentially hazardous situation which could result in minor or moderate injury. It may also be used to alert against unsafe practices.*

## Equipment Alert

*Indicates a potentially hazardous situation which could result in equipment damage.*

## Hot Surface

*Indicates the surface can be hot and should not be touched with out taking care.*

# Shipping Symbols

Proper Shipping Orientation

Keep Dry

Fragile

Max. Stacking Height  
*(Refer to number on package)*

Humidity Limit

2 Person Lift Caution

> 34 kg (75 lb)

Temperature Limit

## Transportation / Storage Conditions

### **Equipment Alert**

Water must be drained from the unit's reservoir before transporting or storing below 0°C (+32°F).

<b>Ambient Temperature Range:</b>	-30°C to +60°C (-22°F to 140°F)
<b>Relative Humidity:</b>	10% to 90% (non-condensing)
<b>Atmospheric Pressure:</b>	49.6 kPa to 106.4 kPa (7.2 psia to 15.4 psia)

## Operating Environment

<b>Ambient Temperature Range:</b>	+20°C to +40°C (68°F to 104°F)
<b>Relative Humidity:</b>	< 80% (non-condensing)
<b>Normal Operating Altitude:</b>	< 3000 m (9842 ft.) above sea level
<b>Device approved for INDOOR USE ONLY.</b>	
<b>Device to be operated in a relatively dust-free environment.</b>	
<b>The M9® and M11® will emit 5275 kJ (5000 BTU) per hour during operation.</b>	

# Electrical Requirements

## **WARNING**

*Connect the sterilizer to a power source that provides rated voltage as specified. Refer to the product label on the back of the sterilizer for the applicable electrical rating. Failure to do so may result in electric shock to personnel and / or damage to sterilizer.*

## **Note**

*For safety, the sterilizer must be connected to a properly polarized and grounded receptacle. Always use a power cord with grounding connections that match the receptacles in your location.*

### **115 VAC unit:**

115 VAC, 50/60 Hz, 12 amp  
Dedicated Supply Circuit: 115 VAC, 50/60 Hz, 15 amp  
Max. Power Consumption: 1425 Watts

### **230 VAC unit:**


230 VAC, 50/60 Hz, 6.4 amp  
Dedicated Supply Circuit \*: 230 VAC, 50/60 Hz, 10 amp  
Max. Power Consumption: 1500 Watts

\* The sterilizer is to be connected to a single-phase (split-phase) three-wire supply with 115 V nominal line to ground and 230 V nominal line to line

## **Note**

*115 VAC models may be operated in a voltage range of 103-127 VAC. 230 VAC models may be operated in a voltage range of 207-253 VAC. Overall cycle time will vary based on conditions such as voltage, starting temperature, and altitude. At the lower end of the voltage range, heat-up times will increase, and it may be necessary to run a Pre-Heat cycle prior to the sterilization cycle.*

# Location Requirements

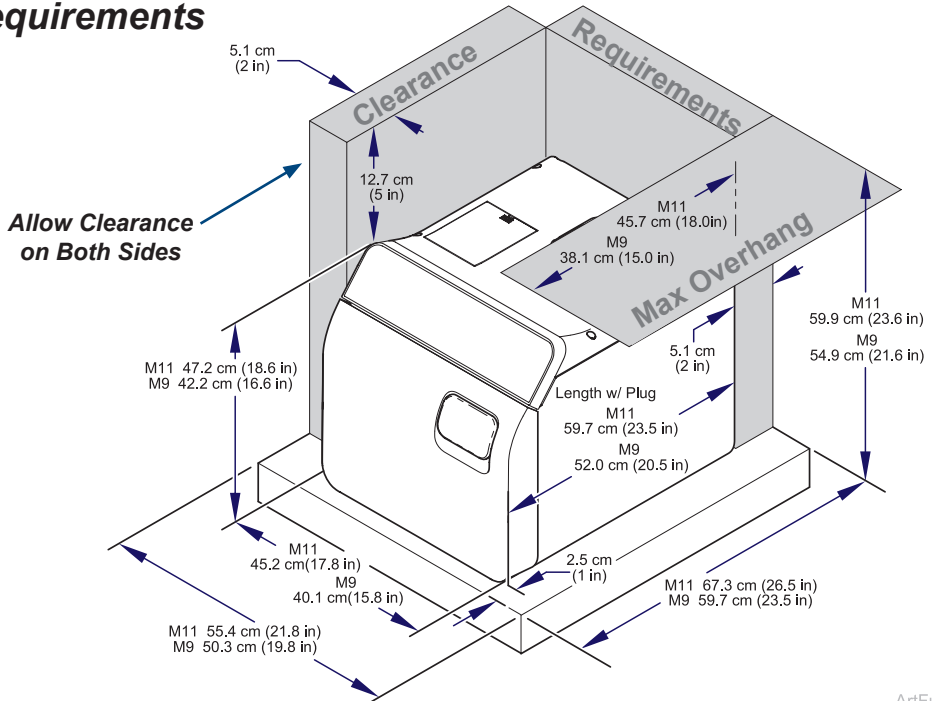ArtFull

## Support Surface

- Surface should be water-resistant material (e.g. laminate, stainless steel, solid surface materials, stone)
- Surface must be level to ensure chamber fills with correct water level. Improper water level in the chamber could cause a sterilizer malfunction.
- Surface including clearances should meet these minimum dimensions:

Depth (front to back w/plug)	M11® - 67.3 cm (26.5 in)	M9® - 59.7 cm (23.5 in)
Width (side to side)	M11 - 55.4 cm (21.8 in)	M9 - 50.3 cm (19.8 in)

## Clearance Requirements


To ensure proper air circulation and to allow access to the power cord, adhere to the minimum clearance requirements listed below. If the sterilizer will frequently be operated without pauses between cycles, position unit where steam will not damage materials or equipment in the surrounding area. If the sterilizer will be used in an enclosed area and will frequently be operated without pauses between cycles, a direct to drain solution is recommended.

- Back of Unit - Back Wall: 5.1 cm (2 in)
- Front Support Surface - Front Sterilizer: 2.5 cm (1 in)
- Sides of Unit - Side Wall: 5.1 cm (2 in)
- Distance above Unit: 12.7 cm (5 in)
- Maximum Upper Overhang: M11 - 45.7 cm (18 in) M9 - 38.1 cm (15 in) (from back of sterilizer)
- Under Cabinet or Shelf: M11 - 59.9 cm (23.6 in) M9 - 54.9 cm (21.6 in)


## Relocation Requirements for Sterilizer

- Disconnect power cord from electrical outlet and allow sterilizer to cool.
- Drain water from reservoir to prevent water from spilling. Refer to **Weekly Care** for process
- Perform Qualification Testing. Refer to **Qualification Testing** for details.
- It is recommended to use two persons when moving or lifting sterilizer.

# Unpacking M11<sup>®</sup> Steam Sterilizer
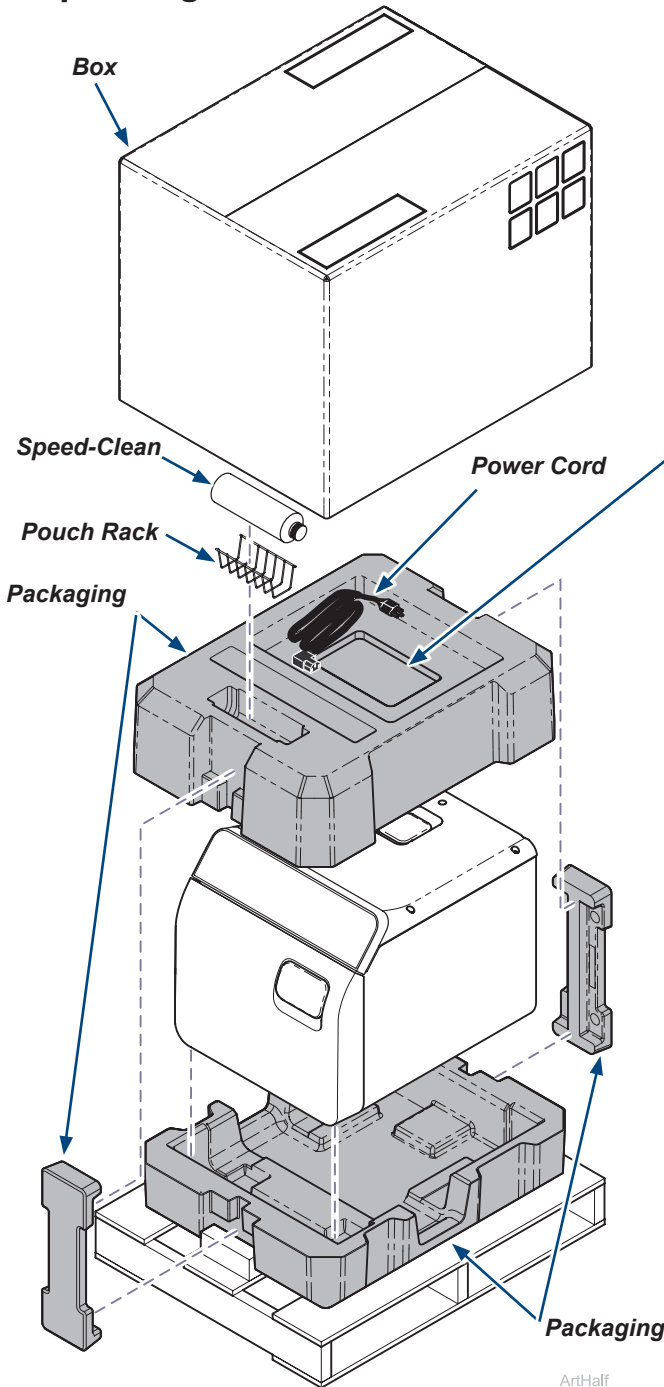
## 2 Small Trays + 2 Large Trays
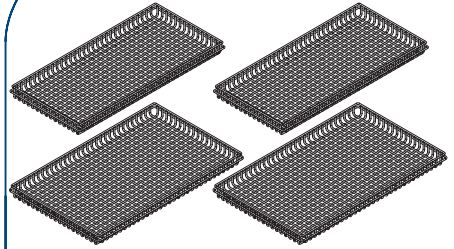
## Speed-Clean SDS Sheet

Midmark Corporation  
07240 L40 Rev. P  
Doc Created: June 2007  
Revised: February 2015  
Page 1 of 4

SAFETY DATA SHEET

**Section 1: Identification**

1.1. Product Identifier  
Speed-Clean

1.2. Recommended Use and Restrictions on Use  
Sterilizer Cleaner

1.3. Supplier Details  
Midmark Corporation  
60 Wood Drive, P.O. Box 286, Versailles, OH 45383, USA  
Address: Telephone: FAX: 614/291-1000  
L40W4MIDMARK

1.4. Emergency Telephone Number  
(973) 258-6400

**Section 2: Hazard Identification**

2.1. Hazard Classification  
Skin Irritation Category 2  
Eye Irritation Category 2A

2.2. Label Elements  
Signal Word  
Warning

**Symbol**  
Corrosion Mark

**Hazard Statements**  
Prolonged or repeated skin contact causes skin irritation  
Contact with eyes causes eye irritation

**Precautionary Statements**  
**Prevention:**  
Wear suitable protective clothing  
Wear eye protection  
Wash hands after handling this material

**Response:**  
IF ON EYES: Immediately flush eyes with water for 15 minutes, intermittently lifting the lower and upper eyelids.  
IF ON SKIN: Flush skin with running water for 15 minutes.

2.3. Hazards Not Otherwise Classified  
None.

## Care / Operation Card

midmark M9<sup>®</sup> and M11<sup>®</sup> Steam Sterilizer

For Models: M9 (400-953) M11 (400-015-053)

Serial Number: 001 Data Matrix

Follow the QR code to register your sterilizer.

Follow the QR code to the User Guide.

Cycle Type	Temperature	Exposure Time	Pressure	Wiping Time	Wiping Solution	Wiping Frequency	Wiping Area	Wiping Method
Wiped	121°C (250°F)	15 min	103 kPa (3.0 psi)	15 sec	70% Isopropyl Alcohol	10 times	100% of all surfaces	Use a clean, lint-free cloth
Wiped	121°C (250°F)	5 min	103 kPa (3.0 psi)	15 sec	70% Isopropyl Alcohol	10 times	100% of all surfaces	Use a clean, lint-free cloth
Wiped	121°C (250°F)	15 min	103 kPa (3.0 psi)	15 sec	70% Isopropyl Alcohol	10 times	100% of all surfaces	Use a clean, lint-free cloth
Wiped	121°C (250°F)	5 min	103 kPa (3.0 psi)	15 sec	70% Isopropyl Alcohol	10 times	100% of all surfaces	Use a clean, lint-free cloth

**Routine Care**

**Daily Routine Care:**

- Wipe down the exterior surfaces with 70% Isopropyl Alcohol.
- Wipe down the interior surfaces with 70% Isopropyl Alcohol.
- Wipe down the door gasket with 70% Isopropyl Alcohol.

**Weekly Routine Care:**


- Wipe down the exterior surfaces with 70% Isopropyl Alcohol.
- Wipe down the interior surfaces with 70% Isopropyl Alcohol.
- Wipe down the door gasket with 70% Isopropyl Alcohol.

**Monthly Routine Care:**


- Wipe down the exterior surfaces with 70% Isopropyl Alcohol.
- Wipe down the interior surfaces with 70% Isopropyl Alcohol.
- Wipe down the door gasket with 70% Isopropyl Alcohol.

For complete Routine Care refer to the User Interface or User Guide.

# Unpacking Steam M9® Sterilizer
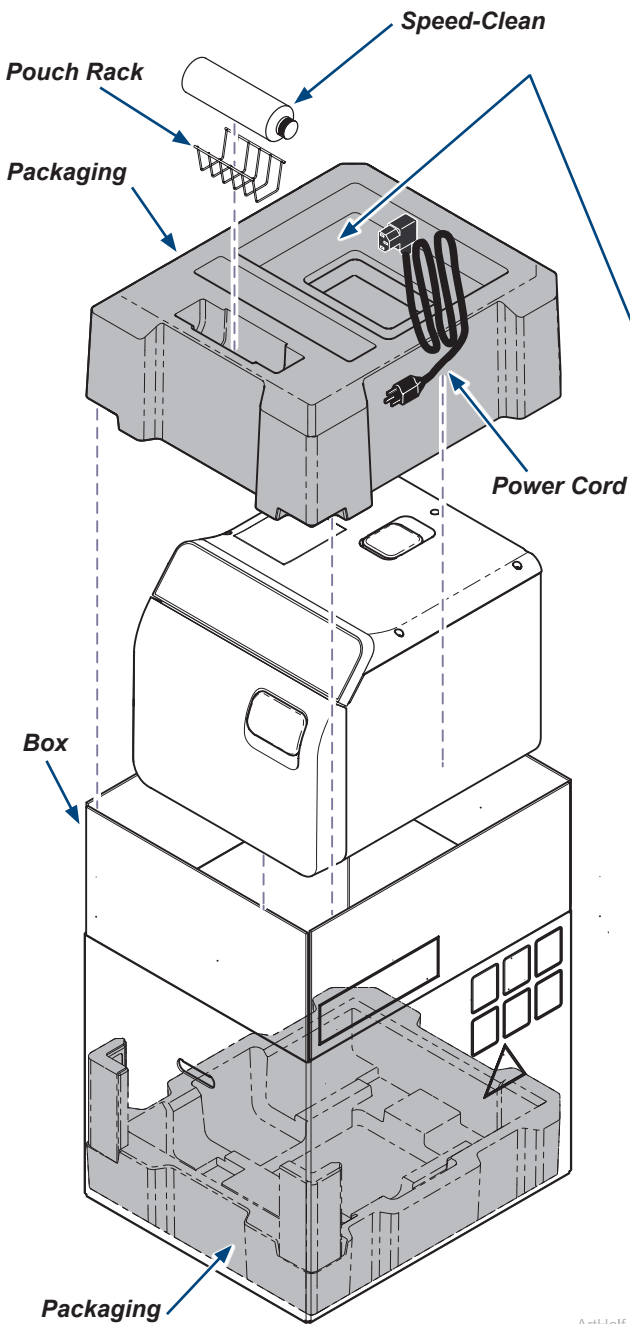
## 2 Small Trays + 2 Large Trays
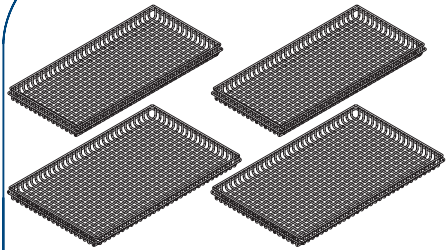
## Speed-Clean SDS Sheet

Midmark Corporation  
07243 L40 Rev. P  
Data Created June 2007  
Revised February 2015  
Page 1 of 4

**SAFETY DATA SHEET**

**Section 1: Identification**

1.1. Product Identifier  
Speed-Clean

1.2. Recommended Use and Restrictions on Use  
Sterilizer Cleaner

1.3. Supplier Details  
Midmark Corporation  
40 West View, P.O. Box 286, Versailles, OH 45383, USA  
Address: LA50MIDMARK  
Telephone:

1.4. Emergency Telephone Number  
(973) 256-6600

**Section 2: Hazard Identification**

2.1. Hazard Classification  
Skin Irritation Category 2  
Eye Irritation Category 2A

2.2. Label Elements  
Signal Word  
Warning

Symbols  
Corrosion Mark

GHS07

**Hazard Statements**  
P201 or repeated skin contact causes skin irritation  
Contact with eyes causes eye irritation

**Precautionary Statements**  
P201  
P273  
P301+P312  
P302+P352  
P303+P361+P531  
P305+P351+P338  
P501

**Exposure Routes**  
Inhalation, skin contact, eye contact

**Acute Toxicity**  
None

**Environment**  
None

**Other Information**  
None

## Care / Operation Card

midmark M9® and M11® Steam Sterilizer

For Models:  
M9 (450-055)  
M11 (450-055-055)

Follow the QR code to register your sterilizer.

Follow the QR code to the User Guide.

Cycle Type	Temperature	Time	Pressure	Notes	min. Cycle Time	min. Cycle Time
Wetpack	121°C (250°F)	15 min	15 PSI	Preheated or pre-soaked items recommended. Do not load items that are not preheated or pre-soaked. Do not load items that are not preheated or pre-soaked.	15 min (M9) / 15 min (M11)	15 min (M9) / 15 min (M11)
Wetpack	134°C (273°F)	3 min	15 PSI	Preheated or pre-soaked items recommended. Do not load items that are not preheated or pre-soaked. Do not load items that are not preheated or pre-soaked.	3 min (M9) / 3 min (M11)	3 min (M9) / 3 min (M11)
Wetpack	134°C (273°F)	3 min	15 PSI	Preheated or pre-soaked items recommended. Do not load items that are not preheated or pre-soaked. Do not load items that are not preheated or pre-soaked.	3 min (M9) / 3 min (M11)	3 min (M9) / 3 min (M11)
Wetpack	134°C (273°F)	3 min	15 PSI	Preheated or pre-soaked items recommended. Do not load items that are not preheated or pre-soaked. Do not load items that are not preheated or pre-soaked.	3 min (M9) / 3 min (M11)	3 min (M9) / 3 min (M11)

**Routine Care**

**Daily Routine Care:**

- Clean sterilizer exterior with mild soap and water.
- Wipe down and clean interior with a damp cloth or Speed-Clean.
- Clean trays after each use with Speed-Clean.

**Weekly Routine Care:**

- Clean Speed-Clean in clean container with hot water.
- Wipe down sterilizer exterior with Speed-Clean.
- Clean trays after each use with Speed-Clean.

**Monthly Routine Care:**

- Clean Speed-Clean in clean container with hot water.
- Wipe down sterilizer exterior with Speed-Clean.
- Clean trays after each use with Speed-Clean.

**Equipment Alert**

• If you notice, allow and call before cleaning packets and internal surfaces or removing trays, they may, in part and time.

• If you notice, allow and call before cleaning packets and internal surfaces or removing trays, they may, in part and time.

• If you notice, allow and call before cleaning packets and internal surfaces or removing trays, they may, in part and time.



# Leveling Procedure

## Step 1: Leveling Your Sterilizer

- A) Place a level on top of support surface to ensure it is level front to back and side to side. If the support surface is level, verify the leveling feet are tightened securely against the shim. There is no need to adjust the sterilizer in this case.
- B) If support surface is not level: Use pliers to remove shims from between leveling feet and base of sterilizer to raise or lower leveling feet located under front corners of sterilizer.
- C) Place a level on top of sterilizer, thread leveling feet in or out of bottom of sterilizer until the sterilizer is level front to back and side to side.

### Note

Use a 11/16 in wrench to adjust leveling feet. Turning leveling foot CW lowers foot, CCW raises foot.

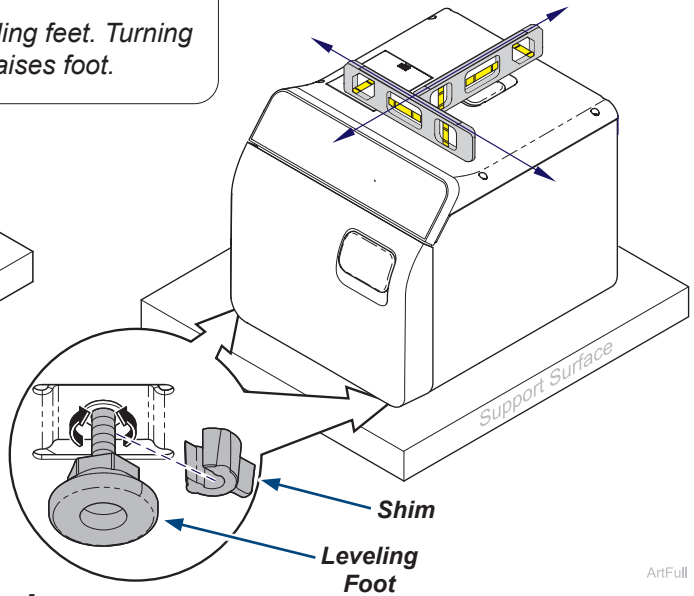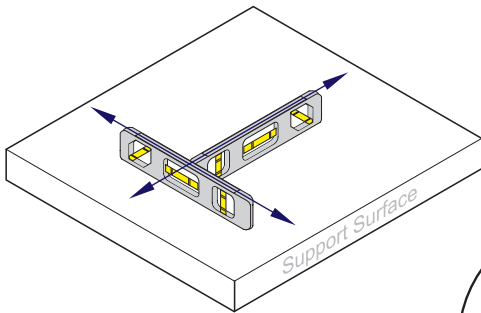ArtFull

## Connecting the Power Cord


### Equipment Alert

For optimal performance, allow sterilizer to reach room temperature before operating.

### Step 1: Connect the power cord


- A) Plug power cord into receptacle on back of sterilizer.
- B) Plug power cord into a properly polarized and grounded receptacle. A dedicated circuit only used for the sterilizer is recommended.
- C) M9® and M11® are not equipped with an on/off switch. The display operates off low power similar to a microwave oven display.
- D) Allow approximately 30-45 seconds for initial start up.
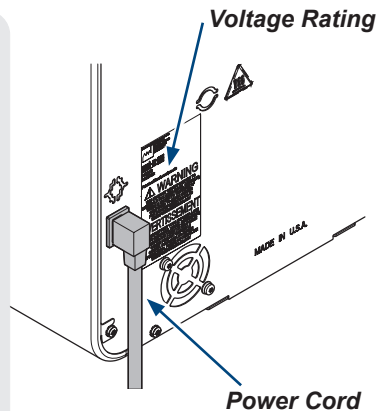
## User Interface Setup

The Midmark M9® and M11® Sterilizers are equipped with a touch screen display. Use the touch screen to perform the user interface setup by following the on-screen prompts and steps.

- Scan the QR code below to register the sterilizer.
- Enter the serial number off of the label located inside the door on the bottom right corner.
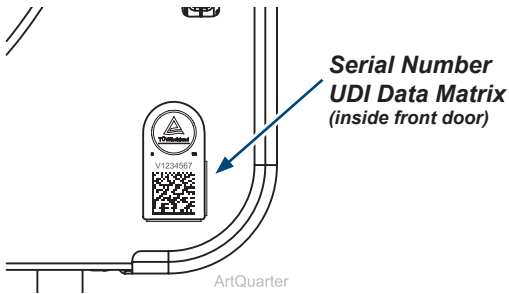
### Note

*These settings can be edited at any time by accessing the **Settings** menu by pressing the Menu icon on the left of the home screen. Refer to the **Settings** section in this guide for details.*

# ***Important Information***

## ***Intended Use***

The Midmark M9® and M11® Steam Sterilizers can be used in medical and dental offices, hospitals, clinics, nursing homes, laboratories and other facilities to sterilize heat and moisture stable reusable items (including dental handpieces) that are compatible with steam sterilization. Refer to ***Standard Cycle Parameters*** in this guide for detailed information.

## ***Sterilization Technology***

The Midmark M9 and M11 Steam Sterilizers utilize a dynamic-air-removal method to effectively remove air in the chamber and load while generating steam. There are two commonly used dynamic-air-removal systems; one type is pre-vacuum which utilizes a vacuum pump, and the other type is steam-flush pressure-pulse (SFPP), which utilizes a series of steam flushes and pressure pulses without the need for additional vacuum components. Midmark sterilizers utilize the SFPP type of air removal.

## ***FCC Compliance Statement***

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the user guide, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Any modifications to this product may violate the rules of the Federal Communications Commission and make operation of the product unlawful.

FCC, IC ID's and wireless transmitter information can be found in the **Settings Menu**.

# Electromagnetic Interference

This sterilizer is designed and built to minimize electromagnetic interference with other devices. However, if interference is noticed between another device and this product:

- Remove interfering device from room
- Plug sterilizer into isolated circuit
- Increase separation between sterilizer and interfering device
- Contact Midmark if interference persists

## **Equipment Alert**

*Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 in) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.*

## **WARNING**

- ***Use of the sterilizer adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, the sterilizer and the other equipment should be observed to verify that they are operating normally.***
- ***The sterilizer is not intended to be used in an environment that contains other equipment that are unique medical emitters such as electrocautery, MRI, electrosurgical units, and diathermy devices.***
- ***Use of accessories, transducers, and cables other than those specified or provided by Midmark could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.***

## **Wi-Fi and Ethernet Connections**

- A remote connection with the myMidmark app is available for 30 minutes after powering on the sterilizer.
- Disconnect and reconnect power to the sterilizer if a connection to the myMidmark app is needed to modify network settings or to re-register the device.
- This sterilizer includes an Ethernet connection and a Wi-Fi module that allows the device to connect to the local network.
- In rare instances, the sterilizer may be unable to connect to a wireless network at power-up due to operating location/environment frequency interference on certain wireless channels. Moving the sterilizer slightly may resolve this issue. Note: Interference conditions may change. Connecting the sterilizer via Ethernet is another option.
- Please consult your IT department/contact to verify the proper Wi-Fi coverage in the area where the device will be used.
- Only products certified to one of the following may be connected to the Ethernet ports.
  - IEC 60950-1 and complying with SELV or TNV-1 limits
  - IEC 62368-1 and complying with ES 0 or ES 1 limits
- Connecting the sterilizer to an IT-network could result in previously unidentified risks to patients, operators or third parties and the system integrator should identify, analyze, evaluate and control these risks.

## ***myMidmark.com Cloud Service Information***

***Loss of Connection*** - If the sterilizer is off-line for:

- Less than a period of 90 days, no action is required to backup or restore authenticated configurations. Midmark's cloud-based services are backed up and will, in the event there is a loss of services, be restored when the device reconnects to the Midmark cloud services.
- A period of 90 to 180 days, the device may need to be re-registered using the myMidmark app to reconnect to Midmark cloud services.
- More than 180 days, the device must be re-registered using the myMidmark app to reconnect to Midmark cloud services.

On reconnection, authenticated configurations stored in the Midmark cloud services will be synchronized to the sterilizer.

## ***Unregistering a Sterilizer***

The sterilizer may be unregistered in the myMidmark.com cloud-based application.

- An unregistered device will be inactivated in the Midmark cloud services. Past activity and data will not be available to any new customer associated to that device. Future activity or data is not associated with a prior customer account.
- Device activity and data is retained in the cloud for a prior customer so that historical data may be retrieved by that customer.
- No data is retained on the device network hardware.

## ***Cybersecurity***

Connecting the sterilizer to a local network could result in previously unidentified cybersecurity risks that the organization responsible for the network should identify, analyze, evaluate, and control. The product captures relevant system and security logs to identify date, time, and source of user activity. A software bill of material (SBOM) is available upon request. If a cybersecurity vulnerability or incident is detected while using the sterilizer, or for any other cybersecurity needs, contact Midmark Technical Service for assistance (see **Calling for Service**).

## ***Software Update***

Within myMidmark.com, a sterilizer software update may be available, and the update may be scheduled to be applied to the sterilizer either immediately or for a specific time.

- If the sterilizer is being used at the scheduled time or up to 30 minutes past the scheduled time, the software update will be canceled and will need to be rescheduled.
- Follow prompts on the sterilizer display as the device may need to be disconnected and reconnected to power.
- If applicable, the device will indicate if qualification cycles should be run after the update (see **Qualification Testing**).

## Authorized Representatives

### **WARNING**

***In the event of any serious incident in relation to the device, please contact Midmark and the appropriate competent authority.***

Customers in Colombia should direct all questions, incidents, and complaints to the Midmark Authorized Representative:

Técnica Electromédica S.A.

Calle 47 # 5-26

Bogota Colombia c.p. 110111

Tel. (+ 57 1) 756 87 87

atencionalcliente@tecnicaelectromedica.com

Customers in Costa Rica should direct all questions, incidents, and complaints to the Midmark Authorized Representative:

Corporacion Almotec SA

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San Jose, Zapote, Costa Rica

Tel: 506 7112 9782

Customers in Panama should direct all questions, incidents, and complaints to the Midmark Authorized Representative:

TECNIDENTAL PANAMA S.A.

VILLA MAGNA, VILLA ZAITA

OFFICE B8

PANAMA

Tel/Fax: 343-4049/ 6614-8029

tecnidentalpanama@gmail.com

Customers in Mexico should direct all questions, incidents, and complaints to the Midmark Authorized Representative:

Sofía Zavala | Asuntos Regulatorios

Carlos Nafarrate SA de CV

Avenida Las Américas 601,

Ladrón de Guevara, 44600,


Guadalajara, Jalisco, México | Tel (33) 3669 9800

## Accessories, Tools and Common Parts List

To maintain reliability and performance use only genuine Midmark replacement parts. Unless noted, items can be used on the M9® and M11® Steam Sterilizers.

<b>Common Service Parts / Tools</b>	<b>Part Number</b>
Auto Fill Kit 115 V	9A694001
Auto Fill Kit 230 V	9A694002
Cassette Rack (Vertical for M11 only)	9A703001
Door and Dam Gasket Kit for M11	002-10880-00
Door and Dam Gasket Kit for M9	002-10879-00
Pouch Rack Kit for M11	002-2108-00
Pouch Rack Kit for M9	002-2108-01
Speed-Clean, 1 (0.47 Liter [ 16 oz.]) Bottle	002-0396-00
Speed-Clean, 1 Case (qty 12 - 0.47 Liter [ 16 oz.]) Bottles	002-0396-05
Tray/Cassette Tool	9A307001
USB Kit	9A698001
VistaCool™ Adapter Kit	9A586007
VistaCool™ Direct-To-Drain Thermal Reduction System	9A586002 (Dual)

# Component Location
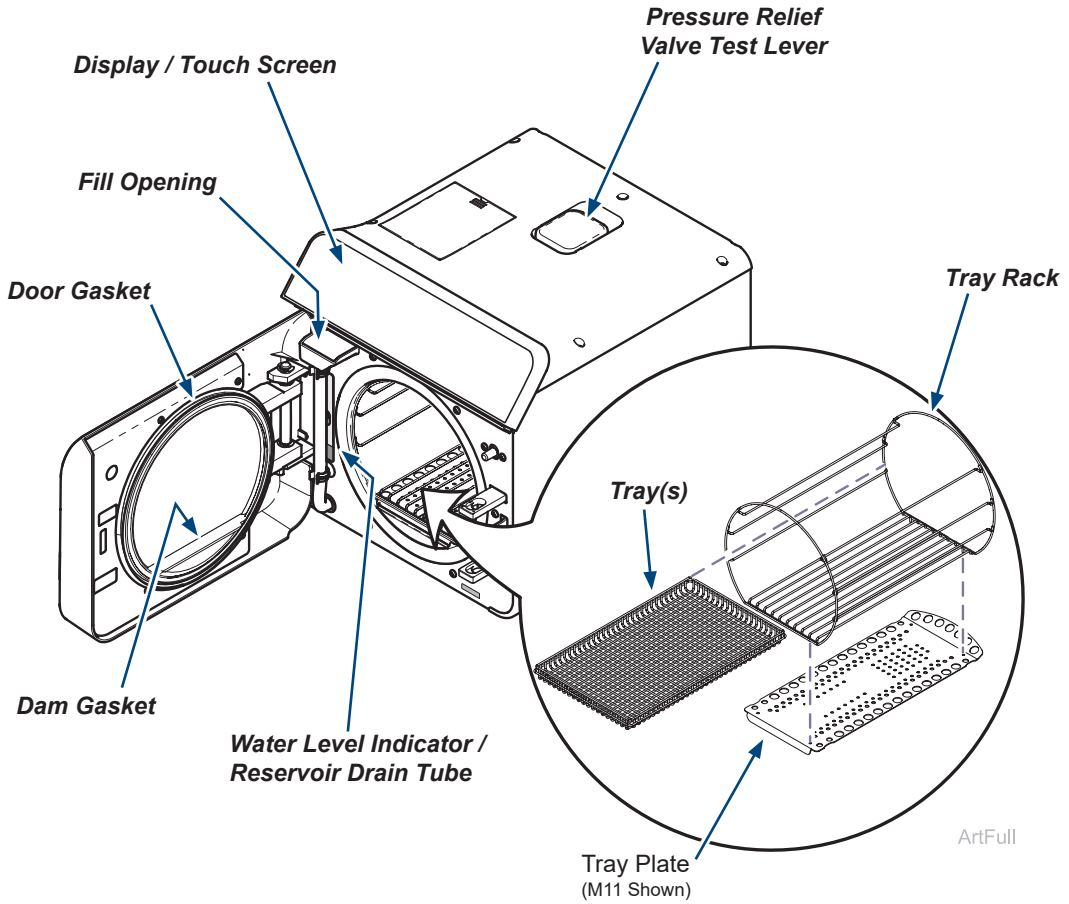
# Sterilization Monitoring Guidelines

## Note

*This information is provided for reference only. Contact appropriate state/local agencies for specific sterilization guidelines for your office. Additional information on infection control is available from the Centers for Disease Control and Prevention (CDC), Organization for Safety and Asepsis Procedures (OSAP), Association for the Advancement of Medical Instrumentation (AAMI) and Association for Professionals in Infection Control and Epidemiology (APIC).*

## Physical Monitors

Temperature and pressure sensing controls are integrated in the sterilizer that monitor the performance of the unit, help detect sterilizer malfunctions and allow suggested data to be verified. The sterilizer uses "Cycle History" to meet these recommendations and will abort the cycle adding a cycle error to "Cycle History" if the sterilizing conditions are outside of the established limits.

## Note

*Use only FDA cleared chemical and biological indicators designed for steam sterilization that are compatible with the particular dynamic-air-removal sterilization cycle temperature and exposure time being monitored (see **Intended Use** for more information). Use chemical indicators with each sterilization load. If a sterilizing cycle is terminated prematurely, reprocess instruments to ensure sterility of the load. Follow manufacturer's instructions for proper disposal of used indicators.*

**Chemical Indicators** Chemical indicators are designed to verify that conditions in the sterilizer chamber were adequate to achieve sterilization. They do not validate that a processed item is sterile. If a chemical indicator shows a failure, items in that load are considered non-sterile. Potential causes for sterilization failure include: improper cleaning, packing, loading, or a sterilizer malfunction. Determine the cause of any sterilization failure and remedy the situation before running the next cycle. Only FDA cleared chemical indicators labeled for use with the steam sterilization cycle parameters (e.g. temperature and exposure time) of the M9® and M11® Steam Sterilizers should be used for monitoring the cycles. Follow the chemical indicator's instructions for proper storage, use, interpretation and disposal.

## Biological Indicators

Biological indicators are microbiological products designed to accompany items being sterilized to monitor adequacy of the sterilization process. If a biological indicator fails, items in that load are considered non-sterile. Potential causes for sterilization failure include: improper cleaning, packing, loading, or a sterilizer malfunction. Determine the cause of any sterilization failure and remedy the situation before running the next cycle. Only FDA cleared biological indicators labeled for use with the steam sterilization cycle parameters, e.g. temperature and exposure time, of the M9® and M11® Steam Sterilizers should be used for monitoring the cycles. Follow the biological indicators instructions for proper storage, use, interpretation and disposal.

A biological indicator and a chemical should be placed in a test pack, which is a representation of a typical load/package processed in the sterilizer. This test pack is to be placed on the bottom tray near the chamber door. In addition to the test pack, add other items to create a fully loaded chamber.

# Operation

## Before Operating the Sterilizer

### **WARNING**

- *Do not use this sterilizer for sterilizing volatile substances or for any purpose other than its intended design. Burns and toxic or explosive conditions could result.*
- *Do not force door handle at any time. Chamber pressure may cause door to open with extreme force. If door handle does not move freely, allow unit to cool and depressurize for 40 minutes before opening door. Failure to adhere could result in serious personal injury.*
- *Do not run the sterilizer without the tray plate in place. If the sterilizer malfunctions, unplug sterilizer and call for service; do not attempt to repair the sterilizer yourself. Doing so could result in serious injury.*


### **WARNING**

*Programmable cycles are provided for those applications requiring sterilization parameters different than the preset cycles. These cycles are not FDA cleared for medical use. All material processed in these cycles must be validated by the user to ensure sterility of the processed load.*

## Controls and Interface

The sterilizer controls and indicators are shown in the illustrations on this and the following pages. The accompanying tables describe the function of each control / indicator.
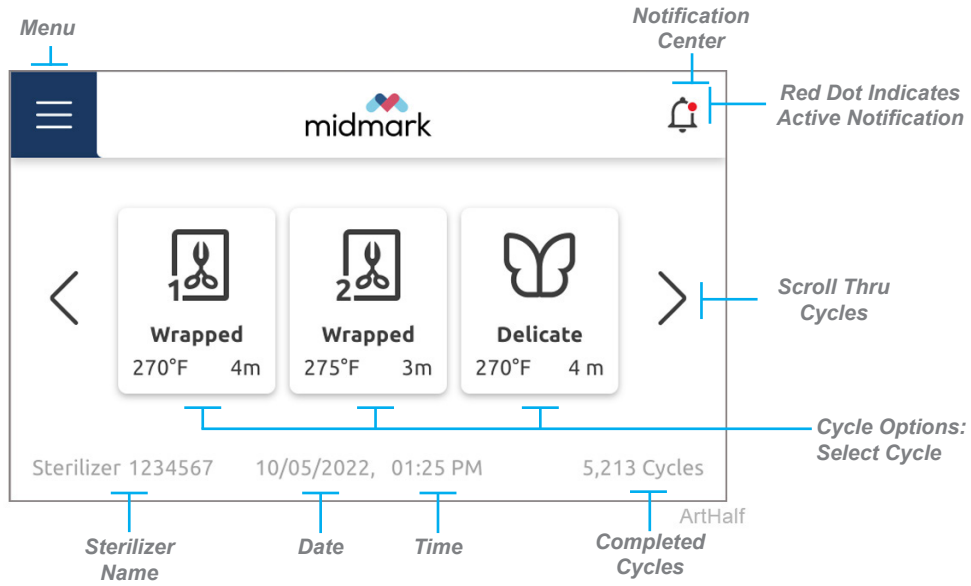
## User Interface Icon Key


Icon	Description	Icon	Description
<	Back to the previous screen Navigate to the left.	✓	Complete
>	Advance to the next screen Navigate to the right	i	Awaiting additional information
✕	Closes current screen / menu	⚠	Status error
→	Shows additional information or options	?	Help menu
↓	Export records	✎	Edit information
⏪	Backspace	≡	Menu items
🏠	Home	⚙	Cycle settings
🔔	Notification Center		

## Additional Menus

Pressing the menu items icon on the left of the home screen will give you access to the following information, menus and settings. Additional menus are shown in more detail within this guide.


**Note**  
All options and menus may not be available for your model sterilizer.
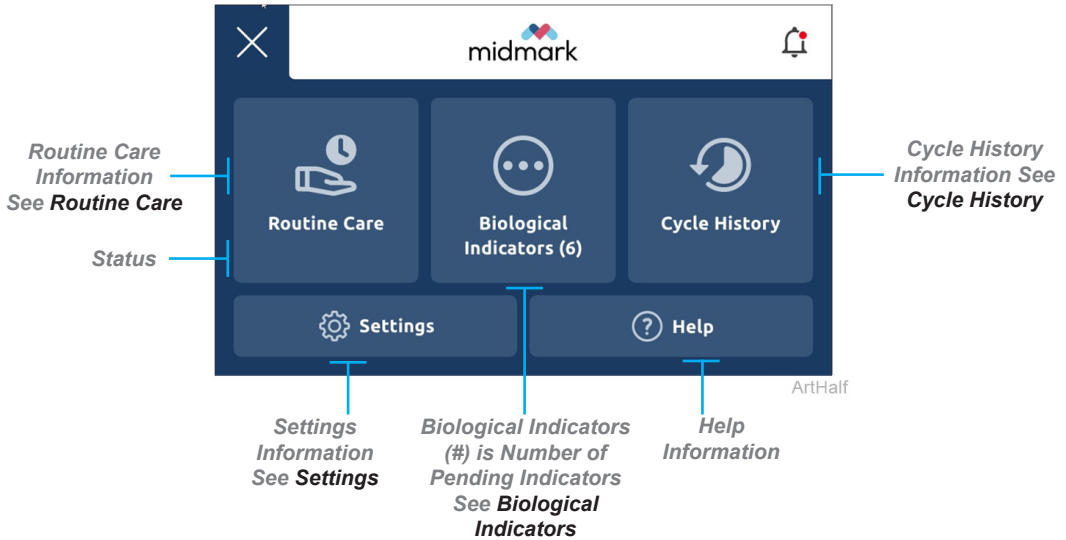
## Biological Indicators (BI)


This screen displays all of the pending BI information as entered by users along with the corresponding cycle information. When the BI results are received they can be entered in the interface by a user as a pass or fail along with any notes. Once the BI results are entered and confirmed for a cycle they are available to view in cycle history.

**Note**  
Notes are for internal use only and not reviewed by Midmark. For additional support contact your Midmark Authorized dealer (see **Calling for Service**).
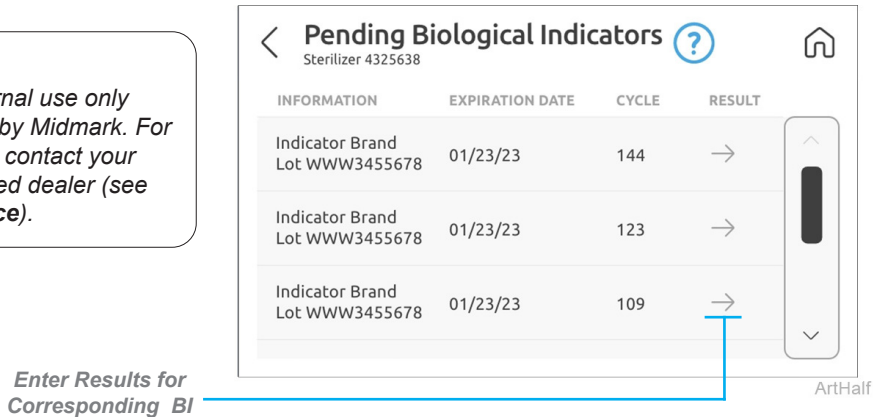
# Cycle History

This screen allows you to view the cycle history with the date, time, cycle type, status and cycle number.

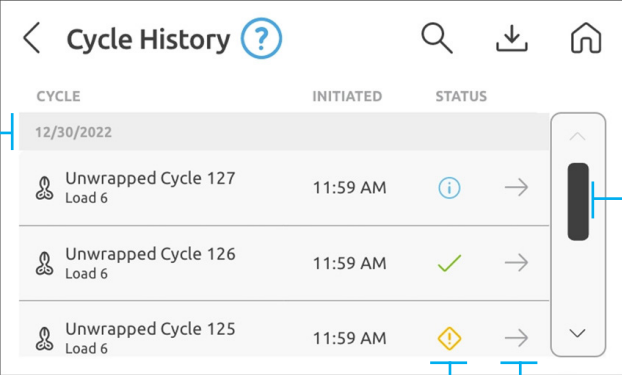*Date of Cycle*

*Scrolls Through Cycle History Information*

*Cycle Status Icon*

*Detailed View of Cycle*


CYCLE STATUS	Description
	Additional information required for cycle such as a pending Biological Indicator (BI) results.
	Cycle completed successfully.
	Indicates the cycle was interrupted; note that symbol color is determined by the severity of the interruption.
	Accesses detailed information: summary of cycle, chemical indicator (CI), BI, notes and details of cycle.

# Settings

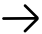Within the Settings menu there are three main tabs across the top **GENERAL**, **CYCLES** and **SYSTEM**. Navigating through these tabs will provide access to information and modifying device settings.

**Note**

*All options and menus may not be available for your model sterilizer.*

General			
Menu Button	Standard	Enhanced	Available Information and Settings
 Service	X	X	Select between the <b>History</b> tab or the <b>Troubleshooting</b> tab.
<b>History</b>			
 Errors	X		This menu shows history of errors including the date when an error occurred, error code, cycle number when error occurred, and description of error.
 Routine Care Events	X	X	This menu shows a history of when routine care occurred and type of routine care completed. Recording who completed routine care is an Enhanced feature.
 Cycles	X		Refer to <b>Cycle History</b> section in this guide for details.
	X		Applicable to all history screens - pressing the right arrow provides detailed information for the selected record.
	X		Applicable to all history screens - this selection requires an optional USB kit that may be used to export records.
<b>Troubleshooting</b>		*** Reference <b>Troubleshooting</b> section in this guide for details *** Ability to run a Diagnostics cycle is an Enhanced feature.	
 About	X	X	View Regulatory, Sterilizer and Software information.
<b>Regulatory</b>		View information such as: the Connectivity Module that contains a radio transmitter, FCC Number, IC number, FCC Compliance Statement and IC Compliance Statements.	
<b>Sterilizer</b>		View information such as: Sterilizer Name, Location, Model Number, Serial Number, Installed and Warranty Thru dates, Midmark Smart View subscription, and Printed circuit board details. Warranty and subscription information are applicable to Enhanced only.	


## Settings - continued...

Menu Button	Standard	Enhanced	Available Information and Settings
 Updates	X	X	View current software information and software update history. Information about scheduled software updates is an Enhanced feature.
 User	X	X	Manage the users by Creating New, Editing or Deleting users. The user list can be sorted from A-Z or from most recent. Disabling the ability to manage users on the sterilizer is an Enhanced feature.
 Water	X		When the water was last changed and when the next change is due. Turn the Auto-Fill optional system on or off.
 Connectivity	X	X	A QR code is available to download the my Midmark mobile app to register the sterilizer. Ability to manage network settings is applicable to Enhanced only.
<b>Cycles</b>			
 Compliance	X	X	This menu allows you to select/deselect these standard features: <ul style="list-style-type: none"> <li>• Require authentication for Cycle Release.</li> <li>• Require authentication for Cycle Interruption.</li> <li>• Require entry of Chemical Indicator (CI) results.</li> <li>• Show Cycle Status in Cycle History and Cycle Details.</li> <li>• Allow Bowie Dick cycle to be ran and require entry of test results.</li> <li>• Require entry of Content Condition for the completion of a cycle.</li> </ul> For Enhanced, these settings are not managed on the sterilizer.
 Instruments		X	Search for, edit or delete existing instruments. Create new instruments. This is an Enhanced feature only.
 Packs		X	Search for, edit or delete existing packs. Create new packs. This is an Enhanced feature only.
 Pre-Heat	X		Toggle Pre-Heat feature on or off. Toggle the Pre-Heat feature on or off for each day of the week. Set the start Pre-Heat time and end Pre-heat time frame for each day.
 Dry Time	X	X	Toggle the ability of the user to decrease dry time on or off. For Enhanced, this setting is not managed on the sterilizer.

# Settings - continued...

Menu Button	Standard	Enhanced	Available Information and Settings
<b>System</b>			
 Date + Time	X		This menu allows editing of the Date + Time and selection of the Date + Time format.
 Language	X		Toggle between English, French and Spanish.
 Notifications	X	X	This menu allows you to select or deselect: <ul style="list-style-type: none"> <li>• Enable Notifications for <b>Overdue Daily Care</b>.</li> <li>• Enable Notifications for <b>Upcoming</b> and <b>Overdue Weekly Care</b>.</li> <li>• Enable Notifications for <b>Upcoming</b> and <b>Overdue Monthly Care</b>.</li> </ul> For Enhanced, these settings are not managed on the sterilizer.
 Units	X		Toggle between Imperial (°F, PSI) or Metric (°C, kPa).
 Appearance	X		Toggle between Light Mode or Dark Mode.
 Sounds + Alerts	X		This menu allows you to turn off non-essential audible alerts.

## Filling the Reservoir
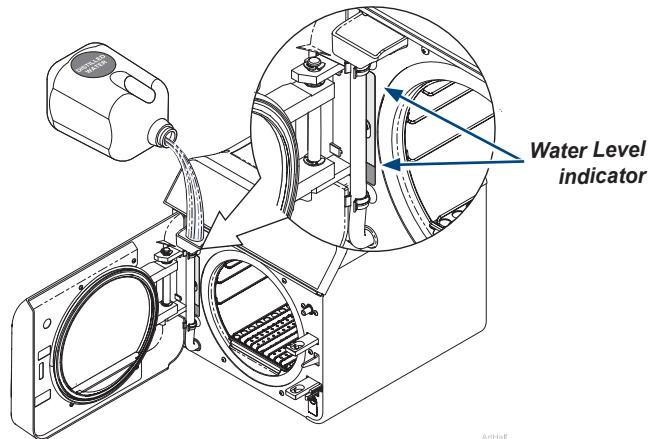
### Step 1: Fill Reservoir

- A) Open door to unit.
- B) Pour distilled water into fill opening until water level reaches the top of the fill level label on the water level indicator tube.

*Note: An Auto-Fill optional accessory is available.*

### Equipment Alert

*Use distilled water or water that meets the referenced water purity specifications. Failure to comply may result in sterilizer malfunction and/or premature failure due to excessive corrosion.*

## Qualification Testing

Qualification testing should be performed prior to placing the sterilizer in service after installation, relocation, major repairs or software updates if indicated in service instructions. Qualification testing should be conducted for the cycle type(s) with the shortest cycle times commonly used in your facility. A test pack that is representative of a load typically processed in the sterilizer should be placed on the bottom tray near the chamber door. Use items that are routinely processed and considered to be the most difficult to sterilize. Also include at least one Biological Indicator (BI) and one Chemical Indicator (CI) in the test pack. Additional items should be placed in the chamber along with the test pack so that the chamber is fully loaded (don't exceed the maximum capacities listed in the **Guidelines for Loading** tables in this guide). Three consecutive test runs, for each cycle type tested, with negative results from the BIs and appropriate readings from all CI's demonstrating complete sterilization, provide verification that the sterilizer will function effectively in the facility in which it is installed. All items processed during qualification testing should be quarantined until the results of the biological testing for all three test runs are available.

## Guidelines for Loading

Items to be sterilized should be arranged to permit free circulation of steam, refer to **General Guidelines** for examples of proper loading.

### **Equipment Alert**

*Do not place items directly on the tray rack or tray plate. Always use a tray or cassette rack - otherwise, serious instrument or equipment damage may occur.*

## Types of Items that can be processed

Before placing any instrument in the M9® and M11® Steam Sterilizer, check with the instrument manufacturer to be sure the materials are compatible with steam sterilization and to verify the acceptability of sterilization parameters. The M9 and M11 Steam Sterilizers are designed to sterilize the following:

- Metal instruments
- Rubber / plastic devices (ex. suction cannulas, impression trays, etc.)
- Steam sterilization packaging for medical devices. (ex. sterilization wrap, paper or plastic pouches, etc.)
- Cassettes (which fit in the sterilizer trays or the cassette rack accessory)
- High / low speed handpieces
- Surgical instruments (ex. Ophthalmology instruments)

### **Equipment Alert**

*Do not sterilize items composed of any of the following materials in the M9 or M11.*

- Corrosion sensitive metals (ex. carbon steel, iron, etc.)
- Fragile items susceptible to breaking under pressure / high temperature
- Biomedical waste
- Plastics that may break down or produce residue when exposed to steam / high temperatures.

#### Examples

*Polyethylene, Styrene, Cellulosics, ABS, PVC, Acrylic (Plexiglass™), PPO (Noryl™), Latex and Neoprene*

## ***Guidelines for Loading - continued...***

### ***Preparing Items for Sterilization***

#### **WARNING**

***Clean and dry instruments thoroughly before placing them in sterilizer. Improper cleaning or drying may result in non-sterile instruments or damage to the unit. Follow instrument manufacturer's guidelines and industry recommendations for handling and cleaning instruments prior to sterilization.***

Instruments must be thoroughly cleaned to remove all residual matter, such as debris, disinfectant residuals, blood, organic tissue, etc. General cleaning guidelines:

- Follow instrument manufacturer's instructions for proper cleaning and preparation.
- Clean instruments immediately after use to avoid drying of residual matter.
- The use of automated cleaning equipment (e.g. ultrasonic cleaner, instrument washer or washer/disinfector) is recommended over manual cleaning for clinician safety and cleaning effectiveness.
- After cleaning, thoroughly rinse instruments with distilled or equivalent water type. After rinsing, instruments should be inspected for damage, debris, detergent residue and then dried before packaging. If necessary repeat the cleaning process.
- If the instrument manufacturer's instructions require lubrication of the instruments after cleaning, wipe off excess lubricant before packing for sterilization or loading into the sterilizer.

### ***Immediate Use Steam Sterilization***

The M9<sup>®</sup> and M11<sup>®</sup> Steam Sterilizers are capable of Immediate Use Steam Sterilization (IUSS)

- Sterilizing instruments for immediate use should only be used in urgent clinical situations.
- Instruments should be thoroughly cleaned and prepared for sterilization.
- Sterilized instruments should be aseptically transferred to the sterile field
- Sterilized instruments should not be stored for later use.
- Sterilization records should indicate that IUSS was utilized.

## **Guidelines for Loading - continued...**

### **Unwrapped**

For an unwrapped cycle, place a surgical cotton towel, paper tray liner or CSR wrap folded to fit on the tray bottom before putting unwrapped items in the tray. Arrange unwrapped items on the towel so they do not touch one another **See Photo 1**.

**Photo 1**

Please consider the following when choosing whether or not to sterilize your instruments unwrapped:

- The sterility of unwrapped instruments is compromised upon exposure to a non-sterile environment. Follow CDC guidelines for using unwrapped, sterilized instruments.
- Due to the sensitive nature of some types of surgery (including, but not limited to ophthalmology), instruments used in such procedures must be wrapped or pouched to maintain sterility and in order to be stored. The water reservoir should also be drained and refilled with fresh distilled water on each day when processing instruments for these procedures on a routine basis.

### **WARNING**

***Do not overload the chamber! Adequate space is required around items in trays for steam circulation and drying. Failure to allow adequate space will compromise sterilization and drying. Items and packaging should be completely dry when removed from the sterilizer to minimize the potential for recontamination.***

### **General Guidelines**

- Use only M9® and M11® trays in their appropriate sterilizer. Using other trays could restrict air / steam flow to items resulting in inadequate sterilization and drying.
- All items must fit within the tray and not extend over the lip of the sterilizer tray. Instruments must not scrape the chamber walls when sliding the tray into the chamber.
- Jointed items must be sterilized in an open position so all surfaces are exposed to the steam (See Photo 2).

## Guidelines for Loading - continued...

### General Guidelines - continued...

#### Note

*Pouches should be loaded paper side down when loaded flat on a tray. This prevents steam from pooling inside the pouch which can prevent proper drying.*


**Photo 2**

- Reusable Instruments, devices and items must be arranged in a single layer on the trays (not piled or stacked), to permit proper steam flow and drying.
- Glassware or utensils capable of holding water (e.g., bottles, basins, beakers) should be positioned on the tray with the open side down so any water condensate drains from the container (See Photo 3). When sterilizing glassware check with the manufacturer to make sure it is compatible with steam sterilization.

**Photo 3**

- Rinse tubing with distilled water or water that meets the referenced water purity specifications (see **Water Purity Specifications**) and do not dry prior to sterilizing. Arrange tubing on the tray so there are no sharp bends and the tubing ends are open and unobstructed (See Photo 4).

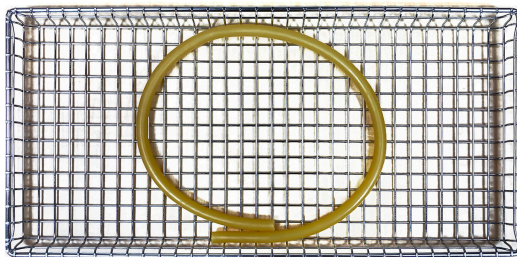**Photo 4**

## Guidelines for Loading - continued...

### General Guidelines - continued...

- Follow the device manufacturer's instructions for disassembly of multi-part instruments prior to packaging and sterilization to assure all parts are adequately exposed to the steam.
- If items are being sterilized and stored for later use they must be packaged (e.g., pouched, wrapped) and completely dry when removed from the sterilizer chamber to avoid potential Recontamination.
- Variations in load configuration, size, wrapping materials and the environment may require the operator to increase the default drying time to assure all packaging and instruments are completely dry. See **Standard Cycle Operation** in this guide for instructions on adjusting drying time. It is also recommended to run a Pre-Heat cycle prior to running sterilization cycles at the start of each day.
- When sterilizing a load with one or more handpieces, utilize the Wrapped 1, Wrapped 2, Unwrapped or Custom cycle depending on the handpiece manufacturer instructions.


### Pouching and Wrapping Items

The M9® and M11® Steam Sterilizers are capable of sterilizing pouched or wrapped items to preserve sterility after processing.

- When pouching or wrapping items, use only sterilization pouches and wraps that have been cleared by the FDA and labeled for use with the steam sterilization cycle being used. Follow the manufacturer's instructions for use.
- Instruments made from different materials (e.g., stainless steel, carbon steel, plastic) should not be mixed in the same pouch or wrapped pack to avoid potential instrument damage.
- Pouches, wraps and their included instruments should not touch the chamber wall. Clearance allows proper steam circulation and avoids potential instrument damage.
- For proper steam circulation and drying the preferred orientation of pouches is resting on their edge, best accomplished using the Midmark Pouch Rack supplied with the M9 and M11 Steam Sterilizers. If additional pouch racks are needed, order the appropriate kit part number. If more than two pouch racks are used in a single load, additional dry time may be required.


**M11 Pouch Rack**  
**(6 slot)**

002-2108-00

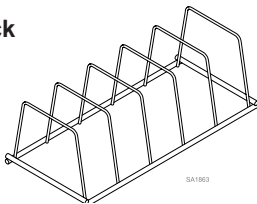**M9 Pouch Rack**  
**(5 slot)**

002-2108-01

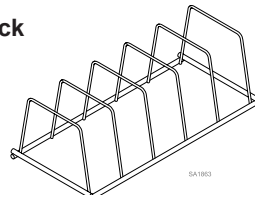- The use of a Pouch Rack requires the removal of some of the sterilizer trays. When using a Pouch Rack, load a single pouch per rack slot. Orient pouches in the rack so the plastic side of one pouch faces the paper side of the adjacent pouch.
- Pouches loaded directly on the trays will dry best if loaded with paper side down.
- Do not layer or overlap pouches in the trays. Pouches should be loosely packed with the handpieces and instruments single height loaded (not piled or stacked), to permit proper steam flow and penetration to the items (See Photo 2).

## Guidelines for Loading - continued...

### Pouching and Wrapping Items - continued...

- If pouches or wrapped packs are labeled using a marking pen the ink should be nontoxic. On paper-plastic pouches only mark on the plastic side of the pouch. On wrapped packs mark the labeling information on the indicator tape or apply a separate label.

Shown here are some recommended pouch load configurations for the M9® (Photo 5 and 6) and M11® (Photo 7 and 8):

#### Equipment Alert

*Do not place items directly on the tray rack or tray plate. Always use a tray or cassette rack - otherwise, serious instrument or equipment damage may occur.*

#### Note

*Heavier loads should be located at the bottom of the sterilizer chamber, unless it will obstruct the transient heat coming from the heater preventing loads above from drying properly.*

M9 Pouches with Pouch Rack  
**Photo 5**

M9 Pouches Mixed  
**Photo 6**

M11 Pouches with Pouch Rack  
**Photo 7**

M11 Pouches Mixed  
**Photo 8**

## Guidelines for Loading - continued...

### Pouching and Wrapping Items - continued...


- When using cassettes in the M9® and M11® Steam Sterilizers, follow the cassette manufacturer's instructions for use.
- Do not wrap items too tightly. If an item is excessively wrapped, sterilization can be compromised and the wrap is more likely to tear if wrapped too tightly.
- Cassettes can be loaded on the trays, but they must fit within the boundaries of the tray and they should not touch each other, the tray above or the chamber wall to allow proper steam circulation. The total instrument load, including the cassettes, should not exceed the limits listed in the maximum capacity tables (see **Load Size**).
- When using the optional Vertical Cassette Rack in the M11, the maximum load including cassettes can be increased to 5.44 kg (12 lb). Additional dry time may be required.

Shown here are some recommended cassette load configurations for the M9 (Photos 9 and 10) and M11 (Photos 11 and 12):

#### **Equipment Alert**

*Do not place items directly on the tray rack or tray plate. Always use a tray or cassette rack - otherwise, serious instrument or equipment damage may occur.*

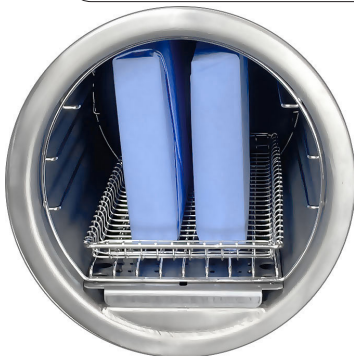
M9 Large Cassette  
**Photo 9**

M9 Mixed Load with Small Cassette and Pouches  
**Photo 10**

M11 Mixed Load with Cassette and Pouches  
**Photo 11**

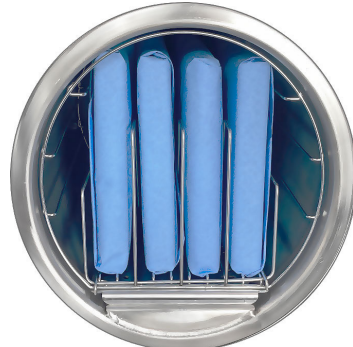M11 with Vertical Cassette Rack  
**Photo 12**

## Guidelines for Loading - continued...


### Textile Loads

- Clean textiles recommended for steam sterilization can be processed in the Delicate Cycle. Verify the acceptability of the Delicate Cycle sterilization parameters with the textile manufacturer to assure the textiles are compatible.
- Textile packs must not exceed the capacity limits specified (see **Standard Cycle Parameters**).
- Textile packs must fit within the boundaries of the sterilizer trays and not touch the chamber wall.
- If multiple packs are placed on a single tray maintain a minimum of 6.4 mm (1/4 in) between packs for proper steam circulation and drying. (Photo 13)

#### **Equipment Alert**

*Do not use towels or packaging containing chlorine bleach residue. Failure to comply may cause rusting or discoloration of the chamber / trays and significantly shorten the life of the sterilizer.*

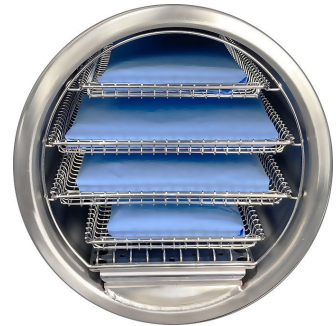Textile load

**Photo 13**

### Load Size

Successful sterilization depends on proper loading of the sterilizer. Do not overload the sterilizer chamber! Adequate space must be maintained around all items placed in the chamber to ensure proper steam circulation and adequate drying. Maintaining proper spacing between all items processed should be the deciding factor in determining the maximum load that can be processed.

# Standard Cycle Parameters

Cycle Type	Cycle Parameters			Drying Time <sup>2</sup>	Items to be Sterilized <i>(Always follow the item manufacturer's recommendations for sterilization.)</i>	M9 Maximum Capacity <sup>4</sup>	M11 Maximum Capacity <sup>4</sup>
	Temperature Minimum	Time	Pressure <sup>1</sup> Reference				
 Wrapped	132°C (270°F)	4 min.	186 kPa (27.1 psi)	30 min.	<ul style="list-style-type: none"> <li>• Pouched or wrapped items manufacturers recommend for exposure at 132°C (270°F) for 4 minutes.</li> <li>• Wrapped cassettes.</li> <li>• Dental Handpieces (wrapped or unwrapped)</li> </ul>	3,629 grams (8.0 lb) or 8 handpieces (2 per tray) with other instruments 3,629 grams (8.0 lb) total	4,082 grams (9.0 lb) or 8 handpieces (2 per tray) with other instruments 4,082 grams (9.0 lb) total
 Wrapped	135°C (275°F)	3 min.	214 kPa (31 psi)	30 min.	<ul style="list-style-type: none"> <li>• Pouched or wrapped items manufacturers recommend for exposure at 135°C (275°F) for 3 minutes.</li> <li>• Wrapped cassettes.</li> <li>• Dental Handpieces (wrapped or unwrapped)</li> </ul>	3,629 grams (8.0 lb) or 8 handpieces (2 per tray) with other instruments 3,629 grams (8.0 lb) total	4,082 grams (9.0 lb) or 8 handpieces (2 per tray) with other instruments 4,082 grams (9.0 lb) total
 Delicate	121°C (250°F)	30 min.	104 kPa (15 psi)	30 min.	<ul style="list-style-type: none"> <li>• Textiles and surgical packs wrapped for sterilization. <sup>3</sup></li> <li>• Items, except liquids, manufacturers recommend for exposure at 121°C (250°F) for 30 minutes.</li> </ul>	3,629 grams (8.0 lb) or 590 grams (1.3 lb) textile load	4,082 grams (9.0 lb) or 2 packs each @ 590 grams (1.3 lb) textile load
 Unwrapped	132°C (270°F)	3 min.	186 kPa (27.1 psi)	30 min.	<ul style="list-style-type: none"> <li>• Instruments loose on a tray.</li> <li>• Open glass or metal canisters.</li> <li>• Tubing not used in surgical procedures. (Max. length - 40 in and Min. inside diameter - .187 in)</li> <li>• Loose items manufacturers recommend for exposure at 132°C (270°F) for 3 minutes.</li> </ul> <p><i>Note: The sterility of unwrapped items is compromised on exposure to a non-sterile environment.</i></p>	3,629 grams (8.0 lb)	4,082 grams (9.0 lb)

- The pressures shown in this table are at sea level and are for reference only. These are the ideal pressure of saturated steam at the sterilization temperature. The pressures on the sterilizer display may be higher.
- Dry time can be changed from 5 to 60 minutes. (USS dry time can be changed from 1 to 5 minutes.) Refer to **Standard Cycle Operation**.
- Allow a minimum of 6.4 mm (1/4 in) space between each pack and from the chamber wall.
- The default dry time may need increased due to variations in load configuration, wrapping materials and the environment to completely dry the chamber contents at these capacities.

**Note**  
Some features may be disabled in Settings, such as ability to reduce dry time or ability to run an IUSS cycle.

# Custom Cycle Parameters

Cycle Type	Cycle Parameters			Drying Time <sup>2</sup>	Items to be Sterilized <sup>3</sup> <i>(Always consult the item manufacturer's recommendations for sterilization.)</i>	M9 Maximum Capacity <sup>4</sup>	M11 Maximum Capacity <sup>4</sup>
	Temperature Minimum	Time	Pressure <sup>1</sup> Reference				
 Custom	121°C (250°F) to 135°C (275°F)	3 min. to 45 min.	104 kPa (15 psi) to 214 kPa (31 psi)	0 min. to 60 min.	Create different programmed cycles for special applications. Sterilization time and temperature, along with drying time and venting procedure can be adjusted or changed to follow instrument manufacturers' instructions that are outside Midmark Standard Cycle Parameters. Refer to <b>Custom Cycle Operation</b> for how to create.	3,629 grams (8.0 lb)	4,082 grams (9.0 lb)

### Caution

These cycles are not FDA cleared and validation of sterility of items processed using a custom cycle is the responsibility of the user.

1. The pressures shown in this table are at sea level and are for reference only. These are the ideal pressure of saturated steam at the sterilization temperature. The pressures on the sterilizer display may be higher.
2. Dry time can be changed from 5 to 60 minutes. (If Slow Vent option is enabled, the dry time will be 0) Refer to **Standard Cycle Operation**.
3. Allow a minimum of 6.4 mm (1/4 in) space between each pack and from the chamber wall.
4. The default dry time may need increased due to variations in load configuration, wrapping materials and the environment to completely dry the chamber contents at these capacities.

**Note**  
The ability to run a Custom Cycle may be disabled in Settings.

# Standard Cycle Operation

## Note

Pre-Heat can be turned on in **Settings** and activated for certain days of the week and for a time range on those days. This feature keeps the chamber pre-heated, potentially shortening the total cycle time.

## Equipment Alert


Using an incorrect sterilization cycle could result in non-sterile items and may damage instruments. Follow instrument manufacturer's recommendations for sterilization.

## Step 1: Select the Correct Cycle

- A) Use right and left arrows to locate correct cycle.
- B) Refer to **Standard Cycle Parameters** chart for time and temperature specifications.

## Note

On units using the metric display, °F will display as °C and psi will display as kPa.
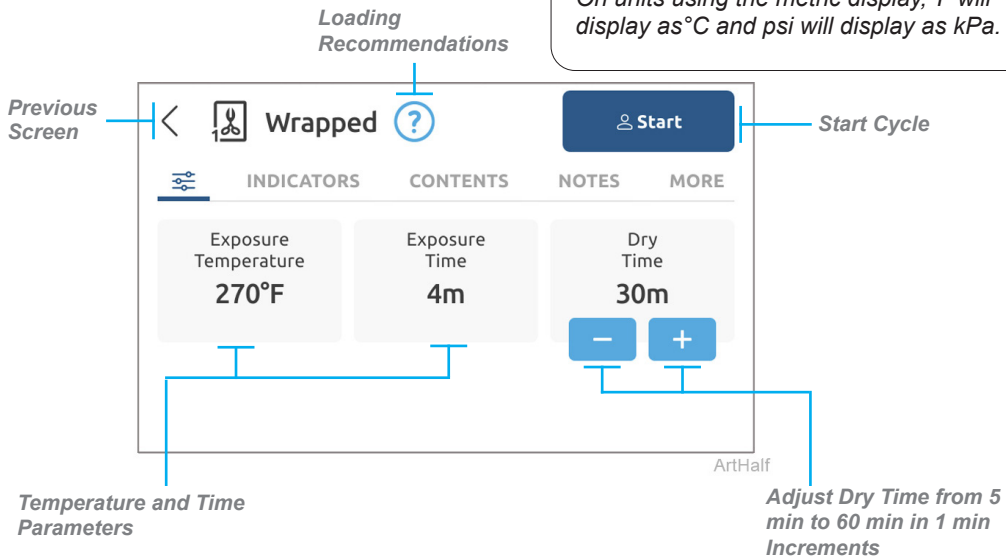
## Step 2: Start Cycle

- A) Press start.
- B) Select user, enter PIN and press start.

## WARNING

When closing the door avoid contacting the bottom of the door handle as it may be hot to touch.

## Standard Cycle Operation - continued...

The sterilizer will progress through the following listed stages. As the cycle progresses the status bar will move to the right. The current temperature and pressure in chamber will be displayed. A timed countdown will start in the sterilization stage and continue until the cycle and drying are complete.

- **Filling** Chamber begins filling with water. When water reaches the proper level the sterilizer proceeds to the next stage.
- **Heating** Heating element turns on heating up the chamber. Status bar progresses, temperature and pressure changes are shown. Air purge routines occur during this phase.
- **Sterilizing** Sterilizing begins when correct temperature and pressure are reached. Time remaining in the cycle counts down while current temperature and pressure in the chamber are continuously updated. Heating element turns on and off, maintaining temperature and pressure. Status bar continues to progress.
- **Venting** When the sterilizing mode completes, the vent valve opens. Steam and water are released back into the reservoir. Status bar progresses, temperature and pressure changes are shown.
- **Door Opening** An audible signal is emitted to indicate that the door is about to open. When pressure in chamber reaches zero, the door actuates to the partially open (drying mode) position.

### **WARNING**

*Keep clear when sterilizer door is ready to open! Failure to do so could result in severe burns from steam being released.*

- **Drying** During the dry cycle the heating element is cycled on and off. Time of Dry Cycle is counted down.

### **WARNING**

*The processed loads may still be wet if the Dry Cycle is aborted prior to completion. To avoid re-contamination, always allow processed loads to dry in the sterilizer before handling.*

- **Cycle Complete** Will appear when the cycle is completed.
- **Cycle Summary** Will appear after the Cycle Complete notification.


### **WARNING**

*STOP button may be pressed at any time to stop or interrupt a cycle. Items must not be considered sterile if this occurs before the Dry Cycle is complete or instruments are wet. Sterilizer will return to SELECT CYCLE mode.*

# Standard Cycle Operation - continued...

## Step 3: Confirm Cycle

- Press next.
- Select PASS/FAIL for the Chemical Indicator (CI).
- Select DRY/WET for the Contents Condition.
- Select confirm, select a user and enter the PIN.

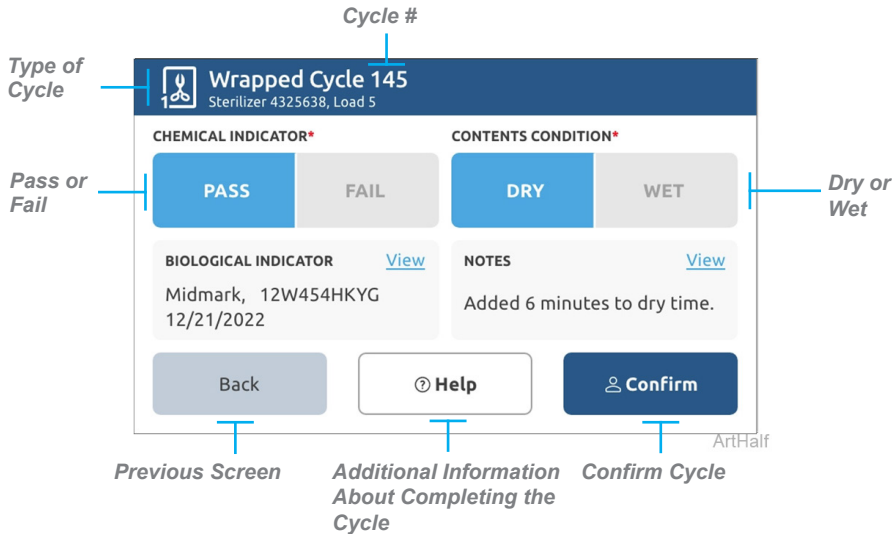- View (BI)** Allows user to see manufacturer, lot number and expiration date of the BI used.
- View Notes** Allows user to view and add notes pertaining to the cycle that just finished.

### Note

Notes are for internal use only and not reviewed by Midmark. For additional support contact your Midmark Authorized dealer (see **Calling for Service**).

### WARNING

- The processed load and inner surfaces will be hot. Avoid contact with hot surfaces. Failure to do so could result in serious burns.
- Do not apply downward pressure to the open sterilizer door. Doing so could cause the sterilizer to tip resulting in serious burns, or injury from trays or instruments sliding out of the chamber.
- If the door does not open due to the chamber not automatically venting, the chamber may be manually vented by pulling up on the pressure relief lever.

# Post-Sterilization Handling

After sterilization is complete, all items must be handled in accordance with accepted and documented standards, such as the Centers for Disease Control and Prevention (CDC) documents, “Guidelines for Infection Control in Dental Healthcare Settings” - 2003, MMWR; 52 (no. RR-17) and “Guidelines for Disinfection and Sterilization in Healthcare Facilities” - 2008, as well as any local requirements that may apply.


Qualified personnel responsible for infection prevention should prepare a protocol for handling sterilized items. This protocol should be followed by all personnel responsible for handling sterilized items.

## Custom Cycle Operation

The Custom Cycle can be used for custom applications that are not covered by standard cycles. These cycles have not been FDA-cleared for medical use. The Custom Cycle allows the cycle to be adjusted or changed to follow instrument manufactures instructions that are outside Midmark Standard Cycle Parameters. Once a Custom Cycle has been saved it will appear in the main cycles menu and can be used by scrolling to it and selecting it. Settings within these Custom Cycles may be changed with options to Save as New Preset or Update This Preset.

**WARNING**  
*These cycles are not FDA cleared. Validation of sterility of items processed using them is the responsibility of the user.*

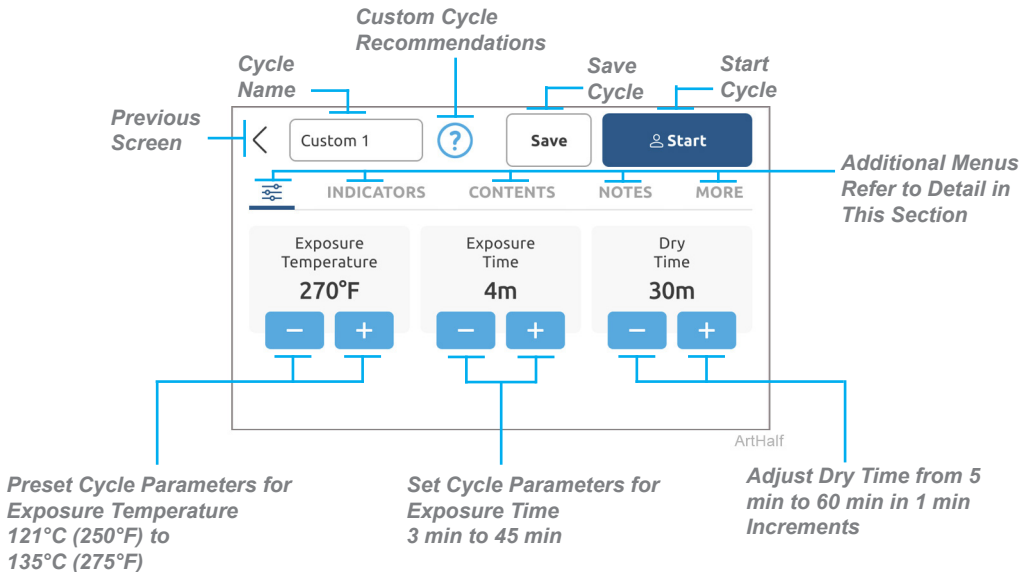**Note**  
When creating a Custom Cycle a custom name can be given to it.

**Note**  
IUSS dry time can be changed from 1 to 5 minutes. If Slow Vent option is enabled the dry time will be 0.

## Custom Cycle Operation - continued...

### **Note:**

*Custom Cycles can be programmed for monthly cleaning cycles with no dry time to expedite monthly maintenance.*

Once a Custom Cycle is started, the sterilizer will progress through the same basic stages of sterilization as a standard cycle. As the cycle progresses, the status bar will continue to move right. The current temperature and pressure in the chamber will be displayed. A timed countdown will start in the sterilization stage and continue until the cycle and drying are complete.

### **Note**

*All options and menus may not be available for your model sterilizer.*

-  Main custom cycle parameters adjustments: Exposure Temperature, Exposure Time and Dry Time.
- **Indicators** Select CI and/or BI indicators to be included in the custom cycle.
- **Contents** Select add past content or add/create a new pack or instrument to be included in the custom cycle.
- **Notes** Add a note to the cycle.

### **Note**

*Notes are for internal use only and not reviewed by Midmark. For additional support contact your Midmark Authorized dealer (see **Calling for Service**).*

- **More** Selection of:
  - Immediate Use Steam Sterilization (IUSS) settings.
  - Designate this cycle as Qualification Cycle.
  - Slow Vent can be activated.

### **Note**

*IUSS dry time can be changed from 1 to 5 minutes.*

# Bowie-Dick Cycle

For the M9® and M11® Steam Sterilizers the Bowie-Dick test is optional and not required. A Bowie-Dick test is a sensitive and rapid means of detecting air leaks and inadequate steam penetration in dynamic-air-removal sterilizers. The test pack should be placed horizontally in the front of the bottom tray with an otherwise empty chamber. Follow the Bowie-Dick test instructions for proper storage, use, interpretation and disposal.

## **Caution**

*This cycle is not intended to be used to sterilize instruments.*

### **Step 1: Start Preconditioning Cycle**

- A) Empty the chamber of any items or test materials.
- B) Press start.

*Note: The preconditioning cycle is performed to raise the sterilizer to operating temperature.*

### **Preconditioning Cycle**

- The Preconditioning Cycle will run a complete cycle (filling, heating, sterilizing, venting, door opening and cycle complete).
- There is no drying operation used for the Preconditioning Cycle.

## **WARNING**

*The sterilizer is at operating temperature - use a tray removal tool or insulated pad to handle the metal trays.*

### **Bowie-Dick Cycle**

#### **Note**

*A test pack must be selected before the Bowie-Dick cycle will run.*

### **Step 2: Start Bowie-Dick Cycle**

- A) Carefully place test pack near front of lower tray.
- B) Select a Test Pack and press start.

### **Step 3: After Cycle is Complete Confirm Cycle**

- A) Press next.
- B) Select PASS or FAIL for the Test Pack.
- C) Add notes if needed.
- D) Select confirm, select a user and enter the PIN.
- E) Press Home.

# Routine Care

To ensure proper operation and maximum sterilizer life, carefully follow all instructions for routine care.


One of the most important ways to prevent problems with your sterilizer is to only use distilled water, or water that meets the referenced water purity specifications (see **Water Purity Specifications**). DO NOT USE TAP WATER. The sterilizer operates with high water temperatures and minerals dissolved in water can form mineral deposits that can prevent valves from opening or closing properly and lead to corrosion in the chamber and tubing.

Routine care reminders will be displayed on the screen at the appropriate intervals to assist the operator.

## Note

Following each cleaning, visually inspect to ensure no material remains. Re-clean if necessary.

This screen has three main tabs: Daily, Weekly and Monthly. Each of these tabs will show it's current status of up to date, due or overdue with the corresponding color bar. Routine care that is due or overdue will have a Run button that steps through the required care.

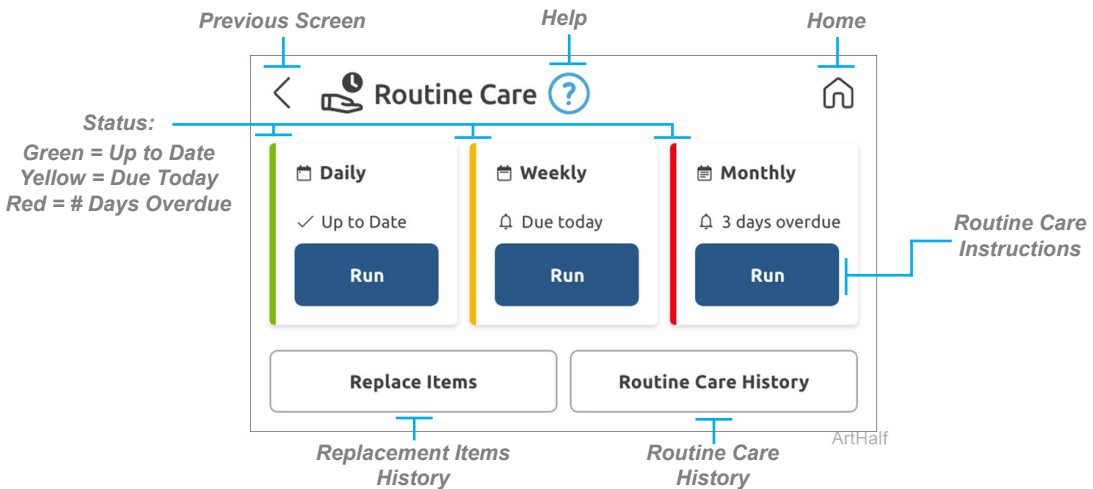Replace Items

The "Replace Items" button allows users to track consumable parts, which are stored in history and can be exported. There are also installation help screens for the user to follow.

Routine Care History

The "Routine Care History" button allows users to view the daily, weekly, and monthly care event history and export these records.

## Daily Care

### **WARNING**

*Wait until the sterilizer is cool to the touch before beginning.*

### **Equipment Alert**

*If you frequently process lubricated instruments, drain and refill the reservoir with distilled water or water that meets water purity specified in the user guide.*


Note: If an Auto-Fill system is enabled proceed to the next step after draining the reservoir.

### **Step 1: Clean and Disinfect External Surfaces**

- A) *Using a soft cloth, clean external surfaces of the sterilizer with soap and warm water.*
- B) *Disinfect the external surfaces of the sterilizer with a quaternary disinfectant according to your facility guidelines for clinical contact surfaces.*
- C) *Follow the instructions provided with the cleaner/disinfectant used regarding rinsing and drying.*

### **Step 2: Inspect and Clean Gaskets**

- A) *Clean exposed surfaces of the door and dam gaskets with a soft cloth moistened with distilled water.*
- B) *If necessary, use Speed-Clean and a soft brush.*
- C) *Inspect door and dam gaskets for damage such as nicks, cuts, shrinkage, or swelling. Replace if damaged.*

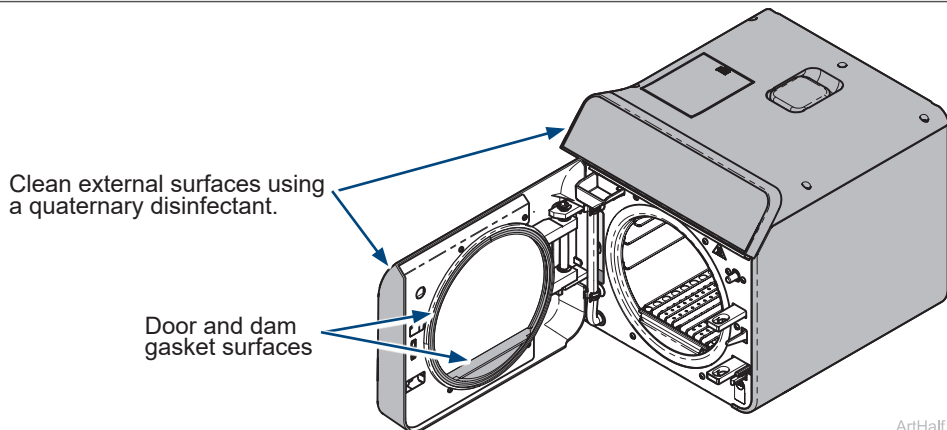Mark Completed

The “Mark Completed” button allows users to confirm they have completed the daily care checklist.

View Guided Walkthrough

The “View Guided Walkthrough” button gives users detailed steps and illustrations to guide them through completing the daily care process.

## Weekly Care

### WARNING

- Wait until the sterilizer is cool to the touch before beginning.
- Handle metal tray components carefully to avoid injury.
- Do not run sterilizer without tray plate in place.


### Equipment Alert

- Failure to change water and clean sterilizer may result in malfunction.
- Use only distilled water or water that meets purity specified in this user guide.
- Do not use bleaching agents or any abrasive materials or substances to clean.
- Use care when wiping the inside of the chamber to avoid damage to sterilizer components.

### Note

On weeks when performing **Monthly Care** the **Weekly Care** is not required.
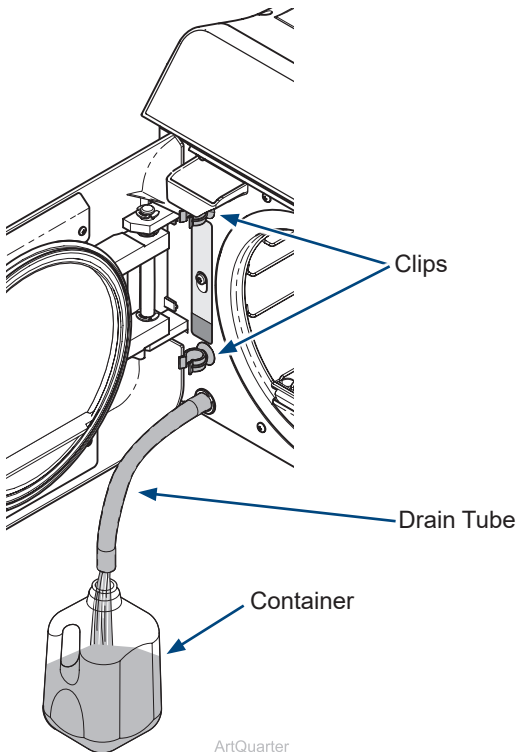
### Step 1: Complete Daily Care

- A) Clean external surfaces and gaskets. See **Daily Care** for details.

### Step 2: Drain Water from Reservoir

- A) Disconnect the drain tube from the clips. Bend tube downward and drain reservoir into a suitable container.
- B) Re-install drain tube into clips.

## Weekly Care - continued...


### Step 3: Remove + Clean Trays, Tray Rack and Tray Plate from Chamber

- Remove trays.
- Pull tray rack and tray plate straight out of chamber then separate them.
- Wash trays, tray rack, tray plate and inside of chamber with Speed-Clean and rinse with distilled water.

#### Equipment Alert

When using Speed-Clean to clean surfaces dilute with 20 parts distilled water to one part Speed-Clean.
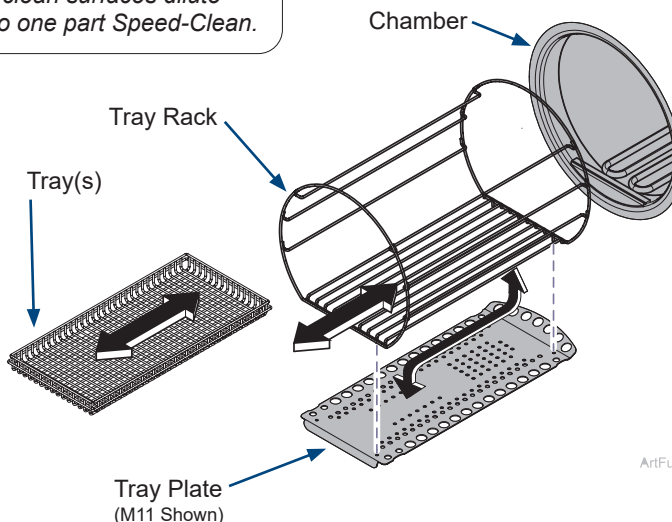
### Step 4: Re-install Trays, Tray Rack and Tray Plate into Chamber

- Insert tray rack back onto the tray plate.
- Install tray rack/plate with curved end of tray plate tipped upward and toward back of the chamber.
- Re-install trays.
- Refill reservoir with water that meets referenced water quality specifications (see **Water Purity Specifications**).

Mark Completed

The "Mark Completed" button allows users to confirm they have completed the daily care checklist.

View Guided Walkthrough

The "View Guided Walkthrough" button gives users detailed steps and illustrations to guide them through completing the daily care process.

## Monthly Care

### WARNING

- Wait until the sterilizer is cool to the touch before beginning.
- Handle metal tray components carefully to avoid injury.
- Do not run sterilizer without tray plate in place.


### Step 1: Complete Daily Care

A) Clean external surfaces and gaskets. See **Daily Care** for details.


### Equipment Alert

- Use only distilled water or water that meets purity specified in this user guide.
- Use only Speed-Clean to flush the system during the cleaning cycle.
- Do not process instruments during cleaning and rinsing cycles.
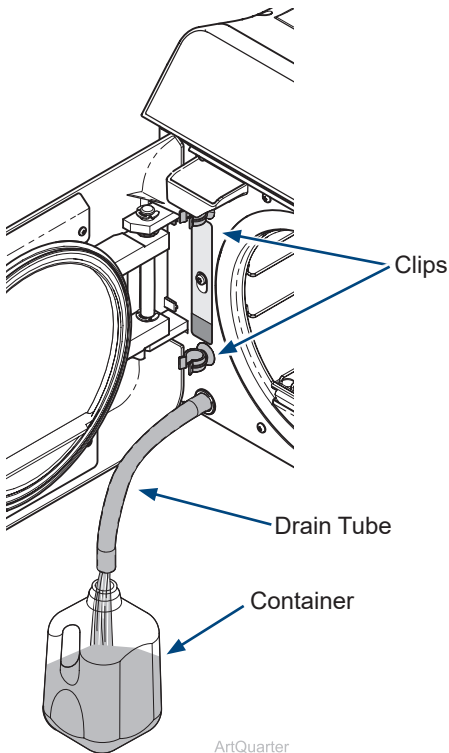
### Step 2: Replace Water in Reservoir


- A) Disconnect the drain tube from the clips. Bend tube downward and drain reservoir into a suitable container.
- B) Re-install drain tube into clips.
- C) If an Auto-Fill system is enabled proceed to the next step. Otherwise, manually refill the reservoir with water that meets referenced water quality specifications.
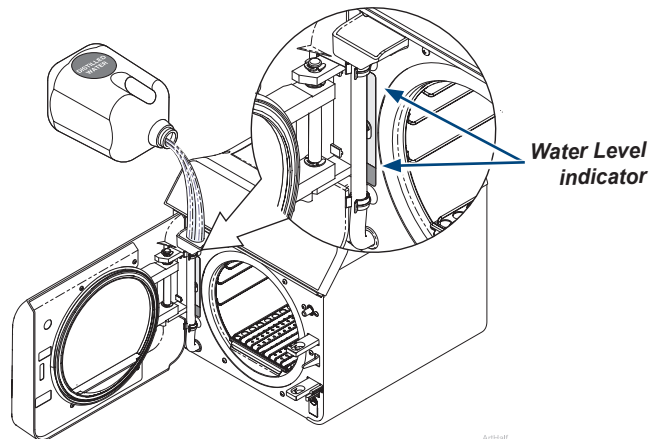
## Monthly Care - continued...

### Step 3: Run Cleaning Cycle

- A) Pour one ounce of Speed-Clean into chamber.
- B) Close the door and start a Cleaning Cycle.
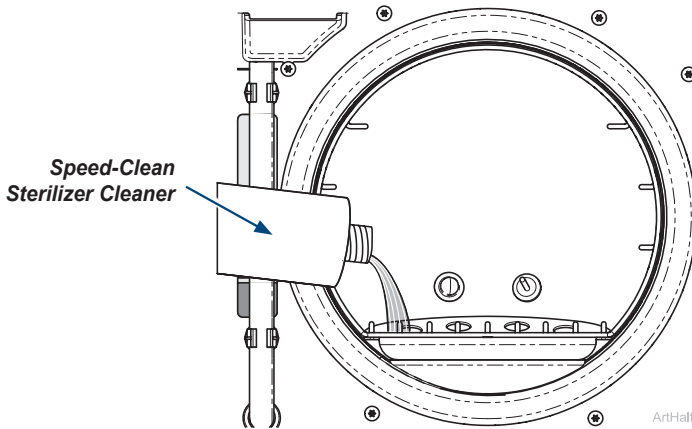
### Step 4: Run Rinsing Cycle

- A) Drain reservoir.
- B) If an Auto-Fill system is enabled proceed to the next step. Otherwise, manually refill the reservoir with water that meets referenced water quality specifications (see **Water Purity Specifications**).
- C) Close the door and start a Rinsing Cycle.


## Monthly Care - continued...

### Note

The user interface will remind the user every 6 months that the pressure relief valve test should be performed.

### Step 5: Test Pressure Relief Valve

- When prompted on the display pull upward firmly on the pressure relief lever for three (3) seconds, then release.
- Steam should discharge freely from beneath rear of the unit when lever is pulled.
- If the valve does not close completely when lever is released, pull lever again and release quickly so that it snaps closed. Repeat until valve reseats properly.
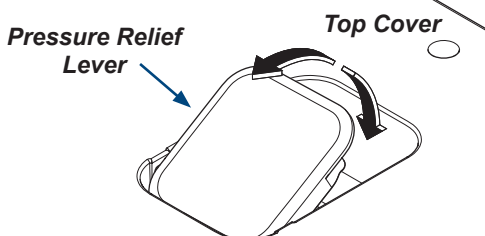
### Equipment Alert

Contact your local supplier/ technician to order a replacement pressure relief valve if: Excessive force is required to open the pressure relief valve or the valve will not reseat properly.


### Step 6: Remove + Wipe Down Trays, Tray Rack, Tray Plate and Chamber

- Drain reservoir.
- If an Auto-Fill system is enabled proceed to the next step. Otherwise, manually refill the reservoir with water that meets referenced water quality specifications (see **Water Purity Specifications**).
- Remove trays.
- Pull tray rack and tray plate straight out of chamber then separate them.
- Wipe any residue off trays, tray rack, tray plate, and chamber with a soft cloth and rinse with distilled water.

### WARNING

**Sterilizer is hot. Wait until the sterilizer is cool to the touch before continuing.**
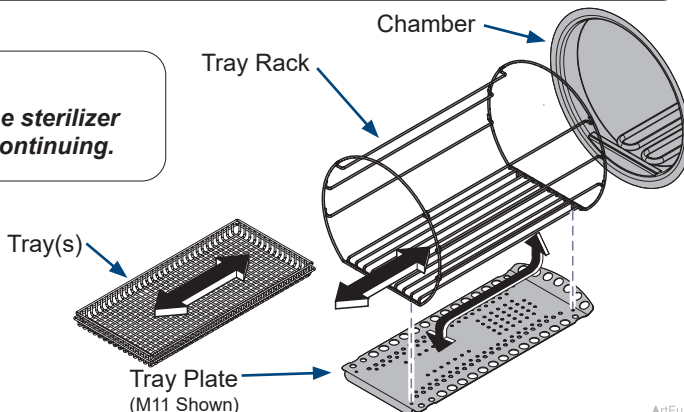
## Monthly Care - continued...

### Equipment Alert

Filters protect sterilizer valves and plumbing. If cycle times become lengthy, filters should be cleaned or replaced.


### Step 7: Clean Filters

- A) Remove filters by grasping filter then gently twisting and pulling away from the chamber wall.
- B) Clean filters using Speed-Clean. A small stiff bristled brush or ultrasonic cleaner may be helpful. Rinse with distilled water.
- C) Replace filters if damaged or if debris cannot be removed by cleaning.
- D) To reinstall filter, firmly press inward toward chamber wall while twisting.

### Equipment Alert

When using Speed-Clean to clean surfaces dilute with 20 parts distilled water to one part Speed-Clean.
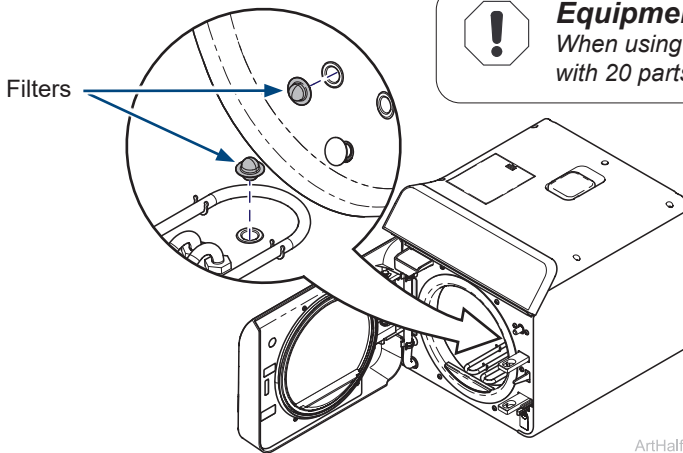
### Step 8: Re-install Trays, Tray Rack and Tray Plate

- A) Insert tray rack back onto the tray plate.
- B) Install tray rack/plate with curved end of tray plate tipped upward and toward back of chamber.
- C) Re-install trays.

## Extended Use Care

The M9® and M11® Steam Sterilizers are designed and tested to provide exceptional reliability throughout their service life. However, like all electro-mechanical devices they are subject to wear and degradation with use. To ensure the integrity, performance and safety of all major components it is the responsibility of the user to have the sterilizer performance / operation verified by a Midmark® Authorized Service Provider at least every 10 years or 25,000 cycles, whichever ever comes first. After 10 years or 25,000 cycles of use, an annual inspection by a Midmark Authorized Service Provider is recommended. Software support may not be available after 10 years.

# Accessory and Tools Details

## Using the Optional Tray and Cassette Tool


### WARNING

- Only use the 9A307001 Tray and Cassette Tool with Midmark® manufactured trays.
- Trays / cassettes may be HOT - Use care when removing or transporting. Hold the tray level and slightly elevated to prevent it from shifting and becoming dislodged. Failure to comply may result in personal injury due to burns.

### Step 1: Remove trays and cassettes

- Slide the proper plate into the handle as shown.
- Hook the top, saw-toothed tab of tool to top center of tray lip / cassette.
- Rotate tool downward until the tool plate is completely under tray / cassette.
- Check to ensure tray / cassette is being held securely and then remove from chamber.

Note: Cassette tool can accommodate cassettes up to 3.8 cm (1 1/2 in) thick.
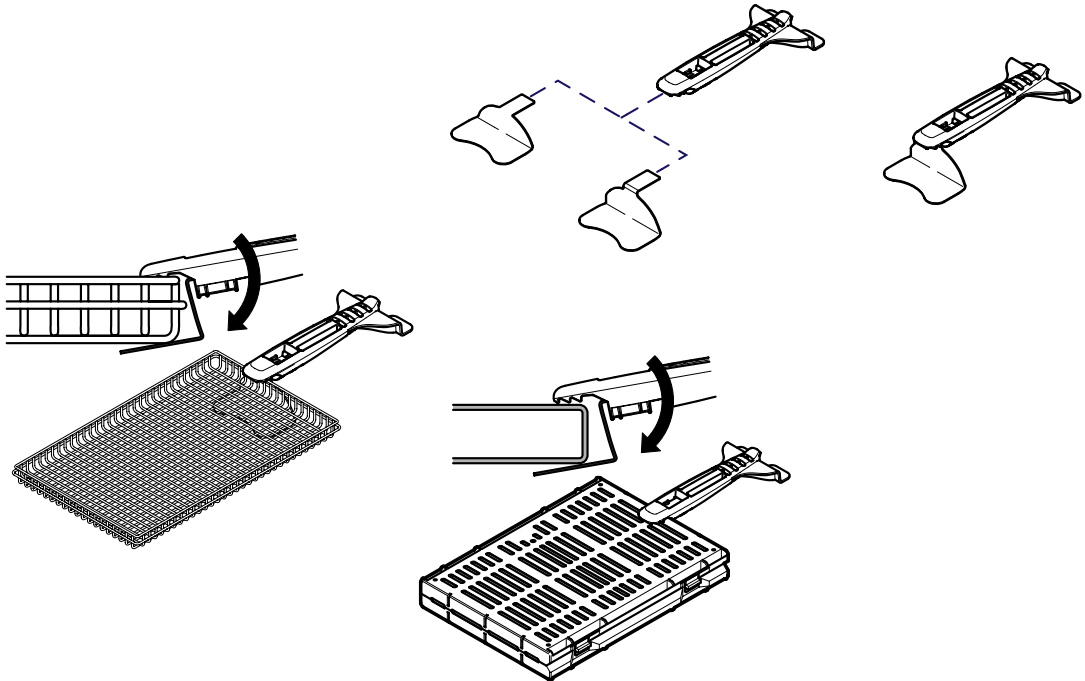
## ***VistaCool™ Direct-to-Drain System***

VistaCool Direct-To-Drain is a system designed to eliminate the need for autoclave condenser/waste bottles and condensation tanks by cooling the exhaust water and automatically sending it directly down the drain.

## ***USB Kit***

The USB kit may be used to electronically transfer sterilization cycle records and routine care records to a USB flash drive. All sterilizer records are stored internally for the life of the device and may be viewed on the sterilizer display. Records will also be transferred to the cloud if the sterilizer is cloud connected. The USB kit is available if the sterilizer is not cloud connected and there is a need to independently store records.

## ***Auto-Fill***

The Auto-Fill system can be connected to a facility water treatment system, see Water Purity Specifications for required water quality. Reservoir water level is monitored and will control a solenoid valve when water level is low to automatically refill the reservoir. The optional VistaCool Direct-to-Drain system may also be used with the Auto-Fill system.

## ***Vertical Cassette Rack (M11® Only)***

The Vertical Cassette Rack allows cassettes to be sterilized in the vertical position. When using the Vertical Cassette Rack in the M11 the maximum load including the cassettes can be increased to 5.44 kg (12 lb). Additional dry time may be required.

# Troubleshooting

Use the following chart to assist in troubleshooting the sterilizer.

## Equipment Alert

*Sterilizers are not intended to be serviced by users. Do not remove covers.*

Problem	Possible Cause	Solution
<b>Sterilizer does not operate (no display).</b>	Sterilizer power cord came loose from supply outlet or back of the sterilizer.	Ensure sterilizer power cord is plugged into outlet and sterilizer.
	No power at sterilizer supply outlet.	Check circuit breaker for supply outlet. If problem recurs, unplug unit power cord and contact a certified electrician.
	Internal issue with sterilizer.	Unplug unit power cord and contact an authorized service technician. See <b>Calling for Service</b> .
<b>Steam escaping from pressure relief valve.</b>	Pressure relief valve is not properly seated.	Reseat pressure relief valve. See <b>Monthly Care</b> .
<b>Sterilization failure evidence from process monitor (e.g., chemical indicator, biological indicator, etc.).</b>	Sterilization conditions were not present at the location of the indicator.	Reload sterilizer in accordance with guidelines for loading the M9® or M11®. If problem recurs, unplug unit power cord and contact an authorized service technician. See <b>Calling for Service</b> .
	Indicator is out-of-date, is inappropriate for sterilization cycle, has otherwise malfunctioned, was incorrectly handled or not completely dry.	Use an indicator, appropriate for the load and cycle selected, that has been stored properly. Contact the indicator manufacturer for additional information on proper selection, use, storage and potential misapplication or malfunction.
	Improper steam conditions.	Use proper water quality per <b>Water Purity Specification Chart</b> .
	Door gasket is leaking.	Clean, remove and reinstall or replace door gasket.
<b>Door gasket leaks.</b>	Door gasket is damaged, dirty or improperly installed.	Clean, remove and reinstall or replace door gasket.
<b>Items are not dry at end of Dry Cycle.</b>	Sterilizer is improperly loaded.	Reload sterilizer in accordance with guidelines for loading the M9 or M11. If the particular load configuration being used is consistently not dry at the end of the dry cycle increase the cycle drying time or reduce the load.
	Filter screen(s) are clogged.	Clean or replace filter screen(s). See <b>Monthly Care</b> .
	Sterilizer is not level.	See <b>Leveling Procedure</b> .
<b>Display goes blank / black / is non legible.</b>	Electromagnetic interference.	Press any key except Stop on the keypad to reset the display without aborting the cycle. Pressing stop will also reset the display but will also abort the cycle in process.


## Troubleshooting - continued...

Problem	Possible Cause	Solution
Discoloration of CSR Wrap, Pouches etc.	Sterilizer is improperly loaded.	Reload sterilizer in accordance with guidelines for loading the M9 or M11.
	Sterilizer is not level.	See <b>Leveling Procedure</b> .
	Heating element not uniformly submerged in water during cycle.   <b>WARNING</b> <i>Wait until the sterilizer is cool to the touch before beginning.</i>	Remove trays, tray rack and tray plate from chamber. Start a sterilization cycle and monitor the display. Press STOP button after the fill mode is complete to interrupt cycle. Heating element should be uniformly submerged in water. If not, adjust front leveling feet so the heating element is uniformly submerged. If issue continues, unplug unit power cord and contact an authorized service representative. See <b>Calling for Service</b> .

## Troubleshooting Menus

### Step 1: Access Troubleshooting

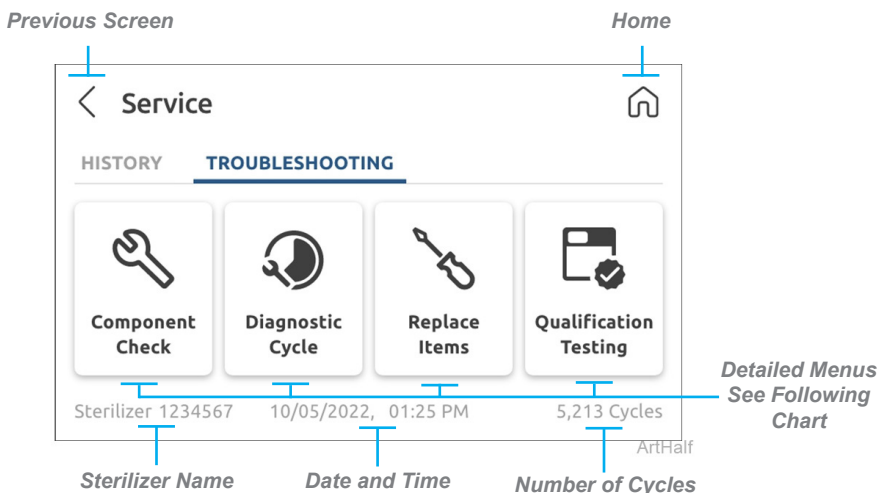- From the home menu press settings.
- Select the Service icon.
- Choose the Troubleshooting Tab.

## Troubleshooting Menus - continued...

### Note

All options and menus may not be available for your model sterilizer.

Settings / General			
Menu Button	Standard	Enhanced	Available Information / Settings
 Service	X	X	Select between the <b>History</b> tab and the <b>Troubleshooting</b> tab
Troubleshooting			
 Component Check	X		Component Check provides digital access to test major components of the sterilizer without running a complete cycle. The status of switches and sensors can be viewed also.
 Diagnostic Cycle		X	A Diagnostic Cycle performs a simulated cycle to verify sterilizer functionality. This cycle may take up to 45 minutes to complete. This is an Enhanced feature only.
 Replace Items	X	X	Replace Items allows for tracking of when commonly replaced parts are replaced. Replacement History can also be viewed. Check boxes to indicate which parts are being replaced. Click on details to learn how to replace items. Recording who completed routine care is an Enhanced feature.
 Qualification Testing		X	Qualification Testing should be run following sterilizer installation, major repairs, and certain software updates. A service technician or a software download may call for the Qualification Testing to be initiated. The user interface will ask the user if the cycle they are attempting to run is a qualification cycle until the three (3) qualification cycles are completed.

## Error Codes

In the event of a malfunction, the sterilizer will generate an error code that will be displayed on the user interface. Additionally, the interface will provide essential information including: the sterility status of the contents, the reason for the error, troubleshooting steps the user can take to potentially resolve the issue, and whether it is necessary to reach out to a authorized service representative.

### WARNING

- **If an error occurs more than once, stop using the sterilizer. Record the message or error code, unplug the unit, then call an authorized service representative. See Calling for Service.**
- **If an error message includes: “Items Not Sterile”, the items in sterilizer shall not be considered sterile and they must be run through a successful sterilization cycle before use.**
- **Items and packaging should be completely dry when removed from the sterilizer to minimize the potential for recontamination.**

# Error Codes - continued...

- Error Code User Interface Examples:**

**Error Status:**  
Red = High  
Yellow = Medium & Low

**Error Title**

 **Low Pressure Error**  
Feb 06, 2023, 09:12 AM

**Error Description** Chamber pressure too low during heating mode.

**Item Status** Instruments not sterile and need to be reprocessed.

**Error Code Number** Error Code 0747-033C

If you continue to receive this message, contact your service provider.

**Troubleshooting**

- Level Sterilizer
- Check Door Gasket
- Check Chamber Filters
- Remove power for 10s

*User Troubleshooting*

**Error Status:**  
Red = High  
Yellow = Medium & Low

**Error Title**

 **Door Sensor Error**  
Feb 06, 2023, 09:12 AM

**Error Description** Sensor not detecting door opening.

**Item Status** Risk of recontaminated instruments.

**Error Code Number** Error Code 0749-1568

If error persists or you need more information, contact your service provider.

**Troubleshooting**

- Check door is opening before drying
- Remove power for 10s if error persists

*User Troubleshooting*

Ok

*Acknowledge Error Code*

## Error Codes - continued...

- Error Status** When an error occurs, the error status is located at the top left of the user interface. It will inform the user of the severity of the error code encountered and if the load is sterile or not.

	<b>RED</b>	High severity error occurred requiring sterilizer power to be removed for a minimum of 10 sec. then reapplied.
		<ul style="list-style-type: none"> <li>When a high severity error message includes: "Instruments not sterile and need to be reprocessed," the items in sterilizer shall not be considered sterile and must be ran through a successful sterilization cycle before use.</li> </ul>
		<ul style="list-style-type: none"> <li>When a high severity error message includes: "Risk of recontaminated instruments". The items in sterilizer shall be considered at risk of recontamination, due to contents not being dry. Items and packaging should be completely dry when removed from the sterilizer to minimize the potential for recontamination.</li> </ul>
		<ul style="list-style-type: none"> <li>Light bar is illuminated red in addition to audible tone.</li> </ul>
	<b>Yellow</b>	Medium and Low severity errors do not require sterilizer power to be cycled off. The error code screen can be exited by pressing the OK button.
		<ul style="list-style-type: none"> <li>When a medium or low severity error message includes: "Instruments not sterile and need to be reprocessed," the items in sterilizer shall not be considered sterile and must be ran through a successful sterilization cycle before use.</li> </ul>
		<ul style="list-style-type: none"> <li>When a medium or low severity error message includes: "Risk of recontaminated instruments," the items in sterilizer shall be considered at risk of recontamination due to contents not being dry. Items and packaging should be completely dry when removed from the sterilizer to minimize the potential for recontamination.</li> </ul>
		<ul style="list-style-type: none"> <li>Light bar is illuminated yellow in addition to optional audible tone turned on or off in settings. Pressing the OK button stops the audible tone and exits the error code screen.</li> </ul>
	<b>Blue</b>	This is a notification indicating an error not involving the sterilization process, such as when exporting records, or updating firmware.

- User Troubleshooting** This information is shown on the right side of the user interface. The user can take these steps to potentially resolve the issue.

Some examples are:

- Check the reservoir water level.
- Ensure door is closed and latched.
- Restart the cycle.
- Level sterilizer.
- Remove power for 10 seconds.
- Check door gasket.

## Calling for Service

### Note

*Model and Serial Number information is required when calling for service. Record any displayed code(s) and relay this information to the service technician.*

If service is required, contact Midmark directly.

Medical US 844.856.1230 / Canada 937.526.8585

Dental US 844.856.1231 / Canada 937.526.8585

Animal Health US 844.856.1232 / Canada 937.526.8585

Monday - Thursday 8:00 AM to 6:00 PM

Friday 8:00 AM to 5:00 PM

[midmark.com](http://midmark.com)

# Specifications Chart: M9<sup>®</sup> Steam Sterilizer

## Physical Dimensions:

Overall Length with Plug	52.0 cm (20.5 in)
Overall Width	40.1 cm (15.8 in)
Overall Height	42.2 cm (16.6 in)
Support Surface	50.3 cm x 59.7 cm (19.8 in x 23.5 in)
Chamber	22.9 cm dia. x 38.1 cm deep (9 in x 15 in)
Standard Tray, Large	18.6 cm x 30.5 cm x 2.2 cm (7 5/16 in x 12 in x 7/8 in) (Usable tray area)
Standard Tray, Small	14.3 cm x 30.5 cm x 2.2 cm (5 5/8 in x 12 in. x 7/8 in) (Usable tray area)

## Weights:

Weight with Empty Reservoir	37.2 kg (82 lb)
Weight with Shipping Carton	41.7 kg (92 lb)

## Electrical Rating:

*Note: A separate (dedicated) circuit is recommended for this sterilizer. Sterilizer should not be connected into an electrical circuit with other appliances or equipment unless circuit is rated for the additional load.*

115 V models	115 VAC, 12 Amp, 50/60 Hz
--------------	---------------------------

## Fuse Ratings:

115 VAC	F1...15 Amp, 250 V, Fast Acting, 6.35 mm x 31.75 mm (0.250 in Dia x 1.250 in L)
	F2...15 Amp, 250 V, Fast Acting, 6.35 mm x 31.75 mm (0.250 in Dia x 1.250 in L)
	F3...2.5 Amp, 250 V, Time Delay/Slow-Blow, 5.20 mm x 20.00 mm (0.205 in Dia x 0.787 in L)

<b>Safety Valve Setting</b>	296 kPa (43 psi)
-----------------------------	------------------

<b>Heat Emission</b>	5275 kJ (5000 BTU) / hr during operation
----------------------	--

<b>Water Reservoir Capacity</b>	4.1 liters (1.1 gallons) to Full Mark. Usable volume is 2.3 liters (0.6 gallons)
---------------------------------	--

<b>Pollution Degree</b>	II (IEC 60664-1:2020)
-------------------------	-----------------------

<b>Overvoltage Category</b>	II (IEC 60664-1:2020)
-----------------------------	-----------------------

# Specifications Chart: M11<sup>®</sup> Steam Sterilizer

<b>Physical Dimensions:</b>	
Overall Length with Plug	59.7 cm (23.5 in)
Overall Width	45.2 cm (17.8 in)
Overall Height	47.2 cm (18.6 in)
Support Surface	55.4 cm x 67.3 cm (21.8 in x 26.5 in)
Chamber	27.9 cm dia. x 45.7 cm deep (11 in x 18 in)
Standard Tray, Large (Usable tray area)	23 cm x 38.7 cm x 2.2 cm (9-1/8 in x 15-1/4 in x 7/8 in)
Standard Tray, Small (Usable tray area)	17 cm x 38.7 cm x 2.2 cm (6-3/4 in x 15-1/4 in x 7/8 in)
<b>Weights:</b>	
Weight with Empty Reservoir	44.9 kg (99 lb)
Weight with Shipping Carton	59.4 kg (131 lb)
<b>Electrical Rating:</b>	
<i>Note: A separate (dedicated) circuit is recommended for this sterilizer. Sterilizer should not be connected into an electrical circuit with other appliances or equipment unless circuit is rated for the additional load.</i>	
115 V models	115 VAC, 12 Amp, 50/60 Hz
230 V models	230 VAC, 6.4 Amp, 50/60 Hz
<b>Fuse Ratings:</b>	
115 VAC	F1...15 Amp, 250 V, Fast Acting, 6.35 mm x 31.75 mm (0.250 in Dia x 1.250 in L)
	F2...15 Amp, 250 V, Fast Acting, 6.35 mm x 31.75 mm (0.250 in Dia x 1.250 in L)
	F3...2.5 Amp, 250 V, Time Delay/Slow-Blow, 5.20 mm x 20.00 mm (0.205 in Dia x 0.787 in L)
230 VAC	F1...8 Amp, 250 V, Fast Acting, 5.20 mm x 20.00 mm (0.205 in Dia x 0.787 in L)
	F2...8 Amp, 250 V, Fast Acting, 5.20 mm x 20.00 mm (0.205 in Dia x 0.787 in L)
	F3...2.5 Amp, 250 V, Time Delay/Slow-Blow, 6.35 mm x 31.75 mm (0.250 in Dia x 1.250 in L)
<b>Safety Valve Setting</b>	296 kPa (43 psi)
<b>Heat Emission</b>	<b>5275 kJ (5000 BTU)</b> / hr during operation
<b>Water Reservoir Capacity</b>	5.3 liters (1.4 gallons) to Full Mark. Usable volume is 3.0 liters (0.8 gallons)
<b>Pollution Degree</b>	II (IEC 60664-1:2020)
<b>Overvoltage Category</b>	II (IEC 60664-1:2020)

## Certification + Standards Compliance

<b>Certifications</b>	<ul style="list-style-type: none"> <li>UL 61010-1:2023</li> <li>CAN/CSA-C22.2 NO. 61010-1-12:2019</li> <li>UL 61010-2-040:2021</li> <li>CSA C22.2 No. 61010-2-040:21</li> <li>ASME Boiler and Pressure Vessel Code, Section VIII, Division 1</li> <li>Pressure Vessel Canadian Registration</li> <li>ISO 13485:2016</li> <li>ISO 9001:2015</li> </ul>
<b>Standards Compliance</b>	<ul style="list-style-type: none"> <li>ANSI/AAMI ST55:2016</li> <li>IEC 61010-1:2017 Edition 3.1</li> <li>IEC 61010-2-040:2020 Edition 3.0</li> <li>IEC 60601-1-2:2020 Edition 4.1</li> <li>ANSI/AAMI/IEC 60601-1-2:2021</li> <li>CAN/CSA C22.2 NO. 60601-1-2:2022</li> <li>IEC62304:2015 Edition 1.1</li> <li>ISO14971:2019 Edition 3.0</li> </ul>

## Electromagnetic Disturbances Guidance

Standard	Test	Level
<b>Emissions</b>		
CISPR 11 edition 6.2 (2019)	Conducted and Radiated RF	Group 1, Class A
<b>Immunity</b>		
IEC 61000-4-2 edition 2.0 (2008)	Electrostatic Discharge	± 8 kV Contact; ± 15kV Air
IEC 61000-4-3 edition 3.2 (2010)	Radiated RF EM Fields	3 V/m 80 MHz – 2.7 GHz, 80% AM at 1 kHz
IEC 61000-4-3 edition 3.2 (2010)	Proximity Fields from RF Wireless Communications Equipment	IEC 60601-1-2, Table 9
IEC 61000-4-4 edition 3.0 (2012)	Electrical Fast Transients	±2 kV at 100kHz repetition frequency
IEC 61000-4-5 edition 3.1 (2017)	Surges	Line-to-Line: ±1 kV
IEC 61000-4-5 edition 3.1 (2017)	Surges	Line-to-Ground: ±2 kV
IEC 61000-4-6 edition 4.0 (2013)	Conducted RF	3V: 0.15MHz – 80MHz 6V in ISM Bands between 0.15MHz – 80MHz
IEC 61000-4-8 edition 2.0 (2009)	Rated Power Frequency Magnetic Fields	30 A/m
IEC 61000-4-11 edition 2.1 (2017)	Voltage Dips	0% U <sub>T</sub> ; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315° 0% U <sub>T</sub> ; 1 cycle
IEC 61000-4-11 edition 2.1 (2017)	Voltage Dips	70% U <sub>T</sub> ; 30 cycles Single phase: at 0°
IEC 61000-4-11 edition 2.1 (2017)	Voltage Interruptions	0% U <sub>T</sub> ; 300 cycles
IEC 61000-4-39 edition 1.0 (2017)	Proximity Magnetic Fields	IEC 60601-1-2, Table 11

### Note

*The emissions characteristics of this equipment make it suitable for use in industrial areas and hospitals (CISPR 11 class A). If it is used in a residential environment (for which CISPR11 class B is normally required) this equipment might not offer adequate protection to radio-frequency communication services. The user might need to take mitigation measures, such as relocating or re-orienting the equipment.*

## Connectivity Options

Connectivity Options Chart		
<b>Network Protocols</b>		Link Layer: IEEE 802.3, IEEE 802.11 Internet Layer: IP Transport: TCP, UDP App Layer: TLS, MQTT, NTP, HTTP, etc.
<b>Network Data</b>	LAN Physical Connection	RJ-45
	LAN Cabling	Ethernet (min. CAT 5)
	LAN Max Speeds	10/100 Mbps
	LAN Protocols	IEEE 802.3u
	Wi-Fi Connection	The device supports connectivity to Wi-Fi using the following standards: <ul style="list-style-type: none"> <li>- 802.11a</li> <li>- 802.11b</li> <li>- 802.11g</li> <li>- 802.11n</li> <li>- 802.11ac</li> </ul>
	Operational Frequency	Connects to a 2.4 GHz or 5.0 GHz
	Wi-Fi Security Standards	When connecting to Wi-Fi, the device uses the following security standards: <ul style="list-style-type: none"> <li>- WPA2-PSK</li> <li>- WPA2 PEAP/MSCHAPv2</li> <li>- WPA2 EAP-TLS</li> </ul>
	Network Interface Card Parameters	Not Applicable

## Water Purity Specifications

Water Quality Specification Chart		
AAMI ST108 (ref.)		
	Units	Specification
Hardness	mg/L	< 1
Conductivity	µS/cm	< 10
pH @ 25C	pH	5 - 7.5
Chlorides	mg/L	< 1
Bacteria	CFU/mL	< 10
Endotoxin	EU/mL	< 10

## Disposal of Equipment

Consult local codes and ordinances for proper disposal of equipment, packaging, accessories and consumable goods.

# Midmark Sterilizer Warranty

Scan the following QR code for access to the sterilizer warranty.

Scan the following QR code for access to the sterilizer warranty registration.

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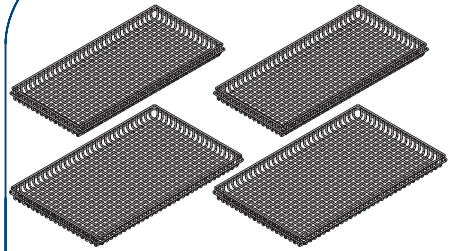
## Hoja SDS de Speed-Clean

Midmark Corporation		07641040 Rev P Doc Control: Jan 2007 Revised: February 2015 Page 1 of 6
SAFETY DATA SHEET		
<b>Section 1: Identification</b>		
1.1. Product Identifier Speed-Clean		
1.2. Recommended Use and Restrictions on Use Sterilizer Cleaner		
1.3. Supplier's Details Developed For: Midmark Corporation Address: 60 Vista Drive, P.O. Box 286, Versailles, OH 45380-0286, USA Telephone: 440-953-8000 1.4. Emergency Telephone Number: (937) 524-4500		
<b>Section 2: Hazard Identification</b>		
2.1. Hazard Classification Skin Irritation Category 2 Eye Irritation Category 2A		
2.2. Label Elements Signal Word Warning		
Symbols Exclamation Mark 		
Hazard Statements Irritating to respiratory skin contact causes skin irritation Contact with eyes causes eye irritation		
Precautionary Statements Prevention: Wear eye face protection. Wear appropriate gloves. Wash hands after handling this material.		
Response: IF IN EYES: Immediately flush eyes with water for 15 minutes, intermittently lifting the lower and upper eyelids. IF ON SKIN: Flush skin with running water for 15 minutes.		
2.3. Hazards Not Otherwise Classified None		

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Cycle Type	Temperature	Exposure Time	Pressure	Velocity	Pre-sterilization Pre-sterilization	min	min
Washed	121°C (250°F)	15	103 kPa (30 inHg)	1.5 m/s (5 ft/s)	Pre-sterilization: 10 min at 121°C (250°F) at 103 kPa (30 inHg) and 1.5 m/s (5 ft/s). Post-sterilization: 10 min at 121°C (250°F) at 103 kPa (30 inHg) and 1.5 m/s (5 ft/s).	10 min (121°C)	10 min (121°C)
Washed	132°C (270°F)	5	138 kPa (40 inHg)	1.5 m/s (5 ft/s)	Pre-sterilization: 10 min at 121°C (250°F) at 103 kPa (30 inHg) and 1.5 m/s (5 ft/s). Post-sterilization: 5 min at 132°C (270°F) at 138 kPa (40 inHg) and 1.5 m/s (5 ft/s).	5 min (132°C)	5 min (132°C)
Spore	121°C (250°F)	15	103 kPa (30 inHg)	1.5 m/s (5 ft/s)	Pre-sterilization: 10 min at 121°C (250°F) at 103 kPa (30 inHg) and 1.5 m/s (5 ft/s). Post-sterilization: 15 min at 121°C (250°F) at 103 kPa (30 inHg) and 1.5 m/s (5 ft/s).	15 min (121°C)	15 min (121°C)
Unwrapped	132°C (270°F)	3	138 kPa (40 inHg)	1.5 m/s (5 ft/s)	Pre-sterilization: 10 min at 121°C (250°F) at 103 kPa (30 inHg) and 1.5 m/s (5 ft/s). Post-sterilization: 3 min at 132°C (270°F) at 138 kPa (40 inHg) and 1.5 m/s (5 ft/s).	3 min (132°C)	3 min (132°C)

**Routine Care**

**Daily Routine Care:**

- Clean control panel and top and bottom covers.
- Wipe and clean gaskets with a damp cloth or Speed-Clean.
- Clean the top and bottom covers with Speed-Clean.

**Weekly Routine Care:**

- Clean the control panel and top and bottom covers.
- Clean the top and bottom covers with Speed-Clean.
- Clean the top and bottom covers with Speed-Clean.


**Monthly Routine Care:**

- Clean the control panel and top and bottom covers.
- Clean the top and bottom covers with Speed-Clean.
- Clean the top and bottom covers with Speed-Clean.


**Equipment Alert**

Use only distilled water or water that meets the water purity specified in the User Guide.

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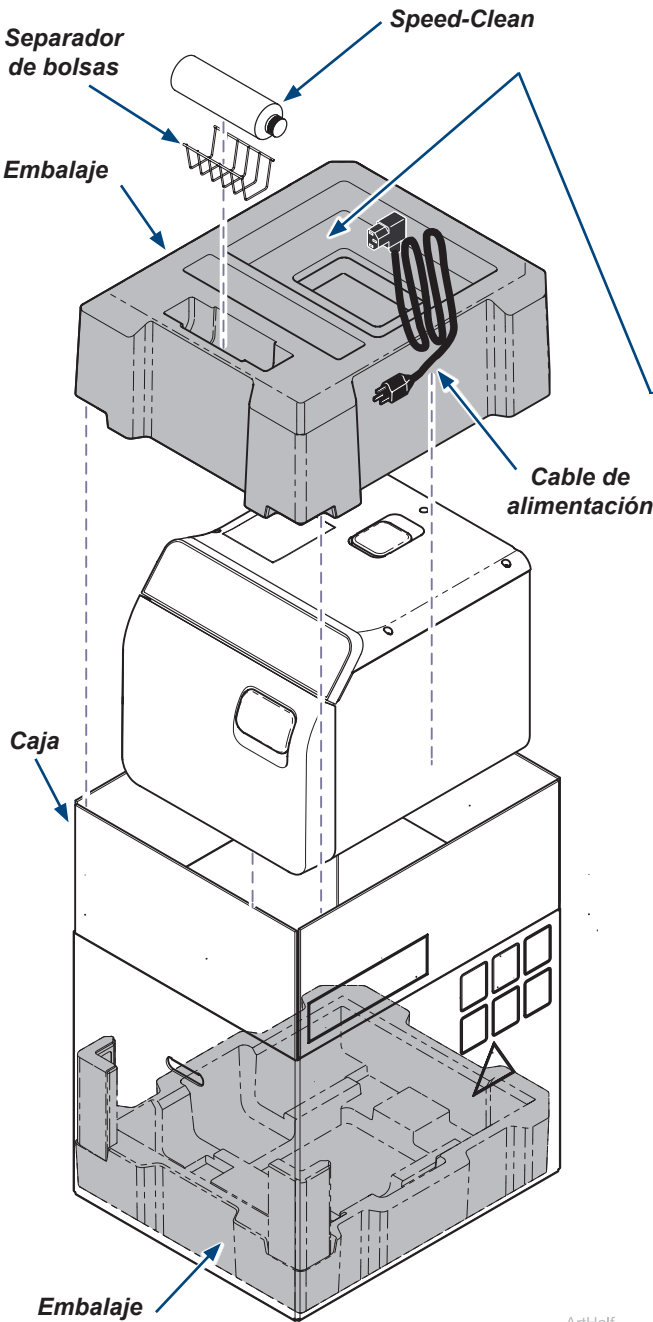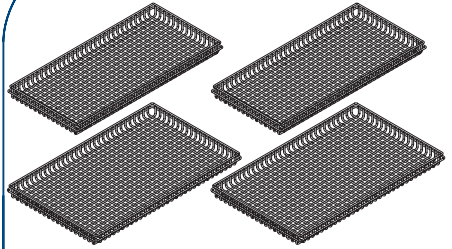## Hoja SDS de Speed-Clean

**Midmark Corporation**  
107481540 Rev. P  
Date Created: June 2007  
Revised February 2015  
Page 1 of 6

**SAFETY DATA SHEET**

**Section 1: Identification**  
1.1. Product Identifier  
Speed-Clean

1.2. Recommended Use and Restriction on Use  
Sterilization

1.3. Supplier's Details  
Developed For: Midmark Corporation  
Address: 600 Vista Drive, P.O. Box 286, Versailles, OH 45384-0286, USA  
Telephone: 1-400-MIDMARK

1.4. Emergency Telephone Number  
(937) 236-6100

**Section 2: Hazard Identification**  
2.1. Hazard Classification  
Skin Irritant Category 2,  
Eye Irritant Category 2A

2.2. Label Elements  
Signal Word  
Warning

**Stichwort**  
Exclamation Mark

**Hazard Statements**  
P303+P361+P353: Washed or repeated skin contact causes skin irritation.  
P305+P351+P338: Contact with eyes causes eye irritation.

**Precautionary Statements**  
P201: Wear protective gloves.  
P202: Wear protective eyewear.  
P273: Avoid release into the environment.

**Response**  
If IN EYES: Immediately flush eyes with water for 15 minutes, intermittently lifting the lower and upper eyelids.  
If ON SKIN: Flush skin with running water for 15 minutes.

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**Routine Care**  
Daily Routine Care:  
• Clean exterior surfaces with soap and water.  
• Wipe control panel with a damp cloth or Speed-Clean.  
• Clean the water reservoir with distilled water.

**Weekly Routine Care:**  
• Clean Speed-Clean trays with distilled water and Speed-Clean.  
• Clean the water reservoir with distilled water.  
• Clean the exterior of Speed-Clean with distilled water.  
• Clean the interior of Speed-Clean with distilled water.

**Monthly Routine Care:**  
• Clean Speed-Clean trays with distilled water and Speed-Clean.  
• Clean the water reservoir with distilled water.  
• Clean the exterior of Speed-Clean with distilled water.  
• Clean the interior of Speed-Clean with distilled water.

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
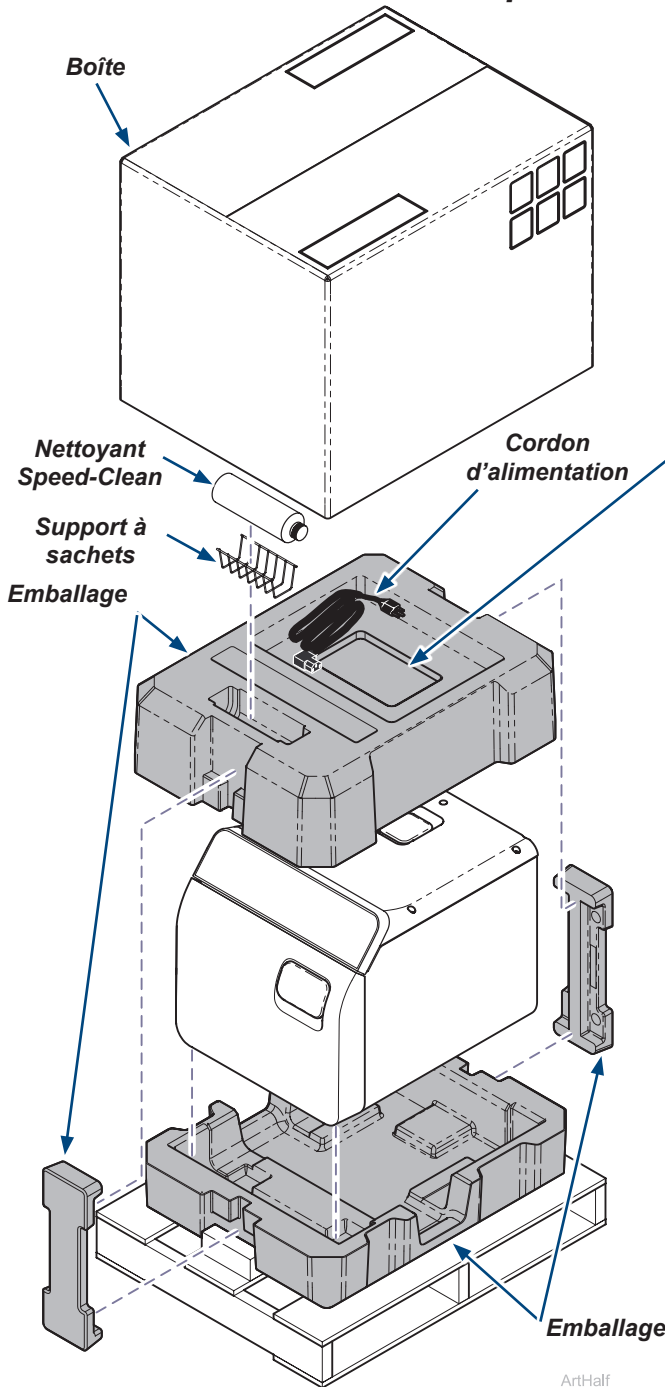
2 petits plateaux + 2 grands plateaux

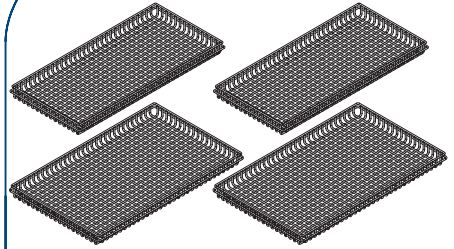FDS du nettoyant Speed-Clean

Midmark Corporation	
072401-040 Rev. P	Doc. Created: June 2007
	Revised: February 2015
	Page 1 of 4
SAFETY DATA SHEET	
Section 1: Identification	
1.1. Product Identifier Speed-Clean	
1.2. Recommended Use and Restrictions on Use Sterilizer Cleaner	
1.3. Supplier Details Midmark Corporation 40 West Dixie, P.O. Box 286, Versailles, OH 45383, USA Telephone: 614-948-9553 L6W4M1DMARK	
1.4. Emergency Telephone Number (977)258480	
Section 2: Hazard Identification	
2.1. Hazard Classification Skin Irritation Category 2 Eye Irritation Category 2A	
2.2. Label Elements Signal Word Warning	
Symbol Corrosion Mark	
Hazard Statements P303+P361+P531: Irritation or repeated skin contact causes skin irritation C314: Contact with eyes causes eye irritation	
Precautionary Statements P201: Hazardous P202: Hazardous P273: Avoid release to the environment P280: Wear eye protection P305+P351+P338: In case of contact with eyes, rinse immediately with plenty of water for several minutes. Remove contact lenses if you wear them and continue rinsing. P501: Dispose of contents and container as instructed.	
Exposure P273: Avoid release to the environment P280: Wear eye protection P305+P351+P338: In case of contact with eyes, rinse immediately with plenty of water for several minutes. Remove contact lenses if you wear them and continue rinsing. P501: Dispose of contents and container as instructed.	
2.3. Hazards Not Otherwise Classified None.	

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Serial Number  
QR Code  
Scan here first

Follow the QR code to the User Guide.

Example: 614-948-9553 ext. 6000  
Example: 614-948-9553 ext. 6000

Cycle Type	Temperature	Pressure	Time	Wiping	Wiping Solution	Wiping Frequency	Wiping Area	Wiping Method
Washed	121°C (250°F)	15 psi (103 kPa)	15 min	Wipe	70% isopropyl alcohol or 70% ethanol	After each cycle	Interior surfaces of the chamber and door	Use a clean, lint-free cloth. Wipe in a circular motion.
Washed	121°C (250°F)	15 psi (103 kPa)	15 min	Wipe	70% isopropyl alcohol or 70% ethanol	After each cycle	Interior surfaces of the chamber and door	Use a clean, lint-free cloth. Wipe in a circular motion.
Washed	121°C (250°F)	15 psi (103 kPa)	15 min	Wipe	70% isopropyl alcohol or 70% ethanol	After each cycle	Interior surfaces of the chamber and door	Use a clean, lint-free cloth. Wipe in a circular motion.
Washed	121°C (250°F)	15 psi (103 kPa)	15 min	Wipe	70% isopropyl alcohol or 70% ethanol	After each cycle	Interior surfaces of the chamber and door	Use a clean, lint-free cloth. Wipe in a circular motion.

**Routine Care**  
To prevent injury, allow unit to cool before cleaning gaskets and interior surfaces or removing trays. Dry rack, tray plate and gasket.


**Daily Routine Care:**  
• Clean control panel with soap and water.  
• Wipe and clean gasket with a damp cloth or Speed-Clean.  
• Clean tray and rack with a damp cloth or Speed-Clean.

**Weekly Routine Care:**  
• Clean control panel with soap and water.  
• Wipe and clean gasket with a damp cloth or Speed-Clean.  
• Clean tray and rack with a damp cloth or Speed-Clean.  
• Clean inside of chamber using Speed-Clean and flush with distilled water.  
• Clean tray plate with a damp cloth.


**Monthly Routine Care:**  
• Clean control panel with soap and water.  
• Wipe and clean gasket with a damp cloth or Speed-Clean.  
• Clean tray and rack with a damp cloth or Speed-Clean.  
• Clean inside of chamber using Speed-Clean and flush with distilled water.  
• Clean tray plate with a damp cloth.

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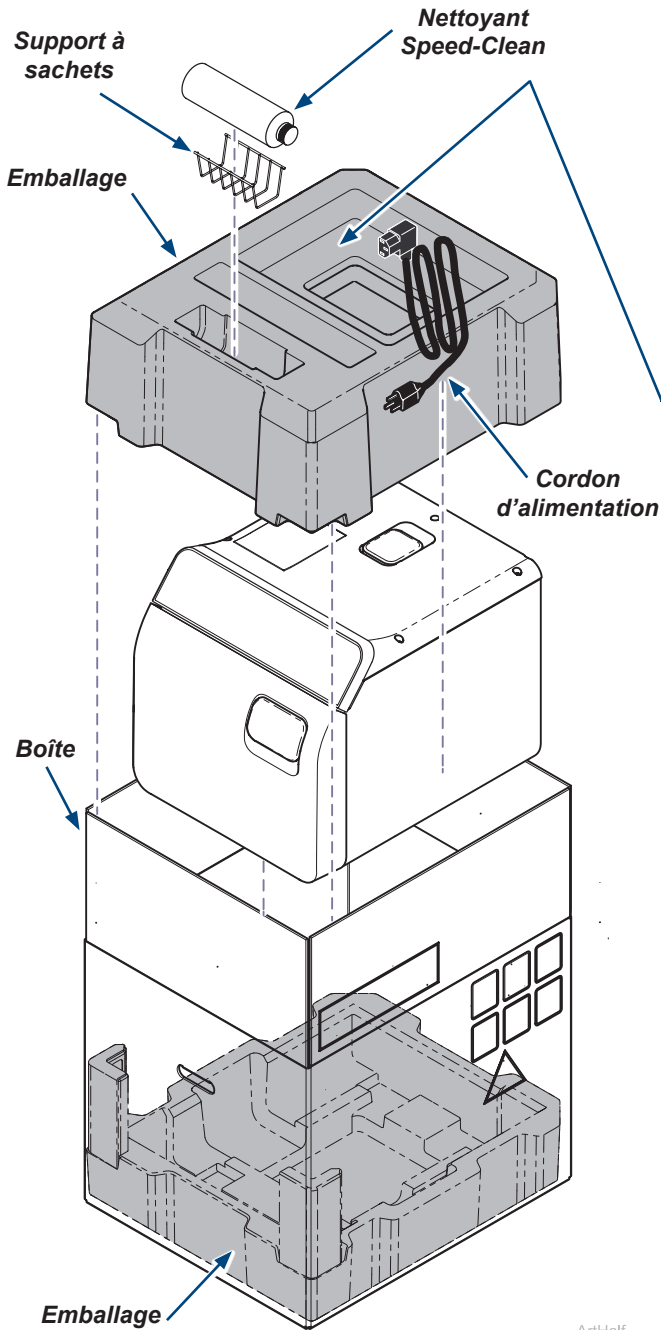2 petits plateaux + 2 grands plateaux
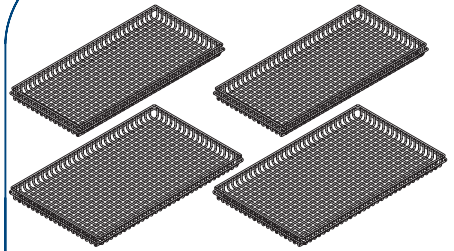
## FDS du nettoyant Speed-Clean

Midmark Corporation  
07240 L40 Rev. P  
Dated Original Issue 2007  
Revised: February 2015  
Page 1 of 6

**SAFETY DATA SHEET**

**Section 1: Identification**  
1.1. Product Identifier  
Speed-Clean  
1.2. Recommended Use and Restrictions on Use  
Sterilizer Cleaner  
1.3. Supplier Details  
Midmark Corporation  
40 Wood Drive, P.O. Box 286, Versailles, OH 45583, USA  
Address: Telephone: 1-800-MIDMARK  
1.4. Emergency Telephone Number  
(937) 256-6600

**Section 2: Hazard Identification**  
2.1. Hazard Classification  
Skin Irritation Category 2  
Eye Irritation Category 2A  
2.2. Label Elements  
Signal Word  
Warning  
2.3. Symbol  
Corrosion Mark  
GHS07  
GHS07

**Hazard Statements**  
P261 or repeated skin contact causes skin irritation  
Contact with eyes causes eye irritation

**Precautionary Statements**  
P201  
P273  
P301+P312  
P302+P352  
P303+P361+P353  
P305+P351+P338  
P501  
Wear eye protection.  
Wear protective gloves.  
Wash hands after handling this material.

**Response**  
IF IN EYES: Immediately flush eyes with water for 15 minutes, intermittently lifting the lower and upper eyelids.  
IF ON SKIN: Flush skin with running water for 15 minutes.

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midmark M9<sup>®</sup> and M11<sup>®</sup> Steam Sterilizer

For Models:  
M9 (450-0505)  
M11 (450-051-050)

Scan the QR code to register your sterilizer.  
Scan the QR code to the User Guide.

Cycle Type	Temperature	Time	Pressure	Notes	min. Cycle	min. Cycle
Wetpack	121°C (250°F)	30 min	150 psig (10.3 bar)	Preheated and/or pre-sterilized contents are required to prevent thermal shock. Check temperature constant or overtemp. Do not use for 100% ethanol.	120 min (M9) / 120 min (M11)	120 min (M9) / 120 min (M11)
Wetpack	134°C (273°F)	3 min	150 psig (10.3 bar)	Preheated and/or pre-sterilized contents are required to prevent thermal shock. Do not use for 100% ethanol.	120 min (M9) / 120 min (M11)	120 min (M9) / 120 min (M11)
Wetpack	134°C (273°F)	3 min	150 psig (10.3 bar)	Preheated and/or pre-sterilized contents are required to prevent thermal shock. Do not use for 100% ethanol.	120 min (M9) / 120 min (M11)	120 min (M9) / 120 min (M11)
Wetpack	134°C (273°F)	3 min	150 psig (10.3 bar)	Preheated and/or pre-sterilized contents are required to prevent thermal shock. Do not use for 100% ethanol.	120 min (M9) / 120 min (M11)	120 min (M9) / 120 min (M11)

**Routine Care**  
Daily Routine Care:  
• Check water level and add water as needed.  
• Check and clean gaskets with a damp cloth or Speed-Clean.  
• Clean the door and handle with Speed-Clean.

**Weekly Routine Care**  
• Clean the door and handle with Speed-Clean.  
• Clean the door and handle with Speed-Clean.  
• Clean the door and handle with Speed-Clean.

**Monthly Routine Care**  
• Clean the door and handle with Speed-Clean.  
• Clean the door and handle with Speed-Clean.  
• Clean the door and handle with Speed-Clean.

**Equipment Alert**  
• Check the water level and add water as needed.  
• Check the door and handle with Speed-Clean.

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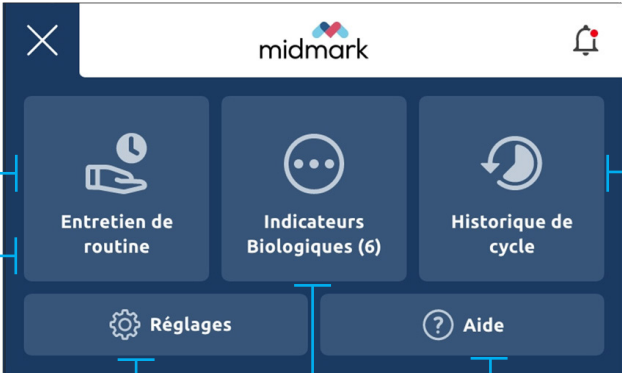
The screenshot shows the Midmark mobile app home screen. At the top, there is a close button (X), the Midmark logo, and a notification bell. Below this are three main menu items: 'Entretien de routine' (with a clock icon), 'Indicateurs Biologiques (6)' (with a speech bubble icon), and 'Historique de cycle' (with a circular arrow icon). At the bottom, there are two more menu items: 'Réglages' (with a gear icon) and 'Aide' (with a question mark icon). Callouts point to each of these items with descriptive text.

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