# STATIM 2000/2000S/2000 G4

CASSETTE AUTOCLAVE™

### Service Manual





Table of Contents 96-108639 Rev 3.0

# STATIM 2000 / 2000S / 2000 G4 Service Manual Table of Contents

About this Service Guide	5
Safety compliance	5
1.Identifying STATIM 2000 Units	7
Identifying software versions	7
Identifying controller board types	
Identifying steam generator types	
2. Tools, Maintenance Schedules,	
Procedures and Testing	
Tools	
The control box	
Calibration cassettes and reference meters	
Calibration cassettes	
Reference meters	
Annual service requirements	31
Routine maintenance schedule and procedures	
Draining the reservoir	
Cleaning the cassette	
Cleaning the water reservoir filter	
Cleaning the water reservoir	
Filters on the STATIM 2000/2000S	36
Bacteria retentive filter	37
Compressor filter	
Filters on the STATIM 5000/5000S	37
Bacteria-retentive filter	37
Compressor filter	38
Replacing the cassette seal	39
Cleaning the cover	40
Shipping the unit	
Upgrading the firmware on STATIM G4 units	41
Using the STATIM G4 remote access function	42
Electrical safety testing	44
Dielectric strength test (hi-pot)	44
Protective bonding impedance test	
(ground continuity)	
Water conductivity testing	45
3. Diagnostics and Troubleshooting Cycle Faults	47
Basic mechanical diagnostics	
Using software as a diagnostic tool	
For units with 2.x/5.x, 6.x controllers	
For units with 7.x controllers (with service menu). Using the service menu on the STATIM G4	48
(Rev. 7 controller board)	64

Troubleshooting cycle faults6	35
Cycle Fault Numbers and descriptions	
for revision 3.x and 4.x controllers	36
Cycle Fault Numbers and descriptions on units	
with revision 2.x/5.x/6.x and 7.x type controller	
boards (1.xx/2.xx/4.xx/5.xx and 6.xx software	37
boards (1.xx2.xx4.xx3.xx and 0.xx software	,,
4. STATIM 2000 Calibration	31
Calibration procedure 1: for Rev. 3.x or 4.x boards 8	31
Calibration procedure 2: for Rev. 2.x or 5.x boards,	
software R1.xx or R2.xx	35
Calibration procedure 3: for Rev. 6.X board,	
software Rev. 4.xx and 5.xx	39
Calibration procedure 4: for Rev. 7.x board,	
software Rev. 6.xx	93
Calibration Procedure 8: for non-S G4 units with	
Rev. 7.x Board, Software Rev. 7.xx, SL 00R1.xx	97
, , , , , , , , , , , , , , , , , , , ,	
5. STATIM 2000S Calibration10	)3
Calibration procedure 5: for Rev. 2.x or 5.x	
boards,software R1.xx and R2.xx10	)3
Calibration procedure 6: for Rev. 6.x boards,	
software R4.xx and R5.xx10	)8
Calibration procedure 7: for Rev. 7.x boards,	
software R6.xx11	13
Calibration Procedure 9: for S-class G4 Units with	
Rev. 7.x Board, Software Rev. 7.xx, SL 00R1.xx 11	18
6. Cover Assembly	23
	20
Removing the cover	
Reinstalling the cover	
Removing and replacing the armature gasket	
Removing and replacing the fascia	
Removing and replacing the membrane keypad 12	
Removing and replacing the LCD12	29
For STATIM 2000 G4	
Removing the cover13	31
Reinstalling the cover13	
Removing and replacing the armature gasket13	
Removing and replacing the LCD overlay	
Removing and replacing the LCD touchscreen 13	
Removing and replacing the Lob touchscreen	
Tiomoving and replacing the lastia	,0

Table of Contents 96-108639 Rev 3.0

# STATIM 2000 / 2000S / 2000 G4 Service Manual Table of Contents

7. Electrical and Electronic Components	141
Controller board	141
Identifying the controller board type	141
Removing the controller board assembly	143
Replacing/refitting the controller board assembly	146
Controller board fuse values	148
STATIM 2000 G4 Electronics	149
Removing and replacing the STATIM 2000 G4	
LCD controller board	149
Removing and replacing the STATIM 2000 G4	
USB port	151
Removing and replacing the STATIM 2000 G4	
speaker	151
Removing and replacing the STATIM 2000 G4	
Ethernet port	
Microprocessor	
Processor types	
Removing the microprocessor	
Replacing/refitting the microprocessor	
Microswitch	
Removing the microswitch	
Replacing the microswitch	
Adjusting the microswitch	
Mains components	
Removing and replacing a fixed power cord	
Removing a detachable power cord	
Replacing/refitting a detachable power cord	
Power (on/off) switch	
Removing the power switch	
Replacing/refitting the A.C. Power switch	
Receptacle / line filter	
Removing the receptacle / line filter	
Replacing/refitting the receptacle / line filter LCD DC Power Source for STATIM 2000 G4	
Removing and replacing the STATIM 2000 G4	100
LCD DC power source	165
LOD DO power source	105
8. Pumps, Reservoir and Compressor	167
Water pumps	
Identifying pump types	
Testing pump flow	
Water pump recovery repair procedure	
Performing filter maintenance	
Performing internal pump filter maintenance	
Performing water reservoir filter maintenance	
Performing inline filter maintenance	
Performing pump tube replacement	
- · · · · · ·	

Removing the Fluid-O-Tech pump assembly	177
Testing the Fluid-O-Tech pump diode	178
Reinstalling the Fluid-O-Tech pump assembly	179
Removing the SciCan or Ulka pump assembly	
Reinstalling the SciCan or Ulka pump assembly	
Reservoir and water sensors	
Removing the reservoir	
Reinstalling the reservoir	
Troubleshooting	
Compressor	
Removing the compressor	
Reinstalling the compressor	
Hemstalling the compressor	13
9. Solenoid Valve and Pressure Transducer	193
Solenoid valve	193
Inspecting and repairing the	
type B solenoid valve	193
Removing the solenoid valve	
Refitting the solenoid valve	
Pressure transducer	
Identifying and checking the	
pressure transducer	197
Removing the pressure transducer	
Reinstalling the pressure transducer	
Removing the pressure interface board	
Reinstalling the pressure interface board	
Removing and replacing the	200
pressure interface board battery	205
10. Armature, Isoplate, and Probe Bracket	
and Chamber Thermocouple	
Armature and isoplate	
Removing armature and isoplate	
Reinstalling armature and isoplate	209
Removing the isoplate from the	
armature-isoplate	210
Replacing and reassembling the isoplate	
into the armature	211
Probe bracket and	
chamber thermocouple	212
Removing the probe bracket and chamber	
thermocouple	212
Replacing and reinstalling the probe bracket	
and chamber thermocouple	215
Aligning the probe bracket and chamber	
thermocouple	218

Table of Contents 96-108639 Rev 3.0

# STATIM 2000 / 2000S / 2000 G4 Service Manual Table of Contents

11. Steam Generator, Check Valve, Thermal Fuse and
Pressure Relief Valve221
Steam generators221
Replacing 1st generation and
stainless steel steam generators (matrix) 222
Removing the aluminum steam generator223
Reinstalling the aluminum steam generator 224
Removing the stainless steel steam generator 226
Reinstalling the stainless steel steam generator 227
Check valves229
Testing the check valve229
Removing and replacing the check valve 230
Steam generator thermal fuses
Removing the thermal fuse
on the aluminum steam generator233
Replacing the thermal fuse
on the aluminum steam generator235
Removing and replacing the thermal fuse
on the stainless steel steam generator236
Pressure relief valve236
Checking the pressure relief valve237
Removing and replacing the pressure relief valve 237
Removing and replacing the pressure relief valve237  12. Printer and Data Logger
12. Printer and Data Logger241
<b>12. Printer and Data Logger241</b> Printer
12. Printer and Data Logger
12. Printer and Data Logger
12. Printer and Data Logger
Printer and Data Logger
Printer and Data Logger
Printer and Data Logger
12. Printer and Data Logger
12. Printer and Data Logger
12. Printer and Data Logger
12. Printer and Data Logger
Printer and Data Logger
Printer and Data Logger
Printer and Data Logger
Printer and Data Logger

13. Installation	259
Environmental considerations2	259
Unit placement2	260
Connecting the waste bottle2	260
Filling the STATIM reservoir2	261
Priming the STATIM pump2	261
Powering up and setting up the unit2	261
Setting the time and date2	262
Language selection2	262
Assigning unit identifier number2	262
Setting/editing serial number2	263
Setting air drying time2	264
Connecting STATIM G4 to a network2	265
Setting up the G4 web portal2	266
Preparing unit for use2	271
Shipping the unit2	271
Unit installation checklist2	271
14. STATIM Operational Diagram	273

# **About this Service Guide**

This STATIM service manual was created to act as reference for the service and repair of the STATIM brand. In the interest of providing one comprehensive global service manual, this book references both the STATIM 2000 and the STATIM 2000S models. Because of this, you will note that there are chapters and sections that do not apply to the STATIM unit that you may have, or are repairing.

If you have a question about the unit you are repairing, please do not hesitate to contact your local SciCan representative for confirmation.



Hazardous voltages are accessible when the cover is removed. Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

If the cover is removed, a dielectric strength test (Hi-Pot) AND a protective bonding impedance test (ground continuity) must be performed on the STATIM when the work is completed and after the cover has been returned to the unit.



The STATIM is heavy. Exercise caution and seek assistance when lifting or carrying the units.



Use only steam-process distilled water in the STATIM.



The STAT*IM* contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.



Ensure that there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.

# Safety Compliance

When a STATIM is serviced, the safety criteria as specified by applicable international safety standards and applicable national, state, provincial, and regional laws and regulations must be observed and maintained.

The following symbols appear in the margins of this book.



A potential hazard to the operator.



A situation which may lead to a mechanical failure.



Important information

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



### **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



## PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# Identifying STATIM 2000 Units

(Does not include 'E' models)

There are two main types of STAT/M 2000 units, the G4 touchscreen version and Statim Classic with keypad. Each of these two versions have a North American specification or "non S" unit, and the European specification or "S" unit. In the Classic series, they are labeled respectively as "STAT/M 2000" and "STAT/M 2000S". The G4 series models do not use the "S" designation in their names, but the G4 series does include S-class models for distribution in markets requiring these cycle types, and non S-class for the North American market.

There are many mechanical identifiers to differentiate between the set types, (e.g., pressure transducer for 'S' class, no pressure transducer for the non 'S' class). The details in the matrix further below will help you identify units based on their factory specifications.


Note that this section does not cover STATIM 2000E models. For more information on these, contact SciCan Technical Service.

### **Identifying Software Versions**

To identify a unit's software version, power it OFF and back ON again. During start up, a number is briefly displayed in the upper right-hand corner of the LCD that indicates which software is installed (See Figure 1). This number is also printed on the label attached to the U14 microprocessor and the U8EEPROM device on the Controller Board (See Figure 2).

Note that units with Revision 3.x or 4.x controller boards WILL NOT display the software revision on the LCD display. These units will display "SELECT A CYCLE".

The software version number is the character to the immediate right of the R as indicated by:

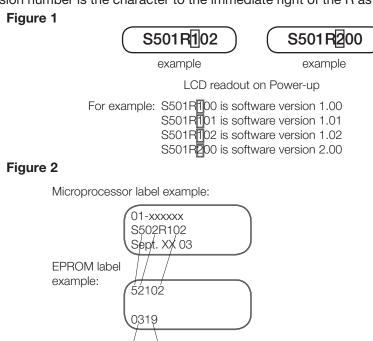

Product Year / week

### Software Versions and Corresponding Calibration Procedures

Model	Approx. Year of Manufacture	Software Revision	Steam generator Type	PCB Revision	Pump	Compressor	Pressure transducer	Validation thermocouple	Calibration procedure
1000	1992 - 1995	N/A	Aluminium Internal	3.x/4.x	<b>√</b>	N/A	N/A	N/A	1
2000	1992 - 1995	N/A	Aluminium Internal	3.x/4.x	<b>√</b>	<b>✓</b>	N/A	N/A	1
2000	1995 - 2005	S2O2R1xx S2O2R2xx	Aluminium Internal	2.x/5.x	<b>√</b>	✓	N/A	N/A	2
2000	2005 - 2007	S2O2R5xx	Aluminium External	6.x	<b>√</b>	<b>✓</b>	N/A	✓	3
2000	2007 →	S2O2R6xx	Aluminium External	7.x	<b>√</b>	<b>✓</b>	N/A	<b>√</b>	4
2000 G4	2010 →	S202R7xx	Aluminium External	7.x	<b>√</b>	✓	N/A	<b>√</b>	8
2000S	1995 - 2004	S2S2R1xx S2S2R2xx	Aluminium Internal	2.x/5.x	<b>√</b>	✓	√ (60 p.s.i.)	N/A	5
2000S	2004 - 2005	S2S2R4xx	Stainless Steel	6.x	<b>√</b>	✓	√ (68p.s.i.)	N/A	6
2000S	2005 - 2007	S2S2R5xx	Aluminium External	6.x	<b>√</b>	✓	√ (68p.s.i.)	N/A	6
2000S	2007 →	S2S2R6xx	Aluminium External	7.x	<b>√</b>	✓	√ (68p.s.i.)	N/A	7
2000S G4	2010 →	S2S2R7xx	Aluminium External	7.x	<b>√</b>	✓	√ (68p.s.i.)	N/A	9

## **Identifying Controller Board Types**

There are four different types of Controller Boards that may be encountered in STAT*IM* 2000/2000S/2000 G4 type units (See 'PCB revision' above). They can be identified as below (Figure 3):

**Figure 3** Revision 3.x/4.x type board (1992 – 1995):



Typical features:

- Revision number bottom right hand side.
- Single, rectangular microprocessor.

### Revision 2.x/5.x type board (1995 – 2004):



### Typical features:

- Revision number bottom right hand side.
- Rectangular microprocessor plus EPROM
- 'W1' jumper for calibration
- Blue 'pressure interface/printer' connector



Revision 6.x type board (2004 - 2007):



Revision 7.x type board (2007 to present):

### Typical features:

- Revision number bottom right hand side.
- Square microprocessor plus EPROM
- 'W1' jumper for calibration
- Blue 'pressure interface/printer' connector

### Typical features:

- Up to revision 7.30, the revision number is on the top right hand side (printed vertically)
- From revision 7.4 onwards, the revision number is on the bottom left hand side (printed horizontally) under the connector J1
- Square microprocessor plus EPROM
- NO 'W1' jumper for calibration
- NO Blue 'pressure interface/printer' connector
- All components integrated on single board
- Surface mount type component.
- 'Push in' yellow thermocouple connectors.

### **Identifying Steam Generator Types**

There are three different types of steam generator (boiler) that may be encountered in STAT*IM* 2000/2000S units. They can be identified as below (Figure 2):

Aluminum steam generator with internal thermocouple (1992 – 2004):





Typical features:

- 70 p.s.i. pressure relief valve (PRV) (no ring pull)
- Thermocouple embedded in side of steam generator
- Software required to drive this steam generator = R1xx/R2xx (revision 3.x/4.x/2.x/5.x PCB)

Stainless steel steam generator (2004 – 2005):



Typical features:

 Software required to drive this steam generator = R4xx (revision 6 PCB only) Aluminum steam generator with external thermocouple, a.k.a. 'ALEX' boiler (2005 – present):





Typical features:

- 43.5 p.s.i. pressure relief valve (PRV) (ring pull on S models, no ring pull on non S models)
- Thermocouple clamped to underside of steam generator
- Blanking plug in the side of steam generator in original thermocouple position
- Software required to drive this steam generator is R5xx for revision 6 PCB and R6xx for revision 7 PCB

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



### **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



## PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# Tools, Maintenance Schedules, Procedures and Testing

### **Tools**

### STATIM specific items from SciCan

Before a STATIM can be serviced, the following special tools are required in addition to your service tool kit. These tools are available from SciCan or your nearest service depot:

1.	Control box	SciCan Part # 01-103141S
2.	Pump tester	SciCan Part # 01-100713S
3.	Water conductivity tester	SciCan Part # 01-103139S
4.	Solenoid plunger tube wrench	SciCan Part # 01-103471S
5.	9/64" ball-end allen-key (hex) screwdriver	SciCan Part # 01103469S
6.	Calibration cassette,	
	SciCan 2000 Non S models only	SciCan Part # 01-103087S
7.	Calibration cassette,	

### Generic reference devices for calibration

- 8. Calibrated digital thermometer with 'K' type probe
- 9. Calibrated digital voltmeter with accuracy of 0.001 Volts
- 10. Calibrated digital pressure meter, 0 400 kPa absolute/0 7 bar absolute

More details on recommended devices can be found in section 'Recommended Reference Meters'

SciCan Part # 01-106366S

### Electrical safety test devices

SciCan 2000 All models

- 11. Hi-Pot tester
- 12. Ground continuity tester
- 13. Static strap
- 14. Static bags


### General tool list

- Phillips screwdriver
- Flat-blade electrician's screwdriver
- Potentiometer trimmer
- Needle nose pliers
- Wire cutters
- Wrench 3/8" A/F
- Wrench 7/16" A/F
- Wrench 9/16" A/F
- Wrench 11/16" A/F
- Wrench 7/8" A/F

**Note:** This tool list is a guide and suggests what is the minimum required to remove and replace components in the STAT*IM*. Other tools may be required as an aid to servicing at the discretion of the individual service engineer.


### Sample equipment set up for calibration of a non-S unit

- Calibration cassette with detachable thermocouple
- 2. Digital voltmeter
- 3. Digital thermometer
- 4. (no pressure meter required)
- 5. Control box
- 2x7 ribbon cable, LCD connector (connect to Controller Board header P3)
- 1x7 cable, keypad connector (connect to Controller Board header P4)
- 2x10 molded socket, test connector (connect to Controller Board header P1)

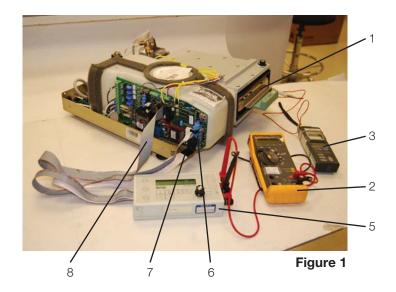

### Sample equipment set up for calibration of S unit

- 1. Calibration cassette.
- 2. Digital voltmeter
- 3. Digital reference thermometer with thermocouple. (Fluke 51 with 80PK-26 probe shown)
- 4. Digital reference pressure meter (S only) (Druck DPI 750 R shown)
- 5. Control box
- 6. LCD connector cable (connect to Controller Board header P3)
- 7. Keypad connector cable (connect to Controller Board header P4)
- 8. Test connector cable (connect to Controller Board leader P1)

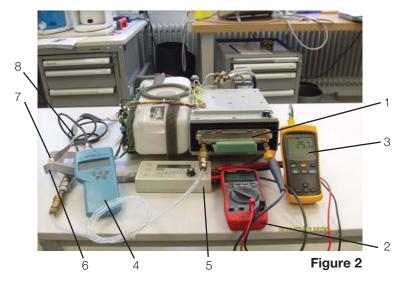


NOTE: The following section describes in greater detail the equipment depicted in Figures 1 and 2.

### The Control Box

NOTE: Late model STAT*IM*s with Rev.7 boards can be calibrated with the cover on using the unit's software and keypad, making use of the Control Box as an option. For older models, use of a control box is required.

The Control Box (5) is a service tool that allows a technician to operate the unit while it is being serviced. There are two kinds of Control Boxes you may encounter: one uses toggle switches, the other uses a membrane keypad arrangement similar to the keypad on the unit.

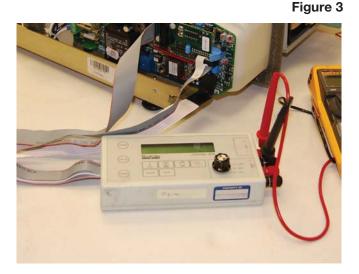

When connecting the Control Box connectors to the Controller Board note the positions of Pin number 1 of the Control Box test connectors and Pin number 1 of the Controller Board headers.

The Control Box has a display (LCD) and a variety of control switches that allow the technician to:

- 1. Operate the unit independently of the keypad mechanism, to manually activate the pump, valve or compressor, or run cycles.
- 2. Make frequently used Controller Board measurements by providing a common output and rotary switch for signal selection during calibration. See the Calibration Instructions in chapters 4 and 5 of this service manual.

If only the display and button functions are being used, connect the LCD cable (6) to Controller Board connector P3, the keypad cable (7) to Controller Board connector P4 and power the STAT*IM* ON.

If the override or measurement functions are being used, connect the test connector cable (8) to connector P1 of the Controller Board in addition to the other cables.

When using a toggle switch model of the Control Box, turn the switch to the ON position to activate the desired device. To turn the device OFF, turn the switch to the AUTO position. To select, start and / or stop a cycle, activate the appropriate push-button switch on top of the Control Box.

Calibration uses the Select out +, Select out - jacks and a setting on the Rotary Switch to set the calibration operation desired for Revision 3.x/4.x Controller Boards (for STATIM 2000 units only). See Calibration Procedure 1 or chapter 4 of this service manual.

When using a keypad model of the Control Box, the keypad switch must be held down to turn the desired device ON. Vref measurements use test leads, the Vref + and Vref - jacks and a voltmeter hooked up in series. Calibration uses the Select out + and Select out - jacks, test leads and a setting on the Rotary Switch to set the calibration operation desired for Rev. 3.x/4.x Controller Boards (for STATIM 2000 units only). The keypad of this Control Box provides the same features as a STATIM keypad. Note: The Control Box is for use with STATIM products ONLY.

### Calibration Cassettes and Reference Meters

Specially designed calibration cassettes allow technicians to take readings from within the cassette while the unit is in operation. In addition, calibration also requires the use of certain reference meters. Use the table below to identify the meters you will require.

Sterilizer type	Part Number	Description	Meter types required
STAT/M 2000 (non S)/ 2000 G4	01-103087S	Calibration cassette STAT <i>IM</i> 2000	Temperature only
STATIM 5000 (non S)/ 5000 G4	01-103088S	Calibration cassette STAT <i>IM</i> 5000	Temperature only
STAT <i>IM</i> 2000S/ 2000S G4	01-103088S	Calibration Cassette 2000S B	Temperature and pressure
STAT <i>IM</i> 5000S/ 5000S G4	01-106367S	Calibration Cassette 5000S C	Temperature and pressure

### **Calibration Cassettes**

The cassette of the STATIM unit is effectively a fixed volume chamber that contains instruments during sterilization. Saturated steam at a specific temperature and pressure is introduced to produce the correct conditions for sterilization. The cassette is part of a sealed system designed to contain the pressurized steam. To calibrate the unit's



Figure 4

temperature and pressure (if appropriate) monitoring devices, a special calibration cassette must be used. The calibration cassette allows a technician to independently monitor the temperature and pressure inside the cassette, ensuring the unit is operating in accordance to its original specifications and to national/international standards.

The relevant cassette types are detailed in the chart above, and the following important information should be noted before choosing the appropriate cassette to use.

- Non S unit cassettes require a 1.6mm (1/16") temperature probe and cannot be used with a 3.2mm (1/8") temperature probe. (See meter details for probe diameter).
- S unit cassettes have temperature fittings (one fitted and one in the cassette accessories supplied) that will allow the use of both 1.6mm (1/16") and 3.2mm (1/8") temperature probe. (See meter details for probe diameter).
- S unit cassettes are designed for both temperature and pressure reference meters and are suitable for all models, S and Non S. If these cassettes are acquired, the non S, single port cassette is not required.
- The non S unit cassettes should not be used to calibrate S class units.
- Please note that when calibrating non-S units with an S unit cassette that the pressure meter will not be required and the test port will require the pressure tube supplied with the cassette to be connected to the pressure fitting to prevent steam leakage.

### Reference Meters

Calibration reference meters are important for accurately setting the STAT*IM* unit so that the correct sterilization conditions (temperature and pressure) occur in accordance with the original specifications of the unit and national/international standards.

When ordering any digital thermometer and temperature probe, ensure that the supplier is aware that the area where the most accuracy is required is between 130°C and 140°C.

Temperature meters and probes should always be calibrated as a matched pair.

Test equipment should be calibrated on a regular basis based on the manufacturer's recommended calibration interval.

Calibration of reference equipment used with autoclaves should ALWAYS be to national or international standards by a certified calibration laboratory.

A number of meters are recommended by SciCan and are referenced as follows:

### For non-S unit calibration

The following reference meters are recommended for use with 01-103087S and 01-103088S test cassettes:

- Omega HH81A single channel multifunction digital thermometer (www.omega.com)
- Omega TJ36-CASS-116G-6-SMP-M Temperature Probe (with SMP-M miniature male connector).

Note: the above probe is 1.6mm (1/16") diameter.

### For S unit calibration

The following options are recommended for use with 01-106366S and 01-106367S test cassettes:

- A set of two independent meters, one pressure and one temperature.
- An integrated hand-held pressure and temperature calibrator.

Note: The thermocouple entry fitting on the S class cassette will need to be changed from the 3.2mm (1/8") diameter fitting to the 1.6mm (1/16") diameter fitting if using the Omega meter.

The two independent meters have been included where the user may already have an Omega meter for non-S use who is upgrading to S class use, when only a pressure meter will be required.

They are also the lower cost options over the integrated device.

# Recommended separate temperature and pressure meters for S unit calibration temperature

- Fluke 51 Series II Digital Thermometer.
- Fluke 80PK-26 SureGrip Tapered Temperature Probe.

#### Notes:

- The flexible thermocouple included with the Fluke digital thermometer will not be used for calibrating SciCan sterilizers. If you wish to use the flexible thermocouple in the future, SciCan advises you to consult the digital thermometer's manual for setting a temporary offset, if required.
- The above probe is 3.2mm (1/8") diameter.
- Technicians already in possession of the Omega HH81A or equivalent, and the associated probe can use this meter as a replacement for the Fluke 51 and 80PK-26 probe. Note that the probes are of different diameter and a cassette compression fitting change may be appropriate.

### **Pressure**

 Druck DPI 705R Absolute Pressure Meter with external 0 - 7 bar absolute pressure transducer with ¼" NPT female thread.

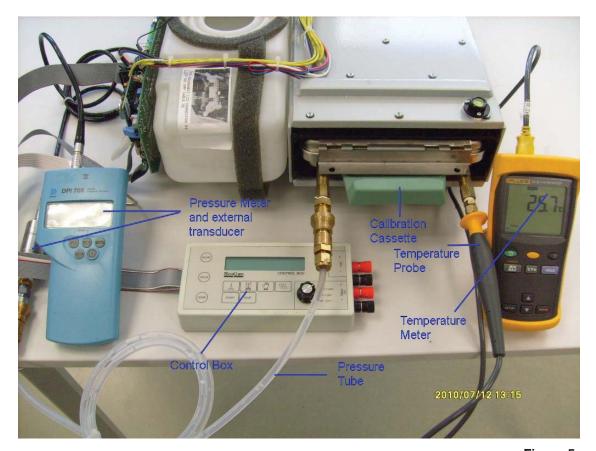


Figure 5

### **Combined Temperature and Pressure Meter**

- Heise PTE-1 Handheld LCD digital calibrator complete with temperature and pressure modules, and PT100 probe as follows:
  - Meter (without data logging capabilities) PTEC = X X 4 4A
  - Pressure Module HQS2 B A A 400 kPa A
  - Temperature Module HQS RT1 PT-100
  - RTD Probe Pt-100 PT5

### Notes:

- These product references are very specific, and the local supplier of Heise equipment (see www.heise.com) should be consulted prior to confirming order codes.
- The meter above is of the non data logging variety. If data logging is required, please consult order code variations from Heise data sheet.
- The above probe is 3.2mm (1/8in) diameter.



Figure 6

### Setting Up your Reference Meters

With certain meter types, some sub-assembly procedures will need to be undertaken to connect the device to your Statim calibration cassette. All cassettes are supplied with the necessary accessories to enable you to convert your meter of choice from those recommended above, to the appropriate fitting on the cassette.

01-103087S and 01-103088S cassettes (non S calibration) with Omega HH81A reference meter and TJ36-CASS-116G-6-SMP-M Temperature Probe.

No sub-assembly required.

To insert the temperature probe (thermocouple) into the compression fittings, loosen the clamp nut attached to the cassette extension fitting and insert the probe into the fitting as far as it will go. A slight resistance will be felt as the probe passes through the seal. Tighten the clamp nut until a steam/air tight seal is achieved. **DO NOT OVERTIGHTEN**.

NOTE: The non S calibration cassette fitting is NOT compatible with the Fluke and Heise 1/8" (3.175mm) temperature probes.

## 01-106366S and 01-106367S cassettes (S-class calibration) with ALL recommended reference meters

S-class calibration cassette kits contain a number of accessories for use with various recommended reference meters and will need to be set up to suit the appropriate meters before use.

Contents of calibration cassette kit:

- Cassette (2000S or 5000S as appropriate) with 'generic' pressure and temperature fittings attached (see Figure 7 below)
- · Accessories as follows: (inside the cassette)
  - 1 x Pressure tube
  - 2 x Male Swagelok B-QC4-S-2PM 'Quick Connect' fittings.
  - 1 x 1/8 NPT female to 1/4 NPT male adaptor
  - 1 x 1/8" (3.175mm) probe fitting.

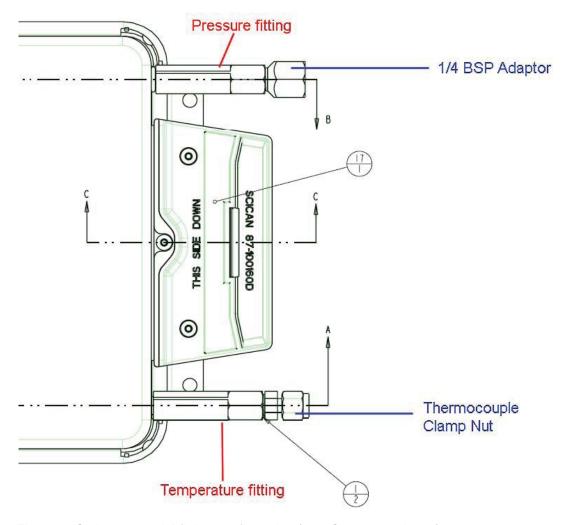


Figure 7. S-class 'generic' fitting configuration (2000S cassette shown)

# Setting up the temperature fittings for use with the recommended temperature reference meter(s)

The temperature fitting (right hand fitting when facing the front of the cassette) is fitted with a compression fitting that consists of a body, clamp nut and conical compression washer. (See Figure 8)

There are two fittings supplied with the cassette, one for a 1/8" (3.175mm) temperature probe and one for a 1/16" (1.6mm) temperature probe.

NOTE: The 1/16" fitting is attached to the cassette when supplied and the 1/8" fitting is in the pack of accessories supplied with the cassette.

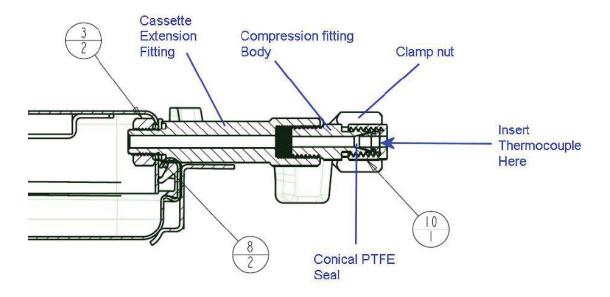


Figure 8. Temperature fitting cut-away

If you have the SciCan supplied **Omega HH81A** type temperature reference meter with the **1/16**" **(1.6mm) temperature probe** (normally used with the Non S-class type calibration cassette) then the compression fitting attached to the S class cassette is correct and does **NOT** have to be removed.

If you have either the **Fluke 51** type temperature reference meter, or **Heise PTE1** type combined temperature and pressure reference meter then the 1/16" (1.6mm) compression fitting (complete) will need to be removed from the extension piece attached to the cassette and replaced with the 1/8" (3.175mm) compression fitting from the accessory kit, as both of these meters have 1/8" (3.175mm) thermocouple probes.

NOTE: PTFE tape should be used when assembling these parts to ensure a steam tight seal.

To insert the relevant temperature probe (thermocouple) into either of the appropriate compression fittings, loosen the clamp nut and insert the probe into the fitting as far as it will go. A slight resistance will be felt as the probe passes through the conical seal. Tighten the clamp nut until a steam/air tight seal is achieved. **DO NOT OVERTIGHTEN**.

# Setting up the pressure fittings for use with the recommended pressure reference meter(s)

NOTE: The pressure fitting (left hand fitting when facing the front of the cassette) is fitted with an adaptor to attach ¼" BSP devices such as hypodermic pressure fittings. (See Figure 9 below). This adaptor is not used with the standard SciCan pressure tube and should be removed from the cassette extension fitting before attaching the pressure tube.

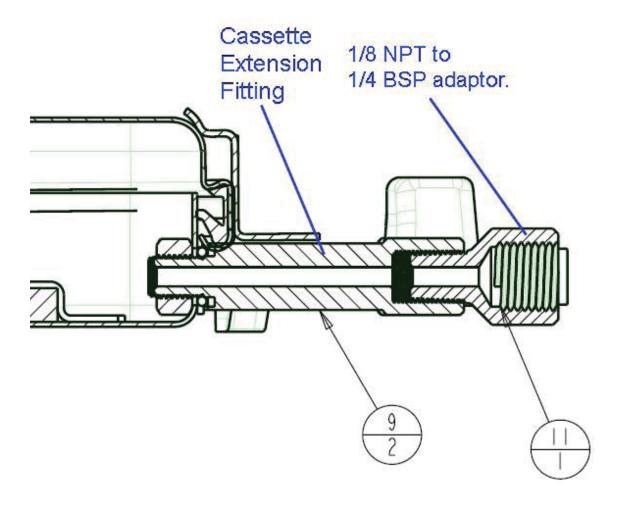


Figure 9. Pressure fitting cut-away

To set up the cassette for use with the recommended **Druck DPI705R** and **Heise PTE1** pressure reference meters, proceed as follows:

Remove the 1/8 NPT to 1/4 BSP adaptor. (The accessories for the cassette include 2 male Swagelok 'Quick connect' fittings (see Figure 10).

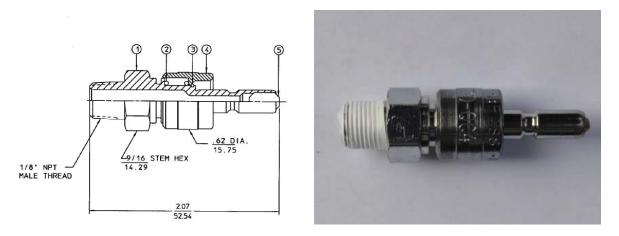


Figure 10. Male Swagelok B-QC4-S-2PM 'Quick Connect' fitting.

Attach and tighten the 'Quick Connect' fitting to the cassette extension fitting in place of the 1/8 NPT to 1/4 BSP adaptor (see Figure 11).

NOTE: PTFE tape should be used when assembling these parts to ensure a steam tight seal.



Figure 11. Male Swagelok 'Quick Connect' fitting and cassette extension fitting assembly.

### Attaching the pressure tube to the cassette

The pressure tube (see Figure 12) supplied with the calibration cassette is designed to insulate the pressure transducer used with the pressure reference meter from the high temperatures experienced during sterilization, which may damage the transducer.



Figure 12

The tube is fitted with a Swagelock, self sealing 'Quick Connect' QC4 female connector (see Figure 13) on either end. These connectors attach to the male 'Quick Connect' fitting shown in Figure 10 above.



Figure 13. Swagelock 'Quick Connect' QC4 female connector

To attach the pressure tube to the cassette (via the male Swagelok connector), press the male and female fitting together until a 'click' is heard and the couplings are firmly attached. (See Figure 14)



Figure 14. Assembled male fitting, female fitting and cassette extension fitting

### Attaching the pressure tube to the recommended pressure reference meter (s)

### Important notes:

- Attaching the pressure tube to the reference meter is achieved using the second male Swagelok b-qc4-s-2pm fitting supplied with the cassette.
- Attaching the tube to the connector is as described in the section above, however, the second male fitting will need to be installed onto the pressure reference meter prior to attempting this.
- How the second male fitting is attached will depend on the type of reference meter used.

Assembling the Swagelok b-qc4-s-2pm 'quick connect' fitting to the Druck dpi 705r absolute pressure meter external 0 - 7 bar absolute pressure transducer with ¼" NPT female thread.

NOTE: To complete this assembly you will need to locate the 1/8 NPT female to 1/4 NPT male adapter supplied with the S-class calibration cassette. (See Figure 15 for assembly).

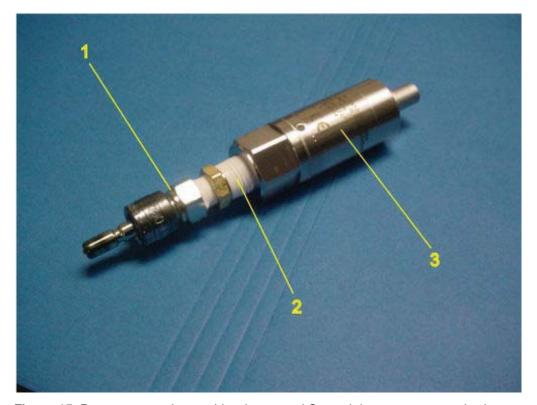


Figure 15. Pressure transducer with adapter and Swagelok connector attached.

Install Swagelok connector (1) onto adaptor fitting (2) DO NOT OVERTIGHTEN FITTING.

NOTE: The use of P.T.F.E. tape as shown is essential to ensure the connection is free from steam leaks. If the joint leaks, steam can enter the transducer and damage the unit internally.

Install assembly (from previous step) onto external pressure sensor (3) supplied with the Druck DPI 705R. **DO NOT OVERTIGHTEN FITTING.** 

NOTE: The use of P.T.F.E. tape as shown is essential to ensure the connection is free from steam leaks. If the joint leaks, steam can enter the transducer and damage the unit internally.

The pressure tube can now be attached to the pressure transducer.

Figure 16, below, shows the meter, pressure transducer and pressure tube attached to a calibration cassette.

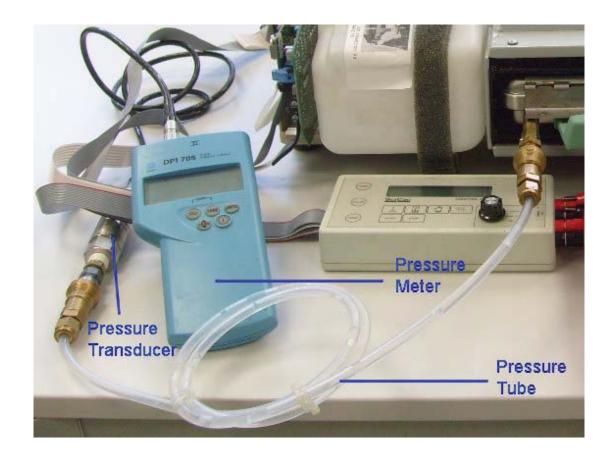


Figure 16

Assembling the Swagelok B-QC4-S-2PM 'Quick Connect' fitting to the internal pressure module supplied with the Heise hand held calibrator.

NOTE: You do NOT need the 1/8 NPT female to  $\frac{1}{4}$  NPT male adaptor for this device. (see Figure 17 for assembly).

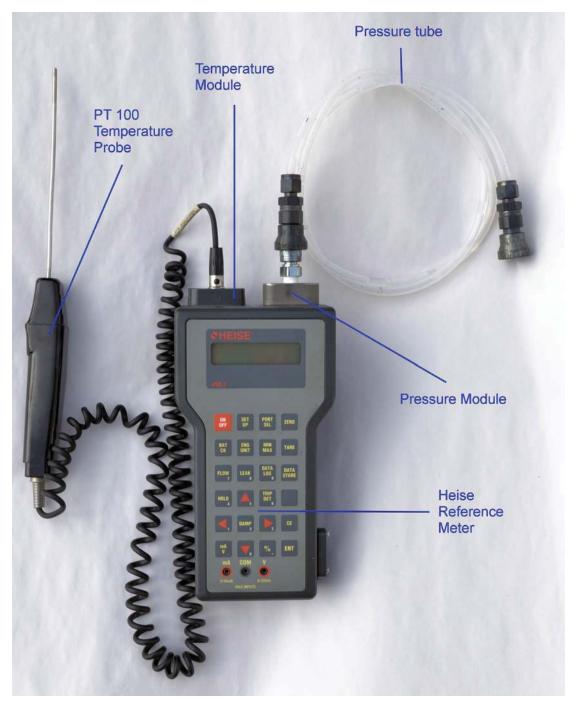


**Figure 17.** Heise meter with Swagelok fitting attached to pressure module.

Insert and tighten the Swagelok fitting into the 1/8 NPT female thread of the pressure module. **DO NOT OVERTIGHTEN FITTING.** NOTE: The use of P.T.F.E. tape as shown is essential to ensure the connection is free from steam leaks. If the joint leaks, steam can enter the module and damage the unit internally.

The pressure tube can now be attached to the pressure module.

Figure 18, below, shows the meter with temperature probe and pressure tube ready to attach to a calibration cassette.



**Figure 18.** Heise hand held calibrator with temperature module, pressure module, PT 100 temperature probe and pressure tube.

# **Annual Service Requirements**

### Important notes:

- The following schedule is designed as a guide for SciCan approved trained service engineers and applies to all STATIM 2000/5000/G4 models.
- The operations contained in this guide are for use during standard service intervals
  on functioning machines. Any operations or components required over and above this
  guide will be at the discretion of the service engineer and the customer and additional
  to the requirements outlined below.
- Components other than the routine annual service components may require changing due to wear and tear or failure but only require changing as and when degradation or failure occurs and should not need to be changed routinely on an annual basis. These may be identified during service and changed as appropriate.
- Consumable items such as seals and filters may require changing at more frequent intervals outside of normal service intervals dependent on the frequency of use of the unit.

SERVICE SCHEDULE			
Standard service parts required for STATIM 2000/2000S and G4 variants:			
Part Number Description			
01-100028S	Seal and Lubricant kit A/B		
01-100207S	Filter for air compressor, STAT/M 2000/2000S (where fitted)		
01-102119S	Filter Biological, B/C (where fitted)		
01-100574S	Check valve, B		
01-100998S	Repair, Solenoid valve (Honeywell) A/B/C/D		

Standard service parts required for STATIM 5000/5000S and G4 variants:	
Part Number	Description
01-101649S	Seal and Lubricant kit C
01-101652S	Filter for air compressor, STAT/M 5000/5000S (where fitted)
01-1021198	Filter Biological, B/C (where fitted)
01-101627S	Check valve, C
01-100998S	Repair, Solenoid valve (Honeywell) A/B/C/D

#### SERVICE PROCEDURES

#### Main unit

#### With cover on

Run cycle on unit to check for leakage or faulting.

Check unit's overall condition, including exhaust tube and bottle. Clean bottle, tighten fittings and check for kinks in exhaust tube as appropriate.

#### With cover off

Check internal condition of unit, pay particular attention to corrosion.

Check water reservoir for contamination. If necessary, disconnect and remove reservoir, clean and rinse with warm distilled water. **Do not use any chemicals. Refit reservoir.** 

Check and clean probes if required.

Remove and fit new bacteriological filter. (where fitted)

Remove and fit new compressor filter. (where fitted)

Remove and fit new solenoid plunger and associated parts.

Remove and fit new check valve.

### Cassette

Remove cassette seal.

Clean process residue from all surfaces of cassette. (Use chlorine-free soap and scrub with cleaning pad designed for Teflon.)

Check cassette for corrosion under seal seat.

Check cassette for damage. Note: pay particular attention to the rear hinge and rear wall of the base as damage in this area may cause steam leakage even if a new seal is fitted.

Fit new cassette seal.

### **Diagnostics and Testing**

### Note:

- Control box, appropriate test cassette (with temperature and pressure port for 'S' units and temperature port only for non 'S' units), appropriate reference meters (temperature and pressure for 'S' units and temperature only for non 'S' units) and voltage meter are required for this section.
- For correct diagnostics and calibration of 'S' units, both temperature AND pressure meters MUST be used.

Check running function of pump, solenoid valve and compressor (control box buttons). Cassette must be disconnected from probes.

Check pump flow in accordance with steam generator type and correct as required.

Check reference voltage as appropriate to board type and adjust as required. **Not required on revision 7 units.** 

Switch off machine and switch on in calibration mode as appropriate to model.

Check steam generator thermocouple calibration as appropriate to board type.

Note: steam generator calibration is not required for software revision numbers XXXXR4XX onwards.

During the above cycle, check pipe work and fittings for steam leaks.

Check and adjust chamber thermocouple calibration as appropriate to board type.

Note: on 'S' class units, *pressure transducer calibration* using the pressure reference meter must be done during this phase. Software revision XXXXR4XX onwards requires this to be done at barometric pressure at the start of the cycle. Previous software variants require this to be done during the sterilization phase.

When set, and during the running cycle, check that the chamber temperature matches the reference meter and that actual and theoretical pressure are within 1 kPa. This is important for the efficient balance of the control and validation circuits. Failure to achieve this may result in unnecessary cycle failures.

Re-check steam generator thermocouple calibration as appropriate to board type.

- Note 1: steam generator calibration is not required for software revision numbers XXXXR4XX onwards.
- Note 2: if the initial steam generator calibration offset and final steam generator calibration offset are > 8 hexadecimal values apart, the steam generator MAY be compromised and further investigation should be undertaken.

After completion of calibration, run a standard cycle to check temperature and pressure values.

Reinstall cover, and check all cycles are available.

Insert customer's cassette and run standard cycle. Check for leaks.

For revision 7 units ONLY: Enter service mode and activate 'Back up NVRAM'.

# Routine Maintenance Schedule and Procedures

To maintain the STATIM in good order between annual services, follow these guidelines:

### **DAILY**

### Water Reservoir

- Replace water as needed.
- For opthalmic use, drain at the end of every workday, leave empty, and refill at the start of the next workday.

### Water Bottle

- Empty the waste bottle every time you refill the reservoir.
- Fill with water, up to MIN line marking. You may also add some chlorine-free disinfectant.

### **WEEKLY**

### Cassette

- Wash the interior of the cassette with dishwashing soap or a mild detergent that does not contain chlorine.
- Scrub the inside with a cleaning pad designed for use with Teflon™-coated surfaces.
- After removing all traces of the detergent, treat interior surfaces of the cassette with the STAT-DRI™ Plus drying agent to enhance the drying process. Order more STAT-DRI™ Plus from SciCan quoting 2OZPLUS, 8OZPLUST, or 32OZPLUS. Please note that STAT-DRI drying agent is not be to used with the U.S. G4 models (G4-121101 and G4-201103).

### Biological and/or Air Filter

Check the filter for dirt and moisture. Replace if dirty. Call for service if wet.

### Water Reservoir Filter

 Check the water reservoir filter every week and clean if necessary. Replace only if necessary.

### **EVERY 6 MONTHS**

### Cassette Seal

Replace every 500 cycles or six months (whichever is first), or whenever necessary.


### Biological and/or Air Filter

Replace every 500 cycles or six months (whichever is first).

### **Draining the Reservoir**

If you have to service the reservoir, ship the STATIM, or move the unit more than a short distance, the reservoir will need to be drained to prevent potential water damage to the internal components (particularly the electronic and electrical components) of the unit.

Also, in some regions of the world, local or national guidelines recommend routine draining and cleaning of the reservoir to reduce the potential of the accumulation of contaminants in the feed water.

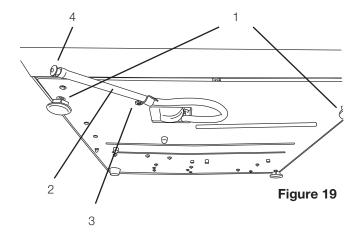

To drain the water from the reservoir, carefully move the STAT*IM* to the edge of the work surface and proceed as follows (see Figure 19):

- The front leveler feet (1) should be approximately 12 mm (1/2 inch) from the edge so the unit remains securely seated on the work surface.
- Lift the front left corner of the STAT*IM* upward and remove the drain tube (2) from the clip (3) located on the underside of the unit. Gently pull the tube out as far as possible so the free end can be positioned over a container when the unit is lowered back to the work surface.
- Remove the stopper (4) from the end of the drain tube and allow the water to drain from the reservoir.
- When water no longer drips from the drain tube, replace the stopper.
- Lift the front left corner of the STATIM upward and re-insert the tube into the clip on the underside of the unit. Push the excess length of tubing back into the unit.

### Cleaning the Cassette

Keeping the STAT*IM* cassette clean is good clinical practice and assists in the proper functioning of the unit. SciCan recommends that the interior surface be cleaned at least once a week. Cleaning the inside of your cassette is very important if you regularly sterilize lubricated instruments.

Use dish washing soap or a mild detergent that does not contain chlorine. Scrub the inside of the cassette with a cleaning pad designed for use with Teflon  $^{\text{TM}}$  coated surfaces. After scouring, rinse thoroughly with water to remove all traces of the detergent.

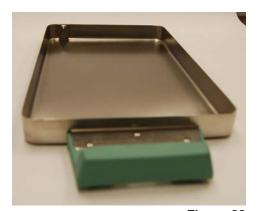


Figure 20

Coating the entire inside surface with STAT-DRI drying agent induces water to form an even coat on the inside surface, without beading. The water in contact with the hot cassette surfaces also evaporates much more efficiently. Spotting is minimized and instruments dry much better.

STAT-DRI should be applied every 10 cycles, and after every cassette cleaning.

STAT-DRI is available from your SciCan Dealer. Please note that STAT-DRI drying agent is not be to used with the U.S. G4 models (G4-121101 and G4-201103).

### Cleaning the Water Reservoir Filter (where fitted)

This water reservoir filter should be cleaned at least once a week or when required. The filter can easily be removed and cleaned by placing the filter upside down under running water to wash away the particles until clean, and then placed back into the reservoir opening. Replacement water reservoir filters are available from you SciCan dealer.



Figure 21

### Cleaning the Water Reservoir

Check the reservoir for dirt or particles. The reservoir may be cleaned by draining followed by cleaning and rinsing with steam process distilled water ONLY. Use of chemicals or cleaning agents is not recommended and could cause damage to the unit.

### Changing the air filter and bacteria retentive filter (where fitted)

Never re-use an old filter. Never run the unit without a clean filter in place.

After each sterilization cycle, the compressor forces air through the cassette to rid it of steam and to dry the instruments. The air is drawn into the compressor through one or two filters located at the back of the unit, depending upon which unit is being serviced.

Note: that some units may not have both or either filters fitted depending on age, unit type, compressor type, or if a compressor is fitted.

### Filters on the STATIM 2000/2000S/2000 G4

The STATIM 2000 type units draw air into the compressor (where fitted) through a circular foam air filter located behind a cover plate on the rear of the compressor. This plate can be accessed from the rear of the unit with the cover still attached. The air is then directed to the steam generator and cassette, by the compressor, through a bacteria-retentive air filter (where fitted) which is held in a bracket attached to the rear cover of the unit. If both filters are present, always change both filters at the same time.

To change the filters, proceed as follows (see Figure 22):

### **Bacteria Retentive Filter**

- Disconnect the inlet tube (1) from the bacteriaretentive filter (3) and remove the filter from the filter bracket (4). Note the orientation of the Arrow mark on the filter and bracket.
- 2. When the filter is free of the bracket, carefully disconnect the outlet tube (2) from the filter. Do not pull on this tube.
- 3. Before installing the replacement bacteria-retentive filter check that the arrow mark on the filter matches the direction of the arrow on the bracket. Push the left hand filter fitting into the outlet tube.
- 4. Gently press the replacement filter into the filter bracket. The arrow mark on the filter is facing out and pointing to the left.
- 5. Re-connect the inlet tube to the right-hand filter fitting.

### **Compressor Filter**

- Remove the screw (5) holding the compressor plate
   (6) to the back of the compressor (7). Remove the plate.
- 2. Remove and discard the old filter.
- Install the new filter (8), and secure the compressor plate to the rear of the compressor using the screw retained from disassembly.




### Filters on the STAT/M 2000/2000S/2000 G4

The STAT*IM* 5000 type unit draws air into the compressor via a tube connected to a cylindrical air filter located on the chassis behind bacteria retentive filter bracket. **Note:** On certain 5000 models (dependent on compressor type) this filter may not be fitted. The air is then directed to the steam generator and cassette, by the compressor, through a bacteria-retentive air filter which is held in a bracket attached to the rear cover of the unit. If both filters are present, always change both filters at the same time (see diagram below).

To change the filters (where fitted), proceed as follows (see diagram above):

### **Bacteria-Retentive Filter**


- Disconnect the inlet tube (1) from the bacteria-retentive filter (3) and remove the filter from the filter bracket (4). Note the orientation of the arrow mark on the filter and bracket.
- 2. When the filter is free of the bracket, carefully disconnect the outlet tube (2) from the filter. **Do not pull on this tube.**

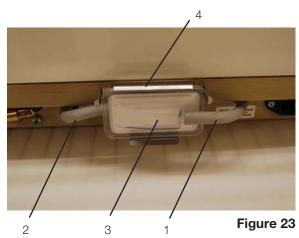

- 3. Before installing the replacement bacteria-retentive filter check that the arrow mark on the filter matches the direction of the arrow on the bracket. Push the left hand filter fitting into the outlet tube.
- 4. Gently press the replacement filter into the filter bracket. The arrow mark on the filter is facing out and pointing to the left.
- 5. Re-connect the inlet tube to the right hand filter fitting.

### **Compressor Filter**

- 1. To enable access to the compressor filter, disconnect the inlet tube (4) from the bacteriaretentive filter (2) and remove the filter from the filter bracket (3). Note the orientation of the arrow mark on the filter and bracket.
- 2. When the filter is free of the bracket, locate and remove compressor filter (5) by unscrewing anti-clockwise using thumb and forefinger. Note: the filter should only be finger tight.
- 3. Install new compressor filter by screwing clockwise using thumb and forefinger. The filter should ONLY be finger tight. Tightening the filter using mechanical means may damage the filter or filter housing.
- 4. Gently press the bacteria retentive filter back into the filter bracket. The arrow mark on the filter is facing out and pointing to the left.
- 5. Re-connect the inlet tube to the right hand filter fitting.

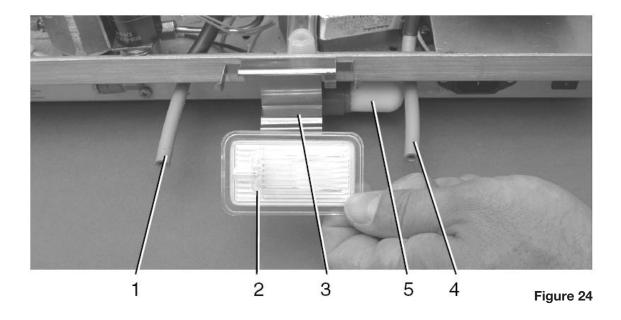


### Replacing the Cassette Seal

To ensure optimum performance of your STATIM cassette autoclave, change the cassette seal every 500 cycles or every six months, whichever comes first. Replacement seals are available from SciCan (order number 01-100028S for STATIM 2000 and 01-101649S for STATIM 5000).


To change the cassette seal, follow these steps:

Place the cassette lid and the new seal on a clean work surface. Examine the position of the old seal in the cassette lid and arrange the new seal in the same orientation, next to the lid.


Remove the old seal by gently easing one corner of the seal out of the seal channel with a small screwdriver or similar device and then pull the remaining seal from the channel by hand. Ensure there are no sharp edges on the removal device that may damage the seal channel. Discard the old seal. Clean any residue out of the seal channel and flush out the channel with distilled water.

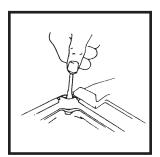


Lubricate the new seal with the liquid seal lubricant provided. Do not use an excess amount of liquid, only use sufficient liquid to lubricate the seal. Note: the seal will appear to be slightly larger in length and width than the seal channel. This is intentional as the seal compresses into the channel when fitted to ensure the correct pressure on the base of the cassette when closed.

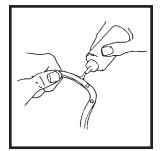

Locate and align the steam inlet holes in the seal by inserting under the lip of the seal channel in the cassette lid and in line with the corresponding holes in the lid. The holes should be concentric to one another and the square tabs adjacent to the holes should be located evenly in the cut-outs in the cassette lid. Ensure seal sits firmly against the back of the channel.

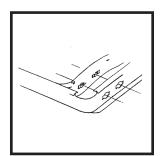

Locate each corner of the seal under the lip of the seal channel in the cassette lid in the corresponding corners of the seal channel. At each corner, the two square tabs should be visible and located on the edges of the corner cut outs in the cassette lid.

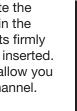



Working from the centre position outwards on each side, locate the sides of the seal by inserting under the lip of the seal channel in the cassette lid and pushing the bottom section in until the seal sits firmly against the back of the channel. Ensure the seal is completely inserted. Running the thumb along the seal applying light pressure will allow you to feel any irregularities and evenly distribute the seal in the channel.

Figure 25

Note: During the first cycle after changing the seal, steam may appear between the lid and the tray as the seal seats. If this persists, remove the cassette and check that the seal is correctly installed.

Be careful. The metal parts will be hot, and the cassette may contain hot steam.

### Cleaning the Cover

Use a soft cloth moistened with a mild cleaning solution or a mild disinfectant to clean all outside surfaces. Do not use solvents or harsh chemicals.

### Shipping the unit

If you must ship your STATIM, follow these directions:

- 1. Drain the reservoir. See Draining the Reservoir.
- 2. Screw in each leveler foot completely.
- 3. Repack your STAT*IM* in the original packing materials.
- 4. Specify heated and insured shipping.

### Upgrading the firmware on STATIM G4 units

To upgrade the firmware, use a USB drive and proceed as follows:


- 1. Download new firmware (available from SciCan upon request).
- 2. Insert a clean USB stick (containing no files) into your computer
- 3. Save the attached zip files (SL00RXXX.zip, SD00RXX.zip) to your desktop and extract the files contained in the zip files directly to the USB stick.
- 4. Power OFF the STATIM G4 unit.
- Remove the USB stick from your computer and insert into your STATIM G4 USB port located on the front left of the unit.
- 6. Power ON the STATIM G4 unit.
- 7. Wait for approximately 10 minutes while the firmware is upgraded. (NOTE: During this process the screen will be blank. To complete the firmware upgrade, the G4 will restart automatically. DO NOT REMOVE THE USB STICK UNTIL THE END OF STEP 13.
- 8. Once you see the standard screen, browse to the User Setup menu.
- 9. Browse to the Themes button. Press the Themes button.
- 10. Select Upgrade Theme.
- 11. Select Upgrade.
- 12. Wait for approximately 24 minutes. (Note: During the first part of this step, text will scroll in a box on the screen. During the second part of this step, the unit will show a screen indicating that a firmware upgrade is in progress and advise you NOT to power OFF the unit. The unit will restart automatically after the theme upgrade is completed.
- 13. Once the LCD displays the standard screen, the process is complete.
- 14. Remove the USB stick.

### Using the STAT*IM* G4 remote access function

Users can allow offsite technicians to remotely access the LCD touchscreens and web portals of STATIM G4 units connected to a network.

For local network remote access, proceed as follows:

- 1. From the TOOLS page, click on the LOCAL CONTROL tab.
- 2. Click on the start button to start a local connection. It will open up a page that mirrors your STAT*IM* unit's touchscreen so that it can be controlled remotely within your local network.

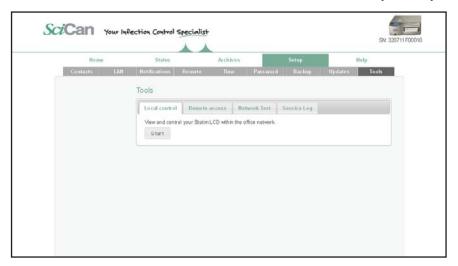


For remote access of a STAT*IM* web portal or touchscreen from outside a local network, proceed as follows:

- 1. Someone onsite with the unit or from within the network must provide access to an outside user by generating a 'token' (or code).
  - To generate a unique token using the web portal, go to the TOOLS page and click on the REMOTE ACCESS tab.




 To generate a unique token using the unit's LCD touchscreen, go to the Settings menu and scroll to Remote access and follow the prompts to enable remote access.

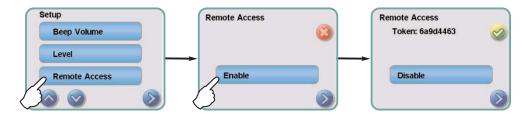

2. The technician attempting to access the G4 from outside the network will need to go to the following URL: http://updates.scican.com and enter their registered email address, password, token and STATIM Serial Number (optional).

To create a new account to enable remote access for a STATIM, click on the CREATE NEW ACCOUNT link, complete the form, and click on the SUBMIT FORM link. The system will send a confirmation email to verify the account. Once confirmed, the account will be ready to use.

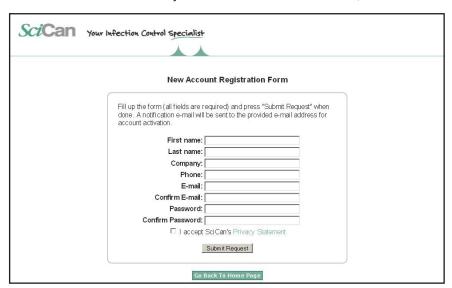

With a valid user name, password and token, a technician can remotely access the STATIM G4 unit's web portal page.

3. Go to SETUP. Another username and password prompt will appear.

### To access the web portal only, use the following default values:

User name: scican

User password: scican (user can change this password)

# To access the web portal and remote access the LCD touchscreen, use the following default values:

User name: scican

User password: s23can173

4. Upon authentication, go to TOOLS and click on REMOTE ACCESS. A page will open that mirrors the STAT*IM* unit's touchscreen so that it can be controlled remotely from outside its local network.

# Electrical safety testing

### Dielectric Strength Test (Hi-Pot)

When mains electrical components are serviced or replaced and when the unit cover is removed and reinstalled following servicing, the dielectric strength of the electrical insulation between the mains and user-accessible conductive parts must be tested to ensure continued compliance of the unit with applicable international safety standards.

Exercise caution while performing this test. Hazardous voltages are present. Do not touch the unit, or allow the unit to touch any conductive surfaces during the test.

Perform the test using a dielectric strength (hi-pot) tester operated in accordance with the manufacturer's written instructions using the following parameters:

STATIM Classic 2000 units

W-Ramp: 2 Seconds

W-Volt: 1500 V

Dwell Time: 2 Seconds

W-High: 3.5mA W-Low: 0.05mA Frequency: 60 HZ Continuity: OFF

STATIM G4 2000 units W-Ramp: 2 Seconds W-Volt: 1500 V

Dwell Time: 2 Seconds

W-High: 6.0mA W-Low: 0.05mA Frequency: 60 Hz Continuity: OFF

Connection Points:

A. Unit chassis

B. Live and neutral terminals of the mains plug connected together.

For the unit being tested to pass, there must be no breakdown of the insulation or any flashover.

### Protective Bonding Impedance Test (Ground Continuity)

Applicable international safety standards require that the impedance between the protective conductor terminal of the power entry connector and any user-accessible conductive parts must not exceed 0.1 ohms. If any components of the protective earthing system are changed or any connections of that system are broken and remade, the impedance of the protective bonding must be tested and verified as being less than 0.1 ohms using test equipment and procedures in compliance with applicable International safety standards and national, state, provincial, and regional laws and regulations.

# **Water Conductivity Testing**

Water conductivity testing determines the amount of dissolved solids in the steam-process distilled water used in the STAT/M unit. Use only steam-process distilled water having 5 ppm or less dissolved solids or a conductivity of less than 10  $\mu\text{S}$  / cm. Follow the manufacturer's instructions to test water conductivity using SciCan conductivity meter, part #01-103139S, or any other equivalent Meter.

Temperature has a significant effect on conductivity readings, therefore the water being tested should be at room temperature.



Water Conductivity Measurement

Figure 26

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



## **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



# PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# **Diagnostics and Troubleshooting**

When a STATIM unit fails to function correctly it is designed to 'fail safe' so that if sterilization conditions are not met, a message or cycle fault will appear on the display to inform you of the failure.

Depending on the age, type and controller variant of the machine, there are a number of things you can do to understand the problem. The following chapter contains guidance on what can be done on the various unit types.

### Basic mechanical diagnostics

If the immediate cause of the problem is unknown and a cycle fault is unavailable, the following functional checks are recommended:

- Inspect the unit for external damage.
- Remove the cover and inspect the unit for internal damage.
- Attach a control box and power ON the unit (if possible) and run a cycle to see if the problem can be replicated.
- While the unit is running, check for steam leaks in the pipe work, cassette or check valve.
- Visually inspect the bacterial filter and compressor filter for signs of contamination.
- Perform component tests to determine the operational status of the solenoid valve, pump and compressor (if fitted).
- Perform a pump flow check.
- Check calibration status of the unit.

### Using software as a diagnostic tool

Depending on the software in the problem unit there are a variety of things you can do to gather information to help with diagnosis. This varies from unit to unit and you will find that the newer the unit, the more options are available.

### For units with revision 2.x/5.x and 6.x controllers (no service menu)

For the above controller boards with version 2.xx and higher software it is possible to print out the last complete or incomplete cycle as a diagnostic tool. If the unit does not have a printer or if the printer was not connected when the fault occurred, a STAT/M printer can still be used as a diagnostic tool, however the time and date printed on the report may not be valid.

To print the last cycle information, proceed as follows:

- 1. Turn the STATIM OFF and remove the cover.
- Install a calibration jumper on controller board header W1 located to the left of the microprocessor. Connect a control box and printer to the controller board and enable the printer.

- 3. Turn the STATIM ON. The unit is now in calibrating mode.
- 4. For version 2 and 3 software without selecting a cycle, press the START button. For other versions of software press the "Wrapped" button.
- 5. The printer will print the last complete or incomplete cycle. This information can be useful in diagnosing intermittent problems when the cycle fault code has not been recorded.

# For units with revision 7.x controllers (with service menu), including G4

All revision units fitted with revision 7.x controllers have a service menu which is designed to:

- Provide information on the setup of the machine.
- Allow the changing or resetting of various machine parameters.
- Operate various aspects of the unit with the cover on as opposed to using a control box.

Certain information provided in this menu may be used to assist in diagnostics.

# Using the service menu on the STATIM Classic (Rev. 7 controller board)

The service menu can be displayed on the LCD screen of any STATIM Classic with a Rev. 7 board by leaving the cover on, or using the LCD screen of the control box if the cover is removed and detached from the unit.

- To enter the Service Menu, power ON the STATIM while holding down the "Unwrapped" and "Wrapped" keys.
- You will then be prompted to enter the service password to allow you to enter the service menu.
- The default password is entered by pressing the "Unwrapped", "Wrapped", "Rubber and Plastics" and "Stop" keys in that order.
- The service technician can change this password. In case the changed password is lost or has been changed by another technician, a backdoor password can be used (which can not be changed) by pressing the "Unwrapped", "Wrapped", "Unwrapped", "Wrapped" keys in this order.

The following menu is now available and the display on the LCD screen of the unit or control box will be a two-line display as highlighted by the **bolded** box in the following table.

Calibration	i
>Calibration Time/Date setup	
Language setup	
Unit ID setup	
Set cycle counter	]
Lock cycles	For S class units only
Conductivity setup	
Water. cnd tmp. comp	]
Last printout	
Stored CF printouts	]
Clear CF printouts	]
Display last CF#	
Devices test On/Off	]
Temperature offset	
Press. atm. offset	For S class units only
Press. ster. offset	For S class units only
Press. sensor type	For S class units only
Validation offset	For non-S class units only
Voltage calibration	]
Voltage setup	1
Repeater mode	]
RS232	1
End of line CR/LF	1
Serial port bitrate	]
Printer user ß char	]
Factory default	
Change password	
Backup NVRAM	
Restore NVRAM	
Save and exit	
Exit	
Water pump type	
Production cycle	

Note that to navigate the menu you can use the keypad buttons that will operate the following

"Unwrapped" button - Selects the next item in the menu (down)
"Wrapped" button - Selects the previous item in the menu (up)

"Rubber and Plastics" button - Enters the indicated sub-menu item (as indicated by the cursor ">")

"Stop" button - Exits the menu and returns the unit to operational mode

selections:


NOTE: For some functions (depending on the function of the sub menu), pressing the "Stop" button may NOT save the changes. When changes made to unit values are NOT saved after pressing the STOP button and the unit returns to the service menu, always select SAVE AND EXIT if you wish to save any new values.

Sub Menu items (Activated by the "Rubber and Plastics" button when the cursor ">" is next to the description):

### >Calibration

This will display the calibration screen relevant to the specific revision 7 unit and allow you to perform calibration operations. For calibration information, refer to calibration procedure 4 in Chapter 4 of this manual (non S class) or calibration procedure 7 in Chapter 5 of this manual (S class).

### >Time/Date Setup

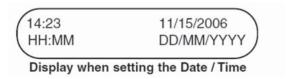


There are 5 selectable items (fields), hours, minutes, date, month and year. After entering the menu, the 'hour' value will be highlighted.

Using the cycle buttons to select and change the selected field's value.

- 1. To increase a field's value, press the UNWRAPPED cycle button. Holding the button down will increase the value until released.
- 2. To decrease the value, press the WRAPPED cycle button.
- 3. To select the next field, press the RUBBER & PLASTICS cycle button.
- 4. To save changes and return to the regular operating mode, press the **STOP** button. The changes will be saved **automatically**.
- 5. To guit without making changes, power the STATIM **OFF**.

### >Language Setup

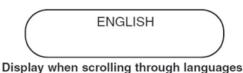


There are many languages available in the STATIM, which will display the operating information on the LCD in the respective language chosen.

The default language in all units is ENGLISH, and on entering the language setup sub-menu, alternative languages can be accessed as follows:

- 1. Press the UNWRAPPED cycle button to scroll to the next language (down).
- 2. Press the WRAPPED cycle button to scroll to the previous language (up).

- 3. To save changes and return to the regular operating mode, press the **STOP** button. The changes will be saved **automatically**.
- 4. To quit without making changes, power the STAT/M OFF.

### >Unit ID Setup

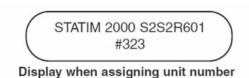

This function allows the user to assign a unique, 3-digit number to a STATIM unit if they have more than one. This allows traceability of recorded data (printed or electronic) to individual units.

On entering the ID sub menu, the first digit will be highlighted and can be changed as follows:

- 1. The UNWRAPPED button will increase the selected value (0, 1, 2, 3, etc.) and the WRAPPED cycle button will decrease it (0, 9, 8, etc.).
- 2. To move to the next digit, press the RUBBER & PLASTICS button.
- 3. To save changes and return to the regular operating mode, press the **STOP** button. The changes will be saved **automatically**.
- 4. To quit without making changes, power the STAT/M OFF.

### >Set Cycle Counter



Screen display when in set cycle counter sub menu.

This function allows the technician to change the cycle count of the unit. This is designed to be used when a complete controller change is necessary and the customer needs to have continuity of the cycle count on their records. A new controller will display 000000, but the 'old' cycle number can be re-introduced.

On entering the Set Cycle Counter sub menu, the first digit will be highlighted and can be changed as follows:

- 1. The UNWRAPPED button will increase the selected value (0, 1, 2, 3, etc.) and the WRAPPED cycle button will decrease it (0, 9, 8, etc).
- 2. To move to the next digit, press the RUBBER & PLASTICS button.
- 3. To save changes and return to the regular operating mode, press the **STOP** button. The changes will be saved **automatically**.
- 4. To guit without making changes, power the STATIM OFF.

### >Lock Cycles (S class units ONLY)

'S' class STATIM units with revision 6.x and 7.x controllers boards have 7 sterilization cycles available as opposed to 3 (pre 2004 2000S units) or 4 (pre 2004 5000S units) sterilization cycles. This means that unlike the earlier (pre 2004) units, each cycle selection button (except extra drying) has 2 or 3 cycles per button as opposed to the single cycle on earlier units. These multiple cycles

per button can be accessed by pressing any specific button repeatedly until the desired cycle is selected. This cycle is automatically the first cycle to show on the display the next time the cycle button is pressed.

Note that units with Revision 6.x controller boards (rev 4.xx/5.xx software) DO NOT have a service menu, so the locking/unlocking procedure is different.

If you do not require some of the extra cycles, each button can have the redundant cycles 'locked out' to simplify cycle selection. Note that in 'locked' mode, each cycle button will only have one cycle available. Single buttons cannot be locked, it must be all or none.

The cycles available on all 'S' class units mentioned above are as follows:

Button	Cycle name	Sterilisation temperature	Sterilisation time (minutes)
	SOLID UNWRAPPED	134°C	3:30
$(\Lambda)$	HOLLOW UNWRAPPED	134°C	3:30
0.0	HOLLOW UNWRAPPED	134°C	18:00
FF	HOLLOW WRAPPED**	134°C	3:30
A	HOLLOW WRAPPED	134°C	18:00
an	RUBBER/PLASTIC	121°C	15:00
	RUBBER/PLASTIC	121°C	30:00
	AIR DRYING ONLY	Not a sterilis	ation cycle

- Solid Unwrapped 134°C/3.5 min (N)
- Hollow Unwrapped 134°C/3.5 min (S)
- Hollow Unwrapped 134°C/18 min (S)
- Hollow Wrapped 134°C/3.5 min (S)
- Hollow Wrapped 134°C/18 min (S)
- Rubber and Plastics 121°C/15min (S)
- Rubber and Plastics 121°C/30min (S)

### Locking procedure for revision 6.x controller boards/software revision R4xx and R5xx

- 1. Turn the unit ON.
- 2. Press "Unwrapped" button repeatedly until the desired cycle is displayed.
- 3. Press "Wrapped" button repeatedly until the desired cycle is displayed.
- 4. Press "Rubber and Plastics" button repeatedly until the desired cycle is displayed.
- 5. Turn the unit OFF.
- 6. Install a calibration jumper on Controller Board pins marked W1.
- 7. Keep button "Rubber and Plastics" pressed and turn the unit **ON**. The unit will enter "Cycle lock/unlock mode"

8. Upon entering this mode the display will show the following for a short time:

```
Cycle UNLOCKED
Available cycles
```

And then, it will display all of the available cycles (7 cycles).

9. Press the "Unwrapped" button and the unit will "lock in" the cycles selected at steps 2, 3 and 4, and "lock out" the unselected cycles. The display will show for a short time:

```
Cycle LOCKED
Available cycles
```

And then, it will display **only** the cycles selected (3 cycles + drying) at steps 2, 3 and 4 plus the drying cycle.

- Please note that if the "Unwrapped" button is pressed again the cycles will be unlocked and so on.
- 11. Turn unit off and remove calibration jumper. The changes will be saved automatically.
- 12. For unlocking the cycles on a unit with locked cycles, follow steps 6 and 7. Upon entering locking/unlocking mode the display will show:

```
Cycle LOCKED
Available cycles
```

And display only the selected cycles.

- 13. Press the "Unwrapped" button to unlock the cycles.
- 14. The display will now show the following, for a short time:

```
Cycle UNLOCKED
Available cycles
```

And then, it will display all of the available cycles (7 cycles).

15. Turn unit off and remove calibration jumper. The changes will be saved automatically.

### Locking procedure for revision 7.x controller boards/software revision R6xx

- 1. Turn the unit ON.
- 2. Press "Unwrapped" button repeatedly until the desired cycle is displayed.
- 3. Press "Wrapped" button repeatedly until the desired cycle is displayed.
- 4. Press "Rubber and Plastics" button repeatedly until the desired cycle is displayed.
- 5. Turn the unit **OFF**.
- 6. Press and hold the "Unwrapped" and "Wrapped" buttons together and switch the unit ON.
- 7. You should now have entered the 'Service Menu' and you will be asked for a password. The default unit password is entered by pressing the "Unwrapped", "Wrapped", "Rubber and Plastics" and "Stop" buttons.
- 8. The display should now show the following:

>Calibration	
Time/Date Setup	

- 9. Using the "Unwrapped" and "Wrapped" buttons, scroll down or up the service menu list until the sub menu title "Lock Cycles" appears next to the ">" cursor.
- 10. Now press the "Rubber and Plastics" button to enter the **Locking** sub menu.
- 11. Upon entering this menu the display will show the following for a short time:

```
Cycle UNLOCKED
Available cycles
```

And then, it will display all of the available cycles (7 cycles).

12. Press the "Unwrapped" button and the unit will "lock in" the cycles selected at steps 2, 3 and 4, and "lock out" the unselected cycles. The display will show for a short time:

```
Cycle LOCKED
Selected cycles
```

And then, it will display **only** the cycles selected (3 cycles + drying) at steps 2,3 and 4 plus the drying cycle.

- 13. Please note that if the "Unwrapped" button is pressed again the cycles will be unlocked and so on.
- 14. To save changes and return to the regular operating mode, press the **STOP** button. The changes will be saved **automatically**.
- 15. To quit without making changes, power the STATIM OFF.
- 16. For unlocking the cycles on a unit with locked cycles, follow steps 6 to 10. Upon entering locking/unlocking mode the display will show:

```
Cycle LOCKED
Selected cycles
```

And display only the selected cycles.

- 17. Press the "Unwrapped" button to unlock the cycles.
- 18. The display will now show the following for a short time:

```
Cycle UNLOCKED
Available cycles
```

And then, it will display all of the available cycles (7 cycles).

- 19. To save changes and return to the regular operating mode, press the STOP button. The changes will be saved **automatically**.
- 20. To quit without making changes, power the STATIM **OFF**.

### >Conductivity Setup



Screen display when in Conductivity setup sub menu.

STATIM units monitor water quality using conductivity. The default settings are based on the minimum level required to detect water in the unit (except units fitted with a float – see below) and the maximum level required to protect the internal components from damage due to mineral deposits (e.g. Calcium) on the internal surfaces.

This function allows the technician to adjust the water conductivity threshold. This may be required after a change of sensor and should be checked after a controller or software change.

Up to 2008, all units used the water sensor to detect both water quality and water level. In 2008, all units were upgraded with a float in the chamber to detect level, which allowed the sensor just to monitor quality.

The default conductivity setup is as follows:

Note that these settings should NOT be changed to suit water used as damage MAY occur to the components of the STATIM.

Units with sensor but WITHOUT float

- LL.L (Lower value threshold / No water threshold. Values lower than this triggers "Refill reservoir empty waste bottle" error) = 0.3uS
- HH.H (Upper value threshold / Bad water quality threshold. Values higher than this triggers "Water quality unacceptable" error) = 10.0uS
- G.GG Water conductivity circuit gain. Default value 1.00

Units with sensor and WITH float

- LL.L (Lower value threshold / No water threshold. Values lower than this triggers "Refill reservoir empty waste bottle" error) = FLOAT
- HH.H (Upper value threshold / Bad water quality threshold. Values higher than this triggers "Water quality unacceptable" error) = 10.0uS
- G.GG Water conductivity circuit gain. Default value 1.00

On entering the Conductivity setup sub menu, the first digit will be highlighted and can be changed as follows:

```
CD=xx.xuS/NNN/y.yppm
L=LL.L H=HH.H G=G.GG
```

The flashing value indicates the currently selected field.

The values shown on the screen are as follows:

- xx.x Actual water conductivity in uS (micro-Siemens) as read by the unit software.
- NNN Water conductivity in ADC (Analog to Digital converter) counts (0...255)
- **y.y** Water quality in ppm (parts per million) as read by the unit software.
- LL.L Lower value threshold / No water threshold/Float status.
- **HH.H** High value threshold / Water quality threshold.
- **G.GG** Water conductivity circuit gain

Note that xx.x, NNN and y.y are for information only as measured by the unit software and can not be changed.

Changes to the threshold value may be changed as follows:

- The UNWRAPPED button will increase the selected value and the WRAPPED cycle button will decrease it.
- 2. To move to the next value, press the RUBBER / PLASTIC button.
- 3. To save changes and return to the regular operating mode, press the **STOP** button. The changes will be saved **automatically**.
- 4. To quit without making changes, power the STAT/M **OFF**.

### >Water.Cnd Tmp. Comp

```
>Water.Cnd Tmp. Comp
On
```

Screen display when in Water.Cnd Tmp. Comp sub menu.

Electrical conductivity can change with temperature and as a result provide a false reading possibly causing the unit to show a water quality/level fault when the water is in fact good. To overcome this, revision 7 STATIM units have water conductivity temperature compensation built in to the software to overcome this potential problem.

The default setting on all STATIM units is **ON**. To change this (NOT recommended) or to reset this condition should it change as a result of a software or controller change, follow these steps:

- The UNWRAPPED or WRAPPED keys will switch the screen display to "ON" or "OFF" depending on the opening display.
- 2. The RUBBER & PLASTICS key will save the change and return the unit to the service menu.
- 3. Pressing the STOP key will exit without saving and return the unit to the operational menu.

### >Stored CF Printouts

When this function is activated, this sub menu will send saved Cycle Faults to a printer or SciCan Data Logger. The unit will then default back to the normal operational menu.

Note that the saved CF printouts will ONLY be sent to the printer or Data Logger when either one is attached and configured.

The following types of errors are saved:

- Cycle Faults (CF)
- Water quality or Water level low errors
- Cycle interrupted due to errors (##)

Pressing the STOP key will exit and return the unit to the operational menu.

### >Clear CF Printouts



Screen display when in Clear CF Printouts sub menu.

This function will reset the cycle fault list. The options are "**YES**" (clear the list) or "**NO**" (do not clear the list.

The default setting on all STATIM units is NO. To clear the list, proceed as follows as follows:

- 1. The UNWRAPPED or WRAPPED keys will switch the screen display to "YES" or "NO" depending on the opening display.
- 2. The RUBBER & PLASTICS key will save your selection and return the unit to the service menu condition.
- 3. Pressing the STOP key will exit without changes and return the unit to the operational menu.

### >Display last CF number (#)

```
>Display last CF#
## (#####)
```

Screen display when in Display last CF# sub menu.

This sub menu is for information only.

This function will allow you to observe the last fault condition of the unit for diagnostic purposes.

- 1. The display will show the CF number (## above) and the cycle number when the fault occurred ((######)).
- 2. Pressing the RUBBER & PLASTICS key will return the unit to the service menu.
- 3. Pressing the STOP key will return the unit to the operational menu.

### >Devices Test ON/OFF

```
>Devices Test On/Off
Pump Off
```

Screen display when in Devices Test ON/OFF sub menu.

This function on revision 7 STAT*IM* units allows the technician to activate a number of components for diagnostic purposes.

The list of components is as follows:

- Pump (water)
- Solenoid Valve (cassette must be disengaged)
- Compressor
- Yellow LED (cycle LED)
- Extra 1L (latched digital output 7000 only not used in 2000/5000)
- Extra 2L (latched digital output 7000 only not used in 2000/5000)
- Drawer relay (STAT/M 900 units only not used on 2000/5000)

When in the service menu (and unit in standby), the default condition of these components is OFF.

Activating (and de-activating) components is achieved as follows:

- 1. Using the UNWRAPPED (down) or WRAPPED (up) buttons will scroll through the sub menu to the next component in the list.
- 2. Pressing the RUBBER & PLASTICS button will activate (switch on) and deactivate (switch off) the chosen component.
- 3. Pressing the STOP key will return the unit to the service menu.

### >Temperature Offset

```
>Temperature Offset
##
```

Screen display when in Temperature Offset sub menu.

This sub menu is for information only.

## = unit chamber thermocouple offset value.

This allows the chamber thermocouple offset value to be read without entering calibration mode. No changes can be made to the offset value in this sub menu.

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Press. Atm.Offset - (S class units ONLY)

```
>Press. Atm. Offset
##
```

Screen display when in Press. Atm. Offset sub menu.

### This sub menu is for information only.

## = unit pressure sensor (transducer) offset value at atmospheric pressure.

This allows the pressure sensor atmospheric offset value to be read without entering calibration mode. **No changes can be made to the offset value in this sub menu.** 

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Press. Ster.Offset - (S class units ONLY)

```
>Press. Ster. Offset
##
```

Screen display when in Press. Ster. Offset sub menu.

### This sub menu is for information only.

## = unit pressure sensor (transducer) offset value for sterilisation conditions.

This allows the pressure sensor sterilisation offset value to be read without entering calibration mode. **No changes can be made to the offset value in this sub menu.** 

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Press. Sensor Type – (S class units ONLY)

```
>Press. Sensor Type
68PSI
```

Screen display when in Press. Sensor Type sub menu.

### This sub menu is for information only.

All S class STATIM units are fitted with a pressure transducer that actively measures the pressure in the chamber. Two types of transducer were fitted (60 p.s.i. and 68 p.s.i.) depending on the age and controller type of the unit.

Both transducers are visually the same but they are **NOT** interchangeable between units as each transducer requires type specific software for the unit to operate correctly.

Notes on how to identify the different types can be found in Chapter 9 of the service manual, but units with revision 7 controllers (including upgraded older units with the correct software) will show the pressure transducer type in this sub menu.

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Validation Offset – (Non S class units ONLY)

```
>Validation Offset
##
```

Screen display when in Validation Offset sub menu.

### This sub menu is for information only.

## = unit validation thermocouple offset value.

This allows the validation thermocouple offset value to be read without entering calibration mode. **No changes can be made to the offset value in this sub menu.** 

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Voltage Calibration.

```
>Validation Calibration
V = VVV VCal. =CCC
```

Screen display when in Validation Offset sub menu.

### This sub menu is for STAT/M 7000 units ONLY and is not used on 2000/5000 units.

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Voltage Setup

```
>Voltage setup
115V
```

Screen display when in voltage setup sub menu.

### This sub menu is for information only.

This allows the unit voltage value to be read. **No changes can be made to the value in this sub menu.** 

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Repeater Mode

```
>Repearer Mode
On
```

Screen display when in Repeater mode sub menu.

Repeater mode allows you to run a selected cycle continuously for bench testing. The cycle selected will run, vent and then dry for 2 minutes before re-setting and repeating.

The default setting on all STATIM units is **OFF**. To activate repeater mode, proceed as follows:

- The UNWRAPPED or WRAPPED keys will both switch the screen display to "ON" (or "OFF" depending on the opening display).
- 2. The RUBBER & PLASTICS key will save the change and return the unit to the service menu condition.
- 3. Pressing the STOP key will exit without saving and return the unit to the operational menu.

>RS232

>End of Line CR/LF

>Serial Port Bitrate

>Printer uses o char

At this point in the service menu you will see a repeat of the items in the user menu used for setting up data recording devices such as printers and Data Logger. For details on setting up these items refer to the user manual or Chapter 12 of the service manual.

### >Factory Default



Screen display when in factory default sub menu.

This function will reset the NVRAM (Non Volatile Random Access Memory) to the factory default settings.

Important Note: Any changes made to the chamber and voltage calibration offsets and conductivity settings will be reset.

The cycle counter will not be reset.

The default setting on all STATIM units is NO.

The alternative is **YES**, **RESET NVRAM**.

To activate a reset, proceed as follows as follows:

- The UNWRAPPED or WRAPPED keys will switch the screen display to "YES, RESET NVRAM".
- 2. The RUBBER & PLASTICS key will action your selection and return the unit to the service menu condition.
- 3. Pressing the STOP key will exit without changes and return the unit to the operational menu.

### >Change Password

This function will allow you to change the default password (Unwrapped, Wrapped, Rubber and Plastics, Stop) to a new key/button sequence.

When entering this sub menu, the unit will prompt you as follows:

• The unit will ask you to enter a new 4-key password.

```
Type New Password
```

- Enter the new password.
- The unit will then ask you to re-enter the new password.

```
Retype New Password
****
```

- Re-enter the password.
- The unit will confirm that the password has been changed, or if changing the password failed, the unit will again ask you to enter a new 4-key password.



### >Backup NVRAM



Screen display when in Backup NVRAM sub menu.

This function will allow the backup of any revised unit settings after calibration for example.

The default setting on all STATIM units is NO.

The alternative is YES.

To activate NVRAM backup, proceed as follows:

- 1. The UNWRAPPED or WRAPPED keys will both switch the screen display to "YES".
- 2. The RUBBER & PLASTICS key will action your selection and return the unit to the service menu condition.
- 3. Pressing the STOP key will exit without changes and return the unit to the operational menu.

### >Restore NVRAM

```
>Restore NVRAM
No
```

Screen display when in Restore NVRAM sub menu.

This function will allow the stored (backed up) unit settings to be restored if operational settings are lost for any reason.

The default setting on all STATIM units is NO.

The alternative is YES.

To activate NVRAM restoration, proceed as follows:

- 1. The UNWRAPPED or WRAPPED keys will both switch the screen display to "YES".
- 2. The RUBBER & PLASTICS key will action your selection and return the unit to the service menu condition.
- 3. Pressing the STOP key will exit without changes and return the unit to the operational menu.

### >Save and Exit

Allows any changes made in the sub menus (where the STOP button does not action save) to be saved and the unit will then exit to the operational menu.

To activate save and exit, press the RUBBER & PLASTICS button.

### >Exit

Allows you to exit the service menu WITHOUT saving changes other than those changes saved by activating the STOP button.

To activate exit, press the RUBBER & PLASTICS button.

### >Water Pump Type



Screen display when in Press. Sensor Type sub menu.

This sub menu is for information only.

This allows you to identify the type of water pump fitted to any revision 7 unit without removing the cover. **No changes can be made to the pump type in this sub menu.** 


There are two specific pump types identified by this sub menu, the ULKA pump (fitted to later machines) and the INVENSYS or SCICAN pump (fitted to earlier units).

The pumps are **NOT** interchangeable between units as the each pump requires type specific software for the unit to operate correctly.

- 1. To return to the service menu, press the RUBBER & PLASTICS button.
- 2. To return to the operational menu, press the STOP key.

### >Production cycle

This function is for manufacturing use only.

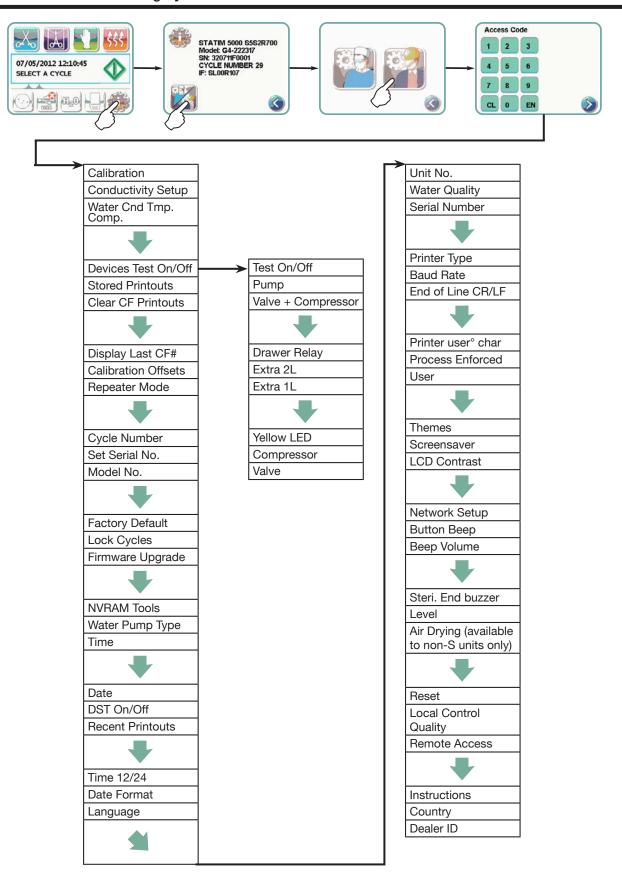

# **Troubleshooting Cycle Faults**

There are three types of cycle fault formats that can be seen on all STATIM 2000 and 5000 Classic and G4 units. They are as follows:

- 1. Cycle fault numbers and descriptions appropriate to **ONLY STAT/M 2000 non S units with revision 3.x and 4.x type controller boards.** 
  - Typically, these will show a number preceded by "Check Cassette".
  - The display on the LCD screen will NOT show any software revision and will display "SELECT A PROGRAM" when switched on.
- 2. Cycle fault numbers and descriptions appropriate to **ALL STAT/M 2000 and 5000 units with revision 2.x/5.x/6.x and 7.x type controller boards** with revision 1.xx/2.xx/4.xx/5.xx/6.xx and 7.xx software as appropriate.
  - These units will show a number preceded by "CF".
  - When the unit is switched on, a software revision number will be displayed for a short time on the top, right hand side of the LCD display (for example S2S2R415), depending on type (2000/5000), model (S or non S) and controller type.
- 3. Cycle fault **descriptions only** with no fault number. These are generic faults that do not necessarily occur as the result of a mechanical or electrical failure of the STAT*IM* unit but of a process error.

After identifying the unit you are trying to diagnose, refer to the following **troubleshooting guides** to assist in the diagnosis.

# Cycle Fault Numbers and descriptions on units with revision 3.x and 4.x type controller boards

CF.	Description of Error Message	Probable Cause of Error Message
number		
Check Cassette #1	The cassette temperature failed to reach 95°C in 3 minutes.	This error message occurs only during warm up. The most probable causes are: 1. Steam generator does not heat up. No power to steam generator. Blown thermal fuse. 2. Check for a large cassette leak or a large load.
Check Cassette #2	The cassette temperature failed to increase from 95°C to 100°C within 1 minute 20 seconds.	Most probable causes are a faulty cassette seal, miscalibration, a damaged cassette, improperly adjusted or failed microswitch allowing unit to function when cassette is not fully inserted. Occasionally, a low pressure weather system in higher altitude sites will alter the boiling temperature significantly and trigger this fault.
Check Cassette #3	The cassette has failed to pressurize and achieve a temperature of 110°C within 70 seconds of pressurization.	This may be caused by a faulty or worn cassette seal, a faulty or dirty solenoid valve, a damaged cassette or a faulty exhaust tube, leaky check valve or pressure relief valve.
Check Cassette #4	The cassette has failed to achieve sterilization conditions within 10 minutes of the chamber first reaching 102°C.	Several conditions may cause this: a damaged or worn cassette seal, a faulty solenoid valve, improperly installed copper tubing, cassette thermocouple misalignment or damage, a faulty steam generator pressure relief valve or a faulty check valve on.
Check Cassette #5	The software causes the pump to activate while between 2 minutes 44 seconds and 2 minutes 24 seconds remaining in the cycle. If a request to pump water occurs outside of 2 minutes 44 seconds, Check Cassette #5 occurs.	This indicates a substantial steam leak from the cassette, a contaminated steam generator or a weak pump.
Check Cassette #6	The steam generator temperature is more than 5°C higher than the chamber temperature.	Most probable causes are a blocked duct, a blocked solenoid valve, miscalibration or a pinched exhaust tube.
Check Cassette #7	The cassette temperature has dropped below 130.5°C during the sterilizing phase of the UNWRAPPED or WRAPPED Cycle, or below 117.5°C during the sterilizing phase of the RUBBER & PLASTIC Cycle.	Most probable cause is a faulty solenoid valve.
Check Cassette #8	The chamber temperature is more than 5°C higher than the steam generator temperature.	Most probable causes are thermocouple miscalibration or very strong pump delivery.
Check Cassette #9	This is a special error message for the U.K. market. The chamber temperature has exceeded 139°C during the UNWRAPPED/WRAPPED Cycle, or has exceeded 126°C during the RUBBER & PLASTIC Cycle.	Most probable causes are a blocked or faulty solenoid valve or a pinched exhaust tube.
Service Needed	The steam generator temperature has exceeded 170°C, chamber temperature has exceeded 147°C, or the thermocouple is faulty.	Most probable causes are a blocked or faulty solenoid valve, a pinched exhaust tube, a broken thermocouple lead, an open, disconnected or defective thermocouple, a weak pump or a contaminated steam generator.
Check Cassette (without a number)	The cassette temperature has failed to drop to 103°C within a timeout period at the end of a cycle.	Most probable causes are a blocked exhaust duct in the cassette, a blocked or faulty solenoid valve or a pinched / kinked exhaust tube.

# Cycle Fault Numbers and descriptions on units with revision 2.x/5.x/6.x and 7.x type controller boards (1.xx/2.xx/4.xx/5.xx/6.xx and 7.xx software)

### Tips for using this guide.

Unlike the previous guide (rev 3.x and 4.x controller boards) where all cycle fault numbers and descriptions are common to those specific controller board revisions irrespective of software, the following guide has been generated to cover all units from the introduction of controller board revision 2.x and have been standardized across the unit range of all STATIM units irrespective of controller board or software revision.

There will be some cases of course where because of the unit type, a particular fault will not be appropriate (e.g. pressure faults on non S units where a pressure transducer is not fitted) so establishing the following criteria before looking up a fault reference will assist in easier diagnosis.

Note also that the 'probable cause(s)' of any fault shown in the following chart are a GUIDE only. There may be multiple faults and the fault displayed may be either the first problem the software has identified or as the result of a different but related issue. For example, a blown thermal fuse is only symptomatic of a boiler overheat, and not necessarily the cause of the failure, which may be through poor water delivery, so it would be useful to undertake an assessment of the unit as described earlier in this chapter under: "Basic mechanical diagnostics."

Before referring to the following guide, establish the following where possible:

- STAT/M type (2000/5000)
- Model type ('S' or non 'S')
- Controller board revision (2.xx/5.xx etc.)
  - TIP: this is usually found on the bottom right hand corner of the controller board except for revision 7 board where it may be located top right or bottom left depending on manufacturing revision. If in doubt, refer to Chapter 1 of this manual.
- Software revision. (\$2\$2R415 for example, the R415 being the important part)
  - TIP: this is usually found on the main microprocessor label, but the easiest method is by powering up the unit. On all of the revisions refered to in this chart, this will appear on the top right hand corner of the LCD display for a few seconds after power up. If in doubt, refer to chapter 1 of this manual.
- Cycle Fault number.

Note that there may be more than one of the same cycle fault which may be the same, but have different causes depending on the derivative of the unit, so when you have established these details, go to the CF number and use the 'unit type' and 'software revision' guide to establish if this particular CF is relevant (/) or not (\*) as the case may be.

ב		1:4	Ü	Carried	Coffice couries	\$		
5 §	5 ₹	type	ร์	וואמום	201210		Description of Fault	Probable Cause (or causes) of Fault
	S	Non	1.xx/	4.xx	5.xx	/xx/9	•	
		တ	2.xx			7.xx		
	<u> </u>	,	,		,	,		Large cassette leak
_	<u>,</u>	<b>&gt;</b>	>	>	>	>	within a time-out period.	Overloaded cassette
								Poor pump delivery
								Blown thermal fuses
								Contaminated steam generator
								<ul> <li>Failed steam generator triac</li> </ul>
								<ul> <li>Steam generator fault (open circuit).</li> </ul>
							The cassette temperature failed to increase from	Major steam leak
7	>	<b>&gt;</b>	>	*	*	*	95°C to 100°C within a time-out period.	Overloaded cassette
I								<ul> <li>Incorrect chamber calibration</li> </ul>
							The cassette has failed to pressurize and	This may be caused by a
~	>	>	>	>	>	>	achieve a temperature of 110°C within a time-	<ul> <li>Faulty cassette seal</li> </ul>
)				,			out period.	Damaged cassette
								Steam circuit leak (tubing)
								<ul> <li>Incorrect probe alignment</li> </ul>
								<ul> <li>Sticking solenoid valve (open)</li> </ul>
								Check valve leak
								PRV leak
								Check Vref is 2.520V± 0.001V (all software revisions
								OTHER THAN revision 6.xx)
							The cassette has failed to achieve sterilization	This may be caused by a
4	<u>`</u>	>	>	>	>	>	conditions within a timeout period.	<ul> <li>Faulty cassette seal</li> </ul>
								Damaged cassette
								<ul> <li>Steam circuit leak (tubing)</li> </ul>
								<ul> <li>Incorrect probe alignment</li> </ul>
								<ul> <li>Sticking solenoid valve (open)</li> </ul>
								Check valve leak
								PRV leak
								Check Vref is 2.520V± 0.001V (all software revisions
								OTHER THAN revision 6.xx)

T.	=	+iul	ď	ftware	Software revision	2		
S	. ≥	type	3	3			Description of Fault	Probable Cause (or causes) of Fault
	S	Non	1 vv/	4 vv	7 × ×	/xx y		
		တ	2.xx	ř	<b>Y</b>	7.xx		
							The software causes the pump to activate at	Leaky cassette
2	>	<b>&gt;</b>	1.xx	*	*	*	predetermined times. If a request to pump water	<ul> <li>Kinked pump tube</li> </ul>
)			7				occurs outside of the predetermined time Cycle	<ul> <li>Leaking solenoid valve</li> </ul>
							Fault #5 occurs.	<ul> <li>Poor pump delivery</li> </ul>
								Overloaded cassette
							The software has detected the filtered steam	Obstruction caused by:
Ç	>	>	2 xx	>	*	*	generator temperature 6°C greater than the	<ul> <li>Seal misalignment</li> </ul>
)			2				chamber during the sterilizing phase of a cycle.	<ul> <li>Foreign material in the seal</li> </ul>
:								channel
Ξ								<ul> <li>Calibration imbalance (boiler</li> </ul>
								thermocouple to chamber
								thermocouple)
								<ul> <li>Contaminated steam generator</li> </ul>
							The software has detected the Validation	Cassette steam leak
9	*	<b>&gt;</b>	*	*	>	>	thermocouple temperature to be 5°C greater	<ul> <li>Sticking solenoid</li> </ul>
							than the chamber during the sterilizing phase of	<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
(11)							a cycle.	<ul> <li>Calibration imbalance (validation</li> </ul>
								thermocouple to chamber
								thermocouple)
								Check Vref is 2.520V± 0.001V (all software
								revisions OTHER THAN revision 6.xx)
							During the sterilization phase of the cycle the	<ul> <li>Sticking solenoid valve (open)</li> </ul>
7	*	<b>&gt;</b>	<b>&gt;</b>	>	>	>	cassette temperature has dropped below a	<ul> <li>Cassette leak</li> </ul>
							threshold value.	<ul> <li>Check valve leak</li> </ul>
								PRV leak

H.		IInit	Ů.	Software revision	revisio	2		
Š.	. ₹	type					Description of Fault	Probable Cause (or causes) of Fault
	S	Non S	1.xx/ 2.xx	4.xx	2.xx	6.xx/ 7.xx		
8	*	*	2.xx	>	*	*	The software has detected the filtered steam generator temperature 6°C less than the chamber during the sterilizing phase of a cycle.	<ul> <li>Contaminated steam generator</li> <li>Failed steam generator</li> <li>Calibration imbalance (hoilar to</li> </ul>
Ξ		_	oniy					chamber)  High pump flow
8	*	>	*	*	>	>	The software has detected the Validation Thermocouple temperature to be 5°C less than	Calibration imbalance (validation thermocouple to chamber
•							the chamber temperature during the sterilizing phase of a cycle.	thermocouple)  High pump flow
<b>E</b>								<ul> <li>Blocked cassette exhaust duct</li> <li>Sticking solenoid valve</li> <li>Check Vref is 2.5200± 0.001V (all software</li> </ul>
6							Not used	TEVISIONS OF THEN THAN TEVISION O.XX)
5	1	,	,	\	`	\	During conditioning (purging), the cassette temperature has failed to drop to115°C during	Blocked cassette exhaust duct     Obstructed exhaust tribing (kink)
2	<b>&gt;</b>	•	•	<b>&gt;</b>	<b>•</b>	<u> </u>	the Unwrapped or Wrapped Cycle OR the temperature has failed to drop to 110°C during	Sticking solenoid valve (closed)     Failed solenoid coil
							the Rubber and Plastics Cycle.	
7	>	`	`	>	>	`	The cassette temperature has failed to drop to 102°C within a timeout period of the end of a	Blocked cassette exhaust duct     Obstructed exhaust tubing (kink)
<b>:</b>	•	•	•	•	•	•	cycle during venting.	Sticking solenoid valve (closed)     Failed solenoid coil
							This indicates a problem with the temperature	Broken or faulty thermocouple
12	>	>	>	>	>	<b>&gt;</b>	measuring system.	Circuit failure on PCB
13							Not used	
4	*	>	>	*	*	*	The steam generator temperature is above 171°C during the sterilization phase of a cycle.	<ul> <li>Poor pump delivery</li> <li>Calibration imbalance (boiler thermocoling to chamber</li> </ul>
								thermocouple)  • Failed water sensor (level)

CF	j	Unit	So	Software revision	revisio	u.		
Š.	tyk	type					Description of Fault	Probable Cause (or causes) of Fault
	S	Non S	1.xx/ 2.xx	4.xx	2.xx	6.xx/ 7.xx		
							The cassette temperature rose above the high	Blocked cassette exhaust duct
15	<u> </u>	<u> </u>	<b>,</b>	>	>	>	threshold during the sterilization phase of a	<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
							cycle or above 138.5°C dufing conditioning or pressurizing phase of the cycle	Sticking solenoid valve (closed)     Ealed colonoid coil
	-						The cheek accordant references meets of	Diologian filters
,	•	•	•	•	•	•	the steam generator temperature went above a	Blocked pump filters
16	<u>,</u>	>	<b>&gt;</b>	>	>	>	unestiola value. (overneat)	Weak pump
			_					<ul> <li>Blocked or kinked pump tube</li> </ul>
			_					<ul> <li>Failed pump</li> </ul>
								<ul> <li>Failed pump triac</li> </ul>
			_					<ul> <li>Failed steam generator triac</li> </ul>
								<ul> <li>Sticking solenoid valve (closed)</li> </ul>
								Failed solenoid coil
								<ul> <li>Faulty seal installation</li> </ul>
17 - 1	18						Not used	
							For software rev. R1xx, R2xx and R4xx the	This occurs when a new controller board
5	5	`	`	1	1	>	steam generator calibration is invalid.	or microprocessor has been installed.
2		<u> </u>	•	<b>&gt;</b>	<b>&gt;</b>	<b>•</b>		This may also happen when the unit has
			_				For software rev. R5xx and R6xx the validation	been subjected to a strong static
							thermocouple calibration is invalid.	discharge corrupting the memory.
			_					In both cases a new steam generator or
								validation thermocouple calibration is
								required.
							The pump has failed to pump water into the	<ul> <li>Blocked pump filters</li> </ul>
20	*	<u> </u>	<b>&gt;</b>	4	4	*	steam generator during a pre-vent pump time-	Weak pump
							out. The steam generator temperature was	<ul> <li>Blocked or kinked pump tube</li> </ul>
			_				greater than 140°C for 3.6 seconds after the	Failed pump
							pump was activated to pump water to cool the	Failed pump triac
							steam generator.	<ul> <li>Failed steam generator triac</li> </ul>
								<ul> <li>Sticking solenoid valve (closed)</li> </ul>
								<ul> <li>Failed solenoid coil</li> </ul>
			_					<ul> <li>Faulty seal installation</li> </ul>
								<ul> <li>Failed water sensor (level)</li> </ul>
71 - 7	24						Not used	

ב	=	1 2 1	ò	0,01179	o io ivio	\$		
S	<b>₹</b>	type	วั				Description of Fault	Probable Cause (or causes) of Fault
	S	Non	1.xx/ 2.xx	4.xx	5.xx	6.xx/ 7.xx		
25	•	,	,	,	_,	_,	The software has failed to detect a need to	±
	>	>	>	>	>	>	purify water within 30 seconds of the start of the cycle.	Purify delivery     Steam generator triac
								0)
							The sterilization phase has failed to start within 3	Calibration imbalance (boiler
26	*	>	>	>	*	*	minutes of the cassette reaching the sterilization	thermocouple to chamber
)							temperature. Note: for rev 4.xx software only,	thermocouple)
:							the unit has to fail on 3 consecutive cycles	<ul> <li>Poor pump delivery</li> </ul>
Ξ							for CF26 to display.	<ul> <li>Sticking solenoid valve</li> </ul>
								<ul> <li>Check Vref is 2.520V± 0.001V</li> </ul>
							The sterilization phase has failed to start within 3	Calibration imbalance (validation)
90	\$	`	5	s	`	`	minutes of the cassette reaching the sterilization	thermocouple to chamber
9	<b>,</b>	•	<b>B</b>	}	•	•	temperature. Note: the unit has to fail on 3	thermocoluple)
							consecutive cycles for CE26 to display	
€							"Cycle Interrinted" displayed for the first 2	
`							Oyde michapied displayed to me mark	Sticking solenoid valve
							occurrences.	<ul> <li>Check Vref is 2.520V± 0.001V</li> </ul>
								(rev 6.xx ONLY)
							The temperature of the steam generator failed to	<ul> <li>Poor pump delivery</li> </ul>
27	>	>	>	>	>	>	drop bellow a set-point temperature (150°C or	<ul> <li>Failed water sensor (level)</li> </ul>
i							165°C) in a timeout period or during a panic	<ul> <li>Contaminated steam generator</li> </ul>
							pump situation.	
							The cassette pressure rose above a ceiling	<ul> <li>Pressure measurement failure.</li> </ul>
28	>	*	>	>	<b>&gt;</b>	>	value.	<ul> <li>Calibration imbalance (pressure</li> </ul>
								transducer to chamber thermocouple)
								<ul> <li>Sticking solenoid valve (closed)</li> </ul>
								<ul> <li>Failed solenoid coil</li> </ul>
								<ul> <li>Blocked cassette exhaust duct</li> </ul>
								<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
								<ul> <li>Blockage or restriction in pressure transducer tube</li> </ul>
29	1						Not used	
2								

H.C	Unit		S.	Software revision	revision	2		
	type	: 4	}	3			Description of Fault	Probable Cause (or causes) of Fault
1	?⊢	2	1					i opanic caase (ci caases) oi i aaii
	2 2	Non S	1.xx/ 2.xx	4.xx	5.XX	6.xx/ 7.xx		
	*	4.5	>	>	>	>	The cassette temperature failed to reach the sterilization temperature within 15 seconds of	<ul> <li>Poor air removal during conditioning (exhaust blockage).</li> </ul>
							the theoretical cassette temperature calculated from the measured cassette pressure reaching	Calibration imbalance (pressure transducer or chamber thermocounts)
							the sterilization temperature.	Blockage or restriction in pressure transducer tube
4	49						Not used	
							For the Rubber and Plastics cycle, the chamber	Poor air removal during conditioning
•	<u>*</u>	4	<b>&gt;</b>	>	>	>	temperature dropped below the sterilization	(exhaust blockage).
							(i.e. Tchm < 121°C).	<ul> <li>Not able to generate steam</li> <li>Steam leak in the system</li> </ul>
							For the Rubber and Plastics cycle, the chamber	Sticking solenoid valve (closed)
	*	٠.	<b>&gt;</b>	>	<b>&gt;</b>	>	temperature rose more than 4 C° above the	<ul> <li>Failed solenoid coil</li> </ul>
							sterilization temperature, allowing for	<ul> <li>Blocked cassette exhaust duct</li> </ul>
							measurement error (i.e. Tchm > 125°C).	<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
							For the Rubber and Plastics cycle, the	<ul> <li>Calibration imbalance (pressure</li> </ul>
-	*	44	<b>&gt;</b>	>	>	>	theoretical chamber temperature as calculated	transducer or chamber thermocouple)
							from the measured chamber pressure was more	<ul> <li>Sticking solenoid valve (closed)</li> </ul>
							than 2°C below the measured chamber	<ul> <li>Failed solenoid coil</li> </ul>
							temperature, allowing for measurement error.	<ul> <li>Blocked cassette exhaust duct</li> </ul>
								<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
								Poor air removal
1							For the Rubber and Plastics cycle, the	Calibration imbalance (pressure
	*	44	<b>&gt;</b>	>	>	>	theoretical chamber temperature as calculated	transducer or chamber thermocouple)
							from the measured chamber pressure was more	<ul> <li>Sticking solenoid valve (closed)</li> </ul>
							than 2°C above the measured chamber	<ul> <li>Failed solenoid coil</li> </ul>
							temperature, allowing for measurement error.	<ul> <li>Blocked cassette exhaust duct</li> </ul>
								<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
								<ul> <li>Poor air removal</li> </ul>
$\dashv$								

CF	_	Unit	So	Software revision	revisio	Ē			
No.	t)	type					Description of Fault	Prob	Probable Cause (or causes) of Fault
	တ	Non S	1.xx/ 2.xx	4.xx	5.xx	6.xx/ 7.xx			
							For the Rubber and Plastics cycle, the	• Ca	Calibration imbalance (pressure
54	>	*	>	>	>	>	theoretical cassette temperature calculated from	tra	transducer or chamber thermocouple)
•	,			,			the measured cassette pressure was below the	• Po	Poor air removal during conditioning
_							sterilization temperature, allowing for	(ex	(exhaust blockage).
_							measurement error (i.e. Pchm < 204.8 kPa).	•	Not able to generate steam
_								• Ste	Steam leak in the system
							For the Rubber and Plastics cycle, the	• Ca	Calibration imbalance (pressure
55	>	*	>	>	>	>	theoretical cassette temperature calculated from	tra	ransducer or chamber thermocouple)
)	,					,	the measured cassette pressure was more than	• Sti	Sticking solenoid valve (closed)
_							4°C above the sterilization temperature, allowing	• Fa	Failed solenoid coil
_							for measurement error (i.e. Pchm > 232 kPa).	• BIC	Blocked cassette exhaust duct
_								8 •	Obstructed exhaust tubing (kink)
_								• Po	Poor air removal
							For a 134°C cycle, the chamber temperature	• Po	Poor air removal during conditioning
9	>	*	>	>	>	>	dropped below the sterilization temperature,	(e)	(exhaust blockage).
)							allowing for measurement error (i.e. Tchm <	• N	Not able to generate steam
							134°C).	• Ste	Steam leak in the system
							For a 134°C cycle, the chamber temperature	• Sti	Sticking solenoid valve (closed)
61	>	*	>	>	>	>	rose more than 4 C° above the sterilization	• Fa	Failed solenoid coil
,							temperature, allowing for measurement error	• BK	Blocked cassette exhaust duct
_							(i.e. Tchm > 138°C).	о •	Obstructed exhaust tubing (kink)
							For a 134°C cycle, the theoretical chamber	• Ca	Calibration imbalance (pressure
62	>	*	>	>	>	>	temperature as calculated from the measured	tra	transducer or chamber thermocouple)
   					,	,	chamber pressure was more than 2°C below the	• Sti	Sticking solenoid valve (closed)
_							measured chamber temperature, allowing for	• Fa	Failed solenoid coil
_							measurement error.	• BIC	Blocked cassette exhaust duct
_								о •	Obstructed exhaust tubing (kink)
								• Po	Poor air removal

S CF	⊃ <b></b>	Unit	So	ftware	Software revision	u	Description of Fault	Probable Cause (seasons) of Eauth
į	<u>`</u> -	20/						riopapie cause (oi causes) oi i auit
	ဟ	Non S	1.xx/ 2.xx	4.xx	5.xx	6.xx/ 7.xx		
							For a 134°C cycle, the theoretical chamber	Calibration imbalance (pressure
63	>	*	>	>	>	>	temperature as calculated from the measured	transducer or chamber thermocouple)
)							chamber pressure was more than 2°C above the	<ul> <li>Sticking solenoid valve (closed)</li> </ul>
							measured chamber temperature, allowing for	<ul> <li>Failed solenoid coil</li> </ul>
							measurement error.	<ul> <li>Blocked cassette exhaust duct</li> </ul>
								<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
								<ul> <li>Poor air removal</li> </ul>
							For a 134°C the theoretical cassette	Calibration imbalance (pressure
64	>	*	<b>&gt;</b>	>	>	>	temperature calculated from the measured	transducer or chamber thermocouple)
,							cassette pressure was below the sterilization	<ul> <li>Poor air removal during conditioning</li> </ul>
							temperature, allowing for measurement error	(exhaust blockage).
							(i.e. Pchm < 304 kPa).	<ul> <li>Not able to generate steam</li> </ul>
								Steam leak in the system
							For a 134°C cycle, the theoretical cassette	<ul> <li>Calibration imbalance (pressure</li> </ul>
65	>	*	<b>&gt;</b>	>	>	>	temperature calculated from the measured	transducer or chamber thermocouple)
)							cassette pressure was more than 4°C above the	<ul> <li>Sticking solenoid valve (closed)</li> </ul>
							sterilization temperature, allowing for	<ul> <li>Failed solenoid coil</li> </ul>
							measurement error (i.e. Pchm > 341.2 kPa).	<ul> <li>Blocked cassette exhaust duct</li> </ul>
								<ul> <li>Obstructed exhaust tubing (kink)</li> </ul>
								<ul> <li>Poor air removal</li> </ul>
- 99	69						Not used	
							The time maintained by the internal timer of the	Transient electromagnetic disturbance
70	>	4	<b>&gt;</b>	>	>	>	processor didn't match the time maintained by	(problem won't repeat).
							the external real-time clock.	<ul> <li>Damaged microprocessor or crystal</li> </ul>
								(STATIM Controller Board).
								<ul> <li>Damaged real-time clock or crystal</li> </ul>
							Pressure reading is outside the possible range	<ul> <li>Misconnected, disconnected or</li> </ul>
71	>	*	>	>	>	>		damaged pressure sensor.
								Disconnected or damaged Pressure     Secondary (25, 25, 25, 25, 25, 25, 25, 25, 25, 25,
	1							Sellsol Iliteriace Boald (Ilot lev 7.XX).

ב	=	121	3	Cofficient Carolina Co	O. O. O.	2		
ş Ş	ξ <u>ς</u>	type	กั	JI Wale	201210		Description of Fault	Probable Cause (or causes) of Fault
	S	Non	1.xx/ 2.xx	4.xx	2.xx	6.xx/ 7.xx		
22	*	>	*	*	>	*	There is a communication error between the microprocessor and the Temperature Adapter Board (TAB).	Check the connection between the square microprocessor and the Printed Circuit Board (PCB).
							(Note: Verify that the capacitor C12 was removed from the main PCB prior to installation	If an adapter board is used check the connection between the adapter board and the PCB. (I ook for hent or broken)
							of the TAB when used in conjunction with the PCB Adapter Controller Rev 2.x, 5.x).	pins.)  Verify that Temperature Adapter PCB (TAB) is properly inserted and secured
08	<b>&gt;</b>	*	*	*	*	*	The boiler calibration offset routine is Being executed at the wrong time. This routine should only run when in Boiler calibration mode.	<ul> <li>This may occur if the microprocessor operates erratically, due to electrical noise.</li> </ul>
20	>	*	>	*	*	*	The state buffer pointer (print data storage buffer) is invalid.	Corruption of pointer is usually due to electrical noise.
82	<b>&gt;</b>	*	*	*	4	4	The calibration buffer pointer (blr-chm data storage buffer) is invalid.	Corruption of pointer is usually due to electrical noise.
83	>	*	<b>&gt;</b>	*	*	4	The EEPROM SIGNATURE string is corrupted. This 16 byte string is located at 0-\$0F in the EEPROM. It should not change unless the EEPROM gets corrupted.	The EEPROM can be corrupted if electrical noise corrupts pointers and counters in RAM, or causes other erratic operation of the microprocessor.
84	>	*	*	*	*	*	The boiler chamber value just calculated is not in the range 8C to FF.	<ul> <li>Bad or miscalibrated temperature measurement circuitry.</li> <li>Thermocouple defective or out of</li> </ul>
								<ul> <li>tolerance.</li> <li>The boiler is old or not conditioned properly.</li> <li>Calibration routine is running at the wrong time due to noise.</li> </ul>
82	>	*	>	*	*	*	Unexpected Serial Interrupt	Interference (electrical noise)

CF	Unit	ıit	So	ftware	Software revisio	on		
<u>.</u>	₽	9e			-	-	Description of Fault	Probable Cause (or causes) of Fault
	<u>~</u> ທ	Non S	1.xx/ 2.xx	4.xx	2.xx	6.xx/ 7.xx		
		)					Unexpected SPI Interrupt.	Interference (electrical noise)
98	<u> </u>	*	<b>&gt;</b>	4	*	*	-	
							Unexpected IRQ Interrupt	Interference (electrical noise)
87	<b>&gt;</b>	*	>	*	*	*		
88	<u> </u>	*	<b>&gt;</b>	4	*	4	The buffer pointers and sample counter are initialized incorrectly.	<ul> <li>This occurs if the boiler calibration routine is executed when it was not expected to, (caused by electrical noise or ESD.)</li> </ul>
68	>	4	>	4	*	*	The software state machine has detected incorrect input parameters.	Usually caused by noise and ESD.
-							Corrupted or not initialized chamber calibration	This occurs when a new controller board
06	<u>`</u>	<u>`</u>	4	*	*	>	value	or microprocessor has been installed.
								<ul> <li>I nis may also happen when the unit has been subjected to a strong static</li> </ul>
								discharge corrupting the memory.  • Re-calibrate the chamber
	-						-	thermocouple.
	<u> </u>	*	4	*	*	>	Corrupted or not initialized pressure calibration	<ul> <li>This occurs when a new controller board or microprocessor has been installed.</li> </ul>
								<ul> <li>This may also happen when the unit has</li> </ul>
								been subjected to a strong static
								<ul><li>discharge corrupting the memory.</li><li>Re-calibrate the pressure sensor.</li></ul>
2	-						Corrupted or not initialized EEPROM	<ul> <li>Failure to read EEPROM messages</li> </ul>
	<u>;</u>	`	*	*	*	>	(G4 units only)	from 24C512 EEPROM (EEPROM
Ļ	+				-		(ylan officer ha) some softonice manage	The TOP contains and the second secon
		•	1	4	4	•	Communication error (G4 units only)	The LCD controller PCB failed to receive     massage from the Doy 7 DCB
	<u>-</u>	<b>&gt;</b>	*	*	*	>		messages from the Rev / PCB.
								<ul> <li>No connection between Colour LCD Controller and Rev 7PCB;</li> </ul>
								Check cable / connection.
								<ul> <li>Verify that MCU is properly inserted into the socket on the Rev 7 PCB.</li> </ul>
								<ul> <li>Verify that 5V is available on the Rev 7 PCB (Green LED ON)</li> </ul>
86							Microcontroller fails to communicate with ADC	Hardware failure.
	`	<b>&gt;</b>	*	<b>&gt;</b>	>	>	(Analog to Digital Converter)	<ul> <li>MCU not inserted properly in the socket</li> <li>Damaged ADC converter, replace PCB.</li> <li>If a PCB adapter is used, verify that PCB adapter is properly secured in the socket and that there are no bent pins</li> </ul>
							To contract of the contract of	County of the society
66	<u> </u>	*	`	4	*	*	timer processor's internal waterladge timer has tripped.	usually noise or ESD related.
l		l						

Cycle Fault	Description and Probable Cause of Fault
"NO CONFIGURATION EEPROM"	<ul> <li>Lack of communication between microcontroller and EEPROM</li> <li>Hardware failure.</li> <li>MCU not inserted properly in the socket, extract MCU, and verify pins are shorted.</li> <li>Damaged wrongly inserted or missing EEPROM, replace software kit.</li> <li>Damaged controller PCB.</li> <li>If a PCB Adapter is used, verify that PCB adapter is properly secured in the socket and that there are no bent pins.</li> </ul>
Message "PRINTER FAULT" (If optional printer is installed)	Printer is not printing.  This may be caused by a paper jam, a failed printer or failed printer electronics. Check that power is ON, check connector attachments. (See also section 12 of this manual)
"CYCLE ABORTED"	This error message is displayed on the printout only, followed by the message "NOT STERILE", as a result of the operator pressing the STOP button to stop the cycle or as a result of any other abnormal cycle termination, including CYCLE FAULT errors.
"CYCLE INTERRUPTED"	This error is displayed on power-up following a power failure occurs during a cycle or whenever the power is turned OFF after an error occurred without pressing the STOP button to reset. (Rev. R200 level code)  This message is displayed when the sterilization phase has failed to start within three minutes of the cassette reaching the sterilization temperature. If it occurs in three consecutive cycles Cycle Fault #26 is displayed. (Software Rev. R4xx, R5xx and R6xx)
"PRESS STOP TO RESET"	This message is displayed for all error faults.  The user MUST press the STOP button on the keypad to reset the unit; otherwise the user will be unable to initiate another cycle.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



## **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

## PERFORM TESTS

• If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



## PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

## STATIM 2000 Calibration

#### Calibration Procedure 1: for Rev 3.x or 4.x Boards

This calibration procedure is for STATIM 2000 non S units (1992-1995) with revision 3.x or 4.x controller boards

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- Do not place the digital reference thermometer in direct sunlight, on the armature or on any other hot surfaces. This may cause improper values on the temperature readout.
- There are specific error messages (cycle faults) for Rev 3.x and Rev 4.x board types. These can be found in the *Troubleshooting Cycle Faults* section.
- Typically, these boards were manufactured with four potentiometers (VREF, DISP, TC-CHM and TC-BLR) and do not have a printer connector.

#### Calibration Procedure


Calibration of these boards should be done in the following order:

- 1. Reference voltage verification / adjustment
- 2. Steam generator thermocouple calibration
- 3. Chamber thermocouple calibration
- 4. Steam generator re-calibration

#### Required equipment:

- Control box 01-103141S
- Digital voltmeter
- Calibration cassette (appropriate to model)
- Potentiometer trimmer
- Digital thermometer


- Before starting calibration, turn the STATIM OFF and remove the cover.
- Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette (6) is correctly engaged in the STAT*IM*.
- Install the thermocouple (7) into the cassette and attach the digital thermometer (8).
- Ensure that the thermocouple and digital thermometer have matched serial numbers.

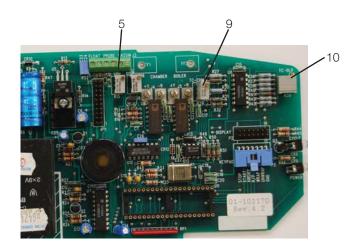

- 1. Control box
- 2. Controller board
- 3. SELECT knob
- 4. Digital voltmeter
- 5. VRH potentiometer
- 6. Calibration cassette
- 7. Thermocouple

- Triefficouple
   Digital thermometer
   TC-CHM potentiometer
   TC-BLR potentiometer

Figure 1

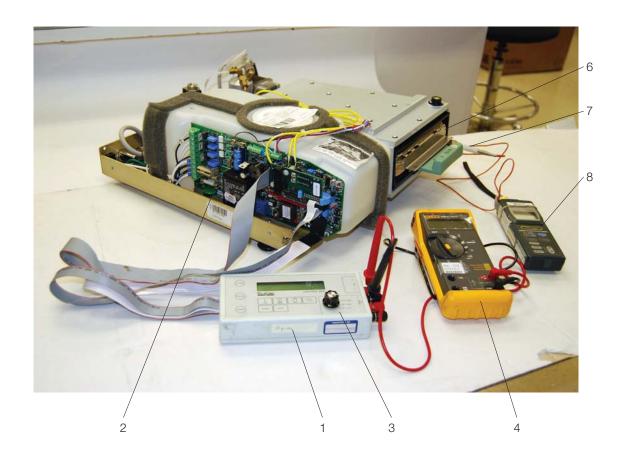

#### 1. Reference Voltage Verification / Adjustment

To verify / adjust the reference voltage:

- Ensure the unit is OFF.
- 2. Connect the Control Box (1) to the Controller Board (2).
- 3. Connect the digital voltmeter (4) to the VREF terminals on the control box.
- 4. Set the voltmeter to read DC Volts with a resolution of 1mV.
- Power the unit ON.
- 6. Adjust the VRH potentiometer (5) until VREF is  $4.000 \text{ V} \pm 0.001$ .
- 7. Once VREF is adjusted, apply a drop of non-conducting lacquer or nail polish to the adjustment screw of the potentiometer.

#### 2. Steam generator thermocouple calibration

To calibrate the steam generator thermocouple:

- 1. Turn the control box SELECT knob to blr-chm. This signal is the difference between the chamber and steam generator thermocouples.
- 2. Set the voltmeter (4) to read DC Volts with a resolution of 1 mV and locate the TC-BLR (10) potentiometer on the controller board.
- 3. Select the UNWRAPPED INSTRUMENT CYCLE button on the control box, then press START.
- 4. Follow the phases of the cycle as they are displayed on the LCD (warming up, purging air, pressurization, sterilization). Calibration occurs during the sterilization phase.
- 5. With 2 minutes and 25 seconds left in the sterilization cycle adjust the TC-BLR potentiometer until the voltmeter reads  $0.000 \text{ V} \pm 0.030 \text{ V}$ .

#### 2. Chamber thermocouple calibration

**Tip:** There are no quick ways to calibrate a unit. Calibration requires a number of readings and adjustments to achieve the required values. Wait for equipment readings to stabilize and record them as quickly and precisely as possible.

To calibrate the chamber thermocouple:

- 1. Install the thermocouple (7) into the calibration cassette (6) and attach the digital thermometer (8). Ensure that the thermocouple and digital thermometer have matched serial numbers.
- 2. Turn the control box (1) SELECT knob (3) to chm X gain.
- 3. Locate the TC-CHM potentiometer (9) on the controller board (2).
- 4. Set the meter (5) to read DC Volts with a resolution of 1 mV.
- 5. Select the WRAPPED INSTRUMENT CYCLE button on the control box and release
- 6. Press START.
- 7. Follow the phases of the cycle as they are displayed on the LCD (warming up, purging air, pressurization, sterilization). Calibration occurs during the sterilization phase.
- 8. After approximately 2 minutes into the sterilization phase look for the digital thermometer (8) reading to remain constant for at least 2 seconds and record the values of the digital thermometer (TCHM) and digital voltmeter (VM).

- 9. Using the TCHM value just recorded, find the corresponding VCHMXGAIN value from the chamber calibration conversion chart below.
- 10. Adjust the TCHM potentiometer until VM matches VCHMXGAIN for the observed TCHM to within  $\pm$  0.005 V.
- 11. Start another cycle and verify that the calibration readings are within acceptable values. If not, repeat steps 4 through 7 until acceptable values are achieved.

Chamber Ca	alibration Con	version - Cha	rt A				
Tchm (°C)	Vchmxgain (VOLTS)	Tchm (°C)	Vchmxgain (VOLTS)	Tchm (°C)	Vchmxgain (VOLTS)	Tchm (°C)	Vchmxgain (VOLTS)
134.0	2.259	135.3	2.300	136.6	2.340	137.9	2.381
134.1	2.262	135.4	2.303	136.7	2.343	138.0	2.384
134.2	2.265	135.5	2.306	136.8	2.347	138.1	2.387
134.3	2.268	135.6	2.309	136.9	2.350	138.2	2.391
134.4	2.271	135.7	2.312	137.0	2.353	138.3	2.394
134.5	2.274	135.8	2.315	137.1	2.356	138.4	2.397
134.6	2.278	135.9	2.318	137.2	2.359	138.5	2.400
134.7	2.281	136.0	2.322	137.3	2.362	138.6	2.403
134.8	2.284	136.1	2.325	137.4	2.365	138.7	2.406
134.9	2.287	136.2	2.328	137.5	2.369	138.8	2.409
135.0	2.290	136.3	2.331	137.6	2.372	138.9	2.413
135.1	2.293	136.4	2.334	137.7	2.375	139.0	2.416
135.2	2.296	136.5	2.337	137.8	2.378		

#### 4. Steam generator thermocouple re-calibration

- 1. Repeat the steam generator thermocouple calibration procedure in 2 above.
- 2. Start another cycle and verify that the calibration readings are within acceptable values. If the readings are not within these values repeat the procedure until the proper voltmeter reading is obtained.
- 3. Once adjusted, apply a drop of non-conducting lacquer or nail polish to the adjustment screw of the TC-ADJ potentiometer.
- 5. Switch off the unit and remove all calibration equipment.

### **Calibration Procedure 2:**

### for Rev 2.x or 5.x Boards, Software R1.xx or R2.xx

This calibration procedure is for STAT*IM* 2000/5000 non S units (1995-2005) with revision 2.x or 5.x controller boards and revision 1.xx or 2.xx software.

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- Do not place the digital reference thermometer in direct sunlight, on the armature or on any other hot surfaces. This may cause improper values on the temperature readout.
- There are specific error messages (cycle faults) for these board/software types, which can be found in the *Troubleshooting Cycle Faults* section.
- In calibration mode, no error messages are displayed.

#### **Calibration Procedure**


#### Calibration of these boards should be done in the following order:

- 1. Reference voltage verification / adjustment
- 2. Steam generator thermocouple calibration (automatic)
- 3. Chamber thermocouple calibration
- 4. Steam generator re-calibration (automatic)

#### Required equipment:


- Control box 01-103141S
- Digital voltmeter
- Calibration cassette (appropriate to model)
- Potentiometer trimmer
- Calibration jumper
- Digital thermometer

- Before starting calibration, turn the STATIM OFF and remove the cover.
- Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette (6) is correctly engaged in the STAT*IM*.
- Install the thermocouple (7) into the cassette and attach the digital thermometer (8).
- Ensure that the thermocouple and digital thermometer have matched serial numbers.
- Install the calibration jumper (11) to controller board pins marked W1 (10).

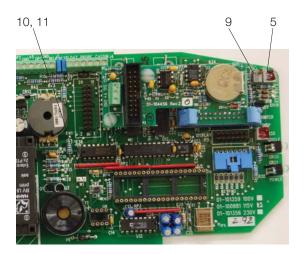

- 1. Control box
- 2. Controller board
- 3. SELECT knob
- 4. Digital voltmeter
- 5. Vref potentiometer
- 6. Calibration cassette
- 7. Thermocouple
- 8. Digital thermometer
- 9. TC-ADJ potentiometer
- 10. Calibration jumper connection
- 11. Calibration jumper (not shown)

Figure 1

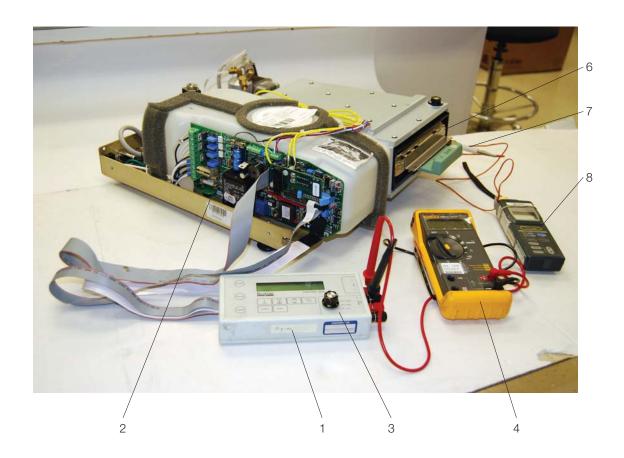

#### 1. Reference voltage verification / adjustment


To verify / adjust the reference voltage:

- 1. Turn the unit OFF.
- 2. Connect the Control Box (1) to the Controller Board (2).
- 3. Connect the digital voltmeter (4) to the Vref terminals on the control box.
- 4. Set the voltmeter to read DC Volts with a resolution of 1mV.
- 5. Power the unit ON.
- 6. Adjust the Vref potentiometer until Vref is 2.520 V  $\pm$  0.001.
- 7. Once Vref is adjusted, apply a drop of non-conducting lacquer or nail polish to the adjustment screw of the potentiometer.

#### 2. Steam generator thermocouple calibration

To calibrate the steam generator thermocouple:

- 1. Turn the unit OFF and install a calibration jumper on Controller Board header W1.
- 2. Turn the unit ON. The LCD should appear and be similar to the diagram below:

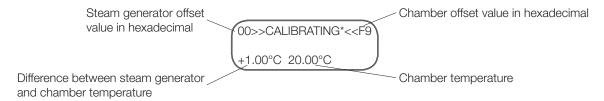

- 3. To start a steam generator (self-calibrating) cycle, **press and hold** the UNWRAPPED cycle button on the control box and press the **START** button. An asterisk \* will appear next to the message 'CALIBRATING'.
- 4. If the asterisk \* does not appear, press the STOP button on the control box twice to reset the unit. Turn the power switch OFF then ON again.
- 5. Hold down the UNWRAPPED cycle button on the control box and press the **START** button to start the self-calibration again.
- 6. During the steam generator warm up (before the cassette temperature reaches 100°C), check that the displayed chamber temperature matches the reading on the digital thermometer to within 2°C.
- 7. If the value is greater than 2°C adjust the TC-ADJ potentiometer (9) so that they match.
- 8. As the cycle proceeds, the chamber temperature reaches the sterilization temperature, drops to 115°C and then regains sterilization temperature.
- 9. After 20-40 seconds, a long beep will sound indicating that steam generator calibration is complete.
- The number in the upper left-hand corner of the LCD has changed to display the new steam generator offset value.



- 11. Press the STOP button on the control box to end the steam generator calibration cycle and to depressurize the cassette.
- 12. Press STOP a second time to reset the unit.

#### 3. Chamber Thermocouple Calibration

To calibrate the chamber thermocouple:

1. Turn the power switch to **OFF**. Ensure the calibration jumper is in place and turn the power switch **ON**. The LCD will read (for example):



- 2. To start a chamber calibration cycle, **press and release** the UNWRAPPED cycle button on the control box and then press the **START** button. The STAT*IM* will run a regular UNWRAPPED cycle, but continue to show calibration information on the LCD.
- 3. Wait for the chamber to reach the sterilization temperature of 134 °C.
- 4. The temperature displayed on the LCD needs to match the temperature displayed on the digital thermometer.
- 5. To change the temperature display on the LCD, adjust the TC-ADJ potentiometer (9) on the Controller Board. Continue adjusting TC-ADJ until the temperatures are within ±0.25°C.
- 6. Observe the temperatures displayed on the LCD and the digital thermometer for 30 seconds. The temperatures registered should remain within  $\pm$  0.25°C of one another. Adjust TC-ADJ if not.
- 7. Press the **STOP** button on the control box when finished. Wait for the beeps to signal that the cassette may be removed.

#### 4. Steam Generator Thermocouple Re-calibration

Repeat the steam generator thermocouple calibration procedure in section 2, *Steam Generator Thermocouple Calibration*, above.

5. Switch off the unit and remove all calibration equipment.

## Calibration Procedure 3: for Rev. 6.x Board, Software Rev. 4.xx and 5.xx

This calibration procedure is for STAT*IM* 2000/5000 non S units (2005-2007) with revision 6.x controller boards and revision 4.xx or 5.xx software.

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- Do not place the digital reference thermometer in direct sunlight, on the armature or on any other hot surfaces. This may cause improper values on the temperature readout.
- There are specific error messages (cycle faults) for these board types. These can be found in the *Troubleshooting Cycle Faults* section.
- IMPORTANT REVISION 4.XX SOFTWARE IS FOR OPERATING UNITS WITH THE STAINLESS STEEL BOILER AND REVISION 5.XX SOFTWARE IS FOR OPERATING THE ALUMINIUM EXTERNAL (ALEX) BOILER. CALIBRATION, HOWEVER, IS THE SAME FOR BOTH.
- Neither boiler requires calibration as these units are fitted with a separate validation thermocouple.
- This calibration procedure may also be used on earlier revision 2.x and 5.x controller boards that have been upgraded to ALEX steam generator (R5.xx software).

#### Calibration Procedure


Calibration of these boards should be done in the following order:

- 1. Reference voltage verification / adjustment
- 2. Chamber thermocouple calibration
- 3. Validation thermocouple calibration (automatic)

#### Required equipment:


- Control box 01-103141S
- Digital voltmeter
- Calibration cassette (appropriate to model)
- Potentiometer trimmer
- Calibration jumper
- Digital thermometer

- 1. Before starting calibration, turn the STATIM OFF and remove the cover.
- 2. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette (6) is correctly engaged in the STAT*IM*.
- 3. Install the thermocouple (7) into the cassette and attach the digital thermometer (8).
- 4. Ensure that the thermocouple and digital thermometer have matched serial numbers.
- 5. Install the calibration jumper (11) to Controller Board pins marked W1 (10).

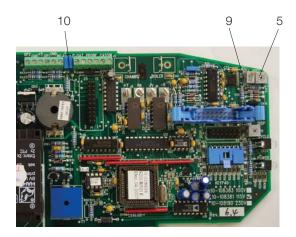

- 1. Control box
- 2. Controller board
- 3. SELECT knob
- 4. Digital voltmeter
- 5. VREF potentiometer
- 6. Calibration cassette
- 7. Thermocouple
- 8. Digital thermometer
- 9. TC-ADJ potentiometer
- 10. Calibration jumper connection
- 11. Calibration jumper (not shown)

Figure 2 (Note: STAT/M 2000 shown)

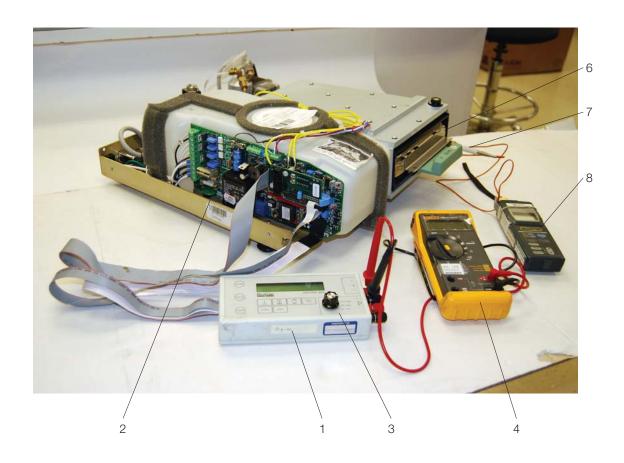

#### 1. Reference voltage verification / adjustment


To verify / adjust the reference voltage:

- Ensure the unit is OFF.
- Connect the control cox (1) to the controller board (2).
- Connect the digital voltmeter (4) to the Vref terminals on the control box.
- Set the voltmeter to read DC Volts with a resolution of 1mV.
- Power the unit ON.
- Adjust the Vref potentiometer until Vref is 2.520 V ± 0.001.
- Once Vref is adjusted, apply a drop of non-conducting lacquer or nail polish to the adjustment screw of the potentiometer.

#### 2. Chamber thermocouple calibration

To calibrate the chamber thermocouple:

 Turn the power switch to OFF. Ensure the calibration jumper is in place and turn the power switch ON. The LCD will read (for example):

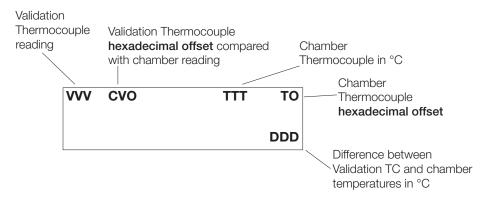

- To start a chamber calibration cycle, press and release the UNWRAPPED cycle button
  on the control box and then press the START button. The STATIM will run a regular
  UNWRAPPED cycle, but continue to show calibration information on the LCD.
- Wait for the chamber to reach the sterilization temperature of 134 °C (TTT).
- The temperature displayed on the LCD needs to match the temperature displayed on the digital thermometer.
- To change the temperature display on the LCD, adjust the TC-ADJ potentiometer (9) on the controller board. Continue adjusting TC-ADJ until the temperatures are within ±0.2°C.
- Observe the temperatures displayed on the LCD and the digital thermometer for 30 seconds. The temperatures registered should remain within ± 0.2°C of one another. Adjust TC-ADJ if not.
- Press the STOP button on the control box when finished.

#### 3. Validation thermocouple calibration

To calibrate the validation thermocouple:

- Turn the power switch to OFF. Ensure the calibration jumper is in place and turn the power switch ON.
- To start a validation thermocouple (self-calibrating) cycle, press and hold the
  UNWRAPPED cycle button on the control box and press the START button. An asterisk
  \* will appear next to the validation thermocouple offset value (CVO) to indicate that a
  validation assembly calibration cycle is running.
- If the asterisk \* does not appear, press the STOP button on the control box twice to reset the unit. Turn the power switch OFF then ON again.
- Hold down the UNWRAPPED cycle button on the control box and press the START button to start the self-calibration again.
- This calibration will take approximately 6 minutes.
- Allow the self-calibration to complete.
- The temperature within the chamber will rise to the sterilization temperature. Wait until sterilization phase of the calibration cycle ends automatically.
- The offset value in the upper left-hand corner of the display (CVO) may have changed to a new offset value.
- Press the STOP button on the control box when finished.
- 4. Switch off the unit and remove all calibration equipment.

## Calibration Procedure 4: for Rev. 7.x Board, Software Rev. 6.xx

This calibration procedure is for STAT*IM* 2000/5000 non S units (2007 on) with revision 7.x controller boards and revision 6.xx software.

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- Do not place the digital reference thermometer in direct sunlight, on the armature or on any other hot surfaces. This may cause improper values on the temperature readout.
- There are specific error messages (cycle faults) for these board types. These can be found in the *Troubleshooting Cycle Faults* section.
- Boiler calibration is not required as these units are fitted with a separate validation thermocouple.
- These units may be calibrated with the cover in place OR using the control box.
- There is NO reference voltage adjustment in these machines.

#### **Calibration Procedure**

Calibration of these boards should be done in the following order:

- 1. Chamber thermocouple calibration
- 2. Validation thermocouple calibration (automatic)

#### Required equipment:


- Control box 01-103141S
- Digital voltmeter
- Calibration cassette (appropriate to model)
- Digital thermometer

- 1. Before starting calibration, turn the STAT*IM* OFF.
- 2. The cover MAY be removed and the control box attached if deemed appropriate, but the unit can be calibrated with the cover on as all adjustments to the calibration offsets are undertaken using the keypad on the unit/control box.
- 3. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette (5) is correctly engaged in the STAT*IM*.
- 4. Install the thermocouple (6) into the cassette and attach the digital thermometer (7).
- 5. Ensure that the thermocouple and digital thermometer have matched serial numbers.



- 1. Control box
- 2. Controller board
- 3. SELECT knob
- 4. Digital voltmeter
- 5. Calibration cassette
- 6. Thermocouple
- 7. Digital thermometer

Figure 3

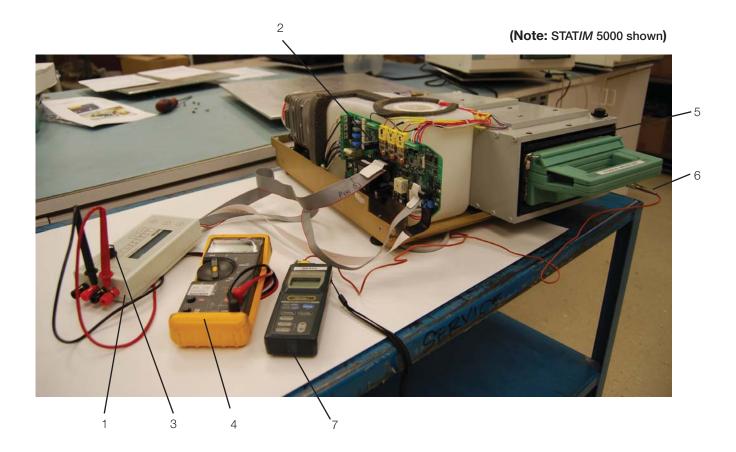


**NOTE:** The cover MAY be removed and the control box attached if deemed appropriate, but the unit can be calibrated with the cover on as all adjustments to the calibration offsets are undertaken using the keypad on the unit.

- 6. Turn ON the unit while keeping Unwrapped and Wrapped button pressed to enter STAT*IM* Service Mode.
- 7. This **Service Mode** is password protected, enter password to continue (default password is: Unwrapped, Wrapped, R&P and Stop keys pressed in this order).

Keypad functions at this time will be:

- Unwrapped Key: Select next item in the menu.
- Wrapped Key: Select previous item in the menu.
- Rubber and Plastics Key: Enter current selection
- Toggle using keypad through the menu selection to reach Calibration option and press R&P key.

The display should be similar to the diagram below:

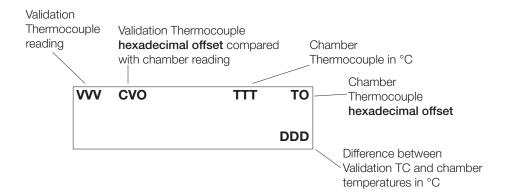

#### 1. Chamber thermocouple calibration

To calibrate the chamber thermocouple:

- 1. To start a chamber calibration cycle, press and release the UNWRAPPED cycle button on the keypad/control box and then press the START button. The STAT/M will run a regular UNWRAPPED cycle, but continue to show calibration information on the LCD.
- 2. Wait for the chamber to reach the sterilization temperature of 134 °C (TTT).
- 3. The temperature displayed on the LCD needs to match the temperature displayed on the digital thermometer.
- 4. To change the temperature display on the LCD, adjust the unwrapped/wrapped keys until the temperatures are within  $\pm 0.2$  °C.

Keypad functions at this time will be:

- Unwrapped Key: increases offset
- Wrapped Key: decreases offset

Observe the temperatures displayed on the LCD and the digital thermometer for 30 seconds. The temperatures registered should remain within  $\pm$  0.2°C of one another. Adjust unwrapped/wrapped keys if not.

5. Press the STOP button on the control box when finished.

#### 2. Validation thermocouple calibration

To calibrate the validation thermocouple:

- 1. Turn the unit off **OFF**.
- Turn ON the unit while keeping Unwrapped and Wrapped button pressed to enter STATIM Service Mode.
- 3. This **Service Mode** is password protected, enter password to continue (default password is: Unwrapped, Wrapped, R&P and Stop keys pressed in this order).
- 4. Toggle using keypad through the menu selection to reach Calibration option and press R&P key.
- 5. The display should show the calibration screen as before
- 6. To start a validation thermocouple (self-calibrating) cycle, press and hold the UNWRAPPED cycle button on the keypad/control box and press the START button. An asterisk "\*" will appear next to the validation thermocouple offset value (CVO) to indicate that a Validation Assy. calibration cycle is running.
- 7. If the asterisk "\*" does not appear, press the STOP button on the control box twice to reset the unit. Turn the power switch OFF then ON again.
- 8. Hold down the UNWRAPPED cycle button on the control box and press the **START** button to start the self-calibration again.
- 9. This calibration will take approximately 6 minutes.
- 10. Allow the self-calibration to complete.
- 11. The temperature within the chamber will rise to the sterilization temperature. Wait until sterilization phase of the calibration cycle ends automatically.
- 12. The offset value in the upper left-hand corner of the display (CVO) may have changed to a new offset value.
- 13. Press the STOP button on the control box when finished.
- 3. Switch off the unit and remove all calibration equipment.

# Calibration Procedure 8: for non-S G4 units with Rev. 7.x Board, Software Rev. 7.xx, SL 00R1.xx

This calibration procedure is for STATIM 2000/5000 G4 non S units with revision 7.x controller boards and revision 7.xx software, SL00R1.xx Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- Do not place the digital reference thermometer in direct sunlight, on the armature or on any other hot surfaces. This may cause improper values on the temperature readout.
- There are specific error messages (cycle faults) for these board types. These can be found in the Troubleshooting Cycle Faults section.
- Boiler calibration is not required as these units are fitted with a separate validation thermocouple.
- These units may be calibrated with the cover in place OR using the control box.
- There is NO reference voltage adjustment in these machines.

#### Calibration Procedure

Calibration of these boards should be done in the following order:

- 1. Chamber thermocouple calibration
- 2. Validation thermocouple calibration (automatic)


**NOTE**: Use of control box is optional. Calibration can be completed using touchscreen. Required equipment:

- Digital voltmeter
- Calibration cassette (appropriate to model)
- Digital thermometer

- 1. Before starting calibration, turn the STATIM OFF.
- 2. The cover MAY be removed and the control box attached if deemed appropriate, but the unit can be calibrated with the cover on as all adjustments to the calibration offsets are undertaken using the touchscreen.
- 3. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette is correctly engaged in the STAT*IM*.
- 4. Install the thermocouple into the cassette and attach the digital thermometer.
- 5. Ensure that the thermocouple and digital thermometer have matched serial numbers.
- 6. Turn ON the unit and access the service menu in the following manner:




- 7. This Service Mode is PIN protected, enter password to continue (default password is: 7919).
- 8. From the service menu, select Calibration.

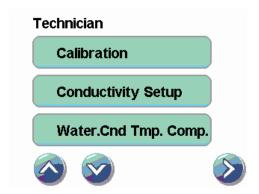


#### 1. Chamber thermocouple calibration

To calibrate the chamber thermocouple:

1. To start a chamber calibration cycle, select Calibration in the Calibration submenu and calibration with start automatically

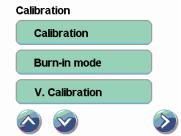

The STATIM will run a regular UNWRAPPED cycle, but continue to show calibration information on the LCD.

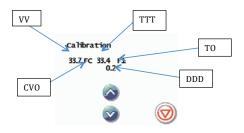

- (VVV) Validation thermocouple reading
- (CVO) Validation thermocouple hexadecimal offset compared with chamber reading
- (TTT) Chamber thermocouple in °C
- (TO) Chamber thermocouple hexadecimal offset
- (DDD) Difference between validation TC and chamber temperatures in °C

2. Wait for the chamber to reach the sterilization temperature of 134 °C (TTT) and for the ADJUST TEMP! message to appear



- 3. The temperature displayed on the LCD needs to match the temperature displayed on the digital thermometer.
- 4. Use the arrow keys to change the temperature display on the LCD until the temperatures are within ±0.2°C. Observe the temperatures displayed on the LCD and the digital thermometer for 30 seconds. The temperatures registered should remain within ± 0.2°C of one another. If not, make adjustments using the arrow keys.
- 5. Press the STOP button on the touchscreen when finished.


#### 2. Validation thermocouple calibration

To calibrate the validation thermocouple:


- 1. Access the service menu as described above.
- 2. From the service menu, select Calibration.



3. To start a validation thermocouple calibration cycle, select V. Calibration in the Calibration submenu. The STATIM will run a Validation Assy. calibration cycle.

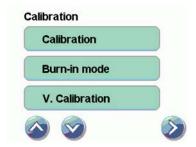

4. This calibration will take approximately 6 minutes.

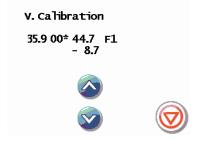

- 5. Allow the self-calibration to complete.
- 6. The temperature within the chamber will rise to the sterilization temperature. Wait until sterilization phase of the calibration cycle ends automatically.



- 7. The offset value in the upper left-hand corner of the display (CVO) may have changed to a new offset value.
- 8. Press the STOP button when finished.
- 3. Switch off the unit and remove all calibration equipment.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



## **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

## PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



## PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

## STATIM 2000S Calibration

## Calibration Procedure 5: for Rev 2.x or 5.x Boards, Software R1.xx and R2.xx

For S-class STAT*IM* 2000 (1995 - 2005) with revision 2.x or revision 5.x type controller boards and software revisions R1.xx or R2.xx.

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- S-Class units contain electronic components that may be damaged or destroyed by electro-static discharge (ESD). Observe appropriate safeguards when calibrating.
- There are specific error messages (cycle faults) for these board/software types. These can be found in the *Troubleshooting Cycle Faults* section.
- In calibration mode, no error messages are displayed.

#### Calibration Procedure

Calibration of these boards should be done in the following order:

- 1. Reference voltage verification / adjustment
- 2. Steam generator thermocouple calibration (automatic)
- 3. Chamber thermocouple and pressure transducer calibration
- 4. Steam generator thermocouple re-calibration (automatic)

#### Required equipment:

- Control box 01-103141S
- Digital voltmeter
- Calibration cassette (appropriate to model)
- Potentiometer trimmer
- Calibration jumper
- Digital thermometer
- Digital pressure meter


- 1. Before starting calibration, turn the STATIM OFF and remove the cover.
- 2. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette (6) is correctly engaged in the STAT*IM*.
- 3. Install the thermocouple (7) into the cassette and attach the digital thermometer (8).

- 4. Ensure that the thermocouple and digital thermometer have matched serial numbers.
- 5. Attach the pressure reference meter to the pressure tube and the pressure tube to the test cassette. (Self locking couplings).
- 6. Install the calibration jumper to Controller Board pins marked W1.
- 7. Locate TC ADJ and Vref potentiometers.



- 1. Control box
- 2. Controller board
- 3. SELECT knob
- 4. VRH potentiometer
- 5. Digital voltmeter
- 6. Calibration cassette
- 7. Thermocouple
- 8. Digital thermometer
- 9. TC-BLR potentiometer
- 10. TC-CHM potentiometer
- 11. Pressure meter

Figure 1

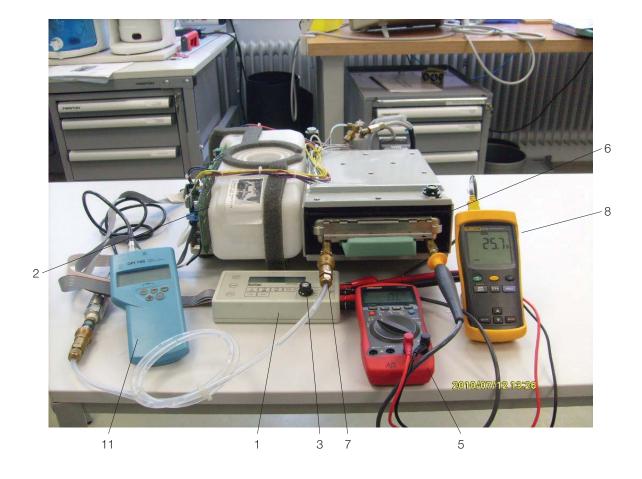

#### 1. Reference voltage verification / adjustment


To verify / adjust the reference voltage:

- 1. Ensure the unit is OFF.
- 2. Connect the control box (1) to the controller board (2).
- 3. Connect the digital voltmeter (5) to the Vref terminals on the control box.
- 4. Set the voltmeter to read DC Volts with a resolution of 1mV.
- 5. Power the unit **ON**.
- 6. Adjust the Vref potentiometer until Vref is  $2.520 \text{ V} \pm 0.001$ .
- 7. Once Vref is adjusted, apply a drop of non-conducting lacquer or nail polish to the adjustment screw of the potentiometer.


#### 2. Steam generator thermocouple calibration

To calibrate the steam generator thermocouple:

- 1. Turn the unit **OFF** and check that there is a calibration jumper on controller board header W1.
- 2. Turn the unit ON. The LCD should appear and be similar to the diagram below:

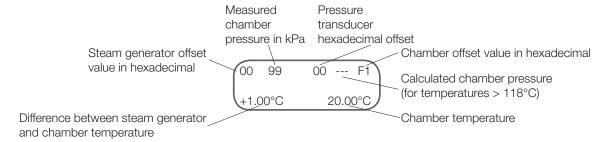

3. To start a steam generator (self-calibrating) cycle, press and hold the UNWRAPPED cycle button on the control box and press the START button. An asterisk "\*" will appear next to the steam generator offset value.

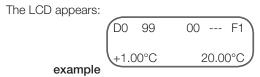

- 4. If the asterisk "\*" does not appear, press the STOP button on the control box twice to reset the unit. Turn the power switch OFF then ON again.
- 5. Hold down the UNWRAPPED cycle button on the control box and press the **START** button to start the self-calibration again.
- 6. If the controller board is new and has never previously been calibrated, adjust the TC-ADJ potentiometer (8) until the chamber thermocouple offset reads F5. If the controller board has previously been calibrated this adjustment is not necessary.
- 7. As the cycle proceeds, the chamber temperature reaches the sterilization temperature, drops to 115°C and then regains sterilization temperature.
- 8. After 20 40 seconds, a long beep will sound indicating that steam generator thermocouple calibration is complete. The number in the upper left-hand corner of the LCD has changed to display the new steam generator offset value.

- 9. Press the **STOP** button on the control box to end the steam generator calibration cycle and to depressurize the cassette.
- 10. Press STOP a second time to reset the unit.

#### 3. Chamber thermocouple and pressure transducer calibration

To calibrate the chamber thermocouple and pressure transducer:

1. Turn the power switch to **OFF**. Ensure the calibration jumper is in place and turn the power switch **ON**. The LCD will read (for example):

- 2. With the calibration cassette removed and the pressure tubing disconnected at both ends, compare the measured chamber pressure reading on the LCD to the pressure displayed on the reference meter. The difference between them is the pressure transducer offset.
- 3. Fully insert the calibration cassette and connect the pressure calibration tubing.
- 4. To start a chamber thermocouple / pressure transducer calibration cycle, press and release the UNWRAPPED cycle button and then press **START**. The STAT*IM* 2000S will run a regular UNWRAPPED cycle, but continue to display calibration information on the LCD.
- 5. Immediately after pressing the START button, use the UNWRAPPED and WRAPPED cycle buttons to make the measured chamber pressure reading displayed on the LCD equal to the reference meter reading. Note that once the boiler begins producing steam the pressure cannot be accurately zero-adjusted to atmospheric pressure. TIP: The offset for the measured chamber pressure is increased or decreased in 0.5 kPa steps. Each press of the UNWRAPPED button adds 0.5 kPa and each press of the WRAPPED button subtracts 0.5 kPa.
- 6. Allow the chamber to reach the sterilization temperature. As pressure builds in the chamber check for leaks in the cassette, associated piping and fittings. A steam leak in the system will introduce errors in the measurement and will result in improper calibration and may compromise sterilization conditions.
- 7. Once the unit has reached sterilization and the temperature and pressure are in equilibrium compare the chamber temperature displayed on the LCD with the temperature displayed on the temperature reference meter. Adjust the TC-ADJ potentiometer on the controller board until the LCD temperature matches the reference temperature to within ± 0.25°C.
- 8. Wait ten seconds after the water pump turns off to allow the pressure and temperature readings to stabilize before making any adjustments.
- 9. IMPORTANT NOTE: When the initial pressure and temperature readings have been adjusted, and to ensure equilibrium of the calibration, compare the measured chamber pressure and the calculated chamber pressure displayed on the LCD. They should be within 1 kPa of each other during equilibrium. If not, adjust the pressure offset using the membrane keypad by ONE UNIT ONLY. If the readings still differ by more than 1 kPa, using the TC-ADJ potentiometer adjust the chamber offset by ONE unit. If the readings still differ by more than 1 kPa abort the cycle and retry the chamber thermocouple / pressure transducer calibration.

- 10. During sterilization and after all necessary adjustments have been made, the measured chamber pressure indicated by the reference meter and the STAT*IM* 2000S LCD should not differ by more than 6 kPa.
- 11. A consecutive failure to achieve the above conditions after adjustment indicates that one or more of the following components may be faulty:
- Chamber thermocouple
- Pressure sensor
- Pressure interface board
- Controller board (Thermocouple amplifier or A to D converter).
- 12. When the calibration is complete, press the STOP button on the control box and wait for the beeps to signal that the cassette may be removed. Press the STOP button again to reset the unit.
- 4. Steam generator thermocouple re-calibration

Repeat the steam generator thermocouple calibration procedure in section 2 described above.

5. Switch off the unit and remove all calibration equipment.

# Calibration Procedure 6: for Rev 6.x Boards, Software R4.xx and R5.xx

For S-class STAT*IM* 2000 (2005 - 2007) with revision 6.x type controller boards and revision 4.xx and 5.xx software.

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- S-Class units contain electronic components which may be damaged or destroyed by electro-static discharge (ESD). Observe appropriate safeguards when calibrating.
- There are specific error messages (cycle faults) for these board/software types. These can be found in the *Troubleshooting Cycle Faults* section.
- In calibration mode no error messages are displayed.
- IMPORTANT REVISION 4.XX SOFTWARE IS FOR OPERATING UNITS WITH THE STAINLESS STEEL BOILER AND REVISION 5.XX SOFTWARE IS FOR OPERATING THE ALUMINIUM EXTERNAL (ALEX) BOILER. CALIBRATION, HOWEVER, IS THE SAME FOR BOTH STATIM UNITS.
- The boiler thermocouples for both boiler types do not require calibration.

### Calibration Procedure

Calibration of these boards should be done in the following order:

- 1. Reference voltage verification / adjustment
- 2. Chamber thermocouple and pressure transducer calibration


### Required equipment:

- Control box 01-103141S
- Digital voltmeter
- Calibration cassette (appropriate to model)
- Potentiometer trimmer
- Calibration jumper
- Digital thermometer
- Digital pressure meter

#### Setting up the unit for calibration:


- 1. Before starting calibration, turn the STAT/M OFF and remove the cover.
- 2. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette (6) is correctly engaged in the STAT*IM*.
- 3. Install the thermocouple (7) into the cassette and attach the digital thermometer (8).
- 4. Ensure that the thermocouple and digital thermometer have matched serial numbers.

- 5. Attach the pressure reference meter to the pressure tube and the pressure tube to the test cassette. (Self locking couplings).
- 6. Install the calibration jumper to controller board pins marked W1.
- 7. Locate TC ADJ and Vref potentiometers.

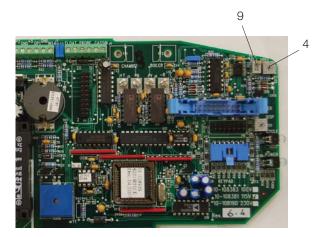

- 1. Control box
- 2. Controller board
- 3. SELECT knob
- 4. VREF potentiometer
- 5. Digital voltmeter
- 6. Calibration cassette
- 7. Thermocouple
- 8. Digital thermometer
- 9. TC-ADJ potentiometer
- 10. Pressure meter

Figure 2

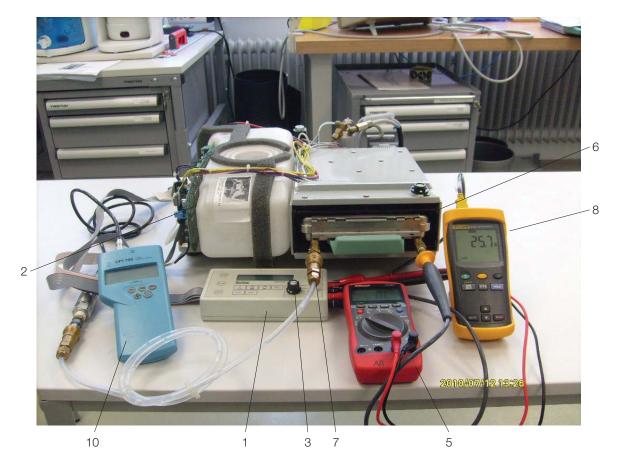

### 1. Reference voltage verification / adjustment

To verify / adjust the reference voltage:

- 1. Ensure the unit is OFF.
- 2. Connect the Control Box (1) to the Controller Board (2).
- 3. Connect the digital voltmeter (5) to the Vref terminals on the control box.
- 4. Set the voltmeter to read DC Volts with a resolution of 1mV.
- 5. Power the unit ON.
- 6. Adjust the VRH potentiometer until Vref is 2.520 V ± 0.001.
- 7. Once Vref is adjusted, apply a drop of non-conducting lacquer or nail polish to the adjustment screw of the potentiometer.

IMPORTANT NOTE: remember – the steam generator thermocouple does NOT require calibration.

#### 2. Chamber thermocouple and pressure transducer calibration

To calibrate the chamber thermocouple and pressure transducer:

1. Turn the power switch to **OFF**. Ensure the calibration jumper is in place and turn the power switch **ON**. The LCD will read (for example):



#### Where

- PPP.P = actual chamber pressure in kPa.
- PO = chamber pressure hexadecimal offset value.
- TTT.T = theoretical pressure calculated from the actual chamber temperature.
- CO = chamber temperature hexadecimal offset value.
- CCC.C = actual chamber temperature in °C.
- 2. To start the calibration cycle, press and release the UNWRAPPED cycle button on the control box and then press the START button. Adjust Press! will appear in the display.



- 3. By using the Unwrapped key (+) and Wrapped Key (-) the pressure offset PO (at atmospheric pressure) can now be adjusted until the reading on the LCD display (as displayed in the PPP.P field) matches the external pressure reference meter's reading. The PO field will then show the new pressure offset. **NOTE**: This has to be done within 30 seconds.
- 4. After 30 seconds, the calibration cycle will continue by entering the heating up phase.
- 5. During calibration the unit will run a normal cycle.

- 6. Allow the chamber to reach the sterilization temperature. As pressure builds in the chamber check for leaks in the cassette, associated piping and fittings. A steam leak in the system will introduce errors in the measurement and will result in improper calibration and potentially non-sterile instruments.
- 7. IMPORTANT NOTE: For units fitted with a Stainless Steel steam generator (revision 4.xx software), at the beginning of the holding phase (sterilization phase) all the devices are turned off for approximately 10 seconds to facilitate the calibration process. This will happen ten times in succession. During these periods chamber temperature could go down to 131°C which is quite normal. After these ten "calm" periods the unit will resume normal operation. This DOES NOT happen on units fitted with aluminium external (ALEX) steam generator (revision 5.xx software).

TIP: To check that the calibration has been done correctly on Stainless Steel/Rev 4.xx units, look at the external temperature reference meter. When it shows 136.5°C you should hear the solenoid valve 'click'. (In this type of unit the software is set to open the solenoid valve at 136.5°C)

8. When Adjust Temp! appears on the display, (see below) and you are confident that the readings are stable, the chamber thermocouple offset CO can now be adjusted using the TC-ADJ potentiometer until the reading on the LCD display (as displayed in the CCC.C



field) matches the external temperature reference meter's reading. The CO field will then show the new chamber temperature offset.

**IMPORTANT NOTE:** When the initial pressure and temperature readings have been adjusted, and to ensure equilibrium of the calibration, compare the measured chamber pressure (PPP.P) and the calculated chamber pressure (TTT.T) displayed on the LCD. They should be within 1 kPa of each other during equilibrium.

- 9. If they are not, then the unit may be fine tuned as follows:
  - These units have a two stage pressure offset and the upper offset can be modified during the sterilization phase of the calibration. This will NOT effect the lower offset value BUT if you subsequently change the lower value after the unit has vented then you will lose the upper value as this will be reset to the same value as the lower offset.
  - If the actual chamber pressure reading on the LCD (PPP.P) is different by more than 1 kPa from the pressure reference meter reading then this can be adjusted using the wrapped and unwrapped keys.
  - If this adjustment corrects the PPP.P reading to within 1 kPa of the TTT.T reading then the unit is now in equilibrium.
  - If these readings are still not within 1 kPa of each other then further adjustment may be undertaken as follows:
    - The theoretical pressure reading (TTT.T) and actual chamber temperature reading (CCC.C) are directly related and when the actual temperature value is increased or decreased, the theoretical pressure value will also increase or decrease respectively.
    - Adjust the theoretical chamber pressure (TTT.T) using the TC-ADJ potentiometer (remember that the theoretical pressure is related to the ACTUAL temperature which is adjusted using the TC-ADJ potentiometer) until PPP.P and TTT.T are within 1 kPa of each other.

- IMPORTANT now re-check the actual temperature reading (CCC.C) with the reading on the temperature reference meter. The readings need to be stable, so allow sufficient time to compare them.
- If you are confident that the temperature readings are still within ±0.25°C, then the unit is in equilibrium and the calibration has finished.
- If the readings are outside of the  $\pm$  0.25°C tolerance then further diagnostics may be required prior to re-calibration.
- 10. Press the STOP button on the control box when finished.
- 3. Switch off the unit and remove all calibration equipment.

# Calibration Procedure 7: for Rev 7.x Boards, Software R6.xx

For S-class STAT*IM* 2000 (2007 - present) with revision 7.x type controller boards and revision 6.xx software.

#### Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- S-class units contain electronic components which may be damaged or destroyed by electro-static discharge (ESD). Observe appropriate safeguards when calibrating.
- There are specific error messages (cycle faults) for these board/software types. These can be found in the *Troubleshooting Cycle Faults* section.
- In calibration mode no error messages are displayed.
- Boiler calibration is not required on these units.
- These units may be calibrated with the cover in place OR using the control box as all adjustments are made with the WRAPPED or UNWRAPPED keys on the relevant keypad.
- There is NO reference voltage adjustment in these machines.

#### **Calibration Procedure**

Calibration of these boards should be done in the following order:

1. Chamber thermocouple/pressure transducer calibration

#### Required equipment and tools:

- Control box 01-103141S
- Calibration cassette (appropriate to model)
- Digital thermometer
- Digital pressure meter
- Digital voltmeter


### Setting up the unit for calibration:

- 1. Before starting calibration, turn the STATIM OFF.
- 2. The cover MAY be removed and the control box attached if deemed appropriate, but the unit can be calibrated with the cover on as all adjustments to the calibration offsets are undertaken using the keypad on the unit.
- 3. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette is correctly engaged in the STAT*IM*.
- 4. Install the thermocouple into the cassette and attach the digital thermometer.
- 5. Ensure that the thermocouple and digital thermometer have matched serial numbers.
- 6. Attach the pressure reference meter to the pressure tube and the pressure tube to the test cassette. (Self locking couplings).



- 1. Control box
- 2. Controller board
- 3. SELECT knob
- Digital voltmeter
   Calibration cassette
- 6. Thermocouple
- 7. Digital thermometer
- 8. Pressure meter

Figure 3

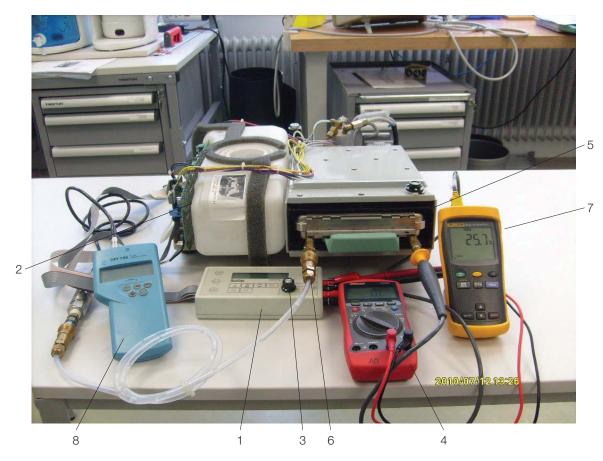

Note: The cover MAY be removed and the control box attached if deemed appropriate, but the unit can be calibrated with the cover on as well as adjustments to the calibration offsets are undertaken using the keypad on the units.

- 7. Turn ON the unit while keeping Unwrapped and Wrapped button pressed to enter STAT*IM* Service Mode.
- 8. This **Service Mode** is password protected, enter password to continue (default password is: Unwrapped, Wrapped, R&P and Stop keys pressed in that order).

Keypad functions (Unit LCD or Control Box) at this time will be:

- Unwrapped Key: Select next item in the menu.
- Wrapped Key: Select previous item in the menu.
- Rubber and Plastics Key: Enter current selection
- Toggle using keypad through the menu selection to reach Calibration option and press R&P key.

The display should be similar to the diagram below.



#### Where

- PPP.P = actual chamber pressure in kPa.
- PO = chamber pressure hexadecimal offset value
- TTT.T = theoretical pressure calculated from the actual chamber temperature.
- CO = chamber temperature hexadecimal offset value
- CCC.C = actual chamber temperature in °C

#### **IMPORTANT NOTES:**

- Remember the steam generator thermocouple does NOT require calibration.
- The VREF does not need adjustment.

1. Chamber Thermocouple and Pressure Transducer Calibration

To calibrate the chamber thermocouple and pressure transducer:

1. To start the calibration cycle, press and release the UNWRAPPED cycle button on the control box / display and then press the START button. Adjust Press! will appear in the display and the PO value will flash. The PO value can only be adjusted while it is flashing.



- 2. By using the Unwrapped key (+) and Wrapped Key (-) the pressure offset PO (at atmospheric pressure) can now be adjusted until the reading on the LCD display (as displayed in the PPP.P field) matches the external pressure reference meter's reading. The PO field will then show the new pressure offset. **NOTE**: This has to be done within 30 seconds.
- After 30 seconds, the calibration cycle will continue by entering the heating up phase and the PO will STOP flashing and the CO value will flash. The CO value can only be adjusted while it is flashing.
- 4. During calibration the unit will run a normal cycle.
- 5. Allow the chamber to reach the sterilization temperature. As pressure builds in the chamber, check for leaks in the cassette, associated piping and fittings. A steam leak in the system will introduce errors in the measurement and will result in improper calibration and may compromise sterilization conditions.
- 6. When Adjust Temp! appears on the display, (see below) and you are confident that the readings are stable, the chamber thermocouple offset CO can now be adjusted using the wrapped and unwrapped keys until the reading on the LCD display (as displayed in the CCC.C field) matches the external temperature reference meter's reading. The CO field will then show the new chamber temperature offset.



**IMPORTANT NOTE:** When the initial pressure and temperature readings have been adjusted, and to ensure equilibrium of the calibration, compare the measured chamber pressure (PPP.P) and the calculated chamber pressure (TTT.T) displayed on the LCD. They should be within 1 kPa of each other during equilibrium.

- 7. If they are not, then the unit may be fine tuned as follows:
  - These units have a two-stage pressure offset and the upper offset can be modified
    during the sterilization phase of the calibration. This will NOT effect the lower offset
    value BUT if you subsequently change the lower value after the unit has vented then
    you will lose the upper value as this will be reset to the same value as the lower offset.
  - If the actual chamber pressure reading on the LCD (PPP.P) is different by more than 1 kPa from the pressure reference meter reading then this can be adjusted using the wrapped and unwrapped keys.
    - Press the Rubber and Plastics key.
    - The CO value will now STOP flashing and the PO value will now flash.
    - Use the wrapped and unwrapped keys to correct the upper PO value.
  - If this adjustment corrects the PPP.P reading to within 1 kPa of the TTT.T reading then the unit is now in equilibrium.

- If these reading are still not within 1 kPa of each other then further adjustment may be undertaken by modifying the TTT.T value.
- The theoretical pressure reading (TTT.T) and actual chamber temperature reading (CCC.C) are directly related and when the actual temperature value is increased or decreased, the theoretical pressure value will also increase or decrease respectively.
- Adjust the theoretical chamber pressure (TTT.T) and actual chamber temperature readings (CCC.C) as follows:
  - Press the Rubber and Plastics key.
  - The PO value will now STOP flashing and the CO value will now flash.
  - Use the wrapped and unwrapped keys to correct the TTT.T value so that it is within 1 kPa of PPP.P.
- IMPORTANT now re-check the actual temperature reading (CCC.C) with the
  reading on the temperature reference meter. The readings need to be stable, so allow
  sufficient time to compare them.
- If you are confident that the temperature readings are still within ±0.25°C, then the unit is in equilibrium and the calibration has finished.
- If the readings are outside of the  $\pm$  0.25°C tolerance then further diagnostics may be required prior to re-calibration.
- Press the STOP button on the control box when finished.
- 2. Switch off the unit and remove all calibration equipment.

# Calibration Procedure 9: for S-class G4 Units with Rev. 7.x Board, Software Rev. 7.xx, SL 00R1.xx

This calibration procedure is for S-class STATIM 2000/5000 G4 units with revision 7.x type controller boards and revision 7.xx software, SL00R1.xx Important:

- Incorrect or inaccurate calibration may cause unsuccessful sterilization of instruments.
- Always calibrate the thermocouples after a thermocouple replacement, thermocouple bend or disconnect, controller board replacement or microprocessor / EEPROM replacement, steam generator replacement or probe bracket replacement.
- Make sure that there is sufficient steam-process distilled water in the unit prior to starting calibration.
- S-class units contain electronic components which may be damaged or destroyed by electrostatic discharge (ESD). Observe appropriate safeguards when calibrating.
- There are specific error messages (cycle faults) for these board/software types. These can be found in the Troubleshooting Cycle Faults section.
- In calibration mode no error messages are displayed.
- Boiler calibration is not required on these units.
- These units may be calibrated with the cover in place OR using the control box as all adjustments are made using the touchscreen.
- There is NO reference voltage adjustment in these machines.

#### Calibration Procedure

Calibration of these boards should be done in the following order:

1. Chamber thermocouple/pressure transducer calibration

Required equipment and tools:

- Calibration cassette (appropriate to model)
- Digital thermometer
- Digital pressure meter
- Digital voltmeter

NOTE: Use of control box is optional. Calibration can be completed using touchscreen.

## Setting up the unit for calibration:

- 1. Before starting calibration, turn the STATIM OFF.
- 2. The cover MAY be removed and the control box attached if deemed appropriate, but the unit can be calibrated with the cover using the touchscreen.
- 3. Check that there is sufficient steam process distilled water in the reservoir, and that the calibration cassette is correctly engaged in the STAT*IM*.
- 4. Install the thermocouple into the cassette and attach the digital thermometer.
- 5. Ensure that the thermocouple and digital thermometer have matched serial numbers.
- 6. Attach the pressure reference meter to the pressure tube and the pressure tube to the test cassette. (Self-locking couplings).

7. Turn ON the unit and access the service menu in the following manner:



- 8. This Service Mode is PIN protected, enter password to continue (default password is: 7919).
- 9. From the service menu, select Calibration.




#### **IMPORTANT NOTES:**

- Remember the steam generator thermocouple does NOT require calibration.
- The VREF does not need adjustment.


## 1. Chamber Thermocouple and Pressure Transducer Calibration

To calibrate the chamber thermocouple and pressure transducer:

1. To start the calibration cycle, select Calibration in the Calibration submenu and calibration will start automatically

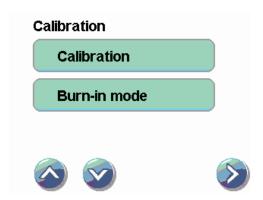

1. Adjust Press! will appear in the display and the PO value will flash. The PO value can only be adjusted while it is flashing.

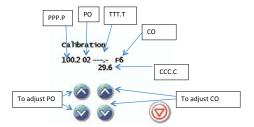

(PPP.P) Actual chamber pressure in kPa.

(PO) Chamber pressure hexadecimal offset value

(TTT.T) Theoretical pressure calculated from the actual chamber temperature.

(CO) Chamber temperature hexadecimal offset value

(CCC.C) Actual chamber temperature in °C

By using the lefthand arrow keys, adjust the pressure offset PO (at atmospheric pressure) until the reading on the LCD display matches the external pressure reference meter's reading. The PO field will then show the new pressure offset.

NOTE: This has to be done within 30 seconds.

- 3. After 30 seconds, the calibration cycle will continue by entering the heating up phase and the PO will STOP flashing and the CO value will flash. The CO value can only be adjusted while it is flashing.
- 4. During calibration, the unit will run a normal cycle.
- Allow the chamber to reach the sterilization temperature. As pressure builds in the chamber, check for leaks in the cassette, associated piping and fittings. A steam leak in the system will introduce errors in the measurement and will result in improper calibration and may compromise sterilization conditions.
- 6. When Adjust Temp! appears on the display, (see below) and you are confident that the readings are stable, the chamber thermocouple offset CO can now be adjusted using the using the righthand arrow keys below it until the reading on the LCD display (as displayed in the CCC.C field) matches the external temperature reference meter's reading. The CO field will then show the new chamber temperature offset.



**IMPORTANT NOTE**: When the initial pressure and temperature readings have been adjusted, and to ensure equilibrium of the calibration, compare the measured chamber pressure (PPP.P) and the calculated chamber pressure (TTT.T) displayed on the LCD. They should be within 1 kPa of each other during equilibrium.

- 7. If they are not, then the unit may be fine-tuned as follows:
  - These units have a two-stage pressure offset and the upper offset can be modified during
    the sterilization phase of the calibration. This will NOT affect the lower offset value BUT if
    you subsequently change the lower value after the unit has vented then you will lose the
    upper value as this will be reset to the same value as the lower offset.
  - If the actual chamber pressure reading on the LCD (PPP.P) is different by more than 1 kPa from the pressure reference meter reading then this can be adjusted using the lefthand arrow keys. If this adjustment corrects the PPP.P reading to within 1 kPa of the TTT.T reading then the unit is now in equilibrium.
  - If these readings are still not within 1 kPa of each other, further adjustment may be undertaken by modifying the TTT.T value.
  - The theoretical pressure reading (TTT.T) and actual chamber temperature reading (CCC.C) are directly related and when the actual temperature value is increased or decreased, the theoretical pressure value will also increase or decrease respectively.
  - Adjust the theoretical chamber pressure (TTT.T) and actual chamber temperature readings (CCC.C) using the righthand arrow keys to correct the TTT.T value so that it is within 1 kPa of PPP.P.
  - IMPORTANT now re-check the actual temperature reading (CCC.C) with the reading on the temperature reference meter. The readings need to be stable, so allow sufficient time to compare them.
  - The unit is in equilibrium and the calibration is complete when you are confident that the temperature readings are still within ±0.25°C.
  - If the readings are outside of the ± 0.25°C tolerance, further diagnostics may be required prior to re-calibration.
- 8. Press the STOP button when finished.
- 2. Switch off the unit and remove all calibration equipment.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



# **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



# PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# Cover Assembly

### For STATIM 2000 Classic


# Removing the Cover

To remove the cover, proceed as follows (see Figure 1):


- 1. Inspect the cover (1) to ensure that it has not been damaged in transit. Position the unit on the bench so the front leveler feet are securely on the work surface and the front cover extends past the work surface edge.
- 2. Check that the LCD (2) and keypad (3) function. This requires that the unit be powered ON. Check for symptoms of the failure to ensure that no further damage occurs.
- 3. Turn the power OFF and unplug the power cord (4) from the wall outlet and remove the cassette from the unit.
- 4. Remove the reservoir cap (5) (and reservoir filter where fitted) from the top of the unit.
- 5. Using a No. 2 Phillips screwdriver, remove seven screws from the cover; four screws from the sides and three at the rear of the unit.
- 6. If a bacteria-retentive air filter is fitted, unclip the filter from the filter bracket before attempting to remove the cover.

Note: Before attempting the next step, be aware that there are cables connecting the cover and the controller. If due care is not taken when removing the cover, damage may occur.

- 7. With the fascia hanging over the edge of the work surface, lift the rear of the cover until it clears the internal components, and then slide the cover forward about 3 cm / 1 inch to clear the front of the armature and the unsupported controller board (8).
- 8. Gently lift the cover up slightly and pivot it around the front left corner of the unit, as shown in Figure 1. When the cover is clear of the unit, rest it on its side. Take care not to strain the ribbon cables connecting the display and keypad.
- 9. Observe the orientation of the ribbon cable connections. Disconnect the LCD and Keypad connectors from the PCB.
  - Tip 1: To disconnect the keypad connector from the controller board, gently release the lock (tab) on the connector housing of the controller board and pull down on the connector.
  - Tip 2: On some units there may be a cable tie securing the LCD cable connector to the PCB. Check for this before attempting to pull off the connector. Remove the cable tie as needed. A new cable tie is not required if the unit is not being shipped by courier.
- 10. Replace the reservoir cap to prevent unwanted materials and debris from entering the reservoir.

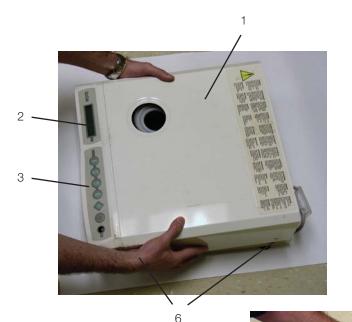

- 1. Cover
- 2. LCD
- 3. Keyboard
- 4. Power cord (not shown)
- 5. Reservoir cap (not shown)
- 6. Screws for side of unit

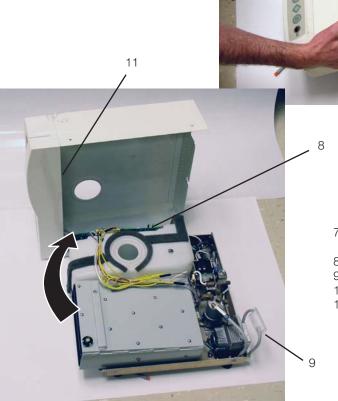

- 7. Screws with lockwashers for rear of unit
- 8. Controller board
- 9. Bacteria-retentive air filter
- 10. Filter bracket (not shown)
- 11. Armature gasket

Figure 1

# Reinstalling the Cover

To reinstall the cover, proceed as follows (see Figure 1).

- 1. Inspect the armature gasket (11) located on the inside front portion of the cover. If it is damaged or fails to adhere to the surface of the cover, it must be replaced. (See Replacing the Armature Gasket.)
- 2. Remove the reservoir cap (5) from the reservoir.
- 3. Rest the cover (1) on the left side of the unit. Reconnect the keypad and LCD connectors to controller board.
- 4. Reconnect the power cord (4) and power the unit ON to check that the LCD and keypad are functioning. If the amber LED on the front cover is flashing, or the solenoid valve is clicking repeatedly, the LCD connector is incorrectly installed. Check that the keypad connector is properly plugged onto the header. Press all of the keypad buttons and check for proper display messages. After checking, turn the power switch OFF and unplug the power cord.
- 5. Carefully lift the cover from the work surface. While rotating the cover slide it forward until the fascia clears the controller board and the front of the armature.
- 6. Tilt the front of the cover down and the rear of the cover upwards. Carefully reposition the front portion of the cover back over the front of the controller board (8) and armature. Ensure that the LEDs at the front left hand corner of the unit fit into the clearance holes in the cover. BE CAREFUL NOT TO PINCH THE RIBBON CABLES.
- 7. When the front of the cover is in place, lower the rear portion and gently push it backwards. Carefully realign the screw holes and reinstall three screws with lockwashers (7) across the rear of the unit. If the unit being serviced has a bacteria-retentive air filter (9), reinstall the filter.
- 8. Reinstall the four remaining screws (6), two on each side. Press firmly on the top of the cover to compress the gaskets and partially reinstall the screws as the holes align. When in place, tighten all the screws. Do not over tighten the screws.
- 9. Replace the reservoir filter (if fitted) and fill the reservoir with steam process distilled water then replace the reservoir cap (5).
- 10. Reconnect the power cord (4). Select a cycle to ensure that the keypad and LCD function, then press the STOP key.

# Removing and Replacing the Armature Gasket

To replace the armature gasket (2), proceed as follows (see Figure 2):

- Remove all traces of the gasket and gasket adhesive from inside the fascia (1). Note that the
  fascia is made of a thermoplastic. DO NOT USE AROMATIC SOLVENTS THEREFORE ON
  THE FASCIA as this will damage the surface. To prepare for installation of the new gasket (2),
  roughen the inside surface of the fascia where the gasket sits, using a fine grit abrasive paper.
- 2. Leave the carrier paper on the adhesive side of the gasket and perform a trial fit. Note that the widest part of the gasket fits at the bottom of the opening and that the gasket does not protrude into the opening.
- 3. Remove the carrier paper from the adhesive side of the gasket. Carefully install the gasket making sure it is seated firmly in place.

Removing and Replacing the Fascia

To remove the fascia (1), proceed as follows (see Figure 2):

- The fascia is attached to the cover (3) using five Plastite<sup>™</sup> screws (4). Remove and retain the screws.
- 2. A bead of silicone was used to provide a water-tight seal between the fascia and cover during manufacturing. Remove the fascia and any silicone residue remaining on the cover.
- The fascia is made of recyclable plastic. If the fascia is to be discarded please recycle wherever possible.

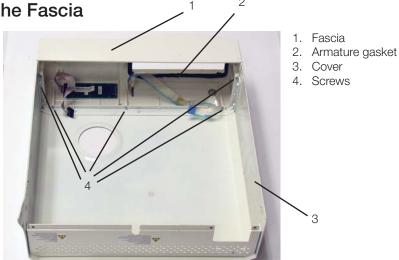


Figure 2

To replace the fascia (1), proceed as follows (see Figure 2):

- 1. Apply a bead of silicone between the facsia and cover (3) to provide a water-tight seal.
- 2. Install the new fascia using the 5 screws (4) retained from the disassembly procedure.

# Removing and Replacing the Membrane Keypad (LCD overlay)

To remove the membrane keypad (1), proceed as follows (see Figure 3):

**Note:** There are two types of membranes. Depending on the age of the unit, it may be a two-part membrane (LCD overlay + Keypad membrane) or a one-piece overlay with integrated keypad membrane. Ensure you have the appropriate new overlay before commencing the replacement.

- 1. The keypad's flexible cable and the LCD ribbon cable may be secured together with double sided adhesive tape. Carefully separate the cables.
- 2. Using a sturdy sharp instrument, carefully lift one corner of the keypad section only or one piece overlay and CAREFULLY peel away from the fascia (2). EXERCISE CAUTION, there may be adhesive on the edges of the LCD. Pass the connector through the slot (5) in the fascia.

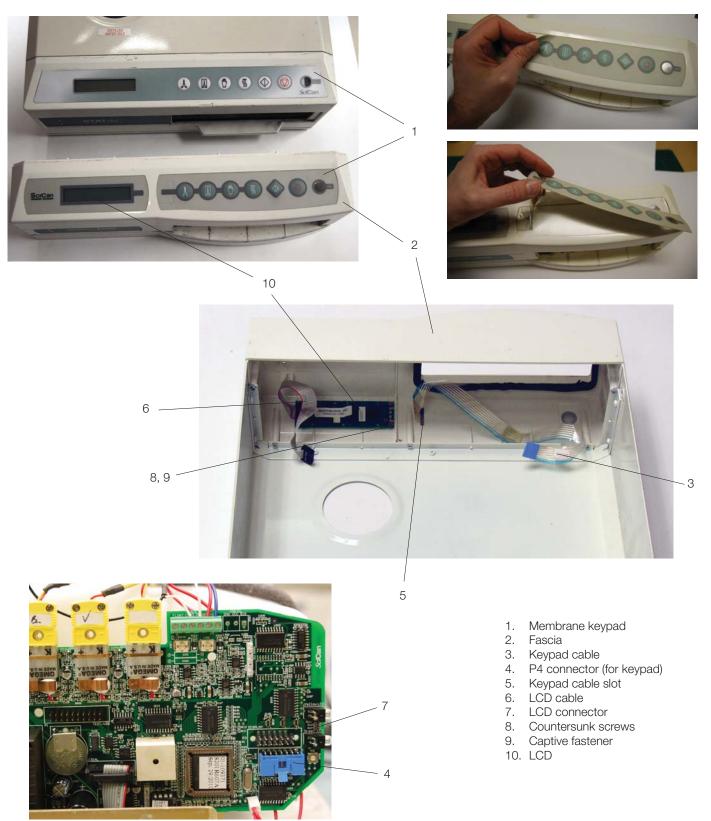


Figure 3

 Carefully remove any residual adhesive or membrane from the recessed areas of the fascia using a sharp blade. DO NOT USE AROMATIC SOLVENTS TO REMOVE ADHESIVE FROM THE FASCIA OR THE MEMBRANE KEYPAD.

To replace the membrane (1), proceed as follows (see Figure 3):

- 1. Before installing the membrane onto the fascia, connect the new keypad cable (3) connector to controller board P4 (4), holding the keypad connector by its tab, insert it into the housing on the controller board and pull up on the tab to lock the connector into place.
- 2. Plug in the power cord and power the unit ON. To test each button on the keypad, ensure that a cassette is not inserted into the armature. Select a cycle and press the START button. The LCD will display "INSERT CASSETTE". Repeat this step for each cycle to ensure the keypad is functioning properly. To test the STOP button, select a cycle. The LCD message will indicate which cycle is selected. Press the STOP button. The LCD message will change to "SELECT A CYCLE".
- 3. Turn off and unplug the unit, and disconnect the keypad cable, by gently releasing the lock on the connector housing of the controller board and pulling down on the connector tab to release it.
- 4. The replacement membrane keypad has carrier paper on the back to protect the adhesive. With the carrier paper in place, feed the keypad cable (4) connector through the slot (5) in the fascia (2).
- 5. Remove the carrier paper and align the bottom edge of the membrane keypad with the bottom edge of the recessed area on the fascia. Carefully lower the membrane into place while continuing to draw the keypad cable connector through the slot. Using finger pressure only, press the membrane into place, eliminating air pockets under the membrane.
- 6. Secure the keypad cable and LCD cable together using double-sided adhesive tape.
- 7. Connect the keypad cable (3) connector to controller board P4 (4) holding the keypad connector by its tab, insert it into the housing on the controller board and pull up on the tab to lock the connector into place.
- 8. Connect the LCD cable (6) connector to controller board (7).
- 9. Test the membrane keypad again by repeating steps 1 and 2.

# Removing and Replacing the LCD

**Note:** There are two types of keypad membranes. Depending on the age of the unit, it may be a two-part membrane (LCD overlay + keypad membrane) or a one-piece overlay with integrated keypad membrane.

To remove the LCD (4), proceed as follows (see Figure 4):

- 1. Carefully remove the LCD overlay or integrated overlay (1) from the fascia. (See Removing and Replacing the Membrane Keypad.) Note that removing the overlay may damage the keypad. Have a replacement available in the event damage occurs.
- Use a No. 1 Phillips screwdriver to remove the four countersunk screws (8) securing the LCD to the fascia. Retain these screws for reassembly. Each screw is held by a captive fastener (5) which is integral to the LCD assembly. In earlier units each screw is held with a nut and washer.
- 3. The keypad cable (3) and the LCD cable (6) may be secured together with double sided adhesive tape. Carefully separate the harnesses and remove the LCD.

To replace the LCD (10), proceed as follows (see Figure 4):

- Align the captive fasteners on the LCD with the four clearance holes in the fascia (2) and reinstall the LCD using the four countersunk screws (8) retained from disassembly. The LCD cable and the keypad cable are on the left side of the cover.
- 2. Carefully reinstall the appropriate overlay (1) on to the fascia. See Membrane Keypad Removal and Replacement. NOTE: When applying pressure to the membrane to ensure the adhesive bonds to the fascia, do not press down in the middle of the LCD assembly. Apply pressure in the areas immediately adjacent the mounting holes only.
- 3. Secure the cables together using double sided adhesive tape.
- 4. Connect the keypad connector to controller board (4).
- 5. Connect the LCD connector to controller board (7).
- 6. Power the unit ON. If the LCD display fails to illuminate, or non alpha-numeric characters appear on the display, the connector is improperly installed or the LCD intensity pot is not adjusted correctly. Other symptoms of an incorrectly installed LCD connector may be a flashing amber LED on the front cover or the solenoid valve activating repeatedly.
- 7. If the LCD intensity or contrast is not suitable, adjust the potentiometer labelled DISP, on the controller board.
- 8. Reinstall the cover.



- 1. LCD overlay
- Fascia
- 3. Keypad cable (not shown)
- 4. LCD
- 5. Captive fastener
- 6. LCD cable
- 7. Controller board
- 8. Countersunk screws



4




Figure 4

## For STATIM 2000 G4

# Removing the STATIM 2000 G4 Cover

To remove the cover, follow these steps (see Figure 5):

- 1. Inspect the cover (1) to ensure that it has not been damaged in transit. Position the unit on the bench so the front leveller feet are securely on the work surface and the front cover extends past the work surface edge.
- 2. Check that the colour LCD functions. This requires that the unit be powered ON Note: Check for symptoms of the failure to ensure that no further damage occurs.
- 3. After checking the LCD, turn the unit OFF, unplug the power cord from the unit's line filter (2) and remove the cassette from the unit.
- 4. Remove the reservoir cap (3) from the top of the unit.
- 5. Remove seven screws from the cover; first remove two on each side (4) and then three with lockwashers (5) at the rear of the unit, using a #2 Phillips screwdriver. The Air Filter (6) must be disconnected from the bracket before the cover can be removed.
- 6. Slowly, lift the rear portion of the cover upwards. When the cover is clear of the internal components, carefully slide the entire cover forward to clear the front of the armature and of the unsupported controller board (8). Slowly pivot the cover the same way a book opens and place it next to the unit. The colour LCD should still be accessible and the internal components exposed.
  - Note: The controller board is unsupported at the front of the unit. Be careful. The LCD cable and LCD DC power source cable are still attached to the unit.
- 7. Disconnect the LCD connector (8) from controller board header P3, and the LCD DC power source cable (9) from the DC power source (10) located at the back of the unit.
- 8. Replace the reservoir cap onto the top of the reservoir.

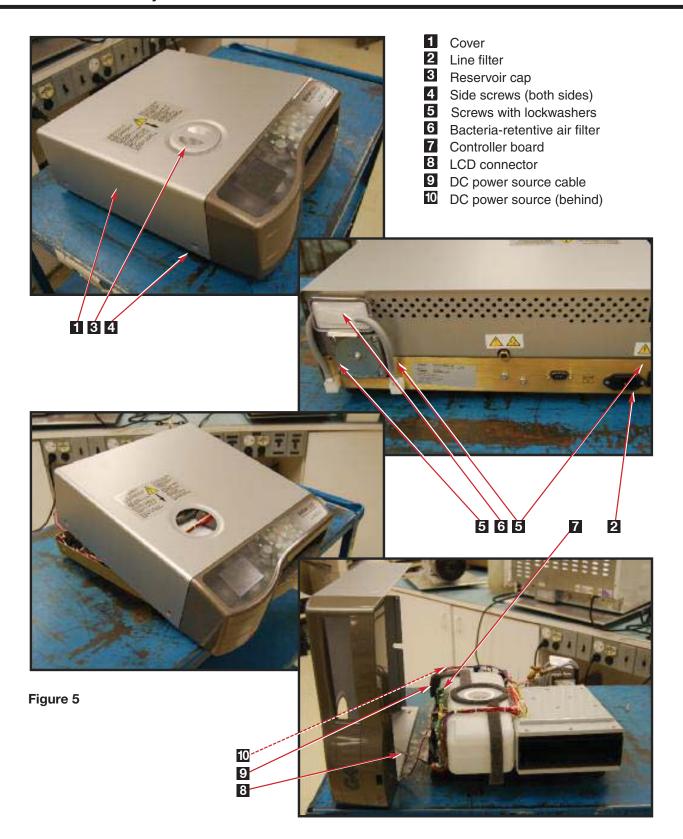

# Reinstalling the STATIM 2000 G4 Cover

To reinstall the cover, follow these steps (see Figure 5).

- 1. Inspect the armature gasket located on the inside front portion of the cover. If it is damaged or fails to adhere to the surface of the cover, it must be replaced. See, *Replacing the Armature Gasket* following.
- 2. Remove the reservoir cap (3) from the top of the reservoir.
- 3. Rest the cover (1) beside the left side of the unit. Reconnect the LCD connector (8) to controller board header P3, and the LCD's DC power source cable (9) to the DC power source (10) located at the back of the unit.
- 4. Reconnect the power cord (2) and power the unit ON to check that the LCD is functioning. After checking, turn the power switch OFF and unplug the power cord.
- 5. Carefully lift the cover from the work surface. While rotating the cover slide it forward until the fascia clears the controller board (7) and the front of the armature.
- 6. Tilt the front of the cover down and the rear of the cover upwards. Carefully reposition the front portion of the cover back over the front of the controller board and armature.

#### BE CAREFUL NOT TO PINCH THE RIBBON CABLE.

- 7. When the front of the cover is in place, lower the rear portion of the cover, and gently push backwards. Ensure the three tabs at the bottom of the fascia are properly lined up with the three holes in the front of the chassis. Carefully realign the screw holes at the back and reinsert three screws with lock washers (5) across the rear of the unit (TIP: start with the rear center screw).
- 8. Re-insert the four remaining screws (4), two on each side. Press firmly on the top of the cover to compress the gaskets and partially re-insert the screws as the holes align. When in place, tighten all the screws. Do not over tighten the screws.
- 9. Place the reservoir cap (3) back onto the reservoir.
- 10. Reconnect the power cord.

# Removing and Replacing the STATIM 2000 G4 Armature Gasket

To remove the armature gasket, remove the cover to access the inside of the fascia and peel the existing armature gasket from the fascia.


To replace the armature gasket follow these steps (see Figure 6).

- 1. Remove all traces of the gasket (1) and gasket adhesive from inside the fascia (2). The fascia is manufactured of polycarbonate plastic. **DO NOT USE AROMATIC SOLVENTS ON THE FASCIA**. This will damage the fascia. To prepare for installation of the new gasket, roughen the inside surface of the fascia using a fine grit sandpaper.
- 2. Leave the carrier paper on the adhesive side of the gasket and perform a trial fit. Note that the widest part of the gasket fits at the bottom of the opening and that the gasket does not protrude into the opening.
- 3. Remove the carrier paper from the adhesive side of the gasket. Carefully install the gasket making sure it is seated firmly in place.

1 Armature gasket

2 Fascia


3 Cover

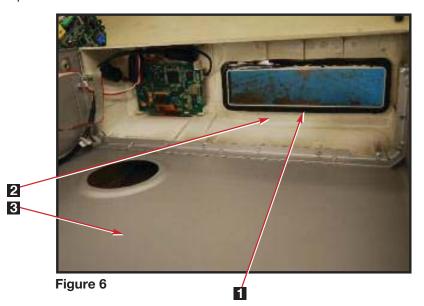

# Removing and Replacing the STATIM 2000 G4 LCD Overlay

To remove the LCD overlay, follow these steps (see Figure 7):

- 1. Using a sturdy sharp instrument, carefully lift one corner of the LCD overlay (1) and CAREFULLY peel it away from the fascia (2). USE CAUTION, there may be adhesive on the edges of the LCD overlay.
- 2. Remove any residual adhesive or membrane from the recessed area of the fascia.

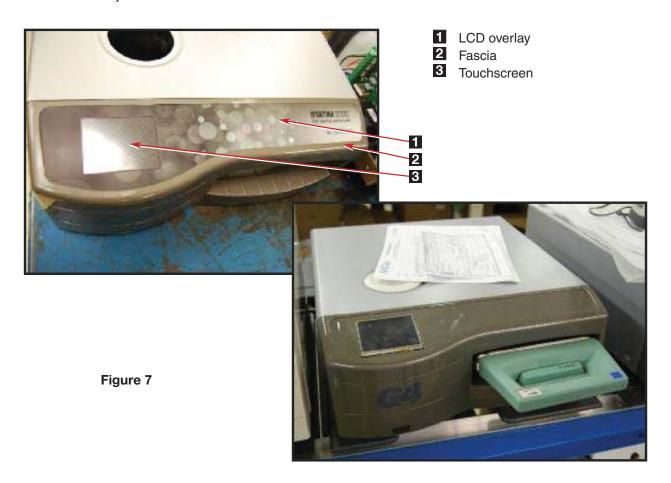

To replace the LCD overlay, follow these steps (see Figure 7):

- 1. The replacement colour LCD overlay has carrier paper on the back to protect the adhesive. Remove the carrier paper and align the bottom edge of the overlay with the bottom edge of the recessed area on the fascia. Carefully lower the overlay into place. With finger pressure only, press it into place, eliminating air pockets under the overlay.
  - **NOTE**: Do not press down in the middle of the colour LCD assembly. Apply pressure around the frame of the LCD.
- 2. Reconnect the power cord. Power the unit ON and test the touchscreen (3).
- 3. Reinstall the cover.

# Removing and Replacing the STAT*IM* 2000 G4 Colour LCD Touchscreen

Note: To replace the colour LCD touchscreen you will also need a new LCD overlay.

To remove the colour LCD touchscreen LCD overlay and remove the cover to access the inside of the fascia (see Figure 8).

- With the cover upside down on a work surface, disconnect the flexible PCB connector (1) cable connecting the colour LCD (2) to the LCD controller board (3). It is connected to the LCD controller board using a hinged catch (4). With your fingernail, pull on the door of the hinged catch to release the flexible PCB connector.
- 2. Turn the cover right side up and carefully remove the LCD overlay (5) from the fascia (6) and discard it.
- 3. The colour LCD touchscreen is held in place with two-sided adhesive tape. If it does not come off easily, put a small screwdriver in the PCB connector cable slit and push upwards.

To replace the colour LCD touchscreen, follow these steps:

- Align the colour LCD touchscreen in the recessed area in the fascia.
   NOTE: Do not press down in the middle of the colour LCD assembly. Apply pressure around the frame of the LCD.
- 2. Replace the LCD overlay with a new overlay.
- 3. Turn the cover upside down to access the inside of the fascia.
- 4. Connect the flexible PCB connector, pushing it down into position and closing the hinged connector. Give the flexible PCB connector a gentle tug to ensure it is properly held in position.
- 5. Connect the LCD controller board power cable to the DC power source and plug in the unit.
- 6. Power ON the unit to test the colour LCD touchscreen.


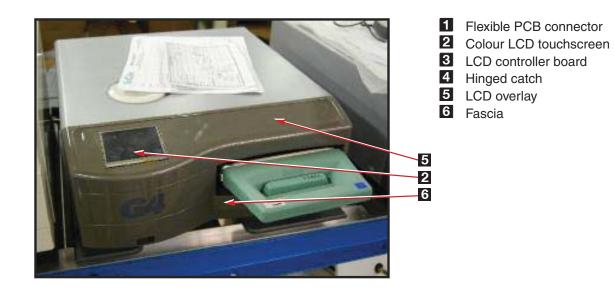


Figure 8

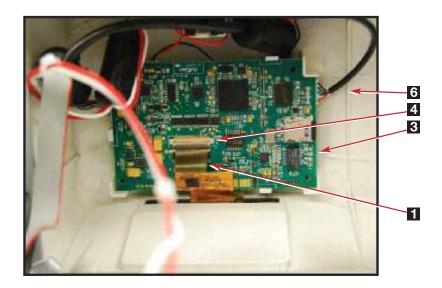

# Removing and Replacing the STATIM 2000 G4 Fascia

To remove the fascia, follow these steps (see Figure 9):

- 1. Place it right side up on a worksurface and remove the LCD overlay and the colour LCD touchscreen (see Removing LCD overlay and removing Colour LCD touchscreen)
- 2. Turn the cover upside down and remove the colour LCD controller board, speaker, Ethernet port and USB port (See removing LCD controller board, speaker, Ethernet port and USB port).
- 3. The fascia is attached to the cover (1) using seven plastite screws (2). Remove and retain the screws.
- 4. A bead of silicone was used to provide a water-tight seal between the fascia and cover during manufacturing. Remove the fascia and any silicon residue remaining on the cover.
- 5. Discard the old fascia. The fascia is made of recyclable polycarbonate. Please recycle wherever possible.

To replace the fascia, follow these steps (see Figure 8):

- 1. Apply a bead of silicon between the fascia and cover (1) to provide a water-tight seal.
- 2. Install the new fascia using the seven screws (2) retained from the disassembly procedure.
- 3. Reinstall the colour LCD touchscreen and LCD overlay.
- 4. Turn it upside down and reinstall the colour LCD controller board, speaker, USB and Ethernet ports.
- 5. Reinstall the cover.
- 6. Verify that the colour LCD functions correctly.
  - 1 Cover
  - 2 Plastite screws
  - 3 LCD controller board

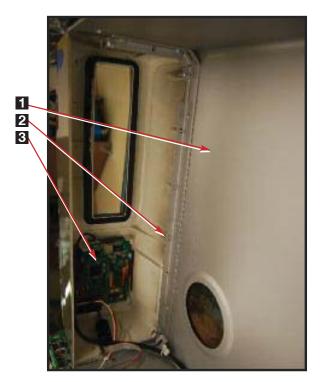


Figure 9

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



# **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

• If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



# PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# **Electrical and Electronic Components**

# **Controller Board**

# Identifying the Controller Board Type

There are different types of controller boards which may be encountered while servicing STAT*IM* 2000/2000S type units (see 'PCB revision' chart in Chapter 1. Identifying Your STAT*IM*). They can be identified as below (Figure 1):

Revision 3.x/4.x type board (1992 - 1995):





Revision 6.x type board (2004 -2007):



Revision 7.x type board (2007 to present):



Figure 1

#### Typical features:

- Revision number bottom right hand side
- Single, rectangular microprocessor

### Typical features:


- Revision number bottom right hand side
- Rectangular microprocessor plus EPROM
- 'W1' jumper for calibration
- Blue 'pressure interface/printer' connector

#### Typical features:

- Revision number bottom right hand side
- Square microprocessor plus EPROM
- 'W1' jumper for calibration
- Blue 'pressure interface/printer' connector

### Typical features:

- Up to revision 7.30 the revision number is on the top right hand side printed vertically
- From revision 7.40 onwards, the revision number is on the bottom left hand side (printed horizontally) under the connector J1
- Square microprocessor plus EPROM
- NO 'W1' jumper for calibration
- NO Blue 'pressure interface/printer' connector
- All components integrated on single board.
- Surface mount type component.
- 'Push In' yellow thermocouple connectors

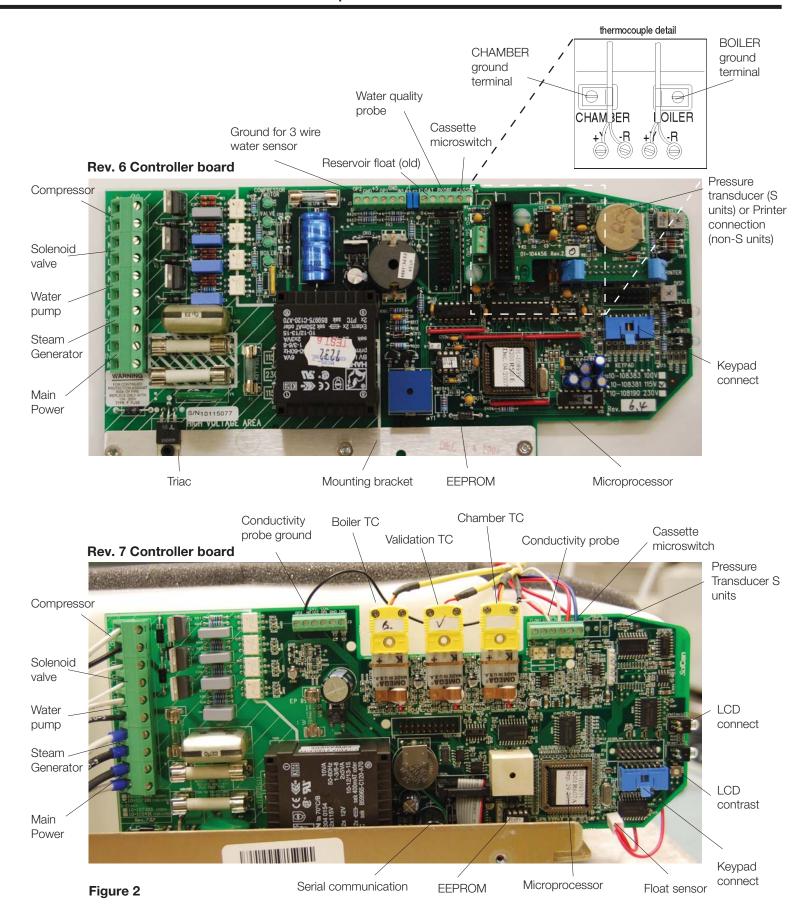


#### Removing the controller board assembly

To remove the controller board assembly, proceed as follows:

NOTE: Exercise extra caution. The thermocouple leads at the board end are very fragile.

1.

- b. For all boards up to and including Rev. 6.x -
  - Detach pressure interface board (1) to access thermocouple wires (S class units only),
    or if a validation thermocouple is present (later or upgraded non S units) detach the
    validation thermocouple interface board. Note that on some non S units where no
    validation thermocouple is present, there may be a printer cable attached in this
    socket. Remove this if present. Rev. 3.x/4.x units DO NOT have this socket.
  - Disconnect the steam generator thermocouple wires (2) from controller board terminal positions BOILER +Y and -R.
  - Disconnect the chamber thermocouple wires (4) from controller board terminal position CHAMBER +Y and -R.
  - Disconnect the ground terminal (3) from the position marked BOILER. Disconnect the ground terminal (5) from the position marked CHAMBER. Retain the ground terminal screws. Leave the screws with contact washers in the terminals. Carefully bend the leads so they do not contact the controller board when it is removed.

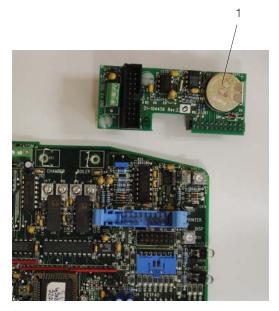



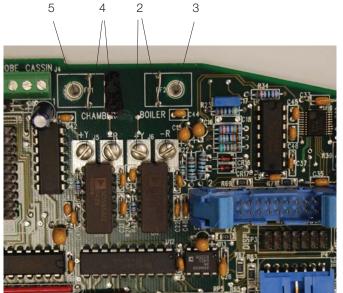


Figure 4

- 1. Pressure interface board
- 2. Steam generator terminals
- 3. Boiler flag terminals
- 4. Chamber thermocouple terminals
- 5. Chamber flag terminals

#### c. For all boards Rev. 7.x onwards -

- Unplug the thermocouple connectors (6). There are three for non S ('boiler', 'chamber' and 'validation') and two for S class ('boiler' and 'chamber').
- Unplug the float connector (7) on units where the water quality and float reservoirs are fitted.

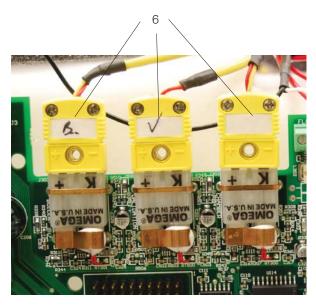




Figure 5

Figure 6

7

- 6. Thermocouple connectors
- 7. Float connector

#### **All Boards**

- 2. Remove the microswitch leads from controller board header terminal positions labeled CASSIN (8).
- Remove the water quality sensor leads (if present) from controller board terminal positions labeled PROBE (9) or the float sensor leads (if present) from terminal positions labeled FLOAT (10). Note that on later units a third, ground wire (11) may be present as well as the two PROBE wires.
- 4. Disconnect all high voltage leads from controller board block terminal J1 (12).
- 5. Remove and retain the three screws with washers that secure the Controller Board assembly to the chassis.
- 8. CASSIN Terminal
- 9. PROBE Terminal
- 10. FLOAT terminal
- 11. Groundwire (some units)
- 12. J1 terminal block



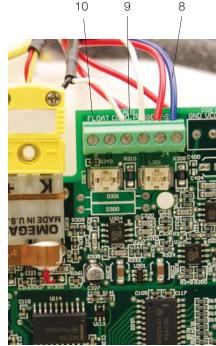




Figure 7


#### Replacing/Refitting the Controller Board Assembly

To replace/refit the controller board assembly, proceed as follows:

- 1. Visually inspect the controller board assembly to ensure that there is sufficient thermal compound on the mounting bracket directly below the steam generator triac.
- 2. Install the controller board assembly using the three screws with washers retained from removal (or fit new if replacing with a new board).
- 3. Connect the high voltage leads from the power switch, steam generator, pump, solenoid valve and compressor to controller board J1 terminal (12). Each lead is numbered with a corresponding J1 designation. See Table B.

TABLE B: CONTROLLER BOARD HIGH VOLTAGE J1 CONNECTOR								
Marking	Fuse Value	Device	Signal	Wire Colour				
COMP	J1-10 J1-9	COMPRESSOR	N L	0.25 x 1.25 inch				
VALVE	J1-8 J1-7	SOLENOID VALVE	N L	WHITE or RED WHITE or RED				
PUMP	J1-6 J1-5	PUMP	N L	WHITE BLACK				
BLR	J1-4 J1-3	STEAM GENERATOR	N L	WHITE BLACK				
INPUT	J1-2 J1-1	POWER SWITCH	L N	BLACK WHITE				

- 4. Connect the microswitch leads to controller board header positions labeled CASSIN (8).
- 5. Connect the water quality sensor leads (if present) to controller board positions labeled PROBE (9) or the float sensor leads (if present) to terminal positions labeled FLOAT (10). Note that on later units a third, ground wire may be present as well as the two PROBE wires.

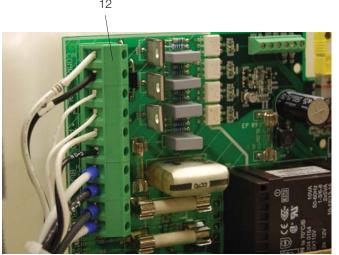

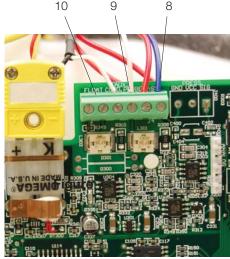


Figure 8

6.

- a. For all boards up to and including Rev. 6.x -
  - Check the pre-bend on both thermocouple leads to ensure that they are the required shape to go under the washers on the connections without touching anything other than the terminal.

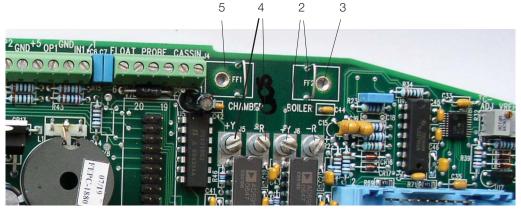


Figure 9

- Note: One lead is colour coded: the unmarked lead is positive, +Y. the red lead is negative, -R. Connect the unmarked lead to the terminal marked +Y on the Controller Board. Connect the red lead to the terminal marked -R on the controller board. Make sure there is extra lead length so that the wires do not break as the screw is tightened. Ensure the wires are seated securely behind the respective washers.
- Reconnect the steam generator thermocouple wires to the controller board terminal positions BOILER +Y and -R (2).
- Reconnect the chamber thermocouple wires to the controller board terminal positions CHAMBER +Y and -R (4).
- Reconnect the ground terminal from the position marked BOILER (3). Reconnect the ground terminal from the position marked CHAMBER (5).

The two leads must not touch one another or any other component. Do not calibrate a thermocouple until it is properly installed and positioned in the unit.

- b. For all boards Rev. 7.x onwards -
  - Plug in the thermocouple connectors.
     There are three for non S ('boiler', 'chamber' and 'validation') and two for S class ('boiler' and 'chamber') (6).
  - Plug in the float connector on units where the water quality + float reservoirs are fitted (7).

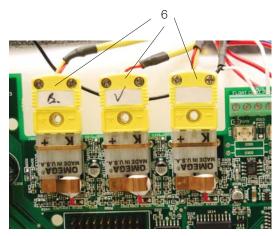


Figure 10

The chamber and steam generator thermocouples MUST be recalibrated. Do not calibrate a thermocouple until it is properly installed and positioned in the unit.

Reattach the pressure interface board, (S class units only), or if a validation thermocouple is present (later or upgraded non S units) the validation thermocouple interface board. Note that on some non S units where no validation thermocouple is present, a printer cable is attached in this socket. Replace this if present. Rev.3.x/4.x units DO NOT have this socket.



7

Figure 11

8. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity should be performed on the STAT*IM* unit at this stage.

**NOTE:** These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.

#### **Controller Board Fuse Values**

Always replace blown fuses with a fuse of the same size and value as indicated directly on the controller board.

## STATIM 2000 G4 Electronics

#### Removing and Replacing the STATIM 2000 G4 LCD Controller Board

**NOTE**: After replacing the LCD controller board, the bubble level must be calibrated. Follow bubble level calibration instructions below.

To remove the LCD controller board, remove the cover to access the inside of the fascia (see Figure 12).

- 1. Disconnect the LCD controller board (1) from the P3 connection on the main controller board and disconnect the LCD's DC power source cable (2) from the DC power source (3) located at the back of the unit.
- 2. Place the cover upside down on a work surface and disconnect the flexible PCB connector (4) cable connecting the colour LCD to the LCD controller board. To do this, slide your fingernail into the centre of the hinged catch (5) and pull away from the board. Once open, the hinge will release the connector.
- 3. To remove the LCD controller board, start by freeing it from the bottom and right clips (6) and lift it out. Remove the cable ferrites from the holding clamps.
- 4. Disconnect the USB cable (7).
- 5. Disconnect the Ethernet port cable (8).
- 6. Disconnect the wires to the speaker (9).

To replace the colour LCD controller board, follow these steps (see Figure 5).


- 1. Re-connect the wires to the speaker (9).
- 2. Re-connect the Ethernet port cable (8).
- 3. Re-connect the USB cable (7).
- 4. Clip the logic board back into position. NOTE: Be certain the board is properly clipped into position. An improperly seated board will affect the functioning of the level.
- 5. Connect the flexible PCB connector, pushing it down into position and closing the hinged connector. Give the flexible PCB connector a gentle tug to ensure it is properly held in position.
- 6. Connect the LCD controller board to the P3 connection on the main controller board, connect the LCD DC power cable to the DC power source and plug in the unit.
- 7. Power ON the unit to test the colour LCD touchscreen

To calibrate the bubble level, proceed as follows:

- 1. Set the unit on a flat surface. Use a bubble level placed on the cover of the unit to ensure the unit is level.
- 2. Access the service menu as described in Chapter 3 in the section titled Using the service menu on the STATIM G4.
- 3. Scroll to the Bubble level calibration screen.



- 4. Press the check mark until the bubble is properly centered in the middle of the circle (wait a few seconds after pressing to allow it to reset).
- When the bubble stays within the screen's center circle, the bubble level is calibrated. Press STOP to exit.

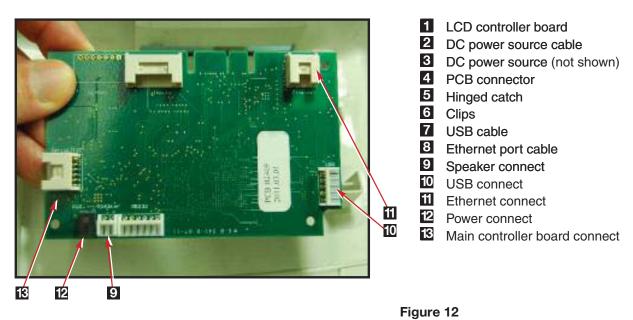


2 8 7

#### Removing and Replacing the STATIM 2000 G4 USB Port

To remove the USB port, first remove the cover to access the inside of the fascia and remove the LCD controller board. Then follow these steps (see Figure 13):

- 1. Disconnect the USB port cable (1) from the LCD controller board (2).
- 2. Lift the tab holding the USB port assembly (3) and pull the assembly out. To replace the speaker, reverse the removal instructions.

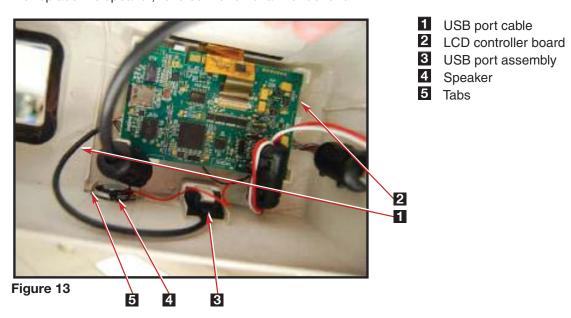

#### Removing and Replacing the Speaker

To remove the speaker, first remove the cover. Remove the LCD controller board and place the cover upside down on a worksurface. Then follow these steps (see Figure 6):


- 1. Disconnect the speaker (4) from the LCD controller board (2).
- 2. The speaker is hot-glued into the tabs. You will have to remove the glue to pull it out. To replace the USB port, reverse the removal instructions (use hot glue to secure speaker).

#### Removing and Replacing the Ethernet port

To remove the Ethernet port, first remove the cover and place it upside down on a worksurface. Remove the LCD controller. Then follow these steps (see Figure 14):


- 1. Disconnect the Ethernet port cable (1) from the LCD controller board (2).
- 2. On the inside rear of the cover, remove the fastening screw (3) holding the Ethernet port assembly in place.
- 3. Unhook the Ethernet port assembly tab (4) from the rear of the cover and remove cable by cutting the cable tie.

To replace, reverse the removal instructions.

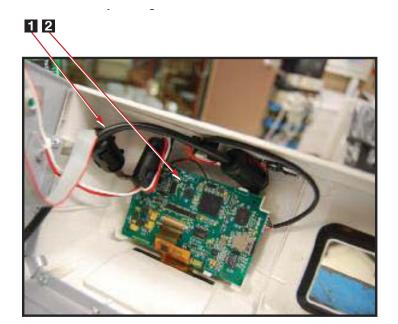

- 1 Ethernet port cable
- 2 LCD controller board
- 3 Screw
- 4 Ethernet port assembly tab

Figure 14

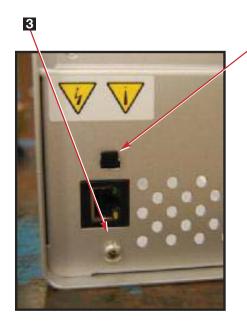




## Microprocessor

Microprocessors are not interchangeable between controller boards of different part numbers and/ or revision, unless otherwise indicated.

#### **Processor Types**

**Note:** a visual guide to the location on the controller board can be found in the section 'Identifying Controller Board Types'

Revision 2.x/5.x, 3.x/4.x controller boards all have the following main processor configuration.

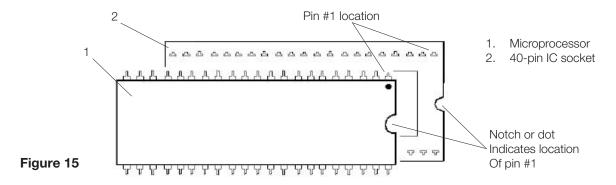

Revision 2.x/5.x controller boards also have a supplementary EEPROM located near to the main processor. (Rev. 3.x/4.x DO NOT have this supplementary device).


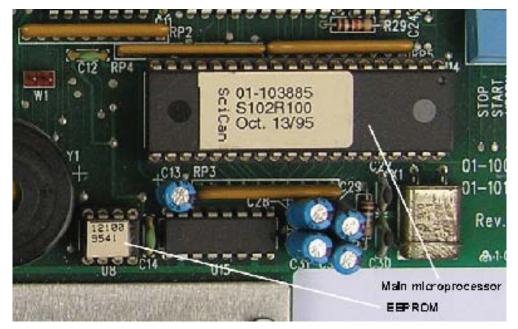



Figure 16

The microprocessor and EEPROM devices are a matched pair set. These devices are supplied together and must be installed together. Any attempt to substitute either device on its own will result in a continuous beeping tone.

Revision 6.x and 7.x controller boards all have the following main processor configuration.

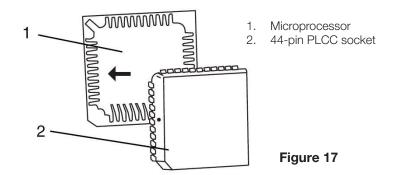

Both Revision 6.x and 7.x of controller boards also have a supplementary EEPROM located near to the main processor.

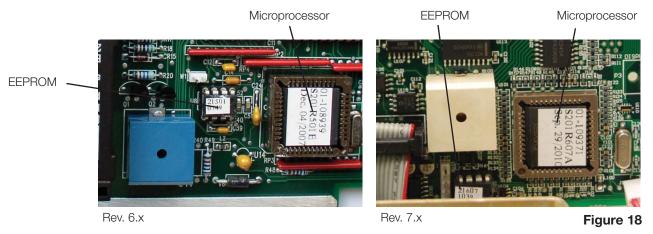

# Removing the Microprocessor (and Pre-Programmed EEPROM where fitted).

Before removing or replacing the devices observe appropriate electrostatic discharge precautions for the work area and technician. Ensure that the unit is powered OFF.

To remove the Microprocessor and pre-programmed EEPROM device if present, proceed as follows:

1.

 For revision 2.x/5.x, 3.x/4.x controller boards — Remove the microprocessor from Controller Board socket using a 40-pin IC puller. Remove the EEPROM device from Controller Board socket using an 8-pin IC puller.

DISCARD the microprocessor and EEPROM devices.

Note: The orientation of Pin 1 of the socket and Pin 1 of the EEPROM

c. For revision 6.x and 7.x controller boards AND units with adaptor boards installed as part of an Alex steam generator upgrade — Remove the microprocessor from Controller Board Socket using a 44-pin PLCC IC puller. Remove the EEPROM device from Controller Board socket using an 8-pin IC puller.

DISCARD the microprocessor and EEPROM devices.

# Replacing/refitting the Microprocessor (and Pre-Programmed EEPROM where fitted).

The microprocessor and EEPROM devices are a matched pair set. These devices must be ordered and installed together. An attempt to substitute either device singly results in a continuous beeping tone.

Note: The orientation of the various main microprocessor types (see Figures 17 and 18 above) to ensure correct alignment. Ensure that the device pins are fully inserted into the socket. Incorrectly installed IC devices may cause damage to the unit.

ALL EEPROM devices are notched in a similar manner to the 40-pin microprocessor.

To replace/retrofit the microprocessor, proceed as follows:

- Determine the part number and revision of the installed controller board. (This number appears
  on the component side of the board in the lower right-hand corner.) Determine the rated voltage
  of the unit by examining the serial number label. Use this information to find the appropriate
  microprocessor replacement kit
- 2.
- a. For revision 2.x/5.x, 3.x/4.x controller boards, use an insertion tool to install the microprocessor into controller board socket. Note the orientation of pin #1 of the socket and microprocessor.
- b. For revision 6.x and 7.x controller boards AND units with adaptor boards installed as part of an Alex steam generator upgrade, the microprocessor can be inserted into controller board socket by hand. Note the orientation of pin # 1 of the microprocessor and align it with the notch in the PLCC socket.
- 3. For all controller boards where an EEPROM is fitted, using an insertion tool, insert the EEPROM device into controller board socket. Note the orientation of pin #1 of the socket and EEPROM.
- 4. Connect the keypad connector to controller board.
- 5. Connect the LCD connector to controller board.
- 6. Connect the printer connector to controller board, if present.
- 7. Power the unit ON. Ensure that the version number displayed briefly when the unit is first powered (not Rev. 3.x/4.x) matches the version number printed on the microprocessor. If the LCD fails to display the "select a cycle" message, review the wiring connector placement and check that the microprocessor and EEPROM are positioned properly in the sockets.
- 8. Calibrate the unit after the installation of a new microprocessor.

## Microswitch


The function of the microswitch is to sense the presence of the cassette when it is inserted in the armature. When the cassette activates the switch, it is possible to start the unit. With the cassette out, the unit will not start.

Failure of the switch will produce the following symptoms:

- Failed 'ON' the unit may start WITHOUT the cassette being engaged.
- Failed 'OFF' the unit will not start even with the cassette fully engaged.

#### Removing the Microswitch

To remove the microswitch, proceed as follows (see Figure 19):

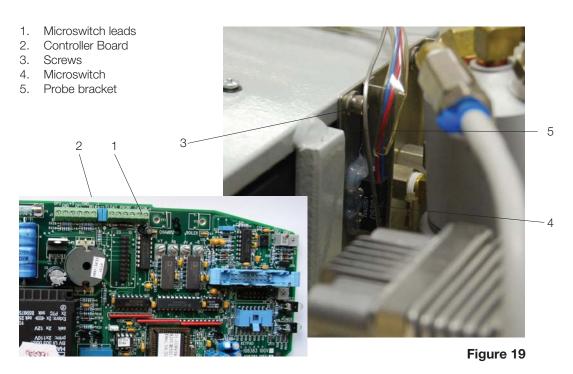

- 1. Turn the power switch OFF, and unplug the unit.
- 2. Disconnect the microswitch leads (1) from the controller board (2) terminal block positions labelled CASSIN J4-1 AND J4-2 (J3-1 and J3-2 on Rev. 3.0 boards) and carefully cut the cable ties securing the microswitch leads to the wiring harness.
- 3. There may be RTV silicone sealant on the heads of the two screws (3) securing the microswitch (4) to the side of the probe bracket (5). Trim away the excess RTV silicone and remove the screws. Retain these for assembling the replacement microswitch to the probe bracket.
- 4. Remove the microswitch.
- 5. Remove all RTV silicone residues from the probe bracket surface.

#### Replacing the Microswitch

To replace and refit the microswitch, proceed as follows (see Figure 19):

- 1. Using the two screws (3) retained from the removal procedure, attach the microswitch to the probe bracket with the hinge end of the switch arm up and facing towards the front of the unit.
- 2. Connect the microswitch leads (1) to controller board (2) block terminal positions labelled CASSIN. The red wire is inserted into CASSIN J4-1 (J3-1) for Rev 3.0 boards), the blue wire is inserted into CASSIN J4-2 (J3-2 for Rev 3.0 boards). Secure the leads.
- 3. Route the microswitch leads back into the wiring harness and secure the complete harness using cable ties every 2-3 inches. Secure these to the top of the armature using the clamps provided.
- 4. There are two methods of testing the microswitch.
  - a. Plug-in the power cord and turn the power switch ON. When the unit is powered ON and the cassette is inserted, the solenoid valve is activated. Carefully insert the cassette. If the solenoid valve "clicks" before the cassette is fully inserted, the microswitch needs to be adjusted. If the solenoid valve does not "click" after the cassette is fully inserted, the microswitch needs to be adjusted or it is defective.
  - b. Using a multimeter set at continuity, touch the meter probes across controller board terminal positions labelled CASSIN. If continuity is achieved before the cassette is fully inserted, the microswitch needs to be adjusted. If continuity is not achieved after the cassette is fully inserted the microswitch needs to be adjusted or it is defective, see Adjusting the microswitch.
- 5. Once installed and functioning, apply a thin bead of RTV silicone sealant along all edges of the microswitch that are in contact with the probe bracket. Allow the silicone to cure per manufacturer's instructions.
- 6. Run a sterilization cycle and observe the area around the microswitch for leaks.

#### Adjusting the Microswitch


1. If adjustment is required, power the unit OFF. Remove the screws securing the microswitch to the probe bracket and remove the microswitch. Using fine needle-nose pliers, adjust the activation point of the switch by carefully bending the microswitch arm. Reinstall and re-test the switch.

## Mains Components

## **Power Cords**

You will encounter two types of power cord on STATIM units:

- Hard wired type, which are attached through the back of the chassis using a bulkhead style strain relief and are wired directly to the power switch using 'fast-on'/'spade' terminals.
- Detachable type, which are modular and plug into a panel mounted A.C. inlet receptacle or receptacle / line filter in the back of the chassis.

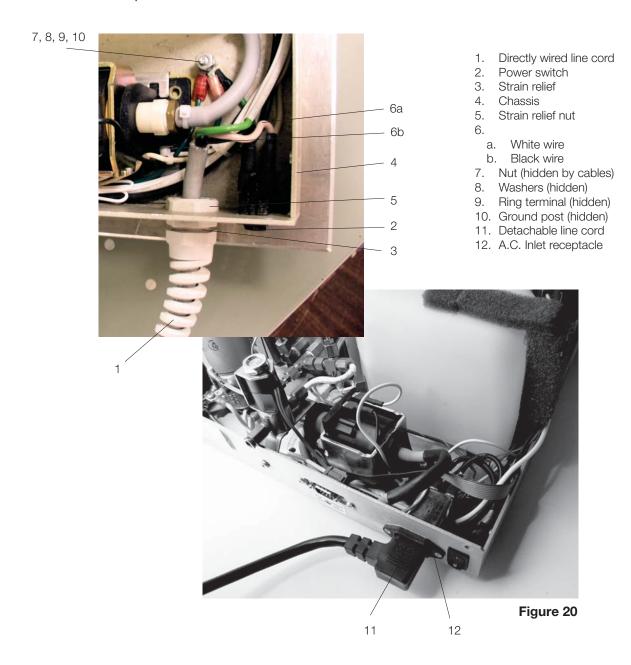

#### Removing a Directly Wired Power Cord

To remove a power cord (1) wired directly to the power switch (2), proceed as follows (see Figure 20):

- 1. Turn the power switch OFF, and unplug the power cord.
- 2. Use a 7/8-inch wrench and a 15/16-inch wrench or two adjustable wrenches to undo the bulkhead strain relief (3). Do not crush any wires between the inside of the chassis (4) and the wrench. Retain the attaching nut (5) from inside the chassis.
- 3. Remove the white wire (6a) from the upper spade terminal, labelled1, and the black wire (6b) from the upper spade terminal, labelled 2, on the power switch.
- 4. Remove the nuts (7) and washers (8) securing the ground wire ring terminal (9) to the ground post (10) using an 11/32-inch wrench.

Several other ground wires are attached at this point. Note the sequence of washers, nuts and ring terminals.

The line-in ground / earth wire should always be in immediate contact with the chassis and secured with a separate washer and nut.

5. Remove the power cord from the unit.

#### Replacing a Directly Wired Power Cord

Use replacement cords with appropriate safety agency ratings and approvals only. The line-in ground / earth wire should always be in immediate contact with the chassis and secured with a separate washer and nut.

Do not attach the input leads from a power cord to any other leads or components in the unit.

To replace/reinstall a power cord, proceed as follows (see Figure 20):

- 1. Ensure the power switch is OFF, and the line cord is unplugged from the wall outlet.
- 2. Pass the end of the power cord with the spade terminals through the opening in the chassis (4) and place the attaching nut (5) over the cord leads.
- 3. Thread the attaching nut to the bulkhead strain relief (3) but do not tighten. The length of cord inside the chassis may require adjustment. A minimum of 5/8-inch / 16 mm of the external cover of cord insulation should project into the unit when the stain relief is positioned correctly. When the cord is adjusted, finger tighten the attaching nut, and then tighten one full turn using a 7/8-inch wrench on the external part of the strain relief and a 15/16-inch wrench on the attaching nut.
- 4. Attach the green ground / earth wire ring terminal (9) to the ground post (10). Note the order of terminals, washers (8) and nuts (7) during disassembly. Use an 11/32-inch wrench to tighten each nut securely.
- 5. Attach the terminal from the black wire (6b) of the power cord to the upper power switch spade terminal labelled 2, nearest the pump.
- 6. Attach the terminal from the white wire (6a) of the line cord to the upper power switch spade terminal labelled 1.
- 7. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STAT*IM* unit at this stage.
  - **NOTE:** These tests must be performed on the STAT*IM* again once the work is completed and the cover has been returned to the unit.
- 8. Plug the line cord into the wall receptacle and turn the power switch ON. Observe the LCD and indicator lights to determine that power is present.

#### Removing a Detachable Power Cord

To remove a detachable power cord (11) connected to the unit through a panel mount A.C. inlet receptacle, proceed as follows (see Figure 20):

- 1. Turn the power switch (2) OFF, and unplug the power cord from the wall outlet.
- 2. Disconnect the power cord from the panel mount A.C. inlet receptacle (12).

#### Replacing/refitting a Detachable Power Cord

Use replacement cords with appropriate safety agency ratings and approvals only.


To replace/reinstall a detachable power cord (11), proceed as follows (see Figure 20):


- 1. Select the appropriate replacement cord from the spare parts list.
- 2. The A.C. inlet receptacle (12) is keyed. Note the shape of the cord connector and the corresponding shape of the receptacle. Plug in the connector.
- 3. Plug the cord into the wall receptacle and turn the power switch (2) ON. Observe the LCD and indicator lights to determine that power is present.

## Power (on/off) Switch

#### Removing the Power Switch

To remove the power switch (1), proceed as follows, (See Figure 21):

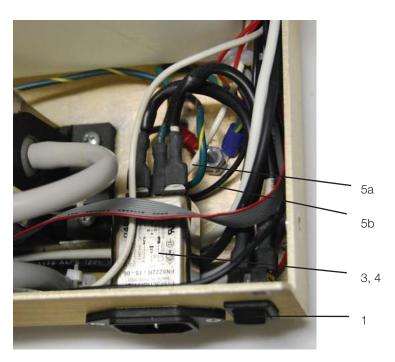

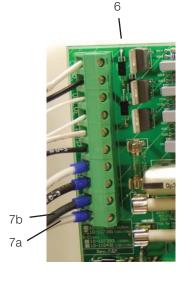

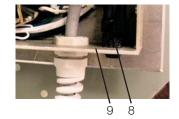


Figure 21


- 1. Power switch
- 2. Line cord (not shown)
- 3. A.C. receptacle
- 4. Line filter
- 5.
  - a. White wire
  - b. Black wire

- 6. Controller Board
- a. White wire (J1-1)
- b. Black wire (J1-2)
- 8. Power switch tabs
- Chassis wall
- 1. Turn the power switch OFF, and unplug the power cord from the wall outlet.
- 2. Remove the cord (2) from the unit. See, Removing a Directly Wired Power Cord or Removing a Detachable Power Cord.
- 3. If an A.C. inlet receptacle (3) or receptacle / line filter (4) is present, disconnect the leads to the power switch. At the switch, disconnect the white wire (5a) from the upper spade terminal, labelled 1, and the black wire (5b) from the upper spade terminal, labelled 2, nearest the pump.
- 4. At the switch, disconnect the white wire (7a) extending from controller board (6) terminal J1-1 to the lower power switch spade terminal labelled 1a and the black wire (7b) extending from controller board terminal J1-2 to the lower power switch spade terminal, labelled 2a, nearest the pump.
- 5. The panel mount style switch is held into the panel with tabs. Compress the tabs (8) on the top and bottom of the switch and push the disconnected power switch out of the chassis (9) wall.

#### Replacing/refitting the A.C. Power Switch

To replace/refit the power switch (1), proceed as follows, (See Figure 22):

- 1. Orient spade terminals 1a and 2a downwards and press the power switch into the clearance hole in the chassis. Apply pressure evenly top and bottom until the bezel rests against the chassis (7) wall.
- 2. Connect the white wire (fig. 21, item 7a) extending from controller board terminal J1-1 to the lower right-hand power switch spade terminal labelled 1a and the black wire (fig. 18, item 7b) extending from controller board terminal J1-2 to the lower left-hand power switch spade terminal, labelled 2a.
- 3. Connect the white wire (4a) from either the attached power cord, the A.C. inlet receptacle or the A.C. inlet receptacle / line filter to the upper right-hand power spade terminal labelled 1, and the black wire (4b) to the upper left-hand spade terminal labelled 2, nearest the pump.
- 4. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STAT/M unit at this stage.
  - **NOTE:** These tests must be performed on the STAT*IM* again once the work is completed and the cover has been returned to the unit.
- 5. Plug the power cord (3) into the wall receptacle and turn the power switch (2) ON. Observe the LCD and indicator lights to determine that power is present.

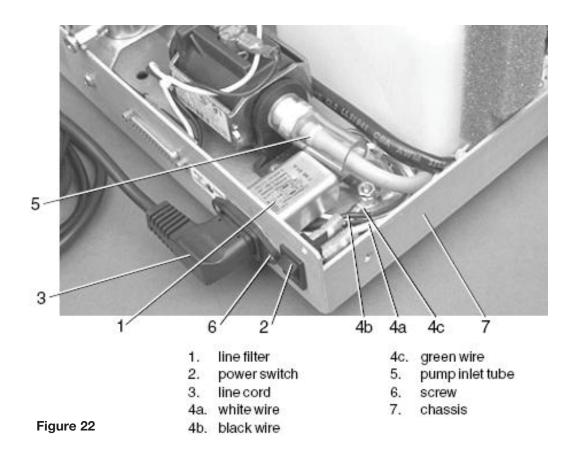


## Receptacle / Line Filter

Some STAT*IM* units contain an A.C. inlet receptacle / line filter. It is difficult to determine whether a line filter has failed or not. If the unit blows mains fuses in the service panel, there may be a short in the line filter. Disconnect all leads from the mains input and output and test for short circuits.

#### Removing the Receptacle / Line Filter

To remove a line filter (1), proceed as follows, (See Figure 23):

- 1. Turn the power switch (2) OFF, and unplug the power cord (3) from the wall outlet.
- 2. Remove the power cord from the unit. See, Removing a Detachable Power Cord.
- 3. Disconnect the white wire (4a) from receptacle line filter position N, the black wire (4b) from line filter position P and the green wire (4c) from line filter position. If access to the back of the line filter is impeded, disconnect the pump inlet tube (5).
- 4. Remove the screws (6) holding the filter to the chassis (7) and remove the filter.

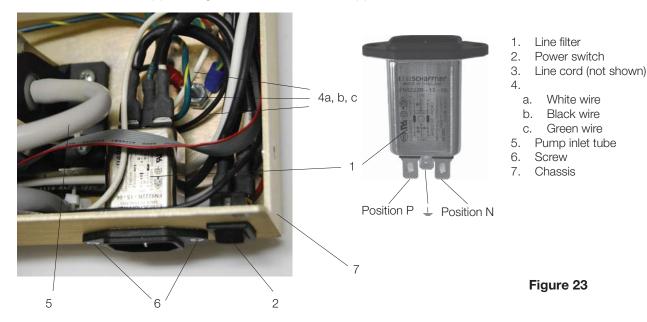

#### Replacing/refitting the Receptacle / Line Filter

To replace/refit a line filter (1), proceed as follows, (See Figure 23):

- 1. Insert the line filter in the opening in the chassis (7). The P and N fast-on spade terminals should face up, and the ground terminal should face down.
- 2. Insert and tighten the screws (6) holding the line filter, using Loctite in the threaded holes in the chassis.
- 3. Connect the white wire (4a) from the power switch to line filter position N and the black wire (4b) from the power switch to line filter position P.
- 4. Connect the green wire (4c) from the ground post to line filter position.
- 5. If the pump inlet tube (5) was moved, reinstall the tube.

- A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.
   NOTE: These tests must be performed on the STATIM again once the work is completed and
  - the cover has been returned to the unit.
- 7. Plug the power cord (3) into the line filter and the wall receptacle and turn the power switch (2) ON. Observe the LCD and indicator lights to determine that power is present.


#### LCD DC Power Source for STATIM 2000 G4

# Removing and Replacing the STAT*IM* 2000 G4 LCD DC Power Source

To remove the LCD DC power source, remove the cover and place it on its side next to the unit. Then follow these steps (see Figure 24):

- 1. Disconnect the DC power source wires (1) from the power source (2) located at the back of the unit.
- 2. The DC power source bracket is held to the chassis by two screws, one from the back of the chassis. Remove the remaining screw (4) and cut the cable tie to lift clear from the unit.

To replace the LCD DC power source, reverse the removal instructions.

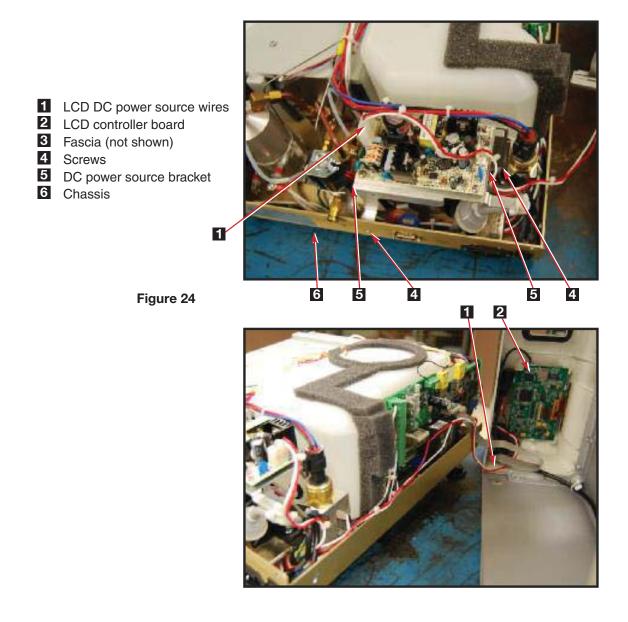

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



## **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



## PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

## Pumps, Reservoir, and Compressor

## **Water Pumps**

STAT*IM* units have been manufactured using three different pumps. Before servicing, you must first determine which model of pump the unit contains (See Figure 1).


**Important Note:** Always make sure that there is sufficient steam-process distilled water in the reservoir prior to testing STAT*IM* pumps.

#### **Identifying Pump Types**



#### Fluid-O-Tech Pump

The Fluid-O-Tech pump is generally factory fitted to older units with revision 3.x/4/x controller boards up to approximately 1994-1995. It is a small rectangular pump with a diode plug attached to the top. The pump rests on a vibration-damping pad on the bottom of the chassis and is held in place with a wire bracket. The pump is secured to the bracket by a cable tie at each end. The black lead is attached to controller board J1-5 (LINE) and the white lead is attached to controller board J1-6 (NEUTRAL).

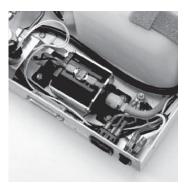

#### SciCan Pump

The SciCan pump is generally factory fitted to units with revision 2.x/5.x/6.x and some early 7.x controller boards from 1995 up to mid 2008. It has a similar footprint to the Fluid-O-Tech pump, however the diode is now integrated into the pump. The pump is mounted on two vibration damping rubber brackets. Each bracket is secured to the chassis with two screws. The black lead is attached to controller board J1-5 (LINE) and the white lead is attached to controller board J1-6 (NEUTRAL).



#### Ulka Pump

The Ulka pump is factory fitted to units with revision 7.x controller boards only from mid 2008 onwards. It is slightly smaller than the SciCan pump and can be identified by the configuration of the wires on top of the pump. The two visible connections are located adjacent and parallel to each other whereas the two visible connections on the SciCan pump are opposing each other and have a thermal switch connecting them. The pump is mounted on two vibration damping rubber brackets. Each bracket is secured to the chassis with two screws. The black lead is attached to controller board J1-5 (LINE) and the white lead is attached to controller board J1-6 (NEUTRAL).

Figure 1

Note: Important information on pump interchangeability:


- The Fluid-O-Tech pump is no longer available, but units with this pump can be retrofitted with the SciCan pump ONLY.
- The SciCan pump and Ulka pump are NOT interchangeable. They are both still available. If replacing the pump during maintenance or servicing, the same pump must be fitted to ensure the unit works correctly.

#### Testing pump flow (for all types)

To test a water pump, proceed as follows (see Figure 2):

The control box should be attached to the unit, or select 'pump' in the device test sub menu of the revision 7.x controller board service menu (if cover is removed but still connected). The unit should be powered on with the cassette removed.

- Disconnect the Teflon™ steam generator inlet tube (1) from the top of the steam generator (2) using a 3/8-inch wrench.
- 2. Connect the disconnected end of the steam generator inlet tube to the pump tester (3). Do not cross thread the fittings. Do not overtighten.
- 3. Using the button on the control box or appropriate keypad button (revision 7.x controllers), activate the pump for 2 seconds to purge any air that is trapped in the fittings. Empty any water that enters the pump tester. Recap the pump test bottle. The pump tester MUST be empty before starting the pump test.
- 4. Be prepared to record the time it takes for the water level displayed on the pump tester to reach the top of the line marked MIN (see below). Activate the pump.
- 5. If the water level reaches the top of the line marked MIN within the time allotted in the Pump Flow Parameters chart (on next page), the pump is performing correctly. If the water level reaches the top of the line marked MIN in less than that time or if the water level does not reach the top of the line marked MIN within the allotted time, follow the Pump Recovery Repair procedures described later in this chapter.

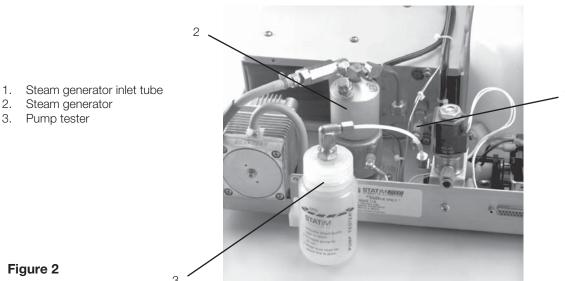


#### Pump flow parameters

The following pump parameters are valid for this procedure.

Unit Voltage Rating	100 V 50Hz		100 V 60 Hz		110 V 60 Hz		230 V 50 Hz	
For units fitted with alluminum steam generators, the following times apply.								
Pump type	min.	max.	min.	max.	min.	max.	min.	max.
Fluid-O-Tech Pump	23	24	25	26	25	26	23	24
SciCan Pump	20	21	23	24	23	24	23	24
Ulka Pump	20	21	23	24	23	24	23	24
For units fitted with Stainless Steel steam generators, the following times apply.								
SciCan Pump	N/A	N/A	N/A	N/A	21.5	23	17.0	18.5

Note: max. / min. numbers represent time in seconds required to fill the bottle to the min. line.

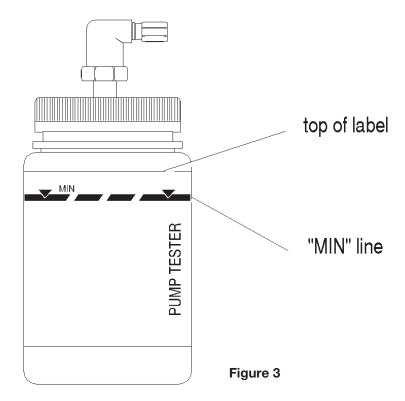


## Water Pump Recovery Repair Procedure

Pump performance may be compromised because of foreign material caught inside the assembly, blockage of an orifice or mechanical wear of internal parts. Before replacing a pump that is performing poorly, attempt to recover the unit by cleaning the filter(s) and/or resizing the pump tube.

#### **Performing Filter Maintenance**

Early STATIM 2000 / 2000S units fitted with SciCan pumps have a coarse mesh filter in the inlet fitting and a fine mesh filter in the outlet elbow fitting located internally in the pump.



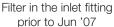




Filter on the outlet fitting prior to Jun '07

Figure 4

In 2007, SciCan removed these two filters due to the implementation of the water reservoir filter in 2006 and subsequent introduction of the in-line filter in 2007.

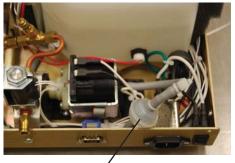

In-line filter


Figure 5

The in-line filter was added between the water reservoir and the water pump inlet to address the issue of customers not using the reservoir filter. The purpose of this in-line filter is to stop debris from entering the STAT*IM* plumbing system when the reservoir filter is not used. The in-line filter kit is available as part # 01-106637S and includes the filter and all necessary tubing. This filter should be replaced every 2000 cycles or every two years.

#### Performing internal pump filter maintenance

For filter maintenance on internal filters located in the pump, proceed as follows (see Figure 6):

Tip: This can be undertaken with the pump in place, but it is easier if the pump is removed from the unit.

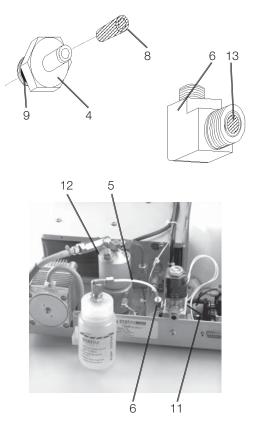

6, 7 (obscured) 10 2

- Reservoir
- Cable tie
- Pump inlet tube
- Pump inlet fitting
- Teflon™ tube
- 6. Right-angle pump outlet fitting
- 7. Pump body fitting
- 8. Mesh filter, integral to item 4
- 9. Rubber seal, integral to item 4
- 10. Pump
- 11. Pump body outlet fitting (obscured)
- 12. Steam generator
- 13. Insert with mesh filter, integral to item 6


Figure 6

If the choice is to have the pump still located in the unit:

- 1. Drain the reservoir or crimp the pump inlet hose to prevent water leaking from the disconnected hose, cut the cable tie (2) holding the hose and remove it from the pump inlet fitting (4).
- 2. Using a 3/8-inch wrench disconnect the white Teflon™ pump tube (5) from the pump outlet fitting (6).
- 3. Using a 9/16-inch wrench on the pump inlet fitting (4) and a 9/16-inch wrench on the pump body inlet fitting (7), unscrew the inlet fitting.
- 4. There is a mesh filter (8) inside the pump inlet fitting and a rubber seal (9) on the outside threaded end of the fitting. From the threaded end of the fitting, insert a blunt instrument and gently push the filter out. If the filter does not come out, soak the fitting until any deposits have been dissolved.

- 5. Inspect the filter and remove any debris. Rinse the filter in clean water.
- 6. Insert the filter and reassemble the inlet fitting to the pump body inlet fitting. Ensure that the rubber seal is in place. Finger tighten the fitting and using a 9/16-inch wrench, tighten one half turn ONLY so as not to damage the internal thread of the plastic pump fitting.
- 7. Using a 7/16-inch wrench, disconnect the right angle fitting (6) and collar sub assembly while holding the pump body outlet fitting (11) with a 9/16-inch wrench and clean the threads of any debris. Note the orientation of the fitting relative to the steam generator (12).

Tip: Draw a black line through the three parts on the output side of the pump to ensure correct orientation when re-assembling.

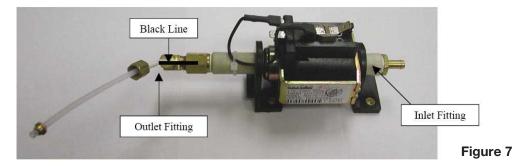

8. Clamp the right angled filter fitting using an adjustable wrench or a vice if available, and unscrew the collar with a 9/16-inch wrench.



Filter fitting and collar assembly (dis-assembly)

Figure 8

- 9. There is an insert (13) containing a fine mesh filter in the end of the right angle fitting. This filter cannot be removed from the fitting. If there are any deposits on the mesh, clean off with a stiff, non-metallic brush. If available, an air line can be used to purge the debris from the filter if appropriate. If there is any scaling on the fitting, it can be removed by soaking the fitting and filter in a mild solution suitable for removing scale or mineral deposits (i.e. vinegar). Rinse thoroughly with clean water.
- 10. If the outlet fitting has become damaged beyond use, then the in-line filter should be removed, and the outlet fitting replaced with replacement fitting 01-111115S, which has no filter. The unit should be upgraded with an in-line filter and reservoir filter, if appropriate.
- 11. If replacing the outlet fitting or re-fitting when the components are clean, apply a small amount of Teflon™ tape or pipe thread compound on to the external thread of the fitting and reassemble the right angle fitting to the collar ensuring the correct orientation of the two components.

- 12. Apply a small amount of Teflon<sup>™</sup> tape or pipe thread compound to the external thread of the collar and insert the fitting/collar sub assembly into the pump body outlet fitting. Finger tighten only at this stage, then using a 9/16-inch wrench, tighten one half turn ONLY. Ensure that the Teflon<sup>™</sup> pump tube will reach and the outlet assembly is orientated correctly.
- 13. Reconnect the pump tube to the pump outlet. **Do not kink the Teflon™ tube**.
- 14. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

#### Performing water reservoir filter maintenance

For maintenance of units where no internal filters are fitted and a water reservoir filter is installed, proceed as follows:



- 1. Remove the reservoir cap.
- 2. Grasp the post in the middle of the filter.
- 3. Remove the reservoir filter.
- 4. Rinse under running water.
- 5. Replace water reservoir filter and reservoir cap.

Figure 9


Where a reservoir filter is not installed, it can be retro-fitted, but note the filter will not work with old reservoir caps and therefore ordering a new reservoir cap is recommended to upgrade existing field units. All replacement caps will come with a new water reservoir filter. Spare Part Kit Numbers are: Reservoir Cap # 01-101783S, Water Reservoir Filter # 01-109300S.

#### Performing inline filter maintenance

For maintenance of units with in-line filters, proceed as follows (See Figure 10):

This filter should be replaced every 2000 cycles or every two years.

- 1. Drain the reservoir or crimp the inlet tube (1) from the reservoir to prevent water leaking from the disconnected tubes.
- 2. Cut and remove the cable ties (2) connecting the filter and inlet/outlet tubes.
- Remove the filter (3) from the tubes.
   Note: the filter is directional, so check the orientation of the filter before removing.

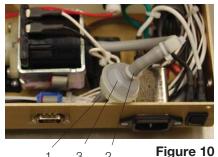

- i iguic iv
- 1. Inlet tube
- Cable tie
   In-line filter
- 4. Insert the new filter in to the inlet and outlet tubes. Check the flow orientation of the filter.
- 5. Secure the tubes to the filter with new cable ties.
- 6. Run a sterilization cycle and observe filter and tubes for leaks. Check LCD read-outs for messages indicating cycle status.
- 7. The in-line filter kit can be retro-fitted and is available as part # 01-106637S and includes the filter, all necessary tubing and installation instructions.


#### Performing Pump Tube Replacement

The pump tube located between the pump and the steam generator is designed to control the flow of water into the steam generator so that it receives the precise amount of water required at specific points in the sterilization cycle. The tube is directional (pump to steam generator) and contains a precision orifice tube inside the main tube that has a precise diameter that controls the volume passing through into the steam generator. The orifice tube is located at the steam generator end of the main pump tube. The diameter of the fitting is represented on the outside of each pump tube by a numeric value. The larger the value, the larger the diameter of the orifice.



Pump tube

Figure 11

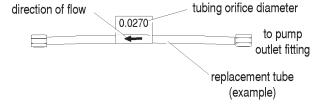

If the filters are clean and the pump flow parameters are not correct, then it may be possible to recover the correct water delivery by the replacement of the pump tube with one of larger or smaller diameter as appropriate.

Single pump tubes are available in a large number of diameters from 0.0210" diameter to 0.0290" diameter in 0.0005" increments.

**Note:** The number as it appears on the tube may only be represented by 3 digits, e.g. a 0.0245" tube will be represented by '245'.

As a guide to the tube size required if changing the tube is considered, the nominal tube sizes for STATIM 2000 models as fitted during production should be as follows:

Unit / Voltage Rating	Nominal Tubing		
*STATIM 100 V	0.0270		
*STAT <i>IM</i> 115 V	0.0255		
*STAT <i>IM</i> 230 V	0.0245		

Note: These are for aluminium boilers only. If a unit is fitted with a stainless steel steam generator, the value will likely be in the region of a difference of 0.0040" larger. For example, a 0.0245" tube fitted to an aluminium steam generator may be 0.0285" on a stainless steel steam generator.

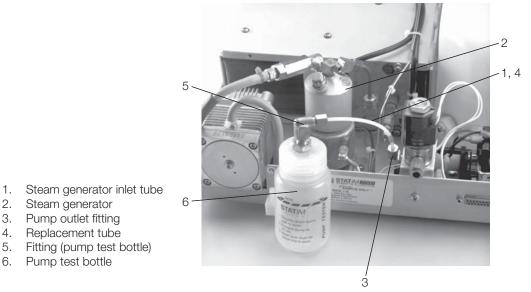


Figure 12


If, as a result of a failed pump flow check after a filter maintenance, the tube may be changed as follows (see Figure 12):

1. Turn the unit OFF.

2. Steam generator 3. Pump outlet fitting 4. Replacement tube 5. Fitting (pump test bottle) 6. Pump test bottle


- 2. Using a 3/8-inch wrench disconnect the pump tube (1) from the top of the steam generator (2) or pump test bottle (6) as appropriate, and disconnect the other end from the pump outlet fitting (3). Record the orifice diameter and discard the tube.
- 3. Select a smaller or larger size tube (depending on whether pump is over or under delivering). Note: There is no prescribed formula over which diameter to select. Selection will need to be based on perception of the under or over delivery amount.
- 4. Connect the pump end of the new tube to the pump outlet fitting. A label on the tube indicates the tube orifice diameter. An arrow on the label indicates the direction in which the tube must be installed. The arrow must point towards the steam generator when installed.
- 5. Thread the other end of the replacement tube (4) to the fitting (5) on the pump test bottle (6) until finger tight. Do not cross thread the fittings. Do not overtighten.

- 6. Using the button on the control box or appropriate keypad button (revision 7.x controllers), activate the pump for 2 seconds to purge any air that is trapped in the fittings. Empty any water that enters the pump tester by removing the bottle from the cap. Recap the pump test bottle. The pump tester MUST be empty before starting the pump test.
- 7. Be prepared to record the time it takes for the water level displayed on the pump tester to reach the top of the line marked MIN (see Figure 13). Activate the pump.
- 8. If the water level reaches the top of the line marked MIN within the time allotted in the Pump Flow Parameters chart, the pump is performing correctly. If the water level reaches the top of the line marked MIN in less than that time or if the water level does not reach the top of the line marked MIN within the allotted time, remove the tube and try a different size, repeating the test procedure until the correct size is established.
- 9. If flow requirements cannot be met using the largest or smallest diameter tubes, pump replacement will be necessary.
- 10. If a test proves successful, disconnect the test bottle. Using a 3/8-inch wrench connect the replacement tube to the top of the steam generator.
- 11. Make sure there is sufficient steam-process distilled water in the reservoir and activate the pump using the control box or appropriate keypad button (if using revision 7.x controllers) for 5 seconds to ensure that the pump is functioning and observe for leaks.
- 12. If the pump is replaced, a dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.
  NOTE: These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.
- 13. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

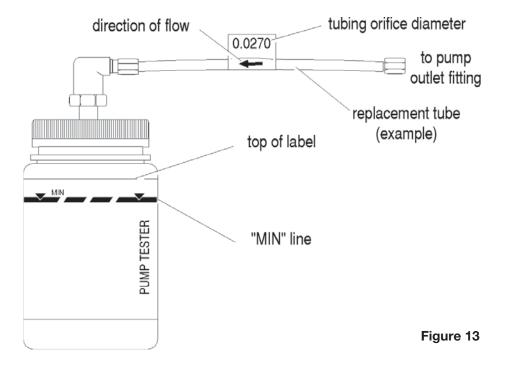

#### Removing the Fluid-O-Tech Pump Assembly

To remove the Fluid-O-Tech pump, proceed as follows (see Figure 14):

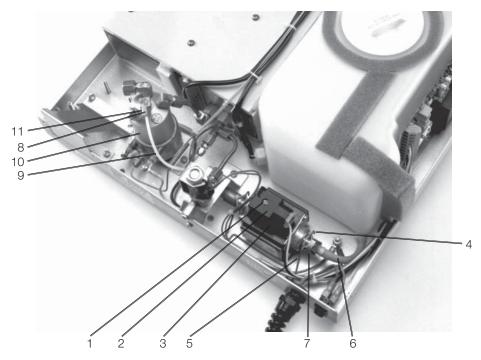


- 1. Diode plug screw
- 2. Diode plug
- 3. Pump
- 4. Cable ties
- 5. Mounting bracket
- 6. Pump inlet tube
- 7. Pump inlet tube fitting
- 8. Compression nut
- 9. Teflon™ tubing
- 10. Steam generator
- 11. Inlet fitting

Figure 14


- 1. Remove the diode plug screw (1) and diode plug (2) from the top surface of the pump (3). If the diode is functioning, retain for reassembly. If not, replace with SciCan Part # 01-104159S. See *Testing the Fluid-O-Tech Pump Diode* in the following section for diode orientation.
- 2. Cut and remove the three cable ties (4): one at each end of the pump mounting bracket (5) and one holding the inlet tube (6) to the pump inlet fitting (7). Do not nick the tubing while cutting the cable ties. Remove the rubber tube from the inlet fitting and clamp or stop the end of the tube securely.
- 3. Using a 3/8-inch wrench, remove the compression nut (8), attached to the Teflon™ tube (9), from the inlet fitting (11) on the top of the steam generator (10).
- 4. The pump assembly is secured to the base with the pump mounting bracket. Lift the pump assembly out of the mounting bracket and clear of the unit.

#### Testing the Fluid-O-Tech Pump Diode

To test the pump diode, proceed as follows (see Figure 15):

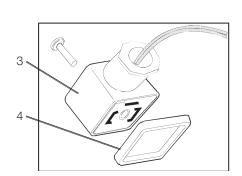

- 1. Center screw
- 2. Diode plug
- 3. Black cap
- 4. Rubber boot

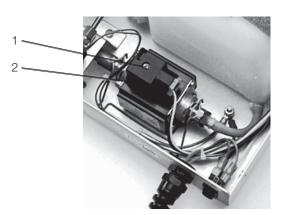

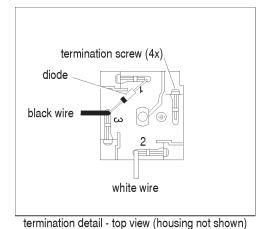


Figure 15

If the pump has not been removed:

- 1. Remove the center screw (1) from the diode plug (2), and retain for reassembly.
- 2. Detach the diode plug from the pump body.

If the pump has been removed, or after detaching the diode plug:

- 1. Remove the rubber boot (4) from the diode plug.
- 2. Using a small screwdriver or similar instrument, pry away the black cap (3) from the plug.
- 3. Test the diode using the diode test setting on a multimeter.


#### Reinstalling the Fluid-O-Tech Pump Assembly

To reinstall the pump assembly, proceed as follows (see Figure 14):

- 1. Slide the outlet side of the pump under the solenoid valve bracket, locate the pump assembly into the mounting bracket and snap it into place. Do not kink the Teflon™ tube (9).
- 2. Push the open end of the rubber inlet tube (6) extending from the water reservoir on to the pump inlet fitting (7) as far as it will go. Secure the tube to the fitting using a cable tie (4).
- 3. Attach the compression nut (8) on the Teflon™ tube to the inlet fitting (11) on top of the steam generator (10) finger tight. Tighten this nut using a 3/8-inch wrench. **Do not overtighten.**
- 4. Re-attach the diode plug (2) using the diode plug screw (1).
- A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.
   NOTE: These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.
- 6. Make sure there is sufficient steam-process distilled water in the reservoir and activate the pump using the control box for 5 seconds to ensure that the pump is functioning.
- 7. Fasten one cable tie to each end of the pump mounting bracket.
- 8. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

#### Removing the SciCan or Ulka Pump Assembly

To remove the pump assembly, proceed as follows (see Figure 16):

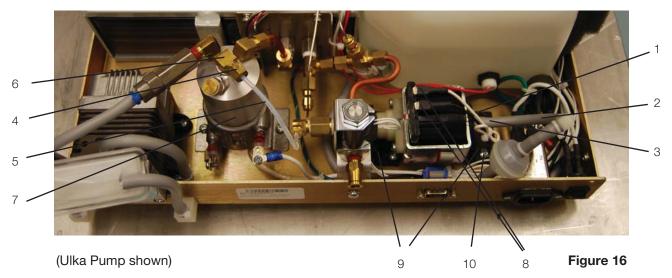

- 1. Cable tie
- 2. Inlet tube
- 3. Pump inlet fitting
- 4. Compression nut
- 5. Teflon tube
- Inlet fitting
- 7. Steam generator
- 8. Terminal connectors
- 9. Rubber brackets
- 10. Screws
- 1. Cut and remove the cable tie (1) holding the inlet tube (2) to the pump inlet fitting (3). Do not nick the tubing while cutting the cable tie. Remove the rubber tube from the inlet fitting and clamp or stop the end of the tube securely.
- 2. Using a 3/8-inch wrench, remove the compression nut (4) attached to the Teflon™ tube (5), from the inlet fitting (6) on the top of the steam generator (7).
- 3. Disconnect the appropriate terminal connectors (8) from the pump. Note the position of each terminal before removal.
- 4. The pump assembly is secured to the chassis by rubber brackets (9) and four shoulder screws (10). Remove the four shoulder screws and retain for re-assembly. Remove the pump assembly from the chassis.

#### Reinstalling the SciCan or Ulka Pump Assembly

To reinstall the pump assembly, proceed as follows (see Figure 5):

- Slide the outlet side of the pump under the solenoid valve bracket and place the pump onto the chassis. Position the rubber brackets (9) so that the mounting holes align with the threaded holes in the chassis. Secure the four shoulder screws (10) using Locktite® Threadlock Permalock compound LM113 or equivalent. Do not kink the Teflon™ tube.
- 2. Push the open end of the rubber pump inlet tube (2) extending from the water reservoir on to the pump inlet fitting (3) as far as it will go. Secure the tube to the fitting using a cable tie (1).
- 3. Thread the compression nut (4) on the Teflon™ tube (5) to the inlet fitting (6) on top of the steam generator (7) until finger tight then tighten this nut using a 3/8-inch wrench. **Do not overtighten.**
- 4. Connect the appropriate terminal connectors (8) to the pump.

- 5. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.
  - **NOTE:** These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.
- 6. Make sure there is sufficient steam-process distilled water in the reservoir and activate the pump using the control box for 5 seconds to ensure that the pump is functioning.
- 7. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

#### **Reservoir and Water Sensors**

STATIM 2000 / 2000S reservoirs have been installed in three different configurations with regard to water sensors. The reservoir moulding however is the same in all cases.

On early units up to 1995 with revision 3.x/4.x controller boards the reservoir contained a black cylindrical float sensor, which only detects the level of water in the reservoir.



(cut away)

On units from 1995 to 2007 with revision 2.x/5.x/6.x controller boards, the reservoir contained a water quality sensor with two exposed probes (2), which detect the water level AND quality in the reservoir by conductivity.



(cut away)

On units from 2007 onwards with revision 7.x controller boards, the reservoir contained a water quality sensor with two exposed probes that ONLY detects the water quality in the reservoir by conductivity, AND a float assembly that ONLY detects the water level (3).

Note: The three sensor types are not interchangeable between units with the different controller boards.



Water quality sensor

Water level sensor (3)

Figure 17


#### Removing the reservoir

To remove the reservoir (1), proceed as follows (see Figure 18):

The unit should be disconnected from the electricity supply before performing this procedure as the unit needs to be turned on its side to detach the reservoir securing screw.

- 1. Drain the reservoir.
- Cut the appropriate cable ties (4) to release the wiring so that it does not restrict the removal of the
  reservoir, and remove the appropriate sensor leads from the controller board terminal positions.
   NOTE: The terminal positions on the controller board so that the sensors may be
  re-connected correctly. This will vary depending on the controller type.
- 3. Cut the cable tie (4) securing the reservoir supply tube (5) to the pump fitting (6) and remove the tube from the fitting.
- 4. Tilt the unit on its side so that the reservoir and PCB are at the top of the unit. Locate and remove three nylon cap nuts with washers (7) on the bottom of the chassis securing the reservoir. Retain the cap nuts and washers.

5. Place the unit back on its feet and carefully remove the reservoir from the chassis. BE CAREFUL NOT TO DAMAGE THERMOCOUPLE LEADS.

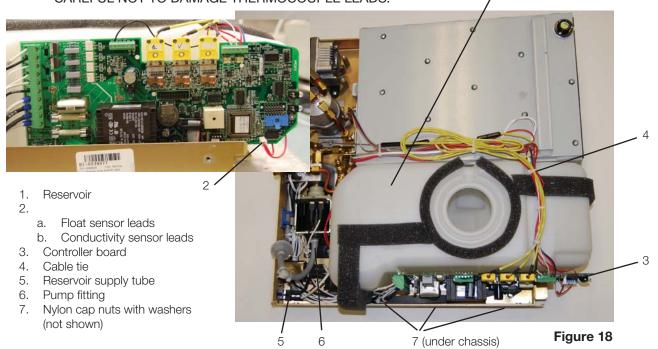

#### Reinstalling the reservoir

To reinstall the reservoir, proceed as follows (see Figure 18):

- 1. Carefully place the reservoir (1) in the chassis taking care not to damage the wiring.
- 2. Tilt the unit on its side so that the reservoir and PCB are at the top of the unit. Locate and install the three nylon cap nuts with washers (7) on the bottom of the chassis securing the reservoir.
- 3. Place the unit back on its feet.
- 4. Connect reservoir supply tube to pump fitting and secure with cable tie.
- 5. Reconnect the appropriate sensor leads from the controller board terminal positions.
- 6. Using cable ties, fasten the leads back in their original positions.
- 7. Fill the reservoir with steam-process distilled water.

**Note:** if filling the reservoir before re-connecting the sensor leads, touch the leads to an earthed point to discharge any static electric charge that may have built up on the reservoir. This reduces the chance of ESD damage to the controller board.

- 8. Reconnect the unit and power the unit ON.
- A dielectric strength test (Hi-Pot) AND a protective bonding impedance test (ground continuity)
  must be performed on the STATIM when the work is completed and after the cover has been
  returned to the unit.

# Troubleshooting reservoirs with water quality sensor only

#### Diagnosing conductivity sensor problems

Before starting diagnostic procedures check that:

- 1. The reservoir is free of any debris and the conductivity sensor is clean.
- 2. The reservoir is filled with good quality steam distilled process water so that the conductivity sensor is completely submerged. If in doubt about water quality, drain the reservoir and refill with steam-process distilled water containing less than 5 ppm total dissolved solids or having conductivity of less than 10 µS/cm.
- 3. The leads of the conductivity sensor are securely connected to controller board. Nothing is connected to controller board terminal positions labeled 'FLOAT'.
- 4. The controller board is clean and dry (both sides). Pay particular attention to the terminal block where the sensor is connected and any components on the upper right hand portion of the board.
- 5. Measure the negative voltage across test connector header positions P1-1 and P1-3. If the voltage reading is not within -8.4V to -9.7V, replace the controller board.

#### Checking the conductivity sensor

To check the conductivity sensor, proceed as follows:

- 1. Power the unit OFF and disconnect the sensor leads from controller board connector positions J4-3 and J4-4.
- 2. Power the unit ON and start an Unwrapped cycle.
- 3. If a "REFILL RESERVOIR, EMPTY WASTE BOTTLE" message appears on the LCD, proceed to step 4. If the cycle starts when the sensor leads are disconnected, the controller board is at fault and should be replaced.
- 4. Short circuit (bridge) controller board connector positions J4-3 and J4-4 together and start a cycle. If the unit displays a "WATER QUALITY NOT ACCEPTABLE" message, proceed to step 5. If no message is displayed then the controller board is at fault and should be replaced.
- 5. Remove the short from J4-3 and J4-4 and short controller board connector positions J4-5 and J4-6 together. If the unit displays a "SELECT A CYCLE" message, it is unlikely that the controller board is damaged.
- 6. If problems persist, leave the unit powered ON to allow internal components to warm up for a period of time (some failures are temperature dependent). Then start any sterilization cycle and repeat steps 1 through 6.
- 7. To test the conductivity sensor, short the sensor posts in the reservoir using a long handle screwdriver. The measured resistance across the unconnected sensor leads should be less than 1.0 ohm. If the reading is other than that the sensor is likely damaged. Replace the sensor.


# Troubleshooting "WATER QUALITY NOT ACCEPTABLE" message

Where water quality is known and measured in the reservoir to be below 5 ppm totally dissolved solids or conductivity of less than 10  $\mu$ S/cm, but the above message is displayed, proceed as follows.

In some STAT*IM* units fitted with water quality sensors, (typically with revision 2.x/5.x controller boards), there may be a case where electro magnetic interference (when the unit/pump gets older) affects the accuracy of the sensor and the "WATER QUALITY NOT ACCEPTABLE" message appears in the LCD display even though nothing is actually wrong with the water. This is rare, but can occur in units prior to those with revision 6.x controller boards onwards. After this, a grounding spring was fitted to the pump and shielding was fitted to the sensor cabling which alleviated the problem.

If this situation occurs, there are three things that can be done:

Option 1. Reverse the polarity of the pump by switching the black and white power leads on the pump side.

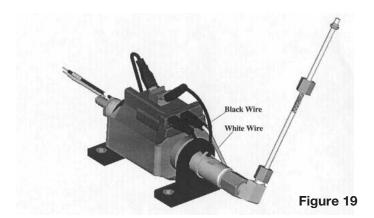

Option 2. Fit a grounding spring to the pump.

To fit a grounding spring, proceed as follows (see Figure 20):







Figure 20


Please Note:
The Spring Pump Ground Assembly is not required for older units that have the metal mounting bracket affixed below

- 1. Remove the pump as described in the section above "Removing the SciCan or Ulka Pump Assembly".
- 2. Insert the grounding spring on the inlet side of the pump as per picture one in figure 10 above.
- 3. Re-install the pump as described in the section above "Installing/reinstalling the SciCan or Ulka Pump Assembly". The spring positioning should appear as per picture two in figure 10 above.

#### Option 3. Upgrade the water quality sensor to the later type (01-103571S).

The original water quality sensor had two wires coming from it whereas the upgraded part has three wires that require attaching as per the configuration below.

- single white wire to J4 pin#3,
- second white wire (sheathed with the black wire) to J4 pin#4, and
- black wire to J3 pin#5.

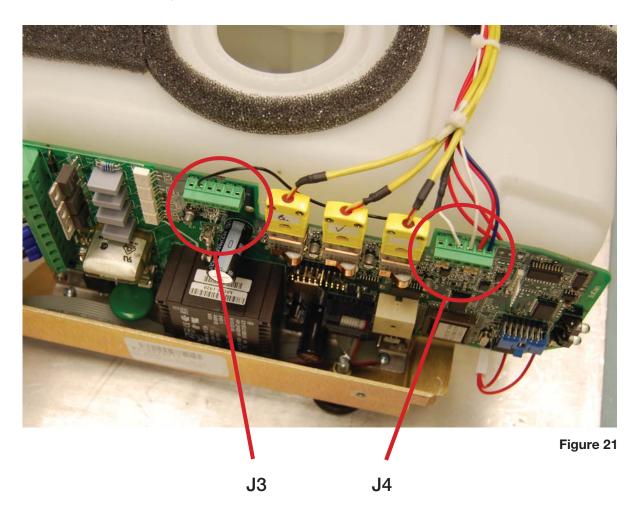


## Troubleshooting reservoirs with float and water quality sensors

In 2008, SciCan introduced a new type of water reservoir fitted with a water quality sensor AND a float switch. In the past, the low water level indication was performed by the water conductivity sensor. If the sensor measured the conductivity to be between 0.0 and 0.3  $\mu$ S/cm, (the minimum default setting in the software is 0.3  $\mu$ S/cm) the unit would notify the user to "REFILL RESERVOIR" even though the reservoir was in fact filled above the sensor.

As the water quality of many water suppliers has improved of late, some Rev. 7.x controller STATIM users may find this "REFILL RESERVOIR" message appear, and as the conductivity sensor in the unit cannot be set to  $0.0~\mu$ S/cm, as the unit would run even without water, the float was introduced to enable extremely pure water to be used without an error message.

To determine if your STAT*IM* unit has a water reservoir float switch, enter the Service Menu and proceed to the "CONDUCTIVITY SETUP" option and observe the L= value.

The display should be similar to the example below.

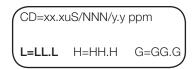

#### Screen Representation

 $xx.x = Water conductivity in \mu S.$ 

NNN = Conductivity measurement in ADC (Analog to Digital) counts (0...255)

y.y = water conductivity in ppm (parts per million).

LL.L = Numerical threshold in uS below which the unit thinks it is out of water OR float identifier (see below).

HH.H = High valve threshold (Bad water threshold - values larger than this trigger

"Water Quality Not Acceptable" error.

G.GG = Water conductivity circuit gain

If the L= value reads "FLOAT", the unit has a water float switch.

If the L= is a numerical value (e.g. 0.3 uS) then the unit has no water float switch.

The float sensor in this reservoir will ONLY work with STAT*IM* 2000/5000 units with revision 7 controller board or higher AND have the following software revisions or higher:

- S201R605, S202R605, S2E2R605, S2S2R605, S2S9R605
- S5S2R605, S501R605, S503R605, S502R605, S5S2R605, S5S9R605

If replacing a reservoir for a unit without a revision 7 controller board, the reservoir without float (with water conductivity sensor) can be ordered.


If replacing a reservoir on a unit with a revision 7 controller board and software lower than as listed above, the reservoir with float switch (and water quality sensor) can only be used if the software is upgraded to the latest revision.

If the unit is being upgraded with the new style reservoir and the choice has been made to upgrade the software to the latest version, the conductivity set-up in the service menu will need to be modified.

#### Modifying the conductivity set up in the service menu

To modify the conductivity set up, proceed as follows:

- 1. Enter Service Menu.
- 2. Select "Conductivity Setup". (Pressing the Unwrapped button will allow you to scroll through the menu).
- 3. When "Conductivity Setup" is on the top line of the LCD, press the Rubber or Plastic to enter the sub menu. The display should be similar to the example below.

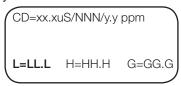

- 4. The L= value at this point SHOULD be a numerical value (e.g. **0.3 uS**). Press the Wrapped button to change the L= value to "FLOAT".
- 5. Press the Stop button. This will save and exit out of the Service Menu.

#### Removing and replacing the float sensor

To remove and replace the float sensor, proceed as follows (See Figure 22)



- 1. Cut the appropriate cable ties holding the float sensor wires in the wiring harness and disconnect the connector from the controller board (see below). Release the wires from the harness so that they are free to be removed from the unit with the sensor.
- 2. Remove nut from float sensor.

3. Position the new float sensor on a flat surface and mark sensor threads as shown in the picture below. The mark indicates the top of the sensor.



Figure 23

- 4. From inside the reservoir, pull the old sensor complete with its wires from the reservoir.
- 5. Insert new sensor (reverse removal instructions).
- 6. Ensuring the mark on the sensor thread is facing up, install the nut that was removed in step 1. Using a torque wrench tighten the nut to 7 in-lb. Inspect and confirm the mark on the threads is still facing up.
- 7. Re-route the float sensor wires in the appropriate position in the wiring harness and reconnect the float sensor. Ensure the wire is run behind the PCB and make the smallest loop possible without causing the wires to kink or bend sharply.
- 8. Cable tie the float sensor wires to the wiring harness.


## Compressor

Most STAT*IM* 2000 type units have been manufactured with a compressor that forces air through the steam generator and cassette to remove residual moisture and to dry the instruments.

Air is drawn into the compressor (where fitted) through a circular foam air filter located behind a cover plate on the rear of the compressor. This plate can be accessed from the rear of the unit with the cover still attached. The air is then directed through a bacteria-retentive air filter (where fitted) which is held in a bracket attached to the rear cover of the unit.

Details on changing these filters can be found in Chapter 2, *Tools, Maintenance Schedules, Procedures and Testing.* 

**Note:** To test the compressor the control box should be attached to the unit or 'compressor' selected in the device test sub menu of the revision 7.x controller board service menu (if cover is removed but still connected).

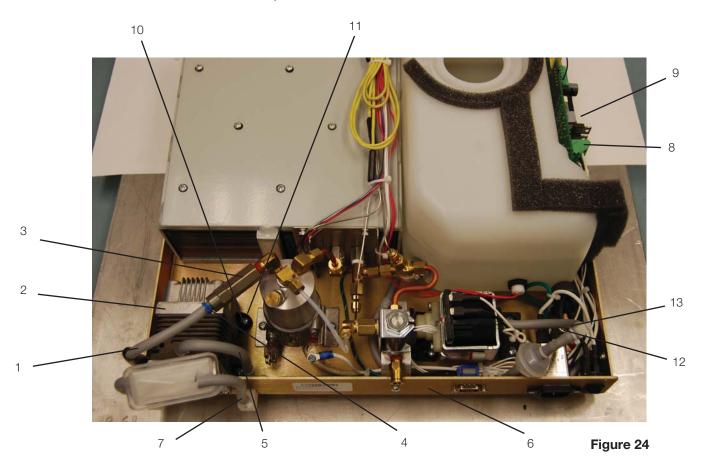

- 1. Silicone tube
- 2. Compressor
- 3. Setam generator check valve
- 4. Cable tie
- 5. Silicone tube
- 6. Chassis
- 7. Filter fitting

- 8. Compressor leads
- 9. Controller board
- 10. Screws (partially obscured)
- 11. Rubber mount
- 12. Ground wire (obscured)
- 13. Ground post (obscured)

#### Removing the Compressor

To remove the compressor, proceed as follows (see Figure 24):

- 1. Turn the power switch OFF, and unplug the unit. Remove the cover.
- 2. For units with the compressor air filter only:
  - a. Locate silicone rubber tube (1) which extends from the compressor (2) to the steam generator check valve (3). Carefully cut the cable tie (4) securing the tube to the compressor and remove the tube from the fitting.

For units with the additional bacteria-retentive filter:

- b. Locate silicone rubber tube (5) which extends from the fitting on top of the compressor through the chassis (6) to the bacteria-retentive input filter fitting (7). Carefully cut the cable tie (4) securing the tube to the compressor and remove the tube from the fitting.
- 3 .Cut the necessary cable ties to release the compressor cables from the wiring harness and disconnect the compressor leads (8) from controller board (9) terminal block positions J1-9 LINE and J1-10 NEUTRAL. Disconnect the ground wire (12) from the chassis ground post (3). Observe the routing of the affected wires for reassembly.
- 4. Remove the four screws (10) that attach the black rubber compressor mounts (11) to the chassis. Retain the screws.
- 5. Remove the compressor from the unit. If removing the compressor from a unit with a bacteria-retentive filter proves difficult, cut the high temperature rated cable tie securing the silicone rubber tube to the check valve and remove the tube from the fitting or attempt to slide the compressor under the tubing.

**Note:** Other than filter changes, the compressor is NOT a serviceable item. If it fails to function it must be replaced.

#### Reinstalling the Compressor

To reinstall the compressor, proceed as follows (see Figure 24):

- 1. Route the ground wire (12) with the wiring harness leading to the controller board and connect it to the ground post (13).
- 2. Reinstall the rubber compressor mounts (11) in the chassis using the four screws (10) retained during disassembly. Apply a small amount of Locktite<sup>™</sup> to the threaded holes in the chassis.
- 3. Route the compressor leads (8) with the wiring harness leading to the controller board (9). Secure the leads and ground wire using nylon cable ties every 2-3 inches. Connect the black wire to controller board terminal block position J1-9 LINE and the white wire to terminal block position J1-10 NEUTRAL.
- 4. For units with the compressor air filter only:
  - a. Connect the silicone rubber tube (1) which extends from the steam generator check valve (3) to the top of the compressor (2). Secure the tube using a cable tie (4).

For units with the additional bacteria-retentive filter (see Figure 24):

- b. Connect the silicone rubber tube (5) which extends from the input side of the bacteriaretentive filter (7) to the top of the compressor (2) and secure the tube using a cable tie (4). If the tube from the check valve fitting was disconnected, re-attach and secure the tube using a high temperature rated cable tie)
- 5. Test the compressor by activating the compressor with the control box or keypad as appropriate.
- 6. A dielectric strength test (Hi-Pot) AND a protective bonding impedance test (ground continuity) must be performed on the STAT/M when the work is completed and after the cover has been returned to the unit.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



#### **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

## PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



#### PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.


## Solenoid Valve and Pressure Transducer


#### Solenoid

STATIM units have been manufactured using one of two types of solenoid valve:

- Type A is a two part square bodied unit (fitted to very early units).
- Type B is a two part cylindrical bodied unit (see Figure 1).

If a Type A unit requires servicing or has failed, replace it with a Type B solenoid valve ONLY.

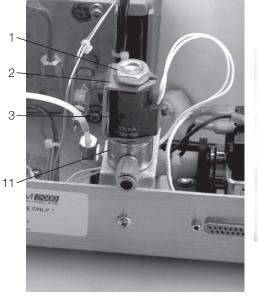

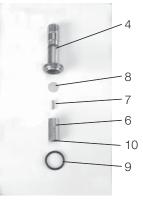

- 1. Nut
- 2. Yoke
- 3. Coil
- 4. Plunger tube
- 5. Plunger tube wrench (not shown)
- 6. Plunger
- 7. Plunger spring
- 8. Bronze disc
- 9. O-ring
- 10. Plunger seal
- 11. Valve body

Type B Solenoid Valve

Figure 1


#### Inspecting and Repairing the Type B Solenoid Valve

The solenoid valve assembly may be inspected and/or repaired while in the unit. To inspect and/or repair the assembly, proceed as follows (see Figures 1 and 2):

 To determine that the solenoid is working, test the solenoid valve by using the solenoid valve switch/button on the Control box. When the switch is activated, the sound of the coil energizing can be heard. Note: the cassette must be removed or unseated so that the microswitch is not active, as the solenoid will be permanently energised and will not be affected by the control box switch/button.

TIP: The solenoid may be also be activated by shorting P1-18 to ground on the PCB.

2. Turn the power switch OFF, and unplug the unit. If the solenoid is activated during disassembly, the coil will be damaged.

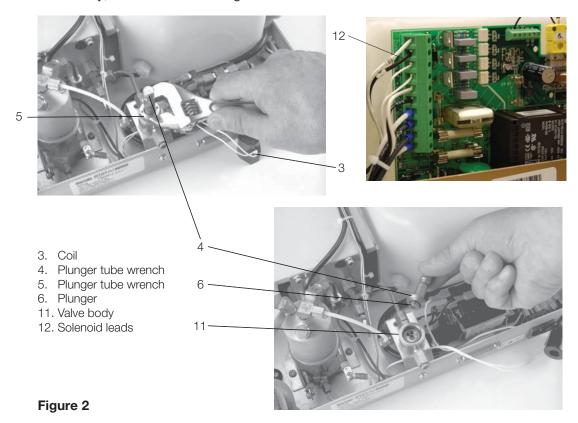

- 3. Disconnect the solenoid leads from the Controller Board terminal block positions labelled J1-7 LINE and J1-8 NEUTRAL. Use a multimeter in diode setting to determine whether or not the coil has failed. If the measured values are approximately equal in each direction, the coil has not failed. If there is no measured value or the values are other than approximately equal, replace the coil assembly.
- 4. Remove the nut (1) from the top of the solenoid valve assembly using an 11/16 inch wrench.
- 5. Pull the yoke (2) and coil (3) off the plunger tube (4) in an upward direction. Note the orientation of the yoke and coil for reassembly. To allow the coil to rest outside the unit, cut any cable ties required.
- 6. Unscrew the plunger tube from the valve body using the solenoid plunger tube wrench (5) and a ¾-inch wrench or adjustable wrench.

TIP: When lifting the plunger tube clear of the valve body, place your finger under the plunger before it drops from the tube and lift the complete sub-assembly off the body to retain the plunger in the tube for checking free movement. Check that the plunger is moving freely in the valve body.

- 7. Remove the plunger (6), plunger spring (7) and bronze disc (8) from the valve plunger tube body.
- 8. Inspect the o-ring (9) and the plunger seal (10) for nicks, compression set or swelling (the seal should not extend more than 0.010 inches from the plunger body).
- 9. Inspect the valve body (11) for debris. If the debris appears to be fragments of rubber, inspect the cassette seal for damage.
- 10. Clean the valve plunger, valve body, and fittings using oil free compressed air to blow debris from the surfaces. If there is residue (e.g. congealed handpiece oil), on the surfaces of the inside of the tube and outside of the plunger they may be cleaned with a moistened cleaning pad designed for use with Teflon™ coated surfaces to remove it. DO NOT USE METALLIC OBJECTS TO REMOVE THE RESIDUE. After cleaning, rinse thoroughly with water to remove all traces of the detergent.
- 11. Once the valve components have been inspected, determine if the assembly needs cleaning only, a new coil, the repair kit for the solenoid valve or an entire valve assembly replacement.
- 12. Reinstall or replace; the bronze disc (8), plunger spring (7), plunger (6) and o-ring (9) to the plunger tube (4). (Order SciCan Kit Part # 01-100998S.)
- 13. Screw the plunger tube into the valve body (11) and tighten using the solenoid tube wrench.
- 14. Place the coil (3) inside the yoke (2) and place onto the plunger tube.
- 15. Reinstall the coil retaining nut (1).
- 16. Connect the leads to controller board terminal block positions labeled J1-7 LINE and J1-8 NEUTRAL. These two wires are interchangeable.
- 17. Test the solenoid by activating the solenoid switch on the control box Note: the cassette must be remove or unseated so that the microswitch is not active as the solenoid will be permanently energized and will not be affected by the control box switch/button.
- 18. If the solenoid valve assembly or the solenoid coil is replaced, a dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) MUST be performed on the unit.
  NOTE: These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.
- 19. Replace any cable ties cut during the procedure.

Figure 3

#### Removing the Solenoid Valve

To remove the solenoid valve proceed as follows (See Figure 3):

- 1. Turn the power switch OFF, and unplug the unit.
- 2. Carefully cut the cable ties holding the solenoid valve leads and disconnect the leads (1) from the Controller Board terminal block positions labelled J1-7 LINE and J1-8 NEUTRAL.
- 3. Disconnect the compression nut (2) holding the outlet probe tube (3) to the outlet fitting (4) on the probe bracket assembly.
- 4. Remove the two screws (5) that attach the solenoid valve bracket (6) to the chassis (7): one on the inside of the chassis, one on the outer rear surface of the chassis.
- 5. Inspect the solenoid valve to determine whether servicing or replacement warranted.



5

- 1. Solenoid leads
- 2. Compression nut
- 3. Outlet probe tube4. Outlet fitting

- 5. Screw (behind)6. Solenoid valve bracket
- 7. Chassis
- 8. Solenoid valve inlet fitting

#### Refitting the Solenoid Valve

To refit the solenoid valve assembly, proceed as follows (see Figure 3):

- 1. Make sure the power switch is OFF, and the unit is unplugged.
- 2. Connect a new outlet probe tube (3) to the solenoid valve inlet fitting (8). Thread the compression nut until finger tight, on the inlet fitting.
- 3. Place the solenoid valve assembly in the unit.
- 4. Connect the other end of the outlet probe tube (3) to the outlet probe fitting on the probe bracket outlet fitting (4). Thread the compression nut until finger tight, on the inlet fitting. Tighten the nuts on both ends of the outlet probe tube using a 7/16 inch wrench. DO NOT OVERTIGHTEN.
- 5. Attach the solenoid valve bracket (6) to the chassis (7) using the two screws (5) retained from disassembly: one on the inside of the chassis, one on the outer rear surface of the chassis.

- 6. Route the solenoid valve leads (1) together with the bundle of wires extending to the controller board and connect the leads to controller board terminal positions labelled J1-7 LINE and J1-8 NEUTRAL. These two wires are interchangeable.
- 7. Bundle the wires together and secure them using cable ties every 2 to 3 inches.
- 8. Test the solenoid by activating the solenoid switch on the Control box. Note: the cassette must be removed or unseated so that the microswitch is not active as the solenoid will be permanently energized and will not be affected by the control box switch/button.
- 9. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) MUST be performed on the unit.
  - **NOTE:** These tests must be performed on the STAT*IM* again once the work is completed and the cover has been returned to the unit.
- Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.


#### **Pressure Transducer**

Note: Only 'S' type (EU) STATIM models have pressure transducers.

STATIM 2000S units use a pressure transducer to measure the steam pressure in the unit. The transducer is attached to a bracket located in the right rear corner of the unit. Transducer wires are routed along the top of the armature to a Pressure Interface Board attached to the P2 printer connector on the controller board, (Revision 2.x, 5.x and 6.x controller boards), or direct to the main PCB (Revision 7.x controller boards).

#### Identifying and checking the Pressure Transducer

Note that there are two types of pressure transducer used on STATIM S models, a 0 - 60 p.s.i. transducer on models up to approximately 2004 and a 0 - 68 p.s.i. transducer on models from 2004 on. The input voltage is the same for both types, but the output voltage is different so the two types are not interchangeable.

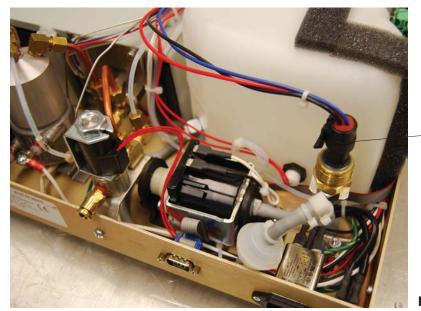

Pressure transducer

Figure 4

#### Identifying the transducer type

The pressure transducer may be identified in two ways:

- by the identification code etched on the side of the upper casing of the transducer
- by the voltage output signal

To identify the pressure transducer by identification code, proceed as follows:

Locate the identification code on the side of the transducer casing. Note that there are several lines to the code. The relevant line of code should contain 9 digits, similar to the following example: 0060AEFY1. The second and third digits in the code identify the pressure range of the transducer and will be 60 or 68.

Tip: the zero in the identifier has a diagonal line through it. 60 may easily be mistaken for 68, so take care when checking the identifier. A magnification device may help.

To identify the pressure transducer by output voltage, proceed as follows:

To measure the voltage output, power up the machine with out a cassette out (so that the machine is measuring atmospheric pressure). Set a multimeter to 'DC volts' and check the reading across the terminals of the BLUE P1-1 'SIG' (signal) and BLACK P1-3 'GND (ground) wires leading from the pressure transducer, at their terminations on the main PCB (Revision 7) or the interface board (all other PCB revisions). Check the atmospheric pressure reading on an independent reference pressure meter. The voltage reading for the appropriate pressure reading should be approximately the value as shown in the table (Figure 5) below for the relevant pressure transducer.

**IMPORTANT:** This identification method is only a guide. Performing a voltage comparison may not be accurate if the pressure transducer is damaged, a high offset value has been applied to it as a result of calibration, the reference voltage is incorrect (n/a revision 7 PCB) or if the unit is installed at high altitude. A more accurate reading can be achieved by finetuning the unit as per the instructions below ( See Checking the transducer).

#### Checking output voltage on pressure transducer

Pressure [kPa]	68PSI sensor reading [V]	60PSI sensor reading [V]
93	1.2934	1.3992
94	1.3020	1.4089
95	1.3105	1.4186
96	1.3190	1.4282
97	1.3276	1.4379
98	1.3361	1.4476
99	1.3446	1.4573
100	1.3532	1.4669
101	1.3617	1.4766
102	1.3702	1.4863
103	1.3788	1.4959
104	1.3873	1.5056
105	1.3958	1.5153

#### **Checking the Transducer**

Figure 5

The pressure transducers used on STATIM units are electro-mechanical in operation. When pressure is applied to the mechanical components located in the base of the transducer, the movements of the components are translated into a voltage signal in the electronics by modifying a fixed input voltage signal to give a value corresponding to the pressure applied. If you suspect that the transducer is faulty then the function of the transducer can be checked. The unit will need to be set up correctly before checking, and the following must be undertaken prior to diagnosis.

- 1. Remove the cover of the unit.
- 2. Attach the control box to the PCB (all models).
- 3. Power up the unit in calibration mode as appropriate to the software revision (see calibration section).
- 4. Connect a multimeter to the control box Vref output and check/adjust the reference voltage as appropriate (n/a on revision 7 units).
- 5. The atmospheric pressure transducer offset value needs to be set at '00' to correct the display pressure reading to perform the following checks. If it is not, then proceed as follows:
  - 5.1. Write down the transducer atmospheric offset reading (e.g. FE).
  - 5.2. Insert a calibration cassette into the unit with the calibrated reference meters (temperature and pressure) attached.
  - 5.3. Start a calibration cycle appropriate to the unit and adjust the pressure transducer offset to '00'
  - 5.4. Stop the calibration and release the cassette from the probe bracket to enable the unit and pressure reference meter to read the atmospheric pressure.
  - 5.5. Using a 3/8-inch wrench, disconnect the compression nut (8) holding the coiled transducer tube (6), from the bottom of the transducer (see Figure 4).

TIP: It is always wise to check the pressure transducer tube for blockages at this time as this may be the source of an incorrect pressure reading or slow transducer response. If the tubing is discoloured, kinked, or of there is debris visible in the tubing, remove the tubing completely for cleaning or fit a new tube.

Now check and record the following values:

Set a multimeter to 'DC volts' and check the reading across the terminals of the BLUE P1-1 'SIG' (signal) and BLACK P1-3 'GND (ground) wires leading from the pressure transducer at their terminations on the main PCB (Revision 7) or the interface board (all other board revisions).

- Compare the reference pressure and the display pressure.
- Compare the voltage reading from the multimeter against the voltage value on the chart above for the relevant pressure reading of the display.

If the transducer is working to manufacturing specifications, the display pressure and reference pressure should be within +/- 4 kPa of each other and the multimeter voltage reading should be close to the value on the chart (Figure 5) taking into account the accuracy/calibration status of the multimeter.

If the voltage readings vary significantly then the electronic components of the transducer may be failing or have failed.

If the pressure readings vary significantly then the mechanical components of the transducer may be failing or have failed.

The pressure transducer is not a serviceable component. It is critical to ensuring the correct conditions occur in the unit during sterilization. If failure has occurred or potential failure is suspected, replacement is recommended.

**NOTE:** Do not forget to reset the pressure transducer offset back to the original setting or, re-calibrate the unit fully if in doubt.


#### Removing the Pressure Transducer


CAUTION: the pressure transducer, steam generator and other components may be hot even if the unit has not been operating. Guard against burns.

Check the atmospheric pressure reading on the independent reference pressure meter AND on the unit display.

- 1. Pressure transducer
- 2. Power switch
- 3. Line cord (not shown)
- 4. Transducer connector
- 5. Compression nut
- 6. Coiled transducer tube
- Solenoid valve inlet "T" fitting
   Compression nut
- 9. Nut and washer
- 10. Transducer bracket
- 11. Screw and washer
- 12. Chassis

Figure 6

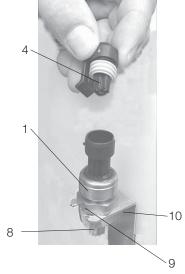

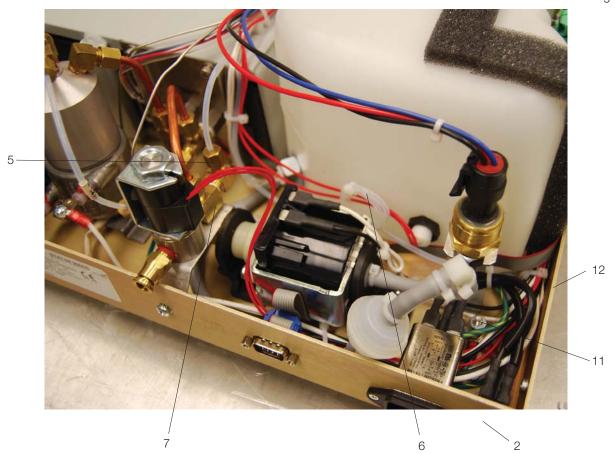

To remove the pressure transducer (1), proceed as follows (see Figure 6):

- 1. Turn the power switch (2) **OFF**, and unplug the mains power cable (3) (not shown). Allow time for the unit to cool.
- 2. Disconnect the transducer connector (4) from the top of the transducer.
- 3. Using a 3/8-inch wrench, disconnect the compression nut (8) holding the coiled transducer tube (6), from the bottom of the transducer.
- 4. Using a 7/8-inch wrench or an adjustable wrench, remove the transducer nut and washer (9) securing the transducer to the bracket (10).
- 5. Should removal of the transducer bracket be necessary, remove the screw and washer (11) holding the bracket to the chassis (12).
- 6. If the transducer is damaged, discard the transducer.

#### Replacing the Pressure Transducer

To replace the pressure transducer (1), proceed as follows (see Figure 6):

- 1. Make sure the power switch (2) is OFF, and the unit is unplugged.
- 2. If the transducer bracket (10) was removed, reinstall the bracket on the chassis (12) using the screw and washer (11) retained from disassembly.
- 3. Install the transducer in the transducer bracket using the nut and washer (9) retained from disassembly. Tighten with a 7/8-inch wrench or an adjustable wrench.
- 4. Reconnect the transducer tube (6) to the bottom of the transducer. Thread the compression nut (8) finger tight, and then tighten the compression nut using a 3/8-inch wrench. **Do not over tighten**.
- 5. If the transducer tube has been removed completely or is being replaced, reconnect the other end of the coiled transducer tube to the solenoid valve inlet "T" fitting (7). Thread the compression nut (5) finger tight, and then tighten the compression nut using a 3/8-inch wrench. **Do not over tighten**.
- 6. Hold the transducer / bracket assembly with one hand and plug the transducer connector (4) onto the transducer. Make sure the tab on the connector snaps in place.
- 7. Check that all fittings have been tightened and that all brackets are secured to the chassis. Reconnect the mains electrical cable (3) (not shown) and turn the power switch (2) **ON**.
- 8. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-out for messages indicating cycle status.

**Important note:** the pressure transducer performs a critical process monitoring function during the sterilization cycle and must be recalibrated prior to the unit being used to process instruments. Refer to the relevent calibration procedure and re-calibrate the unit upon completion of the repair.

#### Removing the Pressure Interface Board

**NOTE:** This is only applicable to Revision 2.x, 5.x and 6.x controller boards. All controller boards from Revision 7 have the pressure circuitry integrated into the main board.

CAUTION: Hazardous voltages are accessible on the controller board when the power is on.

The pressure interface board is a piggyback board assembly, which is attached to the P2 printer connector of the controller board (see note above) using the latch mechanism of the P2 header.

To remove the Pressure Interface Board (1), proceed as follows (see Figure 7):


- 1. Turn the power switch (2) **OFF**, and unplug the power cord (3) (not shown).
- 2. Disconnect the printer cable (4) (if present) from pressure interface board connector P2.
- 3. Disconnect the transducer wires (9) from pressure interface board terminal positions P1-1 'SIG' (blue wire, signal), P1-2 'V+' (red wire V+) and P1-3 'GND' (black wire, ground).
- 4. Press the ejector latches (10) of controller board (11) connector P2 to unseat the pressure interface board and remove the board.

#### Replacing the Pressure Interface Board

To replace the pressure interface board (1), proceed as follows (see Figure 7):

- Orientate the pressure interface board with the component side of the board facing away from
  the controller board (11) and the P2 connector on the left-hand side. Gently seat the pressure
  interface board onto the controller board P2 connector until the ejector latches are in an
  upright position. Do not crush or pinch the thermocouple leads. Support the controller board
  to avoid excessive deflection when pressing the interface board onto the connector.
- 2. Connect the transducer wires (9) to Pressure Interface Board terminal positions P1-1 'SIG' (blue wire, signal), P1-2 'V+' (red wire V+) and P1-3 'GND' (black wire, ground).
- 3. Connect the printer cable (4) (if present) to the pressure interface board connector P2. Ensure that Pin 1 of the keyboard cable aligns with Pin 1 of pressure interface board P2.
- 4. Plug in the power cord (3) (not shown) and turn the power switch (2) ON.
- 5. Run a sterilization cycle and observe LCD read-out for messages indicating cycle status.

**Important:** the pressure transducer performs a critical process monitoring function during the sterilization cycle, and on later units the calibration data is stored on the interface board. The unit <u>must</u> be re-calibrated prior to the being used to process instruments. Refer to the relevant calibration procedure and re-calibrate the unit upon completion of the repair.

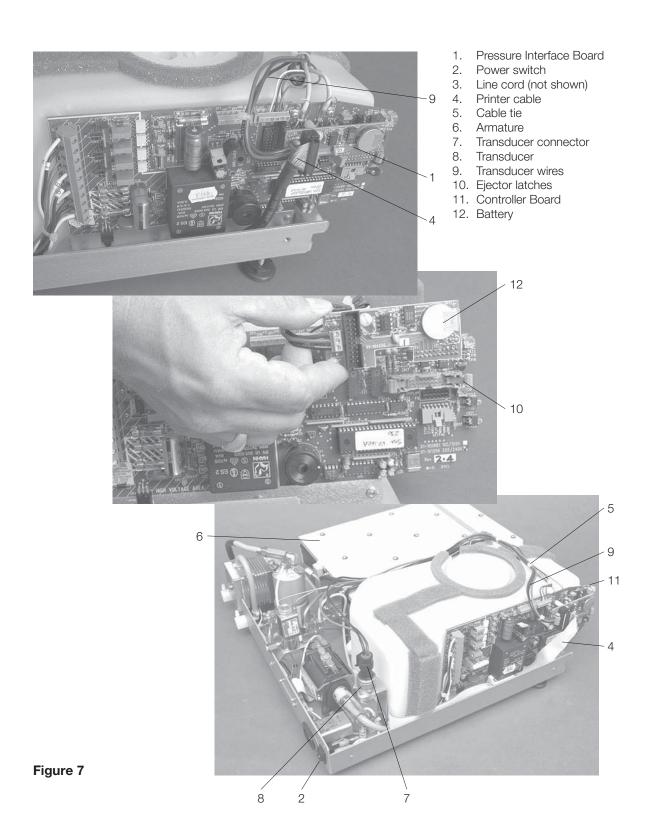

#### Removing and Replacing the Battery

Observe appropriate electrostatic discharge safeguards performing this operation.

When installed and a printer attached, the real time clock and battery on the pressure interface board (or main PCB for revision 7 units) replace the clock and battery found in the printer module.

The battery requires replacement if the date / time are not saved when the power is off.

To remove the battery (12), proceed as follows (See Figure 7 and 8):

- 1. Remove the pressure interface board or main PCB as appropriate to the PCB revision.
- 2. Carefully desolder BAT1 from the solder side of the board. Note the orientation of the anode and cathode. Always replace the battery with a SciCan battery of equal rating and size. Solder the replacement battery in position BAT1.
- 3. Discard the old battery. Observe all applicable environmental laws.
- 4. Replace the Pressure Interface Board or main PCB as appropriate to the PCB revision.

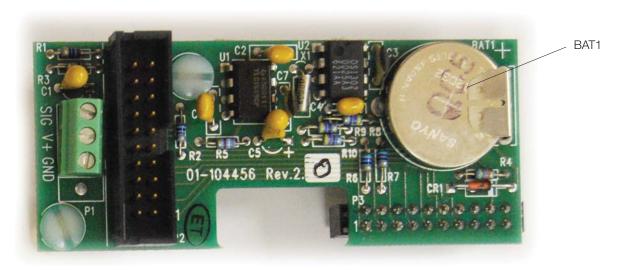


Figure 8

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



#### **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

## PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



#### PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# Armature, Isoplate, and Probe Bracket and Chamber Thermocouple Assembly

### Armature and Isoplate Subassembly


**NOTE:** Replacement armatures are normally supplied complete with isoplates from the factory. The isoplate is available separately should this be the only component requiring replacement in this subassembly.

#### Removing the Armature and Isoplate Subassembly

To remove the armature (1), proceed as follows (see Figure 1):

The unit should be disconnected from the power supply and the cover removed.

- Drain the reservoir.
- 2. Carefully cut the cable ties (3) securing the thermocouple leads (4), the microswitch leads (5) and the reservoir sensor leads (6).
- 3. Remove the microswitch leads from controller board (7) terminal positions labelled CASSIN J4-1 and J4-2.
- 4. Disconnect the compression nut (8) holding the exhaust tube (9) to the inlet probe fitting (10).
- 5. Loosen the two screws (11) holding the steam generator bracket (12) to the chassis (13).
- 6. Disconnect the compression nut (14) holding the exhaust tube (15) to the outlet probe fitting (16).
- 7. Loosen the two screws (17) that attach the solenoid valve bracket assembly (18) to the chassis.
- 8. The probe bracket assembly (19) is held by four cap screws with washers. Remove all screws (20) with washers using a 9/64-inch ball end Allen key. Note the ground strap (21) attached to the bottom right of the probe bracket by the cap screw.
- 9. CAREFULLY bend the thermocouple leads away from the armature.
- 10. Tilt the unit on its side so that the reservoir and PCB are facing upwards. Remove the four cap screws with spring lockwashers (22) securing the armature to the chassis using an Allen key: two at the front on the bottom of the chassis, and two at the rear on the bottom of the chassis.
- 11. Lower the unit to the flat work surface and lift the armature from the chassis. Be careful not to damage the thermocouple in the probe bracket assembly.

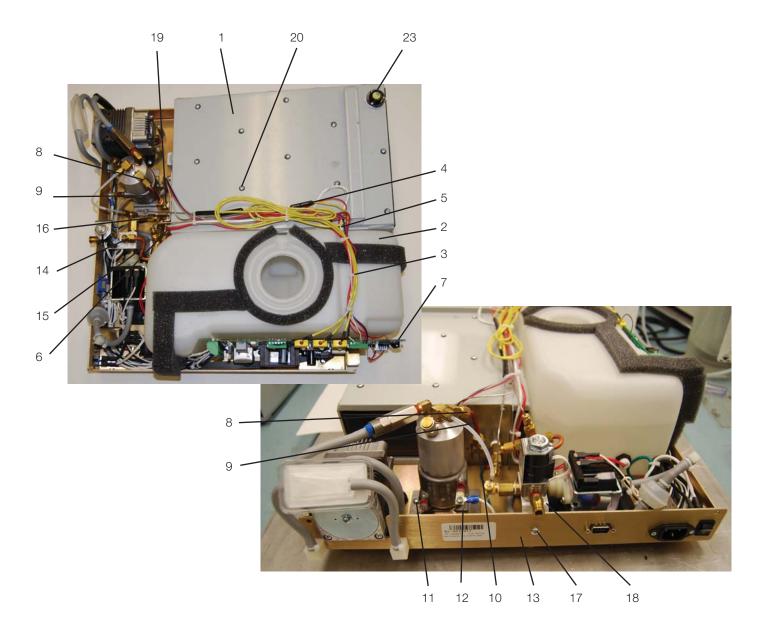

- 1. Armature
- 2. Reservoir
- 3. Cable ties
- 4. Thermocouple leads
- 5. Microswitch leads6. Reservoir sensor leads
- 7. Controller board
- 8. Compression nut for steam generator outlet tube
- 9. Steam generator outlet tube
- 10. Inlet probe fitting
- 11. Screws
- 12. Steam generator bracket assembly
- 13. Chassis
- 14. Compression nut for exhaust tube
- 15. Exhaust tube
- 16. Outlet probe fitting

- 17. Screws
- 18. Solenoid valve bracket assembly
- 19. Probe bracket assembly
- 20. Isoplate screws
- 21. Ground strap (not shown)
- 22. Cap screws with spring lockwashers (located on bottom of chassis)
- 23. Bubble level

Figure 1

#### Reinstalling the armature and isoplate subassembly

To reinstall the armature (1), proceed as follows (see Figure 1):

Align the armature (1) with the chassis (13). Tilt the unit sideways with the reservoir and PCB facing upwards, insert and tighten the four cap screws with spring lockwashers (22) retained from the removal procedure. Ensure that the armature drain tube exits through the hole in the chassis.


IMPORTANT: Be careful not to cross-thread the cap screws in the tapped holes of the armature.

- Be careful not to damage the armature drain tube, the thermocouple in the probe bracket assembly, the ground strap or the bubble level (23).
- 2. Install the probe bracket assembly (19). Partially insert the two bottom cap screws and attach the ground strap (21) to the bottom right hand position of the probe bracket assembly before proceeding.
- 3. CAREFULLY bend the thermocouple leads back into position.
- 4. Complete installation of the probe bracket assembly (19). (See, Probe Bracket and Chamber Thermocouple Subassembly below).
- 5. Secure the solenoid valve assembly (18).
- 6. Tighten the two screws (11) holding the steam generator bracket (12) to the chassis (13).
- 7. Bundle the microswitch leads (5), thermocouple leads (4) and the reservoir sensor leads (6) together using cable ties every 2-3 inches. Secure the wires and thermocouples to the top of the armature using the cable tie anchors provided on top of the armature.
- 8. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.
  - **NOTE:** These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.
- 9. Check that the bubble level (23) on the front right top corner of the armature is not damaged and that the bubble moves freely.
- 10. Refill the reservoir using steam process distilled water.
- 11. Calibrate the thermocouples appropriate to the unit type.
- 12. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

#### Removing the isoplate from the armature-isoplate subassembly

The Isoplate is a plastic box secured inside the armature. It is designed to hold and insulate the cassette in the unit during operation.

To remove the isoplate assembly (1) from inside the armature (2) proceed as follows (see Figure 2):

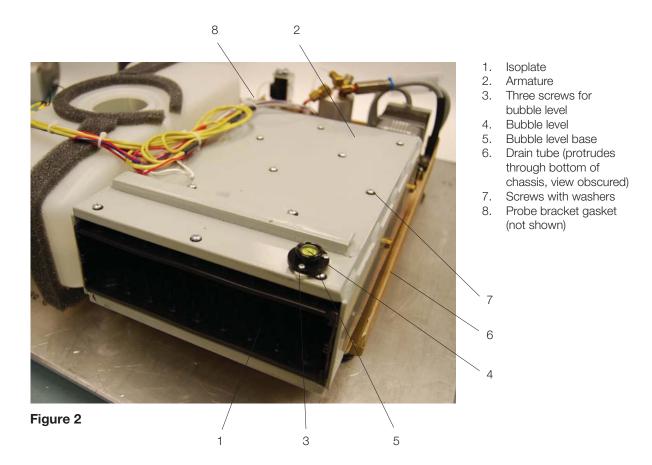

- 1. If not already undertaken, remove the armature from the chassis and place it on a secure work surface (see above).
- 2. Remove the three screws (3) securing the bubble level (4) to the bubble level base (5) on the front right corner of the armature. (The bubble level obscures a fastener that must be removed.)
- 3. The drain tube (6) passes through the bottom of the armature near the probe bracket opening and is attached to the isoplate. Using a pair of pliers, pull the drain tube off and retain.
- 4. The isoplate is secured in the armature using screws with washers (7) on the top and bottom of the armature. Remove and retain all the screws with washers.
- 5. When all fasteners are removed, slide the isoplate out of the armature.

#### Replacing and reassembling the isoplate into the armature

To install the replacement isoplate, proceed as follows (see Figure 2):

- 1. Slide the isoplate (1) into the armature. This may be easier if the cassette is first inserted into the isoplate. Note the orientation of the probe bracket opening in the back wall of the isoplate. Reinstall the screws with washers (7) but do not tighten the screws. When all the screws are reinstalled, slide the cassette into the isoplate if it is not already installed. Now tighten the screws in an "X" pattern, alternating across the isoplate, from the outside to the center.
- 2. Install the drain tube (6) and secure using Locktite™ LM113.
- 3. Reinstall the armature into the chassis (see above).
- 4. Reinstall the bubble level (4) to the bubble level (5) base using the three screws (3) retained from disassembly.
- 5. Before reinstalling the probe bracket assembly to the back of the armature, ensure that all traces of the old probe bracket gasket (8) are removed from the bracket. Fit a new gasket as appropriate. Note that a new gasket comes attached to the replacement isoplate and with a replacement probe bracket. It is also available separately.
- 6. Refit and re-align the probe bracket and chamber thermocouple assembly (see *Aligning the probe bracket and chamber thermocouple assembly below*).
- A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STAT/M unit at this stage.
   NOTE: These tests must be performed on the STAT/M again once the work is completed and the cover has been returned to the unit.
- 8. Refill the reservoir using steam process distilled water.
- 9. Recalibrate the thermocouples.
- 10. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

# Probe Bracket and Chamber Thermocouple Subassembly

#### **IMPORTANT INFORMATION:**

The probe bracket and chamber thermocouple subassembly is a critical part of the STAT*IM* unit. Extreme care should be taken so that the chamber thermocouple is not damaged.


Re-alignment of the assembly in the armature-isoplate subassembly is also critical. It is important that the correct tooling is available to ensure the unit performs as intended.

Note that the probe bracket and chamber thermocouple are NOT available as separate components as the subassembly alignment is set in the factory. They should only ever be replaced as an assembly.

## Removing the probe bracket and chamber thermocouple subassembly

To remove the probe bracket assembly, proceed as follows (see Figures 3, 4 and 5):

For revision 6 controller boards and earlier versions

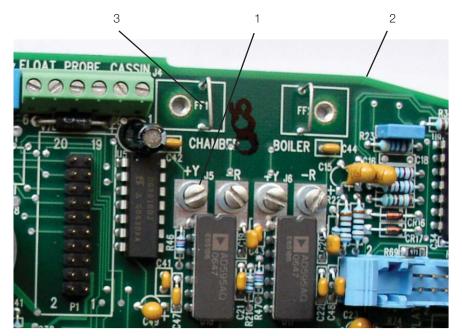


- 1. Chamber terminal positions
- 2. Controller board
- 3. Chamber flag terminal

Figure 3

1.

b. Loosen the screws and disconnect the two fine wires (1) from controller board (2) terminal positions CHAMBER +Y and -R. Disconnect the flag terminal (3) from the position marked CHAMBER directly above the terminals. Leave the screws with contact washers in the terminals.


#### For revision 7.x controller boards

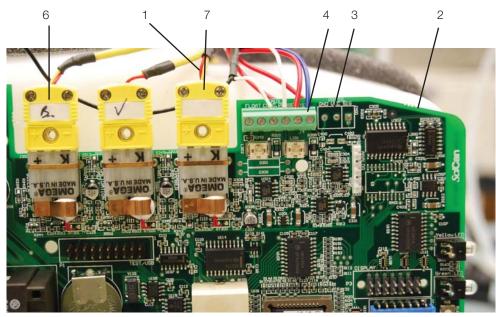

- 1. Chamber thermocouple leads
- 2. Controller board
- 3. Flag terminal
- 4. Microswitch leads
- 5. Cable tie
- 6. Steam generator Thermocouple plug
- 7. Chamber thermocouple plug

Figure 4

- c. Disconnect the chamber thermocouple plug (7) from the socket on the controller board.
- 2. Remove the microswitch (4) leads from controller board block connector positions labelled CASSIN.
- 3. Carefully cut the cable ties (5) securing the microswitch leads and other wires to the top of the armature (26).
- 4. Carefully cut the cable tie (9) holding the compressor tube (10) onto the check valve inlet (11) and pull the tube off the valve.
- 5. Disconnect the compression nut (12) holding the steam generator outlet tube (13) to the steam generator outlet fitting (14).
- 6. Remove the two screws (15) holding the steam generator bracket (16) to the chassis.
- 7. Disconnect the compression nut (17) holding the exhaust tube (18) to the solenoid valve inlet fitting (19).
- 8. Remove the two screws (20) that hold the solenoid valve assembly (21) to the chassis.
- 9. The probe bracket assembly (22) is held by four hex socket cap screws with flat washers (23). The bottom of the probe bracket assembly has two mounting slots. The top of the bracket has two mounting holes. Remove the top two screws with an Allen key. Loosen the bottom cap screws to allow the bracket to slide up and out of the chassis. Note the foam probe bracket gasket (24 not shown) mounted on the plastic end part of the armature (26).
- 10. Disconnect the compression nut (25) holding the steam generator outlet tube (13) to the probe inlet fitting (27). Remove the steam generator outlet tube (13) and the solenoid valve exhaust tube (18) from the probe bracket assembly. (Discard them if fitting a new bracket).

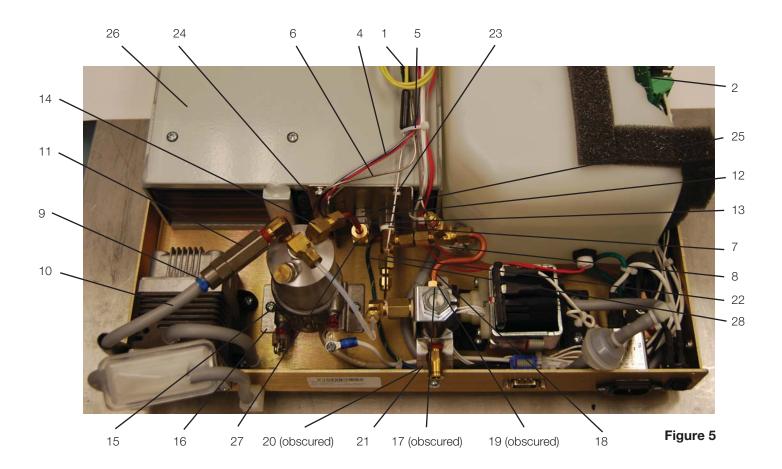

- 1. Chamber thermocouple leads
- 2. Controller board
- 3. Flag terminal (See Figure 3)
- 4. Microswitch leads
- 5. Cable tie
- 6. Steam generator thermocouple
- 7. Chamber thermocouple
- 8. Reservoir sensor leads
- 9. High temperature rated cable tie
- 10. Compressor tube
- 11. Check valve
- 12. Steam generator outlet tube compression nut

- 13. Steam generator outlet tube
- 14. Steam generator outlet fitting
- 15. Screws for item 16
- 16. Steam generator / bracket assembly
- 17. Exhaust tube compression nut (obscured)
- 18. Solenoid valve exhaust tube
- 19. Solenoid valve inlet fitting (obscured)
- 20. Screws for item 21 (one screw obscured)


- 21. Solenoid valve / bracket assembly
- 22. Probe bracket assembly
- 23. Cap screws with flat washers
- 24. Probe bracket gasket (not shown)
- 25. Compression nut for probe inlet fitting
- 26. Armature
- 27. Probe inlet fitting
- 28. Armature ground wire

# Replacing and reinstalling the probe bracket and chamber thermocouple subassembly

To replace the probe bracket assembly, proceed as follows (see Figures 5 and 6):

- 1. Inspect the probe bracket gasket (24) (not shown) on the back of the armature (26) for wear or adhesive failure. If required, remove the damaged gasket and replace it with a new one. NOTE: The gasket must be fitted.
- 2. Inspect the exposed portion of the thermocouple (29) that projects from the outlet probe (30). Note the position of the thermocouple and the orientation of the end of the thermocouple.

NOTE: do not adjust the position of the tip of the thermocouple. The tip is pre-bent during manufacturing to exact specifications.

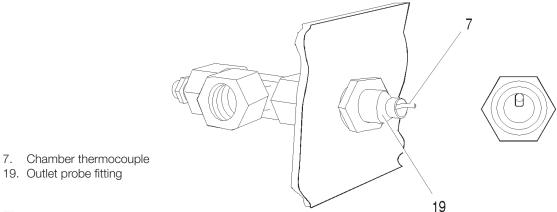

### 19. Outlet probe fitting

### Figure 6

- 3. Verify that the microswitch functions properly. (See Chapter 7. Electrical and Electronic Components: Microswitch).
- 4. Tilt the top of the probe bracket assembly (22) away from the armature and carefully insert the bracket with the bottom two cap screws and flat washers (23) and the armature ground wire terminal (28) in place onto the back of the armature. The washers rest on the outside of the probe bracket assembly. Ensure that the thermocouple leads do not come in contact with the outside of the armature or the probe gasket.
  - NOTE: be careful not to bend the thermocouple leads too sharply (min. Bend radius -3/16 inch / 5mm). Retain the protective sleeve on the end of the thermocouple leads as long as possible.
- 5. Insert two cap screws with flat washers (23) into the upper mounting holes of the probe bracket. Tighten all four cap screws so that the bracket is held in place, but can still be adjusted up, down and sideways.
- 6. Install a new copper tube (13) that extends between the steam generator outlet fitting (14) and the probe inlet fitting (27) on the probe bracket assembly (22). Do not tighten the compression nuts.
- 7. Install a new copper tube (18) that extends between the solenoid valve assembly (21) and the solenoid valve inlet fitting (19) on the probe bracket assembly (22). Do not tighten the compression nuts.

- 8. Attach the steam generator bracket (16) using the two screws (15) previously removed. Do not tighten these screws.
- 9. Tighten the compression nut (12) to the steam generator outlet fitting (14), and the compression nut (17) to the probe inlet fitting (25), using a 7/16 inch wrench.
- 10. Reinstall the solenoid valve assembly (21) using the two screws (20) previously removed. Do not tighten these screws.
- 11. Tighten the compression nuts to the solenoid valve inlet fitting (19), and the outlet probe fitting (31), using a 7/16 inch wrench.
- 12. Push the rubber tube (10) extending from the compressor onto the check valve fitting (11) on the steam generator. Secure the tube to the fitting using a high temperature rated cable tie (9).
- 13. Route the chamber thermocouple leads (7), the microswitch leads (4) and the reservoir sensor leads (8) beside the steam generator thermocouple (6) along the top of the armature over the cable anchors attached to the armature.

For all revision 6 and earlier controller boards...

14.

a. Carefully remove the protective sleeve from the end of the uninsulated chamber thermocouple leads. Pre-bend both leads to achieve the required shape. Reconnect the flag terminal to the threaded lug marked CHAMBER on the Controller Board. NOTE: Be careful, the thermocouple leads at the board end are very fragile.

One lead is colour coded: the unmarked lead is positive, +Y.


the red lead is negative, -R.

Connect the unmarked lead to the terminal marked +Y on the controller board. Connect the red lead to the terminal marked -R on the controller board. Make sure there is extra lead length so that the wires do not break as the screw is tightened. Ensure the wires are seated securely behind the respective washers (see Figure 6).

NOTE: The two leads must not touch one other or any other component. Do not calibrate a thermocouple until it is properly installed and positioned in the unit.

### For Revision 7.x controller boards

b. Connect the wires to the thermocouple connector as follows:

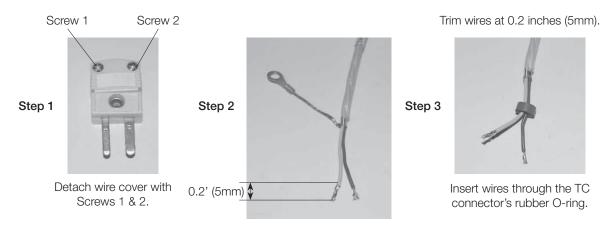

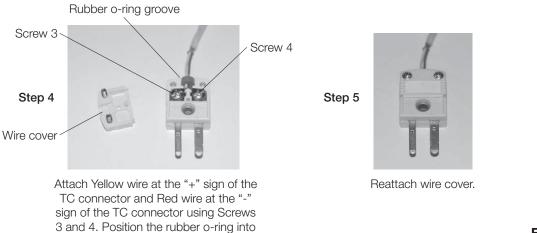


Figure 7

- 15. Connect the microswitch leads (4) to controller board (2) block terminal positions labelled CASSIN.
- 16. Bundle the thermocouple leads, microswitch and reservoir sensor leads together using nylon cable ties and secure the wires to the top of the armature using cable ties and the cable anchors attached to the armature. The leads must not touch the cover when the cover is put back on.
- 17. Align the probe bracket. (See section below).

its groove.

- 18. Ensure that all fittings have been tightened and all tubes are properly secured.
- 19. Use a cable tie to secure the micro-switch leads and reservoir sensor leads to the side of the probe bracket assembly AFTER the probe bracket has been aligned.
- 20. Check the microswitch for proper operation by inserting a cassette with the power ON and listening for the solenoid valve to click.
- 21. Calibrate the unit using the appropriate calibration procedure.

# Aligning the probe bracket and chamber thermocouple subassembly

**IMPORTANT NOTE:** Aligning the probe bracket is critical to the correct functioning of the STAT*IM* unit. It is strongly recommended that the probe bracket alignment jig is used when performing this operation.



Figure 8

Probe Bracket Alignment Jig – Part number 01-106776S (Suitable for both 2000 and 5000 models)

To align the probe bracket assembly, proceed as follows (see Figure 5):

- 1. Be sure that the four cap screws with flat washers (23) are loose enough to allow adjustment of the probe bracket, but there is sufficient friction between the bracket and the armature so that the probe bracket cannot move under its own weight.
- 2.
- a. If using the probe alignment jig:
  - 1. Insert the jig into the armature until it just touches the probe tips.
  - 2. Gently push the jig towards the probe bracket while at the same time adjusting the probe bracket until the probe tips enter the openings in the end of the alignment rods of the jig. PERFORM THIS OPERATION WITH CARE. DO NOT BEND THE TIP OF THE CHAMBER THERMOCOUPLE.
  - 3. Rotate the alignment rods and at the same time finely adjust the position of the probe bracket until the rods rotate, and slide in and out freely.
  - 4. Using the Allen key, tighten the four cap screws with flat washers in the following tightening pattern: upper left hand corner, bottom right hand corner, upper right hand corner and the bottom left hand corner.
  - Recheck the alignment rods to ensure that they still rotate and slide in and out freely. Slacken off the screws and repeat the procedure if necessary if the alignment rods are tight.

### b. If the probe alignment jig is NOT being used:

- 1. Insert a regular cassette (the cassette normally used by the customer with the unit is preferred) or a calibration cassette into the armature until it just touches the probe tips.
- Gently push the cassette towards the probe bracket, while at the same time adjusting the probe bracket until the probe tips enter the openings in the cassette lid. PERFORM THIS OPERATION WITH CARE. DO NOT BEND THE TIP OF THE CHAMBER THERMOCOUPLE.
- 3. Slowly push the cassette to the fully inserted position while continuing to adjust the probe bracket so as to centre the probes in the cassette openings.
- 4. Withdraw the cassette and repeat steps 2 and 3 as often as required to centre the probes correctly in the openings.
- 5. Using the Allen key, tighten the four cap screws with flat washers in the following tightening pattern: upper left hand corner, bottom right hand corner, upper right hand corner and the bottom left hand corner.
- 3. Ensure that all fittings have been tightened and all tubes are properly secured.
- 4. Plug in the power cord and turn the power switch ON. Start an UNWRAPPED cycle and watch the LCD until the PRESSURIZING or CONDITIONING phase of the cycle is displayed.
- 5. At this point, turn the power switch OFF. This will close the solenoid valve and should retain pressure in the machine.

Caution: the steam generator, check valve, probe bracket and associated tubing will be hot. There is steam pressure in the unit so care should be taken when observing the unit. Guard against burns.

- 6. Observe the unit and check for visible or audible signs of steam leakage. Note any leaking joints.
- 7. Turn the unit back on, and allow it to vent.

### Caution: allow the unit to cool sufficiently before attempting adjustments.

- 8. Tighten any leaking joints or fine tune the probe bracket alignment as above if appropriate. Repeat steps 4 to 6 as necessary until no leaks are observed.
- 9. Make sure that all brackets are secured to the chassis and cable ties are installed.
- 10. Calibrate the unit using the appropriate calibration procedure.
- 11. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



# **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



# PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# Steam Generator, Check Valve, Thermal Fuse and Pressure Relief Valve

# **Steam Generators**

There are four types of steam generator you may encounter while servicing a STATIM 1000 / 2000 / 2000S unit.

# 1st Generation aluminium steam generator with internal thermocouple

First generation aluminium steam generator with embedded (internal) thermocouple. These were factory fitted in units up to the end of 2003 and were associated with Rev. 2.x/3.x./4.x and 5.x controller boards and Rev. 1.xx/2.xx software. Later controller boards may have been retrofitted with the appropriate software on field upgraded units.



### Stainless steel steam generator

These were factory fitted from 2004 to 2005 and were associated with Rev. 6.x controller boards and Rev. 4.xx software.



# 2nd Generation aluminium steam generator with external thermocouple (Alex)

These were factory fitted from 2005 and were associated with Rev. 6.x controller boards plus Rev. 5.xx software, or Rev. 7.x controller boards plus Rev. 6.xx software.



Figure 1

# Upgraded 2nd generation Alex steam generator (Boris) NOTE: This version still operates using Rev. 5.x/6.x software.

In mid 2008, SciCan implemented an O-ring seal design to seal the joint between the cap and base of the boiler assembly to improve the manufacturability of the steam generator assembly. The main difference with regard to servicing is that the cap and base can be rotated slightly relative to each other when the boiler assembly is not mounted in the machine. This new seal system does not affect the performance of the sterilizer in any way.

When STATIM steam generator spare part kits are assembled at the factory, they are wrapped and packaged such that the relationship between cap and base will be maintained. If this relationship is inadvertently changed, it can be returned to the correct position by rotating the cap relative to the base. This should be done such that the vertical marks on the boiler base and the boiler cap line up as defined at the factory (see picture) and the gap from the cap to the base is between 1mm (1/32") and 2mm (1/16").



Figure 1 (cont'd)

# Replacing 1st generation and stainless steel steam generators

- These two steam generator types are no longer produced and can ONLY be upgraded with an Alex derived steam generator.
- The upgrade may vary from a complete steam generator and PCB replacement to a steam generator and microprocessor replacement (with or without adaptor boards as appropriate) depending on the model and age of the unit.
- All units factory fitted with, or upgraded to Alex type steam generators will only need a steam generator replacement.
- The following matrix provides the upgrade kit reference numbers and information:


Model	Current Software Revision	Current Steam Generator Type	Current PCB Revision	Alex Upgrade Kit	Kit Contents — Main Components
2000 110 V	-	Aluminium	3.x / 4.x	01-109100S	Alex Boiler, Full R7.x PCB, Vaidation thermocouple, Thermocouple interface PCB
2000 110 V	S201R1xx S201R2xx	Aluminium	2.4-2.9x ONLY	01-108996S	Alex Boiler, Microprocessor adapter board, Vaidation thermocouple, Thermocouple PCB
2000 110 V	S201R4xx	Aluminium	6.x	01-108993S	Alex Boiler, Chipset, Vaidation thermocouple, Thermocouple PCB
2000 230 V	-	Aluminium	3.x / 4.x	01-109717S	Alex Boiler, Full R7.x PCB, Vaidation thermocouple, Thermocouple interface PCB
2000 230 V	S202R1xx S202R2xx	Aluminium	2.4-2.9 / 5.0- 5.3	01-109000S	Alex Boiler, Microprocessor adapter board, Vaidation thermocouple, Thermocouple interface PCB
2000S 230 V	S2S2R1xx S2S2R2xx	Aluminium	2.4-2.9 / 5.0- 5.3	01-108991S	Alex Boiler, Microprocessor adapter board
2000S 230 V	S2S2R4xx	Stainless Steel	6.x	01-108989S	Alex Boiler, R5.x software chipset, Pump tube set for Aluminium boiler

Figure 2

# Removing the Aluminium Steam Generator (for both 1st generation and Alex versions)

CAUTION: The steam generator may be hot if the unit has been operating. Guard against burns.

To remove the steam generator, proceed as follows (see Figure 3):

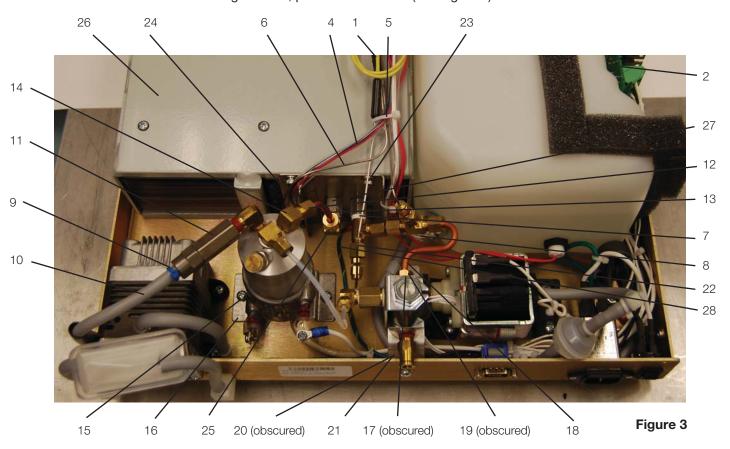

- 1. Chamber thermocouple leads
- 2. Controller board
- 3. Flag terminal (See Figure 3)
- 4. Microswitch leads
- 5. Cable tie
- 6. Steam generator thermocouple
- 7. Chamber thermocouple
- 8. Reservoir sensor leads
- 9. High temperature rated cable tie
- 10. Compressor tube
- 11. Check valve
- 12. Steam generator outlet tube compression nut

- 13. Steam generator outlet tube
- 14. Steam generator outlet fitting
- 15. Screws for item 16
- 16. Steam generator / bracket assembly
- 17. Exhaust tube compression nut (obscured)
- 18. Solenoid valve exhaust tube
- 19. Solenoid valve inlet fitting (obscured)
- 20. Screws for item 21 (one screw obscured)

- 21. Solenoid valve / bracket assembly
- 22. Probe bracket assembly
- 23. Cap screws with flat washers
- 24. Probe bracket gasket (not shown)
- 25. Compression nut for probe inlet fitting
- 26. Armature
- 27. Probe inlet fitting
- 28. Armature ground wire

1. Turn the power switch **OFF**, and unplug the unit. Remove cover.

2.

- a. For all controller boards other than revision 7.x Remove the pressure interface board or printer cable (if fitted) from the blue socket and disconnect the steam generator thermocouple wires (1) from Controller Board (2) terminal positions BOILER +Y and -R and disconnect the ground lead terminal (3) from the position marked BOILER directly above the terminals. Leave the screws with contact washers in the terminals.
- b. For revision 7.x controller boards disconnect the chamber thermocouple plug from the socket on the controller board.
- 3. Disconnect the black thermal fuse (4) wire from controller board connector terminal block J1-3.
- 4. Trace the path of the black wire back to the base of the steam generator (5). Carefully cut all cable ties holding the black wire.
- 5. Carefully cut the cable ties (6) securing the steam generator thermocouple lead (8) and other wires to the armature (11).
- 6. Disconnect the white wire (12) attached to the terminal on the lower half of the steam generator.
- 7. On units where a compressor is fitted, carefully cut the cable tie (13) holding the compressor tube (14) onto the check valve inlet (15) and pull the tube off the valve.
- 8. Disconnect the compression nut (16) holding the Teflon™ tube (17) from the top of the steam generator.
- 9. Disconnect the compression nut (18) holding the steam generator outlet tube (19) to the steam generator outlet fitting (20).
- 10. Disconnect the compression nut (21) holding the steam generator outlet tube to the probe bracket inlet fitting (22).
- 11. Remove the two screws (23) holding the steam generator bracket assembly (24) to the chassis and remove the steam generator.

# Reinstalling the Aluminium Steam Generator (for both 1st generation and Alex versions)

**NOTE:** DO NOT use this section when REPLACING or UPGRADING older 1st generation aluminium steam generators to the newer Alex type. These instructions apply only to RE-INSTALLING an aluminium steam generator (including Alex) that has been removed.

Instructions for UPGRADING to the Alex steam generator vary depending on the status of the unit that is being upgraded (i.e., what controller and software revisions it has) and may also change from time to time. For the most up to date procedures concerning an upgrade, use the matrix in Figure 2 to find the appropriate kit number and follow the instructions included with the kit.

To reinstall the steam generator, proceed as follows (see Figure 3):

- 1. Make sure the power switch is OFF, and the unit is unplugged.
- 2. Place the steam generator in to position.
- 3. Connect the compression nut (18) holding the steam generator outlet tube (19) to the top of the steam generator outlet fitting (20) and the compression nut (21) holding the steam generator outlet tube to the probe bracket inlet fitting (22). Thread the nuts and finger tighten, then tighten with a wrench. **Do not over tighten**.

- 4. Tighten the two screws holding the steam generator bracket assembly to the chassis.
- 5. Connect the compression nut (16) holding the Teflon™ tube (17) to the top of the steam generator. Thread the nut finger tight, then tighten. Do not over tighten.
- 6. Reconnect the white wire (12) to the terminal on the lower half of the steam generator using the retained screw with lockwasher. Tip: If the terminal appears blackened, clean it using fine grit sandpaper. If the terminal threads are stripped, use a nut on the other side of the terminal to hold the screw.
- 7. Route the black thermal fuse wire (4) from the fuse assembly along the back of the chassis along the wiring harness, and connect it to controller board. Re-secure the harness with cable ties.
- 8. Check the bend and route of the steam generator thermocouple lead (8) and place alongside the chamber thermocouple lead to the controller board. The leads must not touch the cover when the cover is assembled to the chassis.
- 9. Reconnect the ground lead terminal (3) to the position marked BOILER directly above the thermocouple terminals on the controller board.

10.

a. For all controller boards other than revision 7.x. Check the pre-bend on both thermocouple leads to ensure that they are the required shape to go under the washers on the connections without touching anything other than the terminal.

### **CAUTION:** The thermocouple leads at the board end are very fragile.

The unmarked lead is positive, +Y, and red lead is negative, -R. Connect the unmarked lead to the terminal marked +Y on the controller board. Connect the red lead to the terminal marked -R on the controller board. Make sure there is extra lead length so that the wires do not break as the screw is tightened. Ensure that the wires are seated securely behind the respective washers (see Chapter 7, Figures 3 and 4).

The two leads must not touch one another or any other component. Do not calibrate a thermocouple until it is properly installed and positioned in the unit.


- b. For revision 7.x controller boards. Reconnect the chamber thermocouple plug to the socket on the controller board (see Chapter 7, Figures 3 and 4).
- 11. Bundle the chamber thermocouple lead (7), the steam generator thermocouple lead, the validation thermocouple lead (where fitted), pressure transducer leads (where fitted), micro switch leads (9) and reservoir sensor leads (10) together using nylon cable ties (6), approximately every 2-3 inches. Secure the wires to the top of the armature using the cable anchors provided.
- 12. Carefully push the compressor tube (14) (where a compressor is fitted) onto the check valve inlet (15) and secure the tube to the valve using a high temperature rated cable tie (13).
- 13. Replace the pressure interface board or printer cable (if fitted) into the blue socket.
- 14. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.

**NOTE:** These tests must be performed on the STATIM again once the work is completed and the cover has been returned to the unit.

- 15. Refill the reservoir using steam process distilled water.
- 16. Recalibrate the unit as appropriate.
- 17. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.
- 18. Reinstall the cover.

# Removing the Stainless Steel Steam Generator

To remove the steam generator, proceed as follows (see Figure 4):

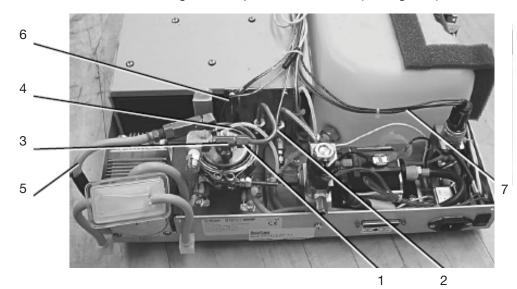

- 1. Compression nut
- Steam generator outlet tube
- 3. T-fitting
- 4. Teflon™ tube
- Compressor tube
- Steam generator thermocouple lead
- 7. Transducer leads

Figure 4

- 1. Turn the power switch **OFF**, and unplug the unit. Remove cover.
- 2. Referring to Chapter 7, Electrical and Electronic components, Figures 3 and 4, remove the pressure interface board or printer cable (if fitted) from the blue socket and disconnect the steam generator thermocouple wires from the controller board terminal positions BOILER +Y and -R and disconnect the ground lead terminal from the position marked BOILER directly above the terminals. Leave the screws with contact washers in the terminals.

### Refer back to Figure 4 (above) for the remainder of the procedure.

- 3. Disconnect the electrical connections of both steam generator power cables at the PC board. Look for the connections on the left side of the board marked "BL". Disconnect wire N and wire L of the "BL" set. Disconnect the black thermal fuse (4) wire from controller board connector terminal block J1-3.
- 4. Trace the path of the wires to the steam generator and cut the cable ties holding them in place. Separate the wires from the wiring harness.
- 5. Carefully cut the cable ties securing the steam generator thermocouple lead (6) and other wires to the armature.
- 6. On units where a compressor is fitted, carefully cut the cable tie holding the compressor tube (5) onto the check valve inlet and pull the tube off the valve.
- 7. Disconnect the compression nut holding the Teflon™ tube (4) from the top of the steam generator.
- 8. Using a 7/16 wrench on the compression nut (1) connected to the steam generator outlet tube and a 3/8 wrench on the T-fitting (3), disconnect steam generator outlet tube (2). When the nut is unscrewed completely, pull the tube gently away from the T-fitting to disengage it from the fitting.
- 9. Loosen the mounting screws that hold the steam generator bracket to the chassis and, using the bracket's key slot, slide the steam generator out of position and remove.

### Reinstalling the Stainless Steel Steam Generator


**NOTE:** This section is NOT for REPLACING or UPGRADING stainless steel steam generators to the newer Alex type. It is only to be used for RE-INSTALLING the stainless steel steam generator that has been removed.

Instructions for UPGRADING to the Alex unit vary depending on the unit that is being upgraded (i.e., what controller and software revisions it has) and may also change from time to time. For the most up to date procedures concerning an upgrade, use the matrix in Figure 2 to find the appropriate kit number and follow the instructions included with the kit.

To reinstall the steam generator, proceed as follows (see Figures 4, 5, 6, 7, and 8):

When servicing a STATIM 2000 / 2000S equipped with a stainless steel steam generator, you may encounter one of two steam generator bracket configurations (see Figures 5-8). Type A has a steam generator bracket and a separate bracket adapter. Type B has an integral modified steam generator bracket.

Both scenarios ensure that, when properly installed, the steam generator terminals are angled slightly downward.

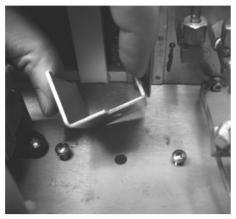

Positioning bracket adapter for Type A Figure 5



Type A configuration





Type A configuration installed

Figure 7



Type B configuration

Figure 8

- Make sure the power switch is OFF, and the unit is unplugged.
- 2.
- a. To refit a Type A configuration, place the bracket into position by aligning the dimple on the bottom of the adapter with the hole in the chassis (See Figure 5). Place the stainless steel steam generator on top of the bracket adapter (See Figure 6). Align the bracket's key slots with the screws on the chassis and tighten just enough to keep the assembly in place (See Figure 7).
- b. To refit a Type B configuration (See Figure 8), align the steam generator and integral bracket key slots with the screws on the chassis. But do not tighten yet.

### Refer back to Figure 4 for the remainder of the procedure.

- 3. Reconnect the compression nut (1) on the steam generator tube to the T-fitting (3) on the steam generator. Only finger tighten at this stage.
- 4. When the nut is fully on the fitting, tighten the two screws holding the steam generator bracket assembly to the chassis.
- 5. Now, using a 3/8 wrench on the T-fitting and a 7/16 wrench on the compression nut, tighten the two together. **Do not overtighten**.
- 6. Reconnect the compression nut of the Teflon™ tube (4) to the top of the steam generator. Thread the nut finger tight, then tighten. **Do not overtighten**.
- 7. Carefully push the compressor tube (5) onto the check valve inlet and secure the tube to the valve using a high temperature rated cable tie.
- 8. Re-route the power cables in the wiring harness along the back of the chassis and connect them to the controller board. The black wire to goes to connection marked L and the white wire to connection marked N. Re-secure the wiring harness with cable ties.
- 9. Check the bend and route of the steam generator thermocouple lead (6) and place alongside the chamber thermocouple lead to the controller board. The leads must not touch the cover when the cover is assembled to the chassis.
- 10. Reconnect the ground lead terminal to the position marked BOILER directly above the thermocouple terminals on the controller board.

Check the pre-bend on both thermocouple leads to ensure that they are the required shape to go under the washers on the connections without touching anything other than the terminal.

### CAUTION: The thermocouple leads at the board end are very fragile.

The unmarked lead is positive, +Y, and red lead is negative, -R. Connect the unmarked lead to the terminal marked +Y on the controller board. Connect the red lead to the terminal marked -R on the controller board. Make sure there is extra lead length so that the wires do not break as the screw is tightened. Ensure that the wires are seated securely behind the respective washers (see Chapter 7, Figures 3 and 4).

The two leads must not touch one another or any other component. Do not calibrate a thermocouple until it is properly installed and positioned in the unit.

- 11. Bundle the chamber thermocouple lead, the steam generator thermocouple lead, the validation thermocouple lead (where fitted), pressure transducer leads (where fitted), micro switch leads and reservoir sensor leads together using cable ties, approximately every 2-3 inches. Secure the wires to the top of the armature using the cable anchors provided.
- 12. Replace the pressure interface board or printer cable (if fitted) into the blue socket.

- 13. A dielectric strength test (Hi-Pot) and a protective bonding impedance test (ground continuity) should be performed on the STATIM unit at this stage.
  - **NOTE:** These tests must be performed on the STAT*IM* again once the work is completed and the cover has been returned to the unit.
- 14. Refill the reservoir using steam process distilled water.
- 15. Recalibrate the unit as appropriate.
- Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.
- 17. Reinstall the cover.

# **Check Valves**

The check valve on STAT*IM* units is basically a non return valve that allows air from the compressor to pass through the steam generator, into the chamber (cassette) to dry the instruments, and exhaust through the solenoid valve to the waste bottle after the unit has vented. It should not allow steam to pass back through the valve and into the compressor. Failure can lead to a drop in pressure in the unit during processing and cause damage to the compressor. Note: As a critical component it is recommended that the check valve is changed annually as part of a routine service.

CAUTION: The steam generator and check valve may be hot. Guard against burns.

The steam generator check valve is only found on STATIM 2000 and 2000S units with compressors.

# Testing the check valve

There are two options of testing the check valve, either with the unit running or manually with the control box.

Option 1: Checking the valve with the unit running.

To perform this test, the STATIM unit must be running with the cover off and the control box attached as if running as test cycle during servicing. For revision 7.x units, the unit can be run with the cover removed but still connected to the controller.

To inspect the check valve, proceed as follows (See Figure 9):

- 1. Carefully pull the compressor tube (2) from the check valve barb fitting.
- 2. Start a normal cycle on the unit and wait for the unit to pressurise.
- 3. Using a reflective, cold surface (e.g. a dental mirror), place it close to the check valve inlet fitting and observe for condensation.
- 4. If condensation is apparent, manually inspect the compressor air filter, the bacteria retentive filter (if present) and the compressor.
- 5. If there is evidence of water in the bacteria retentive filter, replace the check valve and the filter. If there is evidence of water in the compressor or if the compressor filter is wet, replace the steam generator check valve, compressor and filters.


### Option 2: Performing a manual valve inspection.

If it is not possible to run the unit, but it can be powered up, proceed as follows:

- 1. With the power off, remove the cover.
- 2. Attach the control box to the unit and power on the unit, or for revision 7.x controller board units, power the unit on first and select 'compressor' in the device test sub menu of the service menu (if cover is removed but still connected). The unit should be powered on with the cassette removed.
- 3. Turn ON the compressor using the control box or keypad as appropriate, and allow it to run for a few minutes. Allow the compressor to cool until it feels cool to the touch.
- 4. Inspect the compressor air filter, the bacteria retentive filter (if present) and the compressor.
- 5. If there is evidence of water in the bacteria retentive filter, replace the check valve and the filter. If there is evidence of water in the compressor or if the compressor filter is wet, replace the steam generator check valve, compressor and filters.

# Removing and Replacing the Steam Generator Check Valve

To remove and replace the steam generator check valve, proceed as follows (See Figure 9):

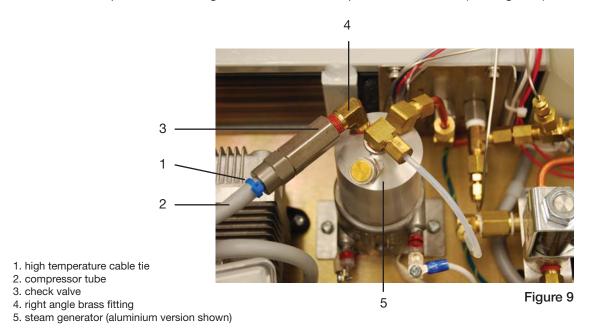

- 1. Turn the power switch OFF, and unplug the unit.
- 2. Carefully cut the cable tie (1) holding the compressor tube (2) onto the check valve (3) inlet and pull the tube off the valve.
- 3. Using a ½" wrench remove the check valve. Apply the wrench to the part of the valve closest to the steam generator (5). Use a wrench to hold the right angle brass fitting so that it does not move. Clean any Teflon™ tape or assembly compound from the threads on the right angle valve fitting.

- 4. Prepare the right angle (4) fitting by wrapping the threads with Teflon™ pipe tape in the rotational direction of the thread. (minimum 4 complete turns).
- 5. Thread the new valve onto the right angle valve fitting on the top of the steam generator (5). Using a wrench on the end of the valve closest to the fitting, tighten the valve until medium resistance is felt. Do not overtighten the new valve. Use a wrench to hold the right angle brass fitting so it does not move.
- 6. Carefully push the compressor tube (2) onto the check valve (3) and secure the tube to the valve using a high temperature application cable tie (1).
- 7. Install replacement filters as required.
- 8. Run a sterilization cycle and observe for leaks.
- 9. Reinstall the cover.


# Steam Generator Thermal Fuses

You may encounter four types of thermal fuse configurations on STAT*IM* units. Three for aluminium steam generators and one for the stainless steel generator. Early 1st generation aluminium steam generators had a single thermal fuse located in series in the main power supply line to the generator (the black wire) which was held in place by a mounting bracket on the underside of the steam generator.

This was replaced with a double fuse, again, in series in the main power supply line and sandwiched between the bracket and the steam generator base. Unlike the single fuse, this protected more of the surface area of the underside of the steam generator.

It should be noted that the steam generator bases/brackets may be different, so the appropriate fuse must be obtained for each steam generator.

The three configurations for aluminium steam generators are as follows (see Figure 10):

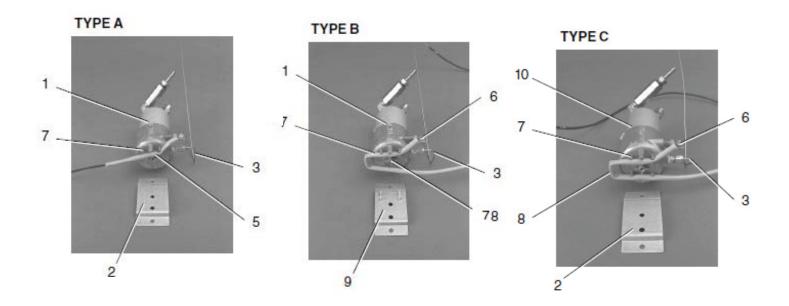

- steam generator
- steam generator bracket
- 3. thermocouple lead
- 4. screw (not shown)
- thermal fuse

- steam generator electrical terminal
- spacing stand-offs
- double thermal fuse
- modified steam generator bracket


Figure 10

Type A is a single thermal fuse wire (5), steam generator bracket (2) and an unmodified steam generator (1).

Type B is a double thermal fuse wire (8) with a modified steam generator bracket (9) and an unmodified steam generator (1).

Type C is a double thermal fuse wire (8), steam generator bracket (2) with a modified steam generator (10).

The fuse configuration for stainless steel steam generators is as follows (see Figure 11):

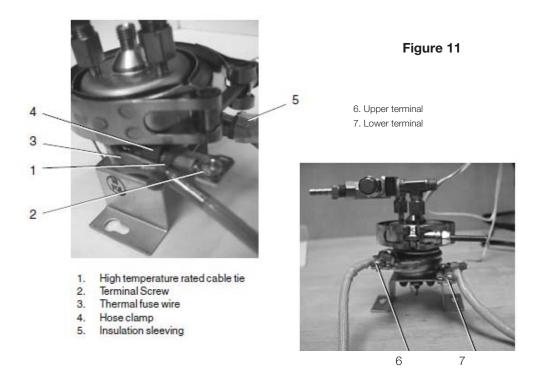


# Removing the Thermal Fuse on the Aluminium Steam Generator (for all models)

To remove the thermal fuse(s) located between the steam generator and the steam generator bracket, proceed as follows (see Figure 12):

- 1. Carefully cut the high temperature rated cable tie (13) holding the compressor tube (14) onto the check valve inlet (15) and pull the tube off the valve.
- 2. Disconnect the black thermal fuse lead wire (4) from controller board connector terminal J1-3 and the white wire (12) from the steam generator.
- 3. Trace the path of the black wire back to the base of the steam generator (5). Cut all the cable ties holding the black wire.
- 4. Using a wrench, disconnect the compression nut (16) holding the Teflon™ tube (17) from the top of the steam generator.
- 5. Using a wrench, disconnect the compression nut (21) holding the steam generator outlet tube (19) to the probe bracket inlet fitting (22) and the compression (18) nut holding the steam generator outlet tube to the steam generator outlet fitting (20).

- 6. Remove the two screws (23) from the steam generator bracket (24). Note: The steam generator will still be attached to the PCB by the thermocouple lead (8).
- 7. Carefully cut the two cable ties closest to the steam generator that hold the thermocouple leads together.
- 8. Gently lift and turn the steam generator assembly onto one side to expose the bottom of the assembly. Be careful not to stress the thermocouple leads. (Min. bend radius 3/16 inch / 5 mm).
- 9. Remove the small screw that attaches the steam generator bracket to the steam generator.
- 10. Disconnect the thermal fuse (4) from the lower power terminal on the steam generator. Note the routing of the thermal fuse assembly between the bracket and spacing lugs on the bottom of the steam generator.
- 11. Remove the thermal fuse.

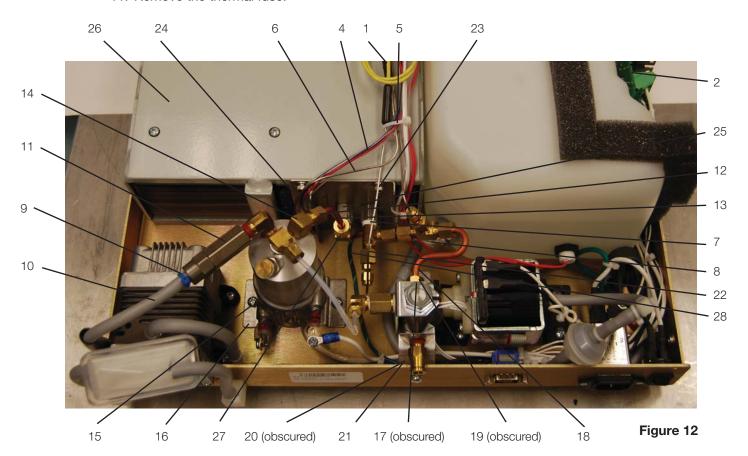

- 1. Chamber thermocouple leads
- 2. Controller board
- 3. Flag terminal (See Figure 3)
- 4. Microswitch leads
- 5. Cable tie
- 6. Steam generator thermocouple
- 7. Chamber thermocouple
- 8. Reservoir sensor leads
- 9. High temperature rated cable tie
- 10. Compressor tube
- 11. Check valve
- 12. Steam generator outlet tube compression nut

- 13. Steam generator outlet tube
- 14. Steam generator outlet fitting
- 15. Screws for item 16
- 16. Steam generator / bracket assembly
- 17. Exhaust tube compression nut (obscured)
- 18. Solenoid valve exhaust tube
- 19. Solenoid valve inlet fitting (obscured)
- 20. Screws for item 21 (one screw obscured)

- 21. Solenoid valve / bracket assembly
- 22. Probe bracket assembly
- 23. Cap screws with flat washers
- 24. Probe bracket gasket (not shown)
- 25. Compression nut for probe inlet fitting
- 26. Armature
- 27. Probe inlet fitting
- 28. Armature ground wire

# Replacing the Thermal Fuse on the Aluminium Steam Generator (for all models)

Determine which configuration is being serviced and use the appropriate parts to service the unit.

To replace the thermal fuse proceed as follows (refer to Figure 10):

1. Connect the terminal of the appropriate thermal fuse (5 for type A - single) or (8 for type B/C - double) to the lower power terminal on the steam generator (1 for type A/B) or (10 for type C).

2.

- a. For Type A configuration, replace the steam generator bracket (2) with a modified steam generator bracket (9). Route the wire from the fuse assembly (5) between the standoffs (7) and the standoffs on the modified steam generator bracket. The fuse must be placed exactly between the standoffs. Attach the steam generator bracket to the steam generator (1) using the small screw (4). Do not pinch or crush the wire between components. Ensure that the screw is tightened securely.
- b. For Type B configuration, route the wire from the fuse assembly (8) between the standoffs (7) on the bottom of the steam generator (1) and the standoffs on the modified steam generator bracket (9). The fuse must be placed exactly between the standoffs. Attach the steam generator bracket to the steam generator (1) using the small screw (4). Do not pinch or crush the wire between components. Ensure that the screw is tightened securely.
- c. For Type C configuration, route the wire from the fuse assembly (8) between the two sets of standoffs (7) on the bottom of the modified steam generator (10).

The fuse must be placed exactly between the standoffs. Attach the steam generator bracket (2) to the steam generator using the small screw (4). Do not pinch or crush the wire between components. Ensure that the screw is tightened securely.

NOTE: For the remainder of the procedure, refer to Figure 12 and proceed as follows.

- 3. Carefully return and attach the steam generator assembly to the chassis using two screws (23). Start the screws into the tapped hole but do not tighten the screws. Be careful not to stress the thermocouple leads. (Min. bend radius 3/16 inch / 5 mm).
- 4. Connect the compression nut (18) holding the steam generator outlet tube (19) to the steam generator outlet fitting (20) and the compression nut (21) holding the steam generator outlet tube (19) to the probe bracket inlet fitting (22). Tighten finger tight.
- 5. Connect the compression nut (16) holding the Teflon™ pump tube (17) to the top of the steam generator. Tighten finger tight, then continue to tighten the nut using a 3/8-inch wrench. Do not overtighten.
- 6. Tighten the two screws (23) that attach the steam generator assembly to the chassis.
- 7. Using a wrench, tighten the compression nut holding the steam generator outlet tube to the steam generator outlet fitting and also the compression nut holding the steam generator outlet tube to the probe bracket fitting.
- 8. Route the black wire (4) from the fuse assembly and connect the wire to Controller Board connector J1-3. Connect the white wire (12) to the steam generator electrical terminal.
- 9. Bundle the loose wires together in the wiring harness and secure using nylon cable ties every 2-3 inches.
- 10. Reinstall the compressor tube (14) onto the check valve inlet (15) and secure the tube to the valve with a high temperature rated cable tie (13).

- 11. A dielectric strength test (Hi-Pot) AND a protective bonding impedance test (ground continuity) must be performed on the STAT/M to confirm that the wire has not been pinched or crushed between the steam generator and the bracket. The tests must be performed again once the work is completed and the cover has been returned to the unit.
- 12. Refill the reservoir using steam process distilled water.
- 13. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.
- 14. Reinstall the cover.

# Removing and Replacing the thermal fuse on the Stainless Steel Steam Generator

Instructions for changing the thermal fuses on the stainless steel steam generator are NOT included in this manual. If the thermal fuses have failed on a stainless steel steam generator, it is strongly recommended that the steam generator be replaced by a modern "Alex" type unit.

# Pressure Relief Valve (PRV)


There are two types of Pressure Relief Valves (PRV) fitted to STATIM units. Note that they are unit specific and the two types SHOULD NOT be interchanged.

The two types are as follow:

### Type 1: 70 p.s.i. PRV without pull ring

This type of PRV are factory fitted on steam generators on the following units:


- ALL STATIM 2000 non S (North American specification) units irrespective of steam generator type.
- STATIM S units (European specification) with 1st generation steam generator from 1996 to 2003.

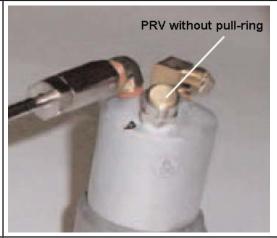

Type 2: 43.5 p.s.i. PRV with pull ring

This type of PRV are factory fitted on steam generators on the following units:

 STATIM S units (European specification) with stainless steel or "Alex" type steam generators (2004 to date).

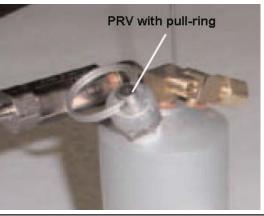

The function of the PRV is to protect the STAT/M unit by allowing steam to vent to atmosphere in the event of an overpressure situation NOT detected by the software.

**Important:** This device is a safety critical item. The checking and replacing of the device may be subject to various pressure regulations in some countries.

It is recommended that you familiarize yourself with the local regulations with regard to pressure systems before attempting to check/replace the PRV.

In the event of a requirement to change the PRV, always replace like with like as the unit will have been designed and validated for that specific PRV.

# Checking the Pressure Relief Valve

Caution: checking of the PRV under running conditions should only be undertaken by visual means. Manually activating the PRV should not be undertaken under any circumstances other than by authorised pressure vessel testing engineers.

To perform this test, the STATIM unit must be running with the cover off and the control box attached as if running as test cycle during servicing. For revision 7.x units the unit can be run with the cover removed but still connected to the controller.

To check the PRV, proceed as follows:

- 1. Start a normal cycle on the unit and wait for the unit to reach sterilization temperature and pressure.
- 2. Using a reflective, cold surface (e.g. a dental mirror), place it close to the PRV and observe for condensation.
- 3. If condensation is apparent when the unit is at optimum operating pressure then the PRV has failed and should be replaced.

Important: Checking of the upper limit of the PRV can only be done with the PRV removed from the steam generator and should only be undertaken in accordance with local regulations and by an authorized pressure vessel testing engineer.

# Removing and Replacing the Pressure Relief Valve

To remove the PRV, proceed as follows:

The unit should be OFF, with the mains cable removed from the socket and the cassette removed from the armature. No pressure or power should be present while performing this operation.

The cover should also be removed.

It is also recommended that the steam generator and associated components are cool to the touch before proceeding.

1. Using a ½" wrench, CAREFULLY unscrew the PRV from the steam generator and remove.

Caution: On aluminium steam generators that have not had the PRV removed or changed for some time, removing the PRV may damage the internal thread on the steam generator.

In certain circumstances the thread may be recovered using a 1/8" NPT taper tap, but the steam genrator will need to be removed and inverted to perform this function so that any debris caused by the re-threading can fall away from the inside of the unit. Under no circumstances should debris be allowed to enter/remain in the steam generator.

If the thread can not be recovered and correct fitting of the PRV can not be safely assured, then the complete steam generator should be replaced.

2. If the internal thread is OK, remove any P.T.F.E. tape CAREFULLY with tweezers so that when the new PRV is inserted the tape debris will not enter the steam generator.

To replace the PRV, proceed as follows:

- Prepare the new PRV by wrapping the threads with Teflon<sup>™</sup> pipe tape in the rotational direction of the thread. (minimum 4 complete turns).
- 2. Thread the new PRV into the steam generator by hand.
- 3. Using the ½" wrench, tighten the PRV until medium resistance is felt. Ideally a MINIMUM of 4 turns should be achieved. Do not over tighten the new valve.
- 4. Run a sterilization cycle and observe all fittings and tubes for leaks. Check LCD read-outs for messages indicating cycle status.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



# **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



# PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# **Printer and Data Logger**


# **Printer**

CAUTION: Never pull printer paper backwards through the printer. This will damage the print head. STATprinter units contain electronic components that may be damaged or destroyed by electro-static discharge (esd). Observe appropriate safeguards when servicing. Always transport electronic components and assembles in static-protected packaging.

The STATPrinter is an optional recording device for use with STATIM 2000/2000S units with revision 2.x/5.x/6.x controller boards and the 25 pin communication port on the rear of the unit. Newer, revision 7 STATIMs with the 9 pin RS 232 communication port CANNOT use the Stat Printer. If a printer is required for a revision 7 unit, consult the list later in this section on page 211.

# Removing and Replacing the STAT Printer Cover

To remove the STATprinter cover (1), proceed as follows, (see Figure 1):

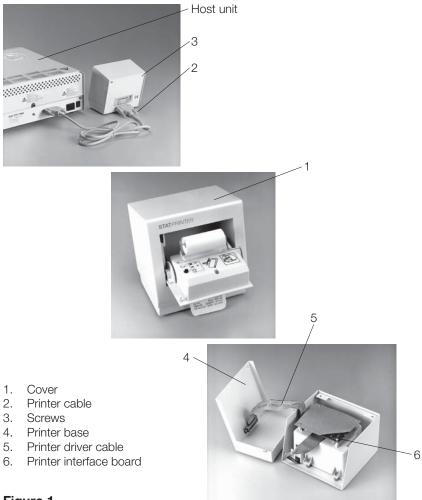

- 1. Turn the unit power OFF and turn Printer OFF.
- 2. Disconnect the printer cable (2) from the STATIM's printer connector port and the STATprinter.
- 3. Remove the four screws (3) securing the cover to the printer base (4) and retain the screws for reassembly. Remove the cover.
- 4. Disconnect the printer driver ribbon cable (5) from the printer interface board (6).


To replace the STATprinter cover, follow these steps, (see Figure 1):

- 1. Connect the printer driver ribbon cable (5) to the printer interface board (6).
- 2. Reinstall the cover (1) on the printer base (4) using the four screws (3) retained from disassembly.
- 3. Connect the printer cable connector (2) to the printer connector port at the rear of the unit and the STATprinter.
- 4. Power the unit ON.
- 5. Power the printer **ON**.

# Removing and Replacing the Printer Module

To remove the printer module assembly (1), proceed as follows (see Figure 2):

Printer interface board shield — with cut-away to show board

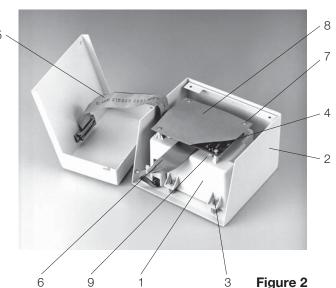

- 1. Printer module assembly
- 2. Printer cover
- 3. Screws (for module)
- 4. Printer interface board
- 5. Printer driver cable
- 6. Printer interface ribbon cable
- 7. Screws (for board)
- 8. Printer interface board shield
- 9. R21 contrast adjustment
- 1. Place the printer cover (2) containing the printer module (1) on a clean work surface to avoid scratching the cover surface.
- 2. Remove four screws (3) securing the printer module assembly to the printer cover. Retain the screws.
- 3. Remove the module from the printer cover.

To replace the printer module (1), follow these steps (see Figure 2):

- 1. Install the printer module assembly (1) in the printer cover (2) using the four screws (3) retained during disassembly.
- 2. If the module contains a new printer interface board, or a new printer, the print quality may require adjustment. See 'Printer Interface Board' and 'Adjusting Print Quality' below.

# Removing and Replacing the Printer Interface Board

To remove the printer interface board (4), proceed as follows (see Figure 2):

- 1. Disconnect the printer driver ribbon cable (5) and printer interface ribbon cable (6).
- 2. Remove the four screws (7) from the printer interface board. Retain the screws.
- 3. Remove and retain the printer interface board and the printer interface board shield (8).

To replace the printer interface board, proceed as follows (see Figure 2):

- 1. Connect the printer interface ribbon cable (6) to printer interface board header P1. Cable assemblies may differ. If the connector is not polarized, note the orientation of Pin 1 of the connector and Pin 1 of the board.
- 2. Connect the printer ribbon cable (5) connector to printer interface board header P2.
- 3. Place the printer interface board (4), component side down, on the module. Replace the printer interface board shield (8) and insert the four screws (7) retained during disassembly.
- 4. If the printer interface board has been repaired, or if a new printer interface board has been installed, the print quality may require adjustment. See 'Adjusting Print Quality' below.

# Removing and Replacing the Battery

You may encounter printers that were manufactured with replaceable batteries or printers manufactured with soldered batteries. Replaceable batteries must be properly stored and handled to avoid discharge. The time and date functions of the printer are battery supported when the unit is not running. If the time or date printout is incorrect, try resetting the time and date as described in the Operator's Manual. Power the printer and STATIM OFF and wait several minutes before powering them back ON. If the time and date are still incorrect, replace the battery on the printer interface board.

To replace the battery, proceed as follows:

- Remove the Printer Interface Board as described in Removing the Printer Interface Board.
- 2.
- a. If there is a replaceable battery, remove and discard the old battery. Install a new battery. Always replace the battery with a fresh battery of equal rating and size.
- b. If the battery is soldered, carefully de-solder BAT1 from the component side of the board. Note the orientation of the anode and cathode. Solder the replacement battery in position BAT1.

# **Adjusting Print Quality**

To adjust print quality, use the R21 contrast adjustment pot (9) located on the Printer Interface Board.


To adjust the print quality, proceed as follows (see Figure 2):


- 1. Remove the cover.
- 2. The printer contrast adjustment pot (R21) is located at the edge of the board. Adjustments may be made using a small adjustment tool while the printer module is still in the cover assembly.

- 3. To adjust and test print quality, connect the printer cable to the STATIM and the STATprinter and turn the host unit power switch **ON**.
- 4. Open the printer module and turn the printer power button ON.
- 5. Start and then quickly stop a cycle. Doing so causes an error message to be printed. While the error message is printing, adjust the pot (R21).
- 6. If further adjustment is required, repeat steps 3 to 5.
- 7. Reinstall the cover.


# Removing the Printer Assembly

To remove the printer assembly (1) from the printer module assembly (2), proceed as follows (see Figure 3):

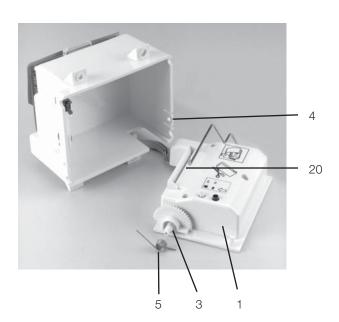

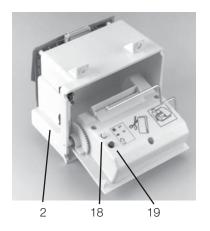



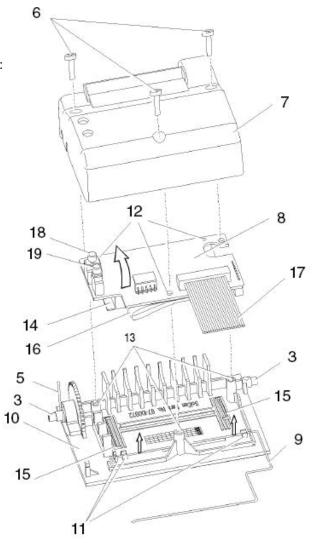

- 1. Cover
- 2. Printer module assembly
- 3. Hinge pin
- 4. Snap mechanism
- 5. Spring
- 6. Screws
- 7. Printer housing
- 8. Printer board
- 9. Paper roll arm
- 10. Printer door

- 11. Paper roll arm retaining clips
- 12. Printer board mounting holes
- 13. Mounting bosses
- 14. Printer body
- 15. Locating ribs
- 16. Print head flexible cable
- 17. Printer ribbon cable
- 18. Power button
- 19. Paper advance button
- 20. Paper feed slot

- 1. Remove the printer module. See 'Removing the Printer Module' above.
- 2. Unlatch the printer assembly from the printer module assembly. The printer assembly is held into the module by two hinge pins (3). The pin on the bottom left of the assembly sits in a recessed slot / retaining hole. The pin on the bottom right hand of the assembly is captured by a snap mechanism. Deflect the snap (4) away from the printer to free the hinge pin, and swing the assembly out of the module housing.
- 3. Remove the printer spring (5) from the left hand hinge pin and retain for re-assembly.
- 4. Using a #1 Phillips screwdriver, remove three screws (6) from the printer housing (7) and set them aside for use in re-assembling the printer.
- 5. Remove the housing. Note the orientation of the printer board (8) and the paper roll arm (9) assembled on the printer door (10).
- 6. Carefully lift the Printer Board upwards and away from the printer door.

Be careful when handling the board. The printer is integral to the wiring board. Do not place strain on the connections of the ribbon cable soldered to the board.

- 7. Remove the paper roll arm from the clips (11).
- 8. Carefully rest the printer board beside the assembly.

# Replacing the Printer

To replace the printer, proceed as follows (see Figure 3):

- 1. Carefully snap the paper roll arm (9), in the position shown, back into the clips (11) on the printer door.
- 2. Place the printer board (8) back into position on the printer door (10). Note the alignment of the printer board mounting holes (12) and the mounting bosses (13) on the printer door. The black plastic printer body (14) rests between the locating ribs (15) on the inside of the printer door.
- 3. Check that the print head flex cable (16) and printer ribbon cable (17) are not pinched between the printer door and the wiring board.
- 4. Place the printer housing (7) on the printer door (10). Check again to be sure that the flexible cables are not pinched between the cover and the door. The power button (18) and the paper advance button (19) must protrude through the openings in the cover and operate freely.
- 5. Secure the printer housing to the printer door with the three screws (6) retained during the disassembly procedure. **Do not overtighten these screws**.
- 6. Place the printer spring (5) on the left hand hinge pin of the printer assembly, with the long arm positioned to align with the long slot on the module housing.
- 7. Place the left hand hinge pin in the recesses slot / retaining hole and align the long spring arm. Swing the right hand hinge pin towards the module housing and push firmly onto the snap mechanism (4).
- 8. Replace the printer interface board. See 'Replacing the Printer Interface Board' above.
- 9. 9. Power the unit ON. Load the thermal paper and adjust the print quality.
- 10. See, 'Installing Thermal Paper into the Statprinter' and 'Adjusting Print Quality'.
- 11. 10. Reinstall the cover. See 'STATprinter Cover Removal and Replacement'.

# Installing Thermal Paper into the STATprinter


NOTE: Use only paper approved for use with the STATprinter. The use of any other paper may damage the printer and will void the warranty. Thermal paper is available from SciCan (SciCan order no. 01-101657S).

Do not operate the printer without paper. If you run out of thermal paper, or if you do not wish to use the printer, turn it OFF.

Never pull the paper backwards through the printer. This will damage the printer mechanism.

To install the paper into the printer, proceed as follows (see Figure 4.1):

- 1. Power the STATIM unit ON.
- 2. Open the printer door (1) by pushing on the top half of the door.
- 3. Power the printer ON.
- 4. Unroll a small amount of paper roll (4) and trim the corners using the paper cutting template included with each box.
- 5. Move the paper roll arm(4) into the loading position. Place the paper roll (3) on the arm so the paper strip feeds from the top of the roll and then carefully insert it into the paper feed slot(5) until it stops. If the paper does not feed from the top, the heat sensitive side of the paper will not be in contact with the print head and the printer will not print.
- 6. With one hand, continue to gently feed the paper strip into the paper feed slot while at the same time pressing the paper advance button with the other hand until the paper feeds by itself. Keep the paper straight when feeding it into the printer or it may jam. Do not force the paper into the slot! If the paper will not feed into the slot, pre-cut the end of the roll again and reload the paper.

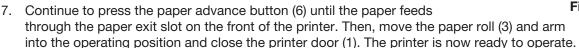

into the operating position and close the printer door (1). The printer is now ready to ope

When you see a red line on one side of the paper, it is time to replace the roll.

If a paper jam occurs, and the paper cannot be removed by pressing the paper advance button (6) do not pull the paper backwards through the printer.

Never put a utensil or tool into the paper exit slot.

For full instructions on how to remove paper jams, see 'Removing Paper Jams' below.

To replace the paper roll (3), proceed as follows (see Figure 4.2):

- 1. With scissors, cut the paper between the roll (3) and the paper feed slot (5)
- 2. Remove the roll from the arm and discard the unused portion.
- 3. Press the paper advance button (6) to feed the paper that remains in the printer out of the slot at the front of the printer.
- 4. Install the new thermal paper roll by following the instructions described above.

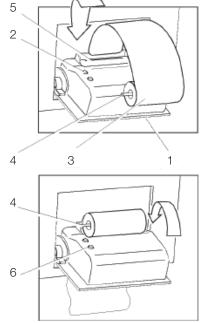


Figure 4.1

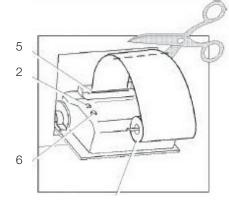


Figure 4.2

# **Removing Paper Jams**

NOTE: If paper jams in the printer and cannot be removed by pressing the paper advance button, the printer must be disassembled. Do not pull the paper backwards through the printer. Never put utensils or tools into the paper exit slot.

When paper is jammed in the printer, proceed as follows (see Figure 3):

- 1. Turn the host unit power OFF.
- 2. Make sure the printer power button (18) is in the **OFF** position and disconnect the printer cable (2) (Figure 1) from the STAT/M printer connector port and the STATprinter.
- 3. Using scissors, cut the paper between the roll and the paper feed slot (20).
- 4. Remove the three screws (6) from the printer housing (7) and set them aside for use in reassembling the printer. Remove the housing.
- 5. Note the orientation of the exposed printed wiring board and the paper roll arm assembled on the printer door (10). Remove the paper roll arm (9) from the clips (11).
- 6. Gently lift the printed wiring board (8) upwards and away from the printer door (10).

Be careful when handling the board. The printer is integral to the wiring board. Do not place strain on the connections of the ribbon cable soldered to the board. Do not remove the connector of the flexible cable (16) from the connector header on the board.

The paper drive mechanism on the underside of the wiring board is now exposed.

- 7. Using a pair of tweezers or fine needle-nosed pliers, **carefully** remove the paper from the mechanism.
- 8. When the paper is removed, proceed as described in Replacing the Printer.

# **External Serial Printers**

(for use on Rev. 7.x controller boards ONLY)

When fitting an external printer to a STATIM unit it is recommended that a 'till roll', 'point of sale' or 'receipt' style printer is used similar to the example shown:

Note that these printers are NOT sold by SciCan and should be purchased in the relevant territory by the user or the dealer for fitting. They are readily available (typically over the internet).

These printers do not normally come with a cable as standard and a Serial Null Modem Cable DB9F – DB25M will be required.



Star SP212FD42-230 'Point of Sale' Printer

# Recommended Printers for STATIMs with Revision 7 Controller Boards

The following printers have been tested by SciCan and are recommended for use as an external printer for STATIM units with Rev. 7 boards. They may be purchased from your local computer store or from an online supplier.

Suggested External Printers by SciCan	End Of Line CR/LF	Serial Port Bitrate	Printer user ° char
Epson TM-U220D (C31C515603)	CR/LF	9600	248 [0xF8]
Citizen IDP-3110-40 RF 120B	CR	9600	N/A
Star Micro SP212FD42-120	CR	9600	210 [0xd2]
Star Micro SP216FD41-120	CR/LF	9600	210 [0xd2]
Star Micro SP512MD42-R	CR/LF	9600	210 [0xd2]

### Important notes:

- The printers referenced above MAY be superseded from time to time. The updated equivalent MAY be able to be used with your STAT/M unit but it is important to ensure that the 'End of Line CR/LF', 'Serial Port Bitrate' and 'Printer uses of char' are appropriate to the new printer.
- Some of the reference numbers above may contain voltage references (e.g. SP212FD42-120). If the voltage in your region is different, then you should ensure that you obtain a printer with the correct voltage (e.g., for a 230V territory, the printer reference should be SP212FD42-230).

# Installing the external printer on your STATIM unit Setting the STATIM unit

Locate the printer port on the rear of the unit and proceed as follows:

- 1. Ensure the STATIM is powered OFF
- 2. Press and hold the STOP button and power ON the STATIM unit to access the USER menu.

The following menu is now available and the display on the LCD screen should be as highlighted.

- 3. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "RS232" and select it by pressing the RUBBER & PLASTICS button.
- 4. You should now have the option of "N/A", "USB/Flash MSD" or "Serial Printer" which can be scrolled through by using the UNWRAPPED and WRAPPED buttons.

# User Mode ➤ Time/Date Setup Language Setup Unit ID Setup Last Printout RS232 End Of Line CR/LF Serial Port Bitrate Printer user ° char Save and Exit Exit

- 5. Scroll to "Serial Printer" and press the RUBBER & PLASTICS button to select and return to the user menu.
- 6. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "End of Line CR/LF" and press the RUBBER & PLASTICS button to access the options.
- 7. Using the UNWRAPPED and WRAPPED buttons, scroll until reaching the desired CR/LF value in the table above (or from the printer's user manual) and press the RUBBER & PLASTICS button to select and return to the user menu.
- 8. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "Serial Port Bitrate" and press the RUBBER & PLASTICS button to access the options.
- Using the UNWRAPPED and WRAPPED buttons, scroll until reaching the desired Bitrate value in the table above (or from the printer's user manual) and press the RUBBER & PLASTICS button to select and return to the user menu.
- 10. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "Printer uses o char" and press the RUBBER & PLASTICS button to access the options.
- 11. Using the UNWRAPPED and WRAPPED buttons, scroll until reaching the desired Hex value in the table above (or from the printer's user manual) and press the RUBBER & PLASTICS button to select and return to the user menu. Important note: the UNWRAPPED button will increase the displayed value by 1, and the WRAPPED button will increase that displayed value by ten.

12. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "Save and Exit" and press the RUBBER & PLASTICS button, which will now save the data, and take you out of the user menu and back to the operating menu.

Note: When printed or electronic records are being kept for future reference, it is important to ensure the correct data is printed, and a check should be undertaken at this stage to ensure that the date, time and unit number (where appropriate) are set before using a data recording device. This may have been undertaken during installation, if not, refer to the user manual or to this service manual's Chapter 13, "Using the service menu" for information on resetting the time, date and unit number.

13. Power OFF the STATIM unit.

# Connecting the Printer

- 1. Connect the printer to the STAT/M unit using a Serial Null Modem Cable DB9F DB25M and secure the cable.
- 2. Connect the printer power supply to the mains power socket as appropriate to the printer being used.
- 3. Power ON the printer.
- 4. Power ON the STATIM unit.
- 5. The unit is now ready for use.

# **Data Logger**

SciCan's Data Logger can record and store cycle information onto a mass storage device (MSD) such as a USB Flash Drive or SD memory card.

The Data Logger is supplied with two cables for connection to

- 25 pin communication ports (STAT/M 2000 units with revision 2.x/5.x/6.x controller boards only) and
- 9 pin communication ports (all 2000 and 5000 units with revision 7.x controller boards)

### Notes:

- STAT/M 2000 units can only be fitted with a data logger if a communication port is present. Early STAT/M 2000 units (with revision 3.x/4.x controller boards) that are NOT fitted with a port can not have a port retrofitted.
- STAT/M 2000 units that have been upgraded to a revision 7 controller board will require the 25 pin port upgrading to a 9 pin port to ensure correct communication.
   A kit is available to enable this which is RS232 Port Kit STAT/M 2000/2000S part number 01-110221S.

#### Installing the SciCan Data Logger

Locate the printer port (9 or 25 pin) on the rear of the unit and proceed as follows:

# Installation on 9 pin revision 7 STATIM units

(2000 or 5000). (Note: these units will have a user and service menu)


For this operation the Null Modem Serial Cable DB9F - DB25M is required

#### Setting the STATIM unit

- 1. Ensure the STATIM is powered OFF
- Press and hold the STOP button and power ON the STATIM unit to access the USER menu.

The following menu is now available and the display on the LCD screen should be as highlighted.

 Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "RS232" and select it by pressing the RUBBER & PLASTICS button.

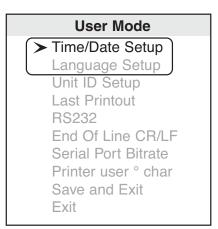

- 4. You should now have the option of "N/A", "USB/Flash MSD" or "Serial Printer" which can be scrolled through by using the UNWRAPPED and WRAPPED buttons.
- 5. Scroll to "USB/Flash MSD" and press the RUBBER & PLASTICS button to select and return to the user menu.
- 6. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "Serial Port Bitrate" and press the RUBBER & PLASTICS button to access the options.
- 7. Using the UNWRAPPED and WRAPPED buttons, scroll (if required) to the Bitrate value of 9600 appears in the display and then press the RUBBER & PLASTICS button to select and return to the user menu.
- 8. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "Printer uses of char" and press the RUBBER & PLASTICS button to access the options.
- 9. Using the UNWRAPPED and WRAPPED buttons, scroll (if required) to the Hex value of 32 [0x20] appears in the display and then press the RUBBER & PLASTICS button to select and return to the user menu. Important note: the UNWRAPPED button will increase the displayed value by 1, and the WRAPPED button will increase that displayed value by ten.

10. Using the UNWRAPPED and WRAPPED buttons, scroll until the cursor is next to "Save and Exit" and press the RUBBER & PLASTICS button, which will now take you out of the user menu and back to the operating menu and will save your changes.

The LCD display should now show the following:

HH:MM DD/MM/YYYY

"MSD NOT DETECTED"/"INSERT MSD/FLASH"/"SELECT A CYCLE"

Messages "MSD NOT DETECTED", "INSERT MSD/FLASH" and "SELECT A CYCLE" will rotate in sequence and not all appear at the same time.

Note: When printed or electronic records are being kept for future reference, it is important to ensure the correct data is printed, and a check should be undertaken at this stage to ensure that the date, time and unit number (where appropriate) are set before using a data recording device. This may have been undertaken during installation, if not, refer to the user manual or to this service manual's Chapter 13, "Using the service menu" for information on resetting the time, date and unit number.

11. Power OFF the STATIM unit.

#### Connecting the Data Logger.

- 1. Connect the Data Logger to the power cable attached to the transformer supplied. Important note: ensure that the plug is fully inserted in the socket on the data logger.
- 2. Attach the adaptor suitable for your regional power supply on to the main body of the transformer. Plug in to the mains supply BUT DO NOT SWITCH ON AT THIS TIME.
- 3. Ensure that the STATIM unit is powered OFF.
- 4. Using the Null Modem Serial Cable DB9F DB25M supplied, connect the Data Logger to the STAT/M unit and secure the cable.
- 5. Power the Data Logger ON.
- Power the STATIM unit ON.

The LCD display should now show the following:

HH:MM DD/MM/YYYY

"MSD NOT DETECTED"/"INSERT MSD/FLASH"/"SELECT A CYCLE"

Messages "MSD NOT DETECTED", "INSERT MSD/FLASH" and "SELECT A CYCLE" will rotate in sequence and not all appear at the same time.

- 7. Insert USB flash drive or SD memory card.
- 8. After a few seconds the LCD display should now show the following:

#### HH:MM DD/MM/YYYY

# USB/FLASH DETECTED/SAFELY REMOVE MSD/SELECT A CYCLE

Messages "USB/FLASH DETECTED", "SAFELY REMOVE MSD" and "SELECT A CYCLE" will rotate in sequence and not all appear at the same time.


9. The unit is now ready for use.


# Installation on 25 pin units (2000 ONLY).

(Note: these units will NOT have a user and service menu)

For this operation the 25 to 25 pin special cable supplied with the data logger is required.

Note that this cable is directional and can only be fitted one way. The cable sockets are marked "STAT/M" and "data logger" to assist in identifying the direction of connection. Failure to connect the cable in the correct direction will cause the LCD screen to go blank, and may cause damage to the unit.

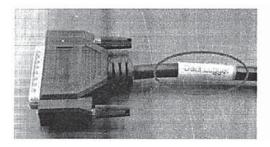

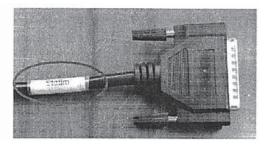

This type of STAT*IM* 2000 unit (no menu facility) will automatically recognize the data logger when connected. To install the data logger, proceed as follows:

- 1. Connect the Data Logger to the power cable attached to the transformer supplied. Important note: ensure that the plug is fully inserted in the socket on the data logger.
- 2. Attach the adaptor suitable for your regional power supply on to the main body of the transformer. Plug in to the mains supply BUT DO NOT SWITCH ON AT THIS TIME.
- 3. Ensure that the STATIM unit is powered OFF.
- 4. Using the 25 pin to 25 pin cable supplied, connect the Data Logger to the STATIM unit and secure the cable. Note the orientation of the cable.

- 5. Insert USB flash drive or SD memory card.
- 6. Power the Data Logger ON.
- 7. Power the STATIM unit ON.
- 8. The display on the unit will now show the time and date (not shown when logger not connected) and the correct time and date will need to be set. To set time, date and unit number, see Chapter 13, Installation.
- 9. When this operation is complete, the unit is now ready for use.

# **Data Logger Troubleshooting**

9 pin revision 7.x units only.

If the display reads: "MSD NOT CONNECTED <> INSERT MSD/FLASH", proceed as follows to ensure everything is setup correctly.

- 1. Check the serial cable connection.
- 2. Check the power connection.
- Ensure the lower red LED is lit. It is located between the Data Logger's serial port and power input.
- 4. Check that the Mass Storage Device (MSD) is properly inserted in the appropriate slot of the Data Logger.
- 5. Repeat the instructions for Installing the Data Logger on your STATIM.

If the display reads: "SAFELY REMOVE MSD <> MSD/FLASH DETECTED", the MSD/Flash can be safely removed without affecting the data.

If the display reads: "MSD/FLASH FULL <> REPLACE MSD", the MSD is full and the data should be downloaded from the MSD.

**If the MSD:** has missing lines of Data, refer to the STAT*IM* screen to confirm successful sterilization.

The Data Logger may be reset by unplugging its power supply, disconnecting the MSD and waiting 10 seconds. Then re-connect the power adapter and insert the MSD into the Data Logger.

If the problem persists, contact the SciCan Service Center.

**If MSD has:** Corrupt or unreadable files and/or directories, refer to the STAT/M screen to confirm successful sterilization.

The MSD may have been unplugged while data was being written to it. The MSD should NOT be unplugged until after "SAFELY REMOVE MSD <> MSD/FLASH DETECTED" is displayed.

The corrupted files or directories may be lost. Reformat the MSD on your computer.

#### 25 pin revision 2.x/5.x/6.x 2000 units only.

Unlike the revision 7.x unit, these units will not have a user or service menu and the LCD will NOT display any messages that can be used for diagnostic purposes. The only way to tell that the logger is not working correctly is either:

- No data will be recorded on the MSD.
- The time and date will not be displayed on the screen of the STATIM unit.

If either of these faults occurs, proceed as follows to ensure everything is setup correctly:

- 1. Check the cable connection, particularly the orientation.
- 2. Check the power connection.
- 3. Ensure the lower red LED is lit. It is located between the Data Logger's serial port and power input.
- 4. Check that the Mass Storage Device (MSD) is properly inserted in the appropriate slot of the Data Logger.
- 5. Repeat the instructions for Installing the Data Logger on your STATIM.

The Data Logger may be reset by unplugging its power supply, disconnecting the MSD and waiting 10 seconds. Then re-connect the power adapter and insert the MSD into the Data Logger.

If the problem persists, contact the SciCan Service Center.

# Transferring data from the USB drive or SD card

9 pin revision 7.x units only

Files are stored on the MSD in txt and PDF format. These can be transferred for storage on a computer hard drive as follows:

- Ensure the LCD screen of the STATIM reads "SAFELY REMOVE MSD <> MSD/FLASH DETECTED" and the STATIM is in standby. THE UNIT MUST NOT BE RUN WHEN REMOVING THE MSD.
- 2. Remove MSD from the Data Logger. The STATIM LCD screen should now show the following:

#### HH:MM DD/MM/YYYY

"MSD NOT DETECTED"/"INSERT MSD/FLASH"/"SELECT A CYCLE"

- 3. Insert the MSD in the relevant port on the computer and wait until the computer has recognised the new device (if transferring for the first time).
- 4. The contents of the MSD should now be available for viewing in the appropriate application (e.g. Windows Explorer).
- 5. Locate, or create the preferred folder for storing the data.
- 6. Data can now be transferred by either:
  - a. Drag and Drop
  - b. Cut and Paste
  - c. Copy and Paste (if you wish to retain the data on the MSD)
- 7. Safely 'eject' the MSD as appropriate to your computer and remove the MSD.
- 8. Re-insert the MSD back in to the Data Logger.
- 9. After a few seconds the STATIM LCD display should now show the following:

#### HH:MM DD/MM/YYYY

USB/FLASH DETECTED/SAFELY REMOVE MSD/SELECT A CYCLE

10. The STATIM and Data Logger are now ready for use.

#### 25 pin revision 2.x/5.x/6.x 2000 units only.

Files are stored on the MSD in txt and PDF format. These can be transferred for storage on a computer hard drive as follows:

- 1. Ensure the STATIM is in standby. THE UNIT MUST NOT BE RUN WHEN REMOVING THE MSD.
- 2. Remove MSD from the logger.
- 3. Insert the MSD in the relevant port on the computer and wait until the computer has recognised the new device (if transferring for the first time).
- 4. The contents of the MSD should now be available for viewing in the appropriate application (e.g. Windows Explorer).
- 5. Locate, or create the preferred folder for storing the data.
- 6. Data can now be transferred by either:
  - a. Drag and Drop
  - b. Cut and Paste
  - c. Copy and Paste (if you wish to retain the data on the MSD)
- 7. Safely 'eject' the MSD as appropriate to your computer and remove the MSD.
- 8. Ensure the STATIM is in standby and re-insert the MSD back in to the Data Logger. THE UNIT MUST NOT BE RUN WHEN REMOVING THE MSD. Tip: Allow a few seconds for the two units to re-establish communication before running the STATIM.
- 9. The STATIM and Data Logger are now ready for use.

# **Data Logger Specifications**

Connection: Serial Cable DB9F-DB25M

Baud Rate: 9600

Electrical Rating: 100V-240V, 0.6A, 50/60 Hz, multiplug

#### ACKNOWLEDGED MSD

#### **USB Flash Drives:**

Cruzer mini 1 GB Memorex Travel drive 1 GB Kingston Travel drive 512 MB PNY Attache 1 GB Lexar Jump Drive 1 GB

#### SD Memory Cards:

Kingston 512 MB Lexar 512 MB PNY 512 MB Ultra 512 MB

Kingston ultimate 512 MB

Other USB Flash Drives and SD Cards may also be compatible with this Data Logger.

# WARNINGS AND PRECAUTIONS

If you have questions about the unit you are repairing, please do not hesitate to contact your local SciCan representative for information. Also, the STAT*IM* is heavy. Exercise caution and seek assistance when lifting or carrying units.



# **EXERCISE CAUTION**

- Hazardous voltages are accessible when the cover is removed.
- Disconnect the power cord before servicing the power mains portion of the controller board and associated devices.

# PERFORM TESTS

 If the cover is removed, a dielectric strength test (Hi-Pot) <u>AND</u> a protective bonding impedance test (ground continuity) must be performed on the STAT*IM* when the work is completed and after the cover has been returned to the unit.



# PROTECT THE UNIT

- Use only steam-process distilled water in the STATIM.
- The STATIM contains electronic circuitry that is static sensitive. Always wear a static strap when working with or near printed wiring boards. In addition, use static footstraps, grounding mats and grounded work surfaces when servicing microprocessor devices. Transport boards and devices in static protected bags.
- Ensure there is sufficient steam-process distilled water in the STAT*IM* before activating the pump.
- In order to ensure adherence to the applicable safety agency approvals, state, provincial, regional and national laws, replace components with SciCan approved parts only.

# Installation

#### **Environmental Considerations**

There are several factors that may affect the performance of your STAT*IM*. Please review these factors, and select a suitable location in which to install the unit.

#### Temperature and Humidity

Avoid installing your STAT/M in direct sunlight or close to a heat source (e.g. vents or radiators). The recommended operating temperatures are 15-25°C (59°F to 77°F) with humidity of 25-70%.

#### Spacing

The vents and openings on the STATIM should remain uncovered and unobstructed. Leave a minimum of 50 mm/2" between the top, sides and back of the unit and any wall or partition.

#### Venting

The STATIM should be operated in a clean, dust free environment.

#### **Work Surface**

The STATIM should be placed on a flat, level, water-resistant surface. Never install and operate the unit on a sloped surface.

#### **Electromagnetic Environment**

The STATIM has been tested and meets applicable standards for electromagnetic emissions. While the unit does not emit any radiation, it may itself be affected by other equipment which does. We recommend that the unit be kept away from potential sources of interference.

#### **Electrical Requirements**

Use properly grounded and fused power sources with the same voltage rating as indicated on the label at the back of your STATIM. Avoid multiple outlet receptacles. If using a surge suppressor power bar, plug in one STATIM only.

#### **Unit Placement**

When placing the unit on a counter top, ensure the following:

The level indicator bubble (1) on the front panel should be balanced in the front right quadrant of the target. This will ensure that the unit drains properly. Adjusting the three leveler feet will help you move the bubble if necessary.

The unit should be stable and all four feet should be securely in contact with the counter surface. This will prevent the unit from moving freely.

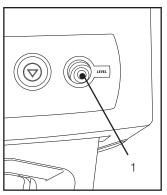



Figure 1

#### Connecting the Waste Bottle

The waste bottle (2) is used to collect the wastewater after it has been converted to steam and then drained from the cassette. To connect

the waste bottle to the STATIM, follow these steps (see Figure 2):

- 1. Insert the exhaust (3) tube into the fitting (4) on the back of the unit and connect tightly.
- 2. Cut the tube to length and slide the waste bottle fitting (5) into place.
- Place the free end of the tube into the hole in the lid of the waste bottle and hand-tighten the fitting. Do not coil the exhaust tube.
- Unscrew the lid and copper condenser coil (6) assembly from the waste bottle. The lid and coil should come out together.
- 5. Fill the waste bottle with water to the MIN line and replace the lid and copper condenser assembly. Empty the waste bottle often to avoid unpleasant odors and discoloration of the contents. (A low-level disinfectant, prepared according to the manufacturer's instructions, may be added to the waste bottle to remedy this situation). As a minimum, empty the waste bottle each time you refill the reservoir.
- 6. Place the waste bottle near the unit. Store the bottle below the unit. The tube can be routed through a hole, (8 mm/0.3" in diameter) in the counter-top and secured with the provided nylon clamps.

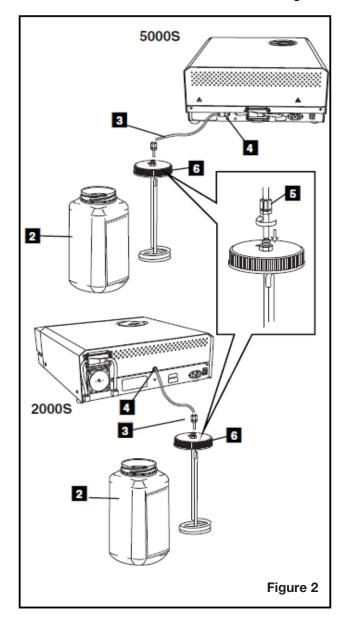

#### Filling the STATIM Reservoir

When filling the reservoir, ensure you only use steam processed distilled water containing less than 5 ppm total dissolved solids (having conductivity of less than 10  $\mu$ S / cm). The impurities and additives in other water sources will cause an error reading on the LCD. If you have a water conductivity meter (available from SciCan, order number 01-103139S) check each new water container before filling the reservoir. If your unit is equipped with a Rev. 7.x controller board or higher (software SxxxR6xx), you can read the water conductivity by referring to the 'Water Quality' sub menu in the main user menu.

To fill the reservoir, follow these steps (see Figure 3):

- 1. Remove the reservoir cap (2)
- 2. Pour steam-process distilled water into the reservoir until almost full (a maximum of 4L/1 U.S. gal). Use a funnel to avoid spillage.
- 3. Replace and secure the cap.

#### Priming the STATIM Pump

To prime the STAT/M pump, follow these steps (see Figure 4):

Move the unit to the edge of the work surface.
 The front leveler feet should be approximately 12 mm/0.5" from the edge.

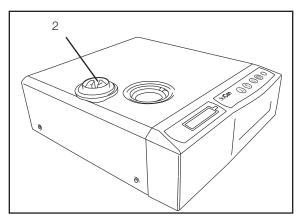


Figure 3

- 2. Lift the front left corner of the unit upward and remove the drain tube (3) from the clip located on the underside of the unit.
- 3. Pull the drain tube outward so the free end can be positioned over a water container.
- 4. Fill the reservoir with steam-process distilled water.
- 5. Remove the plug (4) from the end of the drain tube and allow water to drain from the tube into a container for 30 seconds. When the water flows in a steady stream, replace the stopper.

  4
- 6. Lift the front left corner of the unit upward and reinsert the tube into the clip on the underside of the unit. Push the excess length of tubing back into the space provided.

Make sure the plug in the drain tube is secured.

After installation, and before sterilizing any instruments, run two cycles.

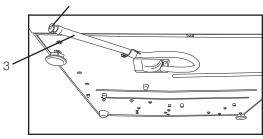


Figure 4

#### Powering up and setting up the unit

- 1. Attach the main power cable to the unit and to the outlet socket.
- 2. Switch on the outlet socket (if appropriate) and then the unit via the power switch on the rear of the unit.
- 3. Check that the display illuminates.

#### Setting the Time and Date

For STATIM G4:

See Using the LCD Service Menu on the STATIM G4 in Chapter 3 to access Time and Date.

 For STAT/M Classic, follow these steps and watch the blinking cursor on the LCD:Power the STAT/M OFF.

- 2. Press and hold the UNWRAPPED cycle button.
- 3. While pressing the UNWRAPPED cycle button, power the STATIM **ON**. The LCD display reads:
- 4. Use the cycle buttons to select and change the selected field's value. To increase a field's value, press the UNWRAPPED cycle button. Hold the button down to increase the value.

18:00 20/09/2005 HH:MM DD/MM/YYYY

Display when setting the Date / Time

- 5. To select the next field, press the START button.
- 6. To save changes and return to the regular operating mode, press the **STOP** button.
- 7. To quit without making changes, power the STAT/M **OFF**.

#### Language Selection

For STATIM G4: See Using the LCD Service Menu on the STATIM G4 in Chapter 3 to access Language Selection.

For STATIM Classic, the messages displayed on the LCD can be changed to a number of different languages.

To change the current language selections, follow these steps:

- N.A. ENGLISH
- Display when scrolling through languages
- 1. Turn the power switch at the back of the unit **OFF**.
- 2. Press and hold the WRAPPED cycle button.
- 3. While depressing the WRAPPED cycle button, turn the power switch at the back of the unit to **ON**.
- 4. Press the UNWRAPPED cycle button to scroll to the next language selection.
- 5. Press the WRAPPED cycle button to scroll to the previous language.
- 6. When the desired language is displayed, press the **STOP** button to save the selection and return to the regular operating mode..

#### **Assigning Unit Identifier Number**

For STATIM G4: See Using the LCD Service Menu on the STATIM G4 in Chapter 3 to access Assigning Unit number

For STATIM Classic, follow these steps:

- 1. Power the STATIM OFF.
- 2. Press and hold the RUBBER / PLASTIC cycle button.
- 3. While pressing the RUBBER / PLASTIC cycle button, power the STATIM ON.
- 4. Using the cycle buttons, select a maximum of 3 digits to be used as the unit's identifier. The UNWRAPPED button will increase the selected value. Use the RUBBER / PLASTIC to move to the next digit.
- 5. To save changes and return to the regular operating mode, press the **STOP** button.

STATIM 2000E S2S2R601 #323

Display when assigning unit number

#### **Setting/Editing Serial Number**

**NOTE:** This function is available on Rev. 7 controller boards with software revision SXXXR612 on STATIM Classic, and software revision SXXXR700 on STATIM G4.

For STATIM Classic, follow these steps::

- 1. Power **ON** the STATIM while holding down the UNWRAPPED and WRAPPED keys.
- 2. Enter the service password to access the service menu.
- 3. The default password is entered by pressing the UNWRAPPED, WRAPPED, RUBBER & PLASTICS and **STOP** keys in that order (a backdoor password can also be used by pressing the UNWRAPPED, WRAPPED, UNWRAPPED, WRAPPED keys in this order).
- 4. Enter the Service Menu and select EDIT SERIAL NUMBER from the list.
- 5. Use the UNWRAPPED button to increase the selected digit and the WRAPPED button to decrease the selected digit.
- 6. Use the RUBBER & PLASTIC button to select the next adjacent digit.
- 7. To save changes and return to the regular operating mode, press the STOP button.

For STATIM G4, follow these steps:

- Access the Service Menu as described in Using the LCD Service Menu on the STATIM G4 in Chapter 3 and scroll to SET SERIAL NUMBER.
- 2. Use the first pair of the arrow buttons to move from one digit to the next. Use the second pair of arrow buttons to change digit values.
- 3. To save changes and return to the regular operating mode, press the **STOP** button

**NOTE:** The serial number can be set also by using a Barcode Scanner connected onto the RS232 port (Settings: 96008N1, CR/LF at the end of string).



#### **Setting Air Drying Time**

**NOTE:** This function is available on Rev. 7 controller boards with software revision SXXXR612 on STAT*IM* Classic, and software revision SXXXR700 on STAT*IM* G4. However, this function is not available on U.S. STAT*IM* models.

For STATIM Classic, follow these steps:

- 1. Power up unit while pressing the **STOP** button to access the User Setup menu.
- 2. Using the UNWRAPPED and WRAPPED buttons, scroll to AIR DRYING and select it by pressing the RUBBER AND PLASTIC button.
- 3. Use the UNWRAPPED and WRAPPED buttons to change the displayed values and the RUBBER AND PLASTIC button to select and move to the next digit. Press the **STOP** button when completed.

For STATIM G4, follow these steps:

- 1. Access the Service Menu as described in Using the LCD Service Menu on the STATIM G4 in Chapter 3 and scroll to DRYING TIME.
- 2. Use the first pair of the arrow buttons to move from one digit to the next. Use the second pair of arrow buttons to change digit values.
- 3. To save changes and return to the regular operating mode, press the **STOP** button.

#### Connecting STATIM 2000 G4 to a Network

The STATIM G4 has a 10/100Base-T Ethernet port located at the back of the unit.

To connect it to a network using a router, follow these steps:


1. Connect the network cable to the Ethernet port at the back of the unit. If the customer's office uses a router, the router should automatically assign the unit an IP address. You will know the IP address has been assigned when the red X over the Network icon disappears.






NOTE: In some circumstances, where there is no router, for example when using Windows Network Sharing, you may have to assign a dedicated or 'static' IP address. To assign a static IP address, contact your local network administrator.

2. From the main screen, press the NETWORKING icon . The Networking screen displays information about the STATIM'S connectivity, including its IP address.

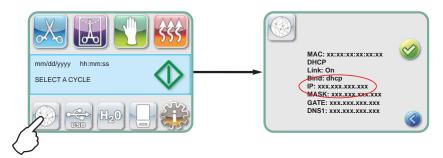

3. Type the IP address displayed on the touchscreen into the browser of any computer connected to the same network as the STAT*IM* to access your unit's web portal. You will know that an IP address has been assigned when the red X over the Network icon disappears.

**NOTE**: Connection time will vary depending on the network speed, and making an initial connection may take longer.

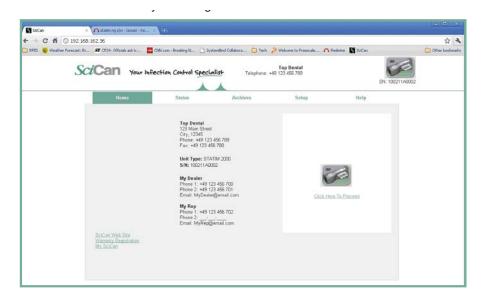

#### Connecting to a wireless network

The STAT*IM* G4 can be configured for wireless use by connecting the Ethernet port to an external wireless bridge / access point. SciCan currently recommends the use of the D-Link® DAP-1522 Xtreme N® Duo Wireless Bridge. Contact the customer's network administrator to learn more about setting up a wireless bridge.

#### Setting up the web portal

This web portal is a direct connection to the STAT*IM* on the customer's network. It is protected by the customer's firewall and not accessible to outside users. To set up the web portal follow these steps:

#### Setting up the CONTACTS page

- 1. Click on the SETUP page and go to the CONTACTS tab.
- 2. Here you will be prompted to provide a user name and password. For both the user name and password, the default is "scican". Use it to enter and thenyou can customize it using the PASSWORD tab.
- 3. The information you use to fill the fields on the CONTACTS page, such as office name, address and uploaded logo, is then used to customize the customer's web portal home page and print outs.



#### Setting up the LAN page

- 1. From the SETUP page, use the LAN tab to check your network settings.
- 2. If the office is using a router, keep the Automatic IP box checked. This will populate the fields for correct network settings automatically.
- 3. If you want to assign different values to these fields, uncheck the box, populate the fields and click Save to preserve the changes.



#### Setting up email


- From the SETUP page, use the REMOTE tab to assign the unit an email address.
- 2. Your STATIM must use an email account set up either with your Internet service provider or with a web-based email service. Populate the email account fields on this page with your provider's information if you are working with an email address from your Internet service provider. To use one of the STATIM web portal's preset email servers (choose from gmail, gmx mail, Yahoo and MSN mail), you must create an account for your STATIM (e.g. drsmithstatim@gmail.com).
- 3. A test email will be sent to verify your settings. Uncheck the box if the customer does not want to receive downloads and promotions from SciCan.



#### Setting up email NOTIFICATIONS

- 1. From the SETUP page, use the NOTIFICATIONS tab to enter the email addresses to which will be sent operating information and/or cycle faults.
- 2. A configurable message box lets you add details to identify the unit, such as its location.
- 3. Uncheck the box if the customer does not do not want SciCan Technical Service to automatically receive the unit's cycle fault emails.

**NOTE**: If you are not receiving emails from the unit, check the spam filter to be certain the STAT*IM* has been identified as an accepted email source.

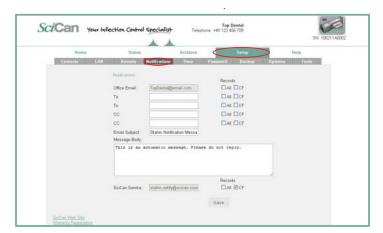

#### Setting the TIME zone and Internet synchronization


1. From the SETUP page, use the TIME tab to set the time and according to the customer's time zone and check the box for the unit to synchronize automatically with dedicated Internet time servers. It will then automatically reset its clock whenever it connects to the network.

NOTE: It will override any manual changes made to time settings directly on the unit.



#### Setting the PASSWORD

- Set a new user name and password to keep unauthorized users from changing the configurations made on ConnectIM as well as the Contacts information the STATIM uses to generate forms.
- 2. From the SETUP page, go to the PASSWORD tab. Note that both the default user name and password is 'scican'.

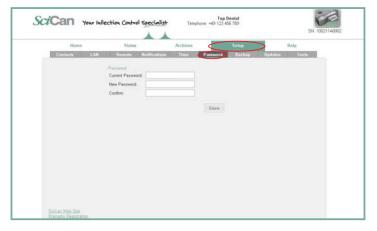


#### Setting up a BACKUP

- 1. From the SETUP page, use the BACKUP tab to back up all the configurable changes and Contacts information made on the STAT*IM*'s web portal by storing a file on the office computer.
  - NOTE: For security reasons, your email password will not be saved.
- 2. Click download settings to automatically store all settings in a file on the computer (TIP: Use the Save function on the computer to choose a specific location, otherwise it will load the file in a default location preset by the browser).
- 3. To restore settings, click Choose file and browse for the file on the computer, select it and click Restore.



#### **Using TESTS**

- 1. From the SETUP page, click on the TOOLS tab.
- 2. From here you can run a simple pass or fail test to check the network connectivity. The results will confirm whether you are accessing the Internet, accessing a router and whether the unit has a proper IP address.

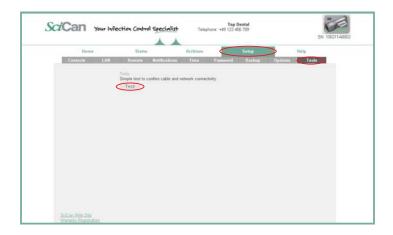



#### **Preparing Unit for Use**

Once the unit is installed and before any instruments are sterilized, run two cycles. Remove the cassette once it has cooled. Clean the top (lid) and bottom (tray) sections using a soft cloth to wipe the inside surfaces and then rinse thoroughly with tap water. Once the cassette is clean and dry, coat the inside surfaces with STAT-DRI. Please note that STAT-DRI drying agent is not be to used with the U.S. models of the G4 models (G4-121101 and G4-201103).

#### Shipping the Unit

Before you move the unit, you will need to drain the reservoir. To do so, follow these steps:

- 1. Place a water container below the unit.
- 2. Using the drain tube (see Priming the STATIM Pump, Figure 4) empty the contents of the reservoir into the water container.
- 3. Remove any remaining water from the reservoir with a non-linting, absorbant towel.
- 4. Screw-in the three leveler feet found underneath the unit.
- 5. Repack the unit in the original packing materials and include all accessories originally included with the unit.
- 6. Specify heated and insured shipping.

#### Unit installation checklist

Listed below you will find a generic check list for use when installing a STATIM unit to ensure installation is complete.

1. INVENTORY CHECK	2. UNIT INSTALLATION
Open box and check contents:	Set up unit as follows:
STAT/M unit	Locate unit and level
Cassette tray and lid	Connect Waste bottle
Unwrapped instrument rack (inside cassette)	Fill Reservoir with steam distilled water
Waste bottle	Prime pump
Waste bottle lid and condenser fitting	Connect mains lead and switch on machine
Tube mounting hardware (clips and screws)	Install printer if appropriate
Power cord	Set time and date
Operators manual	Check and select language as appropriate
Certificate of testing	Assign unit number if required
Exhaust tube	Assemble and insert cassette.
STAT-DRI (not available for U.S. G4 models)	Run unit
Drying plates (5000S only)	
Process Challenge Device (Helix) with test strips (S units only)	

3. TESTING THE UNIT
Power up unit and run as follows:
Check cycle selection on all cycles.
Assemble test helix and integrator as per helix instructions. (S units only).
Remove cassette and place assembled helix on instrument rack.
Insert cassette and run a Hollow Wrapped Cycle (S units)/ Wrapped cycle (non S units).
Check for steam leaks
At end of cycle, (no drying) remove cassette and remove helix.
Check integrator for pass/fail (✓ = pass 🗱 = fail)* (S units only)
Run 2nd Hollow Wrapped/Wrapped cycle

If the integrator shows a fail, do not use the unit. Replace unit and call SciCan representative.

# 14. STAT/// Operational Diagram

